CE-Immundiagnostika GmbH
Karl-Landsteiner-Str. 6
69151 Neckargemünd
Germany

Notified Body Confirmation Letter

Registration no.: D1415300034

To whom it may concern,

Confirmation of the status of a formal application, written agreement, and appropriate surveillance in the framework of Regulation (EU) 2024/1860 amending Regulations (EU) 2017/745 and (EU) 2017/746 as regards the transitional provisions for certain in vitro diagnostic medical devices

This letter confirms that mdc medical device certification GmbH (Kriegerstr. 6, 70191 Stuttgart, Germany), a Notified Body (NB) designated against Regulation (EU) 2017/746 (IVDR) and identified by the number 0483 on NANDO, has received a formal application in accordance with Section 4.3, first subparagraph of Annex VII of IVDR and has signed a written agreement in accordance with Section 4.3, second subparagraph of Annex VII of IVDR with the following manufacturer:

**CE-Immundiagnostika GmbH
Karl-Landsteiner-Str. 6
69151 Neckargemünd
Germany
SRN: DE-MF-000030423**


The devices covered by the formal application and the written agreement mentioned above are identified in the Tables below. Table 1 identifies the devices for which an IVDR application has been received, written agreement concluded and for which the NB is also responsible for appropriate surveillance of the corresponding devices under Directive 98/79/EC. Table 2 identifies the devices for which an IVDR application has been received and a written agreement concluded, but the NB has not yet taken the responsibility for appropriate surveillance of the corresponding devices under Directive 98/79/EC.

In the case of devices covered by certificates issued under Directive 98/79/EC (IVDD) that expired after 26 May 2022 and before 09 July 2024, without having been withdrawn, this letter also confirms that the manufacturer signed the written agreement under IVDR by the date of IVDD certificate expiry; or provided evidence that a competent authority of a Member State had granted a derogation or exemption from the applicable conformity assessment procedure in accordance with Article 54(1) of IVDR or Article 92(1) of the IVDR respectively, by 09 July 2024 for the relevant devices.

The transition timelines that apply to the devices covered by this letter, subject to the manufacturer's continued compliance to the other conditions specified in Article 110.3b of IVDR (as amended by Regulation (EU) 2024/1860), are shown below:

- 31 December 2027 for devices covered by an IVDD certificate regardless of their risk class under the IVDR
- For devices not requiring the involvement of a notified body under the IVDD, but requiring it under the IVDR and for which a declaration of conformity was drawn up prior to 26 May 2022 in accordance with Directive 98/79/EC, the following dates apply:
 - 31 December 2027, for class D devices;
 - 31 December 2028, for class C devices;
 - 31 December 2029, for class B devices and for class A devices placed on the market in sterile condition

Stuttgart, 2025-09-18

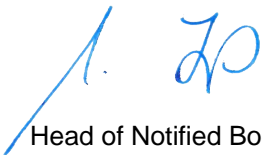Head of Notified Body

Table 1: Devices covered by this letter and for which the NB is also responsible for appropriate surveillance of the corresponding devices under Directive 98/79/EC:

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-A A-11H5 4049863Anti-A_0117J	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300024; NB # 0483 related Supplement D1415300031
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-A Birma-1 4049863Anti-A_0127L	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300024; NB # 0483 related Supplement D1415300031
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-B B-6F9 4049863Anti-B_0217Y	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300024; NB # 0483 related Supplement D1415300031
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-B LB-2 4049863Anti-B_02282	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300024; NB # 0483 related Supplement D1415300031
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-AB A-5E10; B-2D7 4049863Anti-AB_031ZF	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300024; NB # 0483 related Supplement D1415300031
Anti-C MS-24 4049863Anti-C_0418H	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032
Anti-C MS-273 4049863Anti-C_0428K	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032
Anti-c MS-33 4049863Anti-c_051HV	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032
Anti-c MS-35 4049863Anti-c_052HX	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-D MS-201 4049863Anti-D_06192	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032
Anti-D RUM-1 4049863Anti-D_06294	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300032
Anti-D incomplete 4049863Anti-D_06396	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029
Anti-D blend TH-28; MS-26 4049863Anti-D_blend_071GC	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300033
Anti-D blend 175 2; 415 1E4 4049863Anti-D_blend_072GE	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-E MS-80; MS-258 4049863Anti-E_0819K	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-E MS-12; MS-260 4049863Anti-E_0829M	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-e MS-16; MS-21; MS-63 4049863Anti-e_091JX	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-e MS-62; MS-69 4049863Anti-e_092JZ	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300022; NB # 0483 related Supplements D1415300026 D1415300029
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-Kell MS-56 4049863Anti-K_101B2	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483 D1415300023; NB # 0483 related Supplements D1415300027 D1415300030
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-Kell AEK-4 4049863Anti-K_102B4	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300023; NB # 0483 related Supplements D1415300027 D1415300030
Anti- Kell coombsreactive 4049863Anti-K_103B6	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483 D1415300023; NB # 0483 related Supplements D1415300027 D1415300030
Anti-Fy^a coombsreactive 4049863Anti-FyA_123RW	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483
Anti-Fy^b coombsreactive 4049863Anti-FyB_133SC	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483
Anti- Jk^a MS-15 4049863Anti-JkA_141NE	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests <input type="checkbox"/> Class A _{sterile} product	N/A	D1415300025; NB # 0483

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Anti-Jk^a coombsreactive 4049863Anti-JkA_143NJ	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-Jk^b MS-8 4049863Anti-JkB_151NU	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Anti-Jk^b coombsreactive 4049863Anti-JkB_153NY	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
AHG polyspezifisch rabbit 4049863AHG_343QC	<input type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	D1415300025; NB # 0483
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input checked="" type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		

Table 2: Devices covered by this letter and for which the NB is NOT responsible for appropriate surveillance of the corresponding devices under the applicable Directive:

Device name or Basic UDI-DI (under IVDR application)	IVDR Device classification (as proposed by the manufacturer and verified at the pre-application stage)	If the IVDR device is a substitute device, identification of the corresponding IVDD device	IVDD Certificate Reference(s) of the devices under IVDR application, and the NB Identification
Control reagent incomplete 4049863CRI_333W3	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	N/A - Device did not require a Notified Body certificate under Directives
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
Control reagent monoclonal 4049863CRM_331XB	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests	N/A	N/A - Device did not require a Notified Body certificate under Directives
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		
ABO set MP Monoclonal Anti-A Birma-1; Anti-B LB-2 4049863BedSid441HZ	<input checked="" type="checkbox"/> Class D product <input type="checkbox"/> for self-testing <input checked="" type="checkbox"/> for use in near-patient tests	N/A	N/A - Device did not require a Notified Body certificate under Directives
	<input type="checkbox"/> Class C product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class B product <input type="checkbox"/> for self-testing <input type="checkbox"/> for use in near-patient tests		
	<input type="checkbox"/> Class A _{sterile} product		

Confirmation Letter Revision History

Date	NB internal reference traceable to each version of the letter	Action
2025-09-18	D1415300034	Initial

Gamintojas:
CE-IMMUNDIAGNOSTIKA GmbH
Karl-Landsteiner-Str. 6, 69151 Neckargemund, Vokietija
Tel.: +49 6223-80094 00 Faks.: +49 6223-80094 99

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DECLARATION OF CONFORMITY**Manufacturer:**

CE-IMMUNDIAGNOSTIKA GmbH
Karl-Landsteiner-Str. 6, 69151 Neckargemünd, Germany
Tel. +49 6223-80094 00 Fax +49 6223-80094 99

Revision: **2022-05-23**Valid until: **2025-05-26**

As the sole responsible manufacturer, we hereby declare that the products listed below comply with the provisions of Council Directive Annex I 98/79/EC for in-vitro diagnostic medical devices (IVDD). All supporting documents are kept on the premises of the manufacturer.

Supervisory authority :

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Notified body:

Identificationsnumber: 0483

mdc medical device certification GmbH
Kriegerstraße 6
70191 Stuttgart
Germany

Conformity assessment:

For the products named in Lists A and B of Appendix II, the manufacturer has to use the EC declaration of conformity procedure in accordance with Appendix IV (full quality assurance system), for all other products the procedure according to Annex III of the Directive 98/79/EC was applied.

Certificate numbers:

D1415300015_QMS13485
D1415300022_EC Rhesus annex IV.4
D1415300023_EC Kell annex IV.4
D1415300024_EC ABO annex IV.4
D1415300025_EC annex IV excluding section 4 and 6

place, date of issue:

Neckargemünd, 2022-05-23

Sign QMB:

(Angela Grajek)

CE-IMMUNDIAGNOSTIKA GmbH
Karl-Landsteiner-Str. 6
69151 Neckargemünd
Tel. +49 (0)6223 - 80094 00
Fax +49 (0)6223 - 80094 99
www.ce-immundiagnostika.com

Product: List A/List B

REF / Articlenumber	Product designation	Clone	Classification 98/79/EC
01110	Anti-A	A-11H5	Liste A
01210	Anti-A	BIRMA-1	Liste A
02110	Anti-B	B-6F9	Liste A
0221	Anti-B	LB-2	Liste A
03110	Anti-AB	A-5E10+B-2D7	Liste A
04105, 04110	Anti-C	MS-24	Liste A
04205, 04210	Anti-C	MS-273	Liste A
05105, 05110	Anti-c	MS-33	Liste A
05205, 05210	Anti-c	MS-35	Liste A
06110	Anti-D	MS-201	Liste A
06210	Anti-D	RUM-1	Liste A
07110	Anti-D blend	TH-28/MS-26	Liste A
07210	Anti-D blend	D175-2;D415 1E4	Liste A
06310	Anti-D incomplete	-	Liste A
08105, 08110	Anti-E	MS-80/MS-258	Liste A
08205, 08210	Anti-E	MS-12/MS-260	Liste A
09105, 09110	Anti-e	MS-16/MS-21/MS-63	Liste A
09205, 09210	Anti-e	MS-62/MS-69	Liste A
10105, 10110	Anti-Kell	MS-56	Liste A
10205, 10210	Anti-Kell	AEK4	Liste A
10305, 10310	Anti-Kell coombsreactive	-	Liste A
12302, 12305	Anti-Fy ^a coombsreactive	-	Liste B
13302, 13305	Anti-Fy ^b coombsreactive	-	Liste B
14302, 14305	Anti-Jk ^a coombsreactive	-	Liste B
15302, 15305	Anti-Jk ^b coombsreactive	-	Liste B
14102, 14105	Anti-Jk ^a	MS-15	Liste B
15102, 15105	Anti-Jk ^b	MS-8	Liste B
34310	AHG polyspez.-rabbit	-	Liste B

Other Products, self-declaration

REF / Articlenumber	Product designation	Clone	Classification 98/79/EC
30405, 30410	Anti-A1 Lectine	-	Other product
32102, 32105	Anti-C ^w	MS-110	Other product
32302, 32305	Anti-C ^w incomplete	-	Other product
29302	Anti-Di ^a coombsreactive	-	Other product
31105, 31110	Anti-H (A2)	10934C11	Other product
31405, 31410	Anti-H (A2) Lectine	-	Other product
11302, 11305	Anti-k cellano coombsreactive	-	Other product
16302	Anti-Kp ^a coombsreactive	-	Other product
17302	Anti-Kp ^b coombsreactive	-	Other product
33310	Control reagent incomplete	-	Other product
33110	Control reagent monoclonal	-	Other product
18102, 18105	Anti-Le ^a	LEA2	Other product
19102, 19105	Anti-Le ^b	LEB2	Other product
20302	Anti-Lu ^a coombsreactive	-	Other product
21302	Anti-Lu ^b coombsreactive	-	Other product
22102, 22105	Anti-M	11H2	Other product
23102, 23105	Anti-N	1422C7	Other product
26102, 26105	Anti-P1	650	Other product
24102, 24105	Anti-S	MS-94	Other product
24302, 24305	Anti-S coombsreactive	-	Other product
25102, 25105	Anti-s P3BER	P3BER	Other product
25302, 25305	Anti-s coombsreactive	-	Other product
27302	Anti-Wr ^a coombsreactive	-	Other product
28302	Anti-Xg ^a coombsreactive	-	Other product
35310	Anti-IgG monospez. Sheep	-	Other product

41160	Alsever's solution	-	Other product
36110	Albu-LISS reagent	-	Other product
36210, 36250	LISS reagent	-	Other product
37110	Bromelin solution	-	Other product
38110	Bromelin-control reagent	-	Other product
39110	Bovine albumin 22%	-	Other product
40110	Bovine albumin 30%	-	Other product

 CE-Immundiagnostika GmbH Karl-Landsteiner-Str. 6, D-69151 Neckargemünd Phone: +49 6223-80094 00 Telefax: +49 6223-80094 99 www.ce-immundiagnostika.com		
Instructions for Use		
Rev. 002/07-2021		
Description	REF	
Anti-A A-11H5 (monoclonal) 10 ml	01110	
Anti-A BIRMA1 (monoclonal) 10 ml	01210	
Anti-B B-6F9 (monoclonal) 10 ml	02110	
Anti-B LB2 (monoclonal) 10ml	02210	
Anti-AB A-5E10-B-2D7 (monoclonal) 10 ml	03110	

IN VITRO DIAGNOSTIC MEDICAL DEVICE

SUMMARY

In 1900, Karl Landsteiner discovered that the serum of some individuals would agglutinate the erythrocytes of others. Four different blood group phenotypes were recognised: 0, A, B and AB. Subgroups of A and B have since been identified.

Percentage of blood groups AB0 in Caucasians

Antigen		Antibody				AB0 Phenotype	Caucasians %
A	B	A ₁	A ₂	B	0		
+	0	0	0	+	0	A	42
0	+	+	+	0	0	B	10
0	0	+	+	+	0	0	44
+	+	0	0	0	0	AB	4

Table 1: Frequency of antigens

INTENDED USE

Anti-A, anti-B and anti-AB monoclonal murine IgM are designed for the specific and qualitative detection of erythrocytes that carry the corresponding AB0 antigen and are suitable for use by the slide, spot plate, microplate and tube techniques.

The abovementioned techniques are based on the principle of direct haemagglutination. After the addition of erythrocytes to the test reagents, a specific antigen-antibody-reaction occurs if the corresponding antigen is present on the erythrocytes. This reaction is indicated by easily visible agglutination of erythrocytes. Lack of agglutination of the red cells demonstrates – taking into account the limitations of the test methods – the absence of the corresponding antigens (negative test result).

The anti-A A-11H5 and anti-AB A-5E10-B-2D7 products are suitable for agglutinating and, therefore, identifying A_x.

PRODUCT INFORMATION

The monoclonal antibodies used to determine blood groups A, B and 0 are of the IgM isotype and are harvested from murine hybridoma cells. The antibodies are diluted in a buffer containing 0.9% sodium chloride, bovine albumin (no stabiliser), EDTA and reagents that facilitate the resuspension of the cell pellet after centrifugation. Preservative: < 0.1% sodium azide.

The anti-A and anti-B test sera are dyed blue and yellow, respectively, for easy identification.

These antisera have been optimized for use without further dilution or additions.

The lot number and expiry date are detailed on the tube label.

STORAGE

Store the test reagents at 2°C-8°C until the expiry date detailed on the product label. After opening a product for the first time, tightly close again and store at 2°C-8°C.

SPECIMEN COLLECTION AND PREPARATION

Blood specimens for AB0 blood typing should be collected aseptically in EDTA or citrate tubes. The specimen should be tested as soon as

possible following collection. If a delay in testing should occur, store the specimen at 2°C-8°C. Specimens displaying haemolysis or microbial contamination should not be tested with this reagent as this may result in false positive or false negative results.

All blood samples should be washed twice with a 0.9% NaCl solution before being tested by the tube, spot plate or microplate technique. When using the slide technique, prepare a 35-45% suspension of test erythrocytes (whole blood); when using the spot plate technique, use either whole blood or prepare a 10% suspension of test erythrocytes in a 0.9% NaCl solution.

WARNINGS AND PRECAUTIONS

- These reagents are intended for *in vitro* diagnostic laboratory use only.
- These reagents are designed for use by authorized operators trained in serological techniques.
- These reagents are not intended for self-application.
- Do not use these reagents past the expiration date.
- Discard the contents of damaged vials.
- The reagents have been filtered through a 0.2µm capsule to reduce the bio-burden.
- Once a vial has been opened the contents should remain viable up until the expiry date. Discard the contents if turbidity or contamination occur after opening.
- CE-Immundiagnostika GmbH cannot guarantee that products derived from human or animal sources are free from infectious agents. Care must be taken in the use and disposal of each vial and its contents.
- Wear protective clothing such as, for instance, a gown and disposable gloves when using these products

DISPOSAL OF REAGENT AND DEALING WITH SPILLAGES

For information on disposal of the reagent and decontamination of a spillage site, see the Material Safety Data Sheet, available on request from CE-Immundiagnostika GmbH.

CONTROLS AND ADVICE

- Positive and negative control erythrocytes shall be tested in parallel with each batch of tests. Tests must be considered invalid if controls do not show expected reactions.
- Since these reagents do not contain macromolecular potentiators, it is very unlikely that false positive or false negative reactions are caused with IgG coated cells.
- Weak antigens (A_x or B_{weak}) may not be identified. It is advisable to re-test ambiguous results using the tube technique.
- One volume is approximately 35-45µl when using the vial dropper provided.
- The reading and interpretation of results must be carried out by properly trained and qualified personnel in accordance with the requirements of the country where the reagents are in use.
- These reagents should be used only according to these instructions for use.

REQUIRED MATERIAL AND REAGENTS

- 0.9% NaCl solution
- Glass test tubes
- Test tube holder

Gebrauchsinformation

Rev. 002/07-2021

Anti-A, Anti-AB, Anti-B

- Test tube centrifuge
- Volumetric pipettes
- Microplate, plate shaker
- Spot plate
- Glass microscope slide
- Applicator stick
- Positive and negative control erythrocytes
- Chronometer

1. Using whole blood, prepare a 35-45% suspension of test erythrocytes.
2. Place 1 volume of antiserum and one volume of whole blood on a slide.
3. Using a clean applicator stick, thoroughly mix both volumes over an area of about 20x40mm.
4. Slowly move the slide back and forth.
5. Read macroscopically after no more than 2 minutes and record.
6. Incorrect handling or excessive incubation time may lead to drying-induced artefacts, and the test must be considered invalid.

RECOMMENDED TECHNIQUES

A. TUBE TECHNIQUE

1. Prepare a 2-4% suspension of erythrocytes in a 0.9% NaCl solution.
2. Place 1 volume of antiserum and 1 volume of test suspension of test erythrocytes in a labelled test tube.
3. Mix thoroughly and centrifuge straight away at 400g for 1 minute (at 1,500 rcf or for a suitable alternative time and force).
4. Read the result immediately: Gently agitate the tube to dislodge the erythrocyte pellet from the bottom of the test tube and read macroscopically for agglutination; record result.

B. MICROPLATE TECHNIQUE

Preparation of the microplates:

Microplates made by different manufacturers / suppliers have different static properties that may give rise to non-specific red blood cell and protein reactions. It is advisable to pre-treat unused microplates prior to use to reduce the build-up of red blood cells to a minimum. We recommend using "U" wells made from plastic material.

1. Place 1 volume of 22% bovine albumin (BSA) in the appropriate wells.
2. Mix thoroughly, gently agitating or using a microplate shaker, to ensure that wells are covered evenly.
3. Incubate at room temperature (18-25°C) for no less than 10 and no more than 15 minutes.
4. Drain the BSA and discard the content of the wells in a suitable waste disposal container.
5. Flush the microplate at least 10 times with tap water.
6. Then rinse the microplate twice with distilled or deionized water.
7. Tilt and dab the microplate to remove excess water.
8. Allow the microplate to dry prior to use.

Alternative techniques may be used provided that they have been validated by the user.

Procedure:

1. Prepare a 2 - 4% suspension of test erythrocytes in a 0.9% NaCl solution. (Recommendation: 2% suspension)
2. Using the vial dropper, place 30µl of the appropriate antiserum into the marked wells of the microplate.
3. Add 30µl of the previously prepared suspension of test erythrocytes to the microplate.
4. Mix for 30 seconds, either manually or using a shaker.
5. Centrifuge the microplate for 1 minute at 400g (at 1,500 rcf or for a suitable alternative time and force).
6. Briefly agitate the microplate, using a shaker, if necessary.

Record the result and the reaction intensity, testing positive and negative control erythrocytes in parallel. Reading devices, if used, must have been validated. Using additional visual expedients such as, for instance, mirrors or magnifying glasses may facilitate reading the results.

C. SLIDE TECHNIQUE

D. SPOT PLATE TECHNIQUE

1. Using whole blood, prepare a 35-45% suspension of test erythrocytes or a 10% suspension of test erythrocytes in a 0.9% NaCl solution.
2. Place 1 volume of antiserum + 1 volume of suspension of test erythrocytes on a spot plate.
3. Using a clean applicator stick, mix both volumes thoroughly.
4. Incubate at RT for 5-10 minutes.
5. Read macroscopically and record.

Incorrect handling or excessive incubation time may lead to drying-induced artefacts.

INTERPRETATION OF TEST RESULTS

1. **Positive:** Agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the presence of the appropriate AB0 antigen on the test erythrocytes.
2. **Negative:** No agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the absence of the appropriate AB0 antigen on the test erythrocytes.

Discrepancies: If the results obtained for antigens on the test erythrocytes do not correlate with the blood group alloantibodies detected, further investigation is required.

LIMITATIONS

1. AB0 antigens are not fully developed at birth and so weaker reactions may therefore occur with cord or neonatal specimens.
2. When using monoclonal antibodies, blood specimens of weak A or B subgroups may give rise to false negative or weak reactions when tested using slides or microplates. It is advisable to re-test weak subgroups using the tube technique.
3. Stored blood may give weaker reactions than fresh blood.
4. False positive or false negative results may also occur due to:
 - Contamination of test materials
 - Improper storage, erythrocyte concentration, incubation time or temperature
 - Improper or excessive centrifugation
 - Deviation from the recommended techniques
5. When using the spot plate technique and whole blood specimens, rouleaux formation may occasionally occur that resembles weak agglutination and may be interpreted as a false positive reaction. This phenomenon is due to non-immunological causes. Rouleaux formation can be expected in heparin blood and in patients treated with plasma expanders (e. g. dextran) as well as in patients with plasmacytoma (high protein content, altered protein composition), oncological disorders (abnormal hemogram) and impaired coagulation. Blood specimens from these patients must always be tested using the tube technique, as this phenomenon can generally not be verified when using a suspension of erythrocytes. The presence of weak antigens (A₃, A_x, B_{weak}) should be demonstrated using the tube technique, because of greater sensitivity, after an incubation time of 30 minutes, where appropriate. Unusually weak reactions or no reaction at all may possibly be due to the presence of A and B

Gebrauchsinformation

Rev. 002/07-2021

Anti-A, Anti-AB, Anti-B

subgroups. Blood specimens from patients suffering from certain disorders may show false positive/false negative reactions. Cord blood specimens contaminated with Wharton's jelly may show false positive reactions. Do not use the indirect Coombs test (ICT) with monoclonal murine reagents and AHG sera. Reactive autoagglutinins are frequently a potential source of error in AB0 blood typing. The presence of these antibodies is unpredictable. A large enough number of these antibodies may cause unspecific reactions when demonstrating isoagglutinins with A-(1) and B cells as well as unspecific reactions with anti-A, anti-B and anti-AB test sera if unwashed cells or specimens suspended in sera or plasma specimens are tested. AB0 and isoagglutinin assays are therefore indispensable for blood group determination. Any discrepancies between the AB0 and the isoagglutinin assays require further investigation, regardless of the strength of the reaction in the patient's serum or cells. Do not use monoclonal reagents prepared from murine hybridoma cells for direct antiglobulin tests with AHG reagents.

- Always use two different clones for the determination of AB0 and Rh antigens.

STABILITY OF THE REACTIONS

- Read all tube and microplate tests straight after centrifugation.
- Slide tests should be interpreted within 2 minutes to ensure specificity and to avoid the possibility a negative result may be incorrectly interpreted as positive due to drying of the reagent.
- Tests must be considered invalid if they have been performed at temperatures other than those recommended.

SPECIFIC PERFORMANCE CHARACTERISTICS

The performance is assessed in accordance with the Common Technical Specifications (CTS: EU Commission Decision of 03/02/2009).

- The reagents have been tested using all the recommended procedures.
- Each LOT is tested in accordance with the requirements of the Common Technical Specifications by *in vitro* diagnostic medical devices and meets the requirements.
- Specificity of monoclonal antibodies is demonstrated using a panel of antigen-negative erythrocytes.
- The quality control of the reagents was performed using erythrocytes or whole blood that had been washed twice with 0.9% saline solution.
- Tests with more than 1,000 specimens with a sensitivity and specificity of > 99%

DISCLAIMER

- The user is responsible for the performance of the reagents by any method other than those recommended.
- Any deviations from the Recommended Techniques should be validated prior to use.

BIBLIOGRAPHY

- Kohler C. & Milstein C. (1975), Continuous cultures of fused cells secreting antibody of predefined specificity. *Nature*, 256, 495-497.
- Lee H.H., Rouger P., Germain C., Muller A. & Salmon C. (1983). The production and standardisation of monoclonal antibodies as AB blood group typing reagents. Symposium of International Association of Biological Standardisation on monoclonal antibodies.
- Human Blood Groups, by Geoff Daniels, 1st Ed., Blackwell Science, Oxford 1995.

- HMSO, Guidelines for Blood Transfusion Services., 2nd Ed., 1994.

INDEX OF SYMBOLS

	Lot number		For <i>in vitro</i> diagnostic use only
	Catalogue number		Store between +2°C to +8°C
	Expiry date		Manufacturer
	Consult instructions for use (insert)		

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LITERATŪRA

- Kohler C. & Milstein C. (1975), Continuous cultures of fused cells secreting antibody of predefined specificity. Nature, 256, 495-497.
- Lee H.H., Rouger P., Germain C., Muller A. & Salmon C. (1983). The production and standardisation of monoclonal antibodies as AB blood group typing reagents. Symposium of International Association of Biological Standardisation on monoclonal antibodies.
- Human Blood Groups, by Geoff Daniels, 1st Ed., Blackwell Science, Oxford 1995.
- HMSO, Guidelines for Blood Transfusion Services., 2nd Ed., 1994.

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SUMMARY

The rhesus blood group system was first described in 1940. In addition to the ABO antigens, the D antigen is of special clinical importance since a blood transfusion error may cause severe haemolytic transfusion reactions and haemolytic disease of the foetus and newborn (HDFN).

Rhesus positive (D pos.)/rhesus negative (D neg.) mean that a D antigen is present/not present, respectively, on the surface of an individual's erythrocytes.

Distribution of the D antigen

Anti D	Phenotype	Caucasians %	African Americans %
+	Rh D positive	85	72
-	Rh D negative	15	28

Table 1: Frequency of antigens according to Mollison et al.

INTENDED USE

Anti D (monoclonal IgM) is designed for the specific and qualitative detection of erythrocytes that carry the corresponding D (RH1) antigen and is suitable for use by the slide, spot plate, microplate and tube techniques.

The abovementioned techniques are based on the principle of direct haemagglutination. After the addition of erythrocytes to anti D (monoclonal IgM), a specific antigen-antibody-reaction occurs if the corresponding D (RH1) antigen is present on the erythrocytes. This reaction is indicated by easily visible agglutination of erythrocytes. Lack of agglutination of the red cells demonstrates – taking into account the limitations of the test methods – the absence of the corresponding D (RH1) antigen.

Both clones recognize most D_{weak} variants, with the exception of D VI.

PRODUCT INFORMATION

Monoclonal Anti D reagents of the IgM type are harvested from human hybridoma cells.

The dilution medium for this low-protein reagent contains the following potentiators: NaCl, BSA, a buffer and some other select components. Preservative: < 0.1% sodium azide.

When used in accordance with the instructions for use, the reagent agglutinates human erythrocytes if the corresponding D (RH1) antigen is present.

Anti D (monoclonal IgM) will recognize, and directly agglutinate, many red cells of weak D phenotype previously likely to be interpreted as Rh negative (or weak D). This also includes some types of extremely rare partial D cells. The IgM anti D clones show no reactivity with the D Category VI cells that have been tested so far. Weak variants can be detected only using the tube and microplate techniques.

The specificity of each LOT is demonstrated using the tube technique and a panel of D (RH1) antigen-negative erythrocytes.

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The LOT number and expiry date are detailed on the tube label.

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Blood specimens should be collected aseptically in EDTA or citrate tubes. The specimen should be tested as soon as possible following collection. If a delay in testing should occur, store the specimen at 2°C-8°C. Specimens displaying haemolysis or microbial contamination should not be tested with this reagent as this may result in false positive or false negative results.

All blood samples should be washed twice with a 0.9% NaCl solution before being tested by the tube, spot plate and microplate technique. When using the slide technique, prepare a 35-45% suspension of test erythrocytes (whole blood).

WARNINGS AND PRECAUTIONS

- These reagents are intended for *in vitro* diagnostic laboratory use only.
- These reagents are designed for use by authorized operators trained in serological techniques.
- These reagents are not intended for self-application.
- Do not use these reagents past the expiration date.
- Discard the contents of damaged vials.
- Wear protective clothing such as, for instance, a laboratory gown and disposable gloves while handling these products.
- The test sera have been filtered through a 0.2µm capsule to reduce the bio-burden.
- Once a vial has been opened the contents should remain viable up until the expiry date. Discard the contents if turbidity or contamination occur after opening.
- CE-Immundiagnostika GmbH cannot guarantee that products derived from human or animal sources are free from infectious agents. Care must be taken in the use and disposal of each vial and its contents.

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CONTROLS AND ADVICE

- Positive and negative control erythrocytes shall be tested in parallel with each batch of tests. Tests must be considered invalid if controls do not show expected reactions.
- Since these reagents do not contain macromolecular potentiators, it is very unlikely that false positive or false negative reactions are caused with IgG coated cells.
- Weak antigens (partial D) may not be identified.
- 1 volume is approximately 35-45µl when using the vial dropper provided.
- The reading and interpretation of results must be carried out by properly trained and qualified personnel in accordance with the requirements of the country where the reagents are in use.
- These reagents should be used only according to these instructions for use.

REQUIRED MATERIAL AND REAGENTS

- 0.9% NaCl solution
- Glass test tubes
- Test tube holder
- Test tube centrifuge
- Volumetric pipettes
- Microplate, plate shaker
- Spot plate

Instructions for Use

Rev. 002/07-2021

Anti D

- Glass microscope slide
- Applicator stick
- Positive and negative control erythrocytes
- Chronometer

RECOMMENDED TECHNIQUES

A. TUBE TECHNIQUE

1. Prepare a 2-4% suspension of erythrocytes in a 0.9% NaCl solution.
2. Place 1 volume of antiserum and 1 volume of suspension of test erythrocytes in a labelled test tube.
3. Mix thoroughly and centrifuge straight away at 400g for 1 minute (at 1,500 rcf or for a suitable alternative time and force).
4. Read the result immediately: Gently agitate the tube to dislodge the erythrocyte pellet from the bottom of the test tube and read macroscopically for agglutination; record result.
5. Incubate all negative or weak positive results for 15 minutes at room temperature (18-25°C), then repeat 3. and 4. Record the result and the reaction intensity.

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1. Place 1 volume of bovine albumin 22% (BSA) in the appropriate wells.
2. Mix thoroughly, gently agitating or using a microplate shaker, to ensure that wells are covered evenly.
3. Incubate at room temperature (18-25°C) for no less than 10 and no more than 15 minutes.
4. Drain the BSA and discard the content of the wells in a suitable waste disposal container.
5. Flush the microplate at least 10 times with tap water.
6. Then rinse the microplate twice with distilled or deionized water.
7. Tilt and dab the microplate to remove excess water.
8. Allow the microplate to dry prior to use.

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Procedure:

1. Prepare a 2 - 4% suspension of test erythrocytes in a 0.9% NaCl solution.
2. Place 30µl of the appropriate antiserum into the marked wells of the microplate.
3. Add 30µl of the previously prepared suspension of test erythrocytes to the microplate.
4. Mix for 30 seconds, either manually or using a shaker.
5. Centrifuge the microplate for 1 minute at 400g (at 1,500 rcf or for a suitable alternative time and force).
6. Briefly agitate the microplate, using a shaker, if necessary.
7. Record the result and the reaction intensity straight away, testing positive and negative control erythrocytes in parallel. Reading devices, if used, must have been validated. Using additional visual expedients such as, for instance, mirrors or magnifying glasses may facilitate reading the results.

C. SLIDE TECHNIQUE

1. Using whole blood, prepare a 35-45% suspension of test erythrocytes.

2. Place 1 volume of antiserum and one1 volume of whole blood on a slide.
3. Using a clean applicator stick, thoroughly mix both volumes over an area of about 20x40mm.
4. Slowly move the slide back and forth.
5. Read macroscopically after no more than 2 minutes and record.
6. Incorrect handling or excessive incubation time may lead to drying-induced artefacts, and the test must be considered invalid.

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INTERPRETATION OF TEST RESULTS

1. **Positive:** Agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the presence of the D antigen on the test erythrocytes.
2. **Negative:** No agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the absence of the D antigen on the test erythrocytes.

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LIMITATIONS

1. Stored blood may give weaker reactions than fresh blood.
2. False positive or false negative results may also occur due to:
 - Contamination of test materials
 - Improper storage, erythrocyte concentration, incubation time or temperature
 - Improper or excessive centrifugation
 - Deviation from the recommended techniques
3. Rouleaux formation can be expected in heparin blood and in patients treated with plasma expanders (e. g. dextran) as well as in patients with plasmacytoma (high protein content, altered protein composition), oncological disorders (abnormal hemogram) and impaired coagulation. Blood specimens from these patients must always be tested using the tube technique, as this phenomenon can generally not be demonstrated when using a suspension of erythrocytes. The presence of weak antigens should be demonstrated using the tube technique, because of greater sensitivity, after an incubation time of 30 minutes, where appropriate.
4. Blood specimens from patients suffering from certain disorders may show false positive/false negative reactions. Cord blood specimens contaminated with Wharton's jelly may show false positive reactions.
5. Always use two different clones for the determination of D antigens.

STABILITY OF THE REACTIONS

1. Read all tube and microplate tests straight after centrifugation.
2. Slide tests should be interpreted within 2 minutes to ensure specificity and to avoid the possibility a negative result may be incorrectly interpreted as positive due to drying of the reagent.

	CE-Immundiagnostika GmbH Karl-Landsteiner-Str. 6, D-69151 Neckargemünd Phone: +49 6223-80094 00 Telefax: +49 6223-80094 99 www.ce-immundiagnostika.com	
Instructions for Use		
Rev. 002/07-2021	Anti D	

- Tests must be considered invalid if they have been performed at temperatures other than those recommended.

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- Prior to approval, the antisera have been tested using all the recommended procedures.
- Each LOT of anti D (monoclonal IgM) is tested in accordance with the requirements of the Common Technical Specifications for products defined in Annex II, List A of Directive 98/79/EC of the European Parliament and the Council on *in vitro* diagnostic medical devices and meets the requirements.
- The reactivity of each LOT is demonstrated using a panel of erythrocytes in accordance with the techniques recommended in these instructions for use.
- Specificity of monoclonal antibodies is demonstrated using a panel of antigen-negative erythrocytes.
- The quality control of the reagents was performed using erythrocytes or whole blood that had been washed twice with 0.9% saline solution.
- Tests with more than 1,000 specimens with a sensitivity and specificity of > 99%.

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SUMMARY

The rhesus blood group system was first described in 1940. In addition to the ABO antigens, the D antigen is of special clinical importance since a blood transfusion error may cause severe haemolytic transfusion reactions and haemolytic disease of the foetus and newborn (HDFN).

Rhesus positive (D pos.)/rhesus negative (D neg.) mean that a D antigen is present/not present, respectively, on the surface of an individual's erythrocytes.

Distribution of the D antigen

Anti D	Phenotype	Caucasians %	African Americans %
+	Rh D positive	85	72
-	Rh D negative	15	28

Table 1: Frequency of antigens according to Mollison et al.

INTENDED USE

Anti D blend (monoclonal IgM + IgG) is designed for the specific and qualitative detection of erythrocytes that carry the corresponding D (RH1) antigen and is suitable for use by the slide, spot plate, microplate and tube techniques.

The abovementioned techniques are based on the principle of direct haemagglutination. After the addition of erythrocytes to anti D blend (monoclonal IgM + IgG), a specific antigen-antibody-reaction occurs if the corresponding D (RH1) antigen is present on the erythrocytes. This reaction is indicated by easily visible agglutination of erythrocytes. Lack of agglutination of the red cells demonstrates – taking into account the limitations of the test methods – the absence of the corresponding D (RH1) antigen.

The indirect antiglobulin test is used for the detection of D_{weak}.

PRODUCT INFORMATION

Anti D blend (monoclonal IgM + IgG) is a blend of human monoclonal anti D (IgM), clone D175-2 or TH-28, and human monoclonal anti D (IgG), clone D415 1E4 or MS-26.

The dilution medium for this low-protein reagent contains the following potentiators: NaCl, BSA, a buffer and some other select components. Preservative: < 0.1% sodium azide.

When used in accordance with the instructions for use, the reagent agglutinates human erythrocytes if the corresponding D (RH1) antigen is present.

Anti D blend (monoclonal IgM + IgG) will directly agglutinate red cells from most weak D previously likely to be interpreted as Rh negative (or weak D). This also includes some types of extremely rare partial D cells. The IgM anti D clones D175-2 and TH-28 show no reactivity with the D Category VI cells that have been tested so far. IgG clones D415 1E4 and MS-26 are capable of giving a reaction with D Category VI cells.

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Blood specimens should be collected aseptically in EDTA or citrate tubes. The specimen should be tested as soon as possible following collection. If a delay in testing should occur, store the specimen at 2°C-8°C. Specimens displaying haemolysis or microbial contamination should not be tested with this reagent as this may result in false positive or false negative results.

All blood samples should be washed twice with a 0.9% NaCl solution before being tested by the tube, spot plate and microplate technique. When using the slide technique, prepare a 35-45% suspension of test erythrocytes (whole blood); when using the spot plate technique, use either whole blood or prepare a 10% suspension of test erythrocytes in a 0.9% NaCl solution.

WARNINGS AND PRECAUTIONS

- These reagents are intended for *in vitro* diagnostic laboratory use only.
- These reagents are designed for use by authorized operators trained in serological techniques.
- These reagents are not intended for self-application.
- Do not use these reagents past the expiration date.
- Discard the contents of damaged vials.
- Wear protective clothing such as, for instance, a gown and disposable gloves when using these products.
- The test sera have been filtered through a 0.2µm capsule to reduce the bio-burden.
- Once a vial has been opened the contents should remain viable up until the expiry date. Discard the contents if turbidity or contamination occur after opening.
- CE-Immundiagnostika GmbH cannot guarantee that products derived from human or animal sources are free from infectious agents. Care must be taken in the use and disposal of each vial and its contents.

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CONTROLS AND ADVICE

- Positive and negative control erythrocytes shall be tested in parallel with each batch of tests. Tests must be considered invalid if controls do not show expected reactions.
- Since these reagents do not contain macromolecular potentiators, it is very unlikely that false positive or false negative reactions are caused with IgG coated cells.
- Weak antigens (partial D) may not be identified.
- 1 volume is approximately 35-45µl when using the vial dropper provided.
- The reading and interpretation of results must be carried out by properly trained and qualified personnel in accordance with the requirements of the country where the reagents are in use.
- These reagents should be used only according to these instructions for use.

REQUIRED MATERIAL AND REAGENTS

- 0.9% NaCl solution
- Glass test tubes
- Test tube holder
- Test tube centrifuge
- Volumetric pipettes

Instructions for Use

Rev. 003/01-2022

Anti D blend

- Microplate, plate shaker
- Spot plate
- Glass microscope slide
- Applicator stick
- Positive and negative control erythrocytes
- Chronometer

RECOMMENDED TECHNIQUES
A. TUBE TECHNIQUE

1. Prepare a 2-4% suspension of erythrocytes in a 0.9% NaCl solution.
2. Place 1 volume of antiserum and 1 volume of suspension of test erythrocytes in a labelled test tube.
3. Mix thoroughly and centrifuge straight away at 400g for 1 minute (at 1,500 rcf or for a suitable alternative time and force).
4. Read the result immediately: Gently agitate the tube to dislodge the erythrocyte pellet from the bottom of the test tube and read macroscopically for agglutination; record result.
5. Incubate all negative or weak positive results for 5 minutes at room temperature (18-25°C), then repeat 3. and 4. Record the result and the reaction intensity.

INDIRECT ANTIGLOBULIN TEST FOR THE DETECTION OF WEAK D VARIANTS (D_{weak})

1. Prepare a 2-4% suspension of erythrocytes in a 0.9% NaCl solution.
2. Place 1 volume of antiserum and 1 volume of suspension of test erythrocytes in a labelled test tube.
3. Reaction time: 15 minutes at 37°C.
4. Rinse 3 times using isotonic saline solution.
5. Add to volumes of anti-human globulin reagent and centrifuge straight away at 400g for 1 minute (at 1,500 rcf or for a suitable alternative time and force).
6. Holding the tube at an angle, gently agitate it to dislodge the sediment from the bottom and check for agglutination.
7. Record the result and the reaction intensity. The presence of agglutinates indicates a weak expression of the D antigen. All negative reactions should be tested using IgG-sensitized cells (Coombs Control Cells) to test the reactivity of the AHG serum.

B. MICROPLATE TECHNIQUE
Preparation of the microplates:

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Procedure:

1. Prepare a 2 - 4% suspension of test erythrocytes in a 0.9% NaCl solution. (Recommendation: 2% suspension)
2. Place 30µl of the appropriate antiserum into the marked wells of the microplate.
3. Add 30µl of the previously prepared suspension of test erythrocytes to the microplate.
4. Mix for 30 seconds, either manually or using a shaker.
5. Centrifuge the microplate straight away for 1 minute at 400g (at 1,500 rcf or for a suitable alternative time and force).
6. Briefly agitate the microplate, using a shaker, if necessary.

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C. SLIDE TECHNIQUE

1. Using whole blood, prepare a 35-45% suspension of test erythrocytes.
2. Place 1 volume of antiserum and one1 volume of whole blood on a slide.
3. Using a clean applicator stick, thoroughly mix both volumes over an area of about 20x40mm.
4. Slowly move the slide back and forth.
5. Read macroscopically after no more than 2 minutes and record.
6. Incorrect handling or excessive incubation time may lead to drying-induced artefacts, and the test must be considered invalid.

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INTERPRETATION OF TEST RESULTS

1. **Positive:** Agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the presence of the D antigen on the test erythrocytes.
2. **Negative:** No agglutination of the test erythrocytes indicates, within accepted limitations of test procedure (see below), the absence of the D antigen on the test erythrocytes.

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Anti D blend

LIMITATIONS

1. Stored blood may give weaker reactions than fresh blood.
2. False positive or false negative results may also occur due to:
 - Contamination of test materials
 - Improper storage, erythrocyte concentration, incubation time or temperature
 - Improper or excessive centrifugation
 - Deviation from the recommended techniques
3. Rouleaux formation can be expected in heparin blood and in patients treated with plasma expanders (e. g. dextran) as well as in patients with plasmacytoma (high protein content, altered protein composition), oncological disorders (abnormal hemogram) and impaired coagulation. Blood specimens from these patients must always be tested using the tube technique, because of greater sensitivity, after an incubation time of 30 minutes, where appropriate.
4. Blood specimens from patients suffering from certain disorders may show false positive/false negative reactions. Cord blood specimens contaminated with Wharton's jelly may show false positive reactions.
5. Always use two different clones for the determination of D antigens.

STABILITY OF THE REACTIONS

1. Read all tube and microplate tests straight after centrifugation.
2. Slide tests should be interpreted within 2 minutes to ensure specificity and to avoid the possibility a negative result may be incorrectly interpreted as positive due to drying of the reagent.
3. Tests must be considered invalid if they have been performed at temperatures other than those recommended.

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1. The antisera have been tested using all the recommended procedures.
2. Each LOT of anti D blend (monoclonal IgM + IgG) is tested in accordance with the requirements of the Common Technical Specifications for products defined in Annex II, List A of Directive 98/79/EC of the European Parliament and the Council on *in vitro* diagnostic medical devices and meets the requirements.
3. The reactivity of each LOTs was tested using a panel of erythrocytes in accordance with the techniques recommended in these instructions for use.
4. Specificity of monoclonal antibodies is demonstrated using a panel of antigen-negative erythrocytes.
5. The quality control of the reagents was performed using erythrocytes or whole blood that had been washed twice with 0.9% saline solution.
6. Tests with more than 3,000 specimens with a sensitivity and specificity of > 99%.

DISCLAIMER

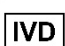
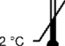
1. The user is responsible for the performance of the reagents by any method other than those recommended.
2. Any deviations from the Recommended Techniques should be validated prior to use.

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IN-VITRO-DIAGNOSTIKUM

SUMMARY

Blood group serological examinations are carried out using reagents with specific antibodies and other components. A control reagent contains all components except the specific antibody and thus produces a negative reaction.

INTENDED USE

Control reagent monoclonal is a negative control for monoclonal reagents.

PRODUCT INFORMATION

The reagents are suspended in a buffered 0.9% NaCl solution that contains bovine albumin (without stabilizer), EDTA and reagents that enable the cell button to be resuspended more easily after centrifugation. Preservative: Na azide.

These reagents are used without further dilution / additives.

The LOT and expiration date are on the vial label.

STORAGE

The test reagents can be used up to the expiry date stated on the label when stored at + 2°C - + 8°C. After opening the test reagents for the first time, store them tightly closed at + 2°C - + 8°C.

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Blood samples should be collected in EDTA or citrate tubes according to standard medical procedures. The evaluation should be carried out as soon as possible after the blood sample has been drawn. If the blood is not to be used immediately, the tubes must be stored at + 2 ° C - + 8 ° C. Blood samples showing haemolysis or microbial contamination should not be used for the test. Such blood tests can give incorrect results.

All blood samples are washed twice in 0.9% NaCl solution for the tube and microtiter plate test before use. Whole blood (35-45% erythrocyte suspension) is used for the slide test, whole blood or a 10% erythrocyte suspension in 0.9% NaCl solution for the spot plate test.

WARNINGS AND PRECAUTIONS

1. The reagents are intended for in vitro diagnostic laboratory use only
2. The reagents may only be used by authorized specialists.
3. The reagents are not intended for personal use.
4. After the expiry date, the reagents may no longer be used
5. Damaged vials must not be used
6. The reagents contain <0.1% sodium azide. In high concentrations, sodium azide can react with lead and copper and form an explosive metal azide compound.
7. Wear protective clothing such as a gown and disposable gloves when using the products
8. The test sera were filtered through a 0.2 µm membrane in order to reduce the bacterial load.
9. Once opened, the contents should be used up to the expiration date. Should it become cloudy or contaminated after opening, the contents should be discarded.
10. CE-Immundiagnostika GmbH cannot guarantee that human and animal raw materials are free from infectious agents, so the products should be used with caution.

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For disposal of the reagents or decontamination in the event of spillage, please request the safety data sheet from CE-Immundiagnostika GmbH.

CONTROLS AND ADVICE

1. Positive and negative control erythrocytes must be included with each experiment. If the controls do not show the expected results, the test batch should be discarded.
2. Since the test reagents do not contain any macromolecular enhancement media, it is very unlikely that false positive or false negative signals will occur, which are caused by IgG-laden cells.
3. 1 drop from the pipette vial corresponds to 35-45µl.
4. Only authorized specialists are allowed to read and evaluate the results.
5. The test reagents may only be used as described here.

REQUIRED MATERIAL AND REAGENTS

- 0.9% NaCl-Solution
- Glas tube
- Tube holder
- Centrifuge
- Microplate, Shaker
- Spot plate
- Glass mikroskope slide
- Applicator stick
- Positive and negative Control test cells
- timer

RECOMMENDED TECHNIQUES**A. METHODE: TUBE TEST**

1. A 2-4% erythrocyte suspension is prepared in 0.9% NaCl solution
2. 1 drop of test reagent and 1 drop of erythrocyte suspension are placed in a labeled tube.
3. Mix well and centrifuge immediately for 1 min. At 400g (1500 rpm, or at an alternative speed with an appropriate time).
4. Read the result immediately: loosen the erythrocyte pellet from the bottom of the tube by gently shaking it and read and record the agglutination strength macroscopically.

B. METHODE: MICROPLATE TECHNIQUE**Preparation of the microplates:**

MTP from different manufacturers / suppliers show different static properties, which can result in non-specific reactions of red blood cells and proteins. It is recommended that unused MTP be pre-treated prior to use to minimize red blood cell adhesion. MTP with a U-profile made of plastic are recommended.

1. Add 1 drop of 22% beef albumin (BSA) to each MTP well.
2. Mix by gentle agitation or on an MTP shaker so that the wells are evenly coated.
3. Let the MTP stand for at least 10 to 15 minutes at room temperature (18-25 ° C).
4. Drain the BSA and transfer the contents of the MTP wells to an appropriate waste container.
5. Rinse the MTP at least 10 times with tap water.
6. Then rinse the MTP twice with distilled or deionized water.

Gebrauchsinformation

Rev. 002/06.2021

Gebrauchsinformation: Control reagent monoclonal REF 33110

7. Tilt and blot the MTP to remove excess water.
 8. Allow the MTP to dry before use.
- Alternative methods are possible, provided they have been appropriately validated by the user.

Procedure:

1. Prepare a 2 - 4% erythrocyte suspension in 0.9% NaCl solution. (Recommendation 2% suspension)
2. Pipette 1 drop of the appropriate anti-serum into the marked MTP wells.
3. Add one drop of the prepared erythrocyte suspension to the MTP.
4. Mix MTP manually or on a shaker for 30 seconds.
5. Centrifuge the MTP for 1 min. At 400g (1500 rpm, or at an alternative speed with an adapted time).
6. Shake the MTP briefly, if necessary with a shaker.

Record the result and reaction strength, positive and negative controls must be included. When using readers, these must be validated. The use of additional visual aids such as test reading mirrors or magnifying glasses can make reading easier.

C. METHODE: Slide Technique

1. Whole blood is used for the slide test (35-45% erythrocyte suspension)
2. Place 1 drop of antiserum and a drop of red cell suspension on the slide.
3. Mix both drops thoroughly with a clean stirring stick over an area of approximately 20x40mm.
4. Slowly move the slide back and forth
5. Read the result macroscopically after a maximum of 2 minutes and record it.
6. In the event of improper handling or an incubation time that is too long, drying artifacts may occur; the test must be discarded.

D. METHODE: Spot Plate Technique

1. Whole blood (35-45% erythrocyte suspension) or a 10% erythrocyte suspension in 0.9% NaCl solution is used for the spot plate test.
 2. Put 1 drop of antiserum + 1 drop of erythrocyte suspension on the spot plate.
 3. Mix both drops thoroughly with a clean stirring rod.
 4. 5-10 min. Incubation time at room temperature.
 5. Read the result macroscopically and record it.
- In the event of improper handling or an incubation time that is too long, drying artifacts can occur.

INTERPRETATION OF TEST RESULT

1. **Positive:** Erythrocyte agglutination. Please note the limits of the test method (see below).
2. **Negative:** no agglutination Please note the limits of the test method (see below).

Discrepancies: A positive reaction rejects the antigen result.

LIMITATION

1. Blood that is not freshly drawn can lead to weaker results.
2. False positive or false negative results can be caused by:
 - Contamination of the material to be tested
 - Incorrect storage, incorrect erythrocyte concentration, incorrect incubation time, incorrect temperature
 - Incorrect centrifugation
 - Deviations from the recommended methods
3. Blood samples can occasionally react in the spotted plate test with whole blood with roll formation, which resembles weak agglutination and can be evaluated as false positive. The phenomenon has non-immunological causes. Rolling of money must be expected in heparin bleeding, patients who have been

treated with plasma expanders (e.g. dextran), as well as in patients with plasmacytomas (high protein content, altered protein composition), oncological diseases (abnormal blood count) and coagulation disorders. A tube test should definitely be used in these patients, as the phenomenon is usually undetectable when using suspension.

STABILITY OF THE REACTIONS

1. Tube and microtiter plate tests must be read immediately after centrifugation.
2. Slide tests must be read within 2 minutes to ensure specificity and to avoid false positive results due to drying artifacts.
3. If a temperature other than the recommended one has been selected, the results must be discarded

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1. The reagents were tested by all recommended methods prior to release.
2. Tested on over 1000 samples with sensitivity and specificity of 100%

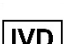
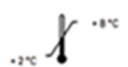
DISCLAIMER

1. The user is liable if methods other than those recommended are used.
2. Any deviations from the recommended test methods must be validated prior to use.

BIBLIOGRAPHY

1. Köhler C. & Milstein C. (1975), Continuous cultures of fused cells secreting antibody of predefined specificity. Nature, 256, 495-497.
2. Mollison P. L. and Engelfriet C. P. (1997), Blood Transfusion in Clinical Medicine, Blackwell Science, 10th Edition.
3. AABB Technical Manual, 13th Edition.
4. Race R. & Sanger R. (1975), Blackwell Science, 6th Edition.

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Anti-AB	
Reference No.:	03110
LOT:	OABM111-3
Clone:	A-5E10+B-2D7
Titer:	<div>A₁</div> <div>A₂</div> <div>B</div> <div>1:512</div> <div>1:128</div> <div>1:512</div>
Manufacture date:	2025-01
Expiry date:	2027-01
Purpose of use:	<i>In-vitro-diagnosticum</i>
Colour:	colourless / clear
PH Value:	7.3
Specificity:	Untested antigens circled: A, B, H, Le ^a , Le ^b , I, K, k, Kp ^a , P ₁ , D, C, E, c, e, C ^w , M, N, Bg ^a , S, s, Lu ^a , Jk ^a , Jk ^b , Di ^a , Fy ^a , Fy ^b , Xg ^a , Do ^a , Do ^b , Yt ^b , Co ^b , Mg, Wr ^a requirements fulfilled <input checked="" type="checkbox"/> not fulfilled <input type="checkbox"/>
Avidity:	tested
Microbial status:	This product was filtered through a filter with a pore size of 0.2 µm and subjected to a bacterial load test on Columbia agar, Tryptic Soy Agar and Sabouraud Dextrose Agar with an incubation time of 2 to 5 days, <1 CFU / ml.
Preservatives:	< 0.1 % Sodium-Azide
Storage:	+2 - +8°C
Origin:	Culture supernatant / monitored animal populations / material of human origin
Virus detection:	Materials of human origin have been tested and are negative for HIV 1 + 2, HCV, HBV (HbsAg). Adequate safety precautions are recommended when handling potentially infectious materials.
Note:	All cattle materials used come from veterinary controlled and certified animal herds. Products of bovine origin are assessed as having a low TSE risk.

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Anti-A	
Referenz No.:	01110
LOT:	OAM168-6
Clone:	A-11H5
Titer: A ₁ : A ₂ B:	1:512 1:128
Manufacture date:	2024-12
Expiry date:	2026-11
Purpose of use:	<i>In-vitro-Diagnosticum</i>
Colour:	blue / clear
PH value:	7.3
Specificity:	<p>Circled antigens not tested:</p> <p>A, B, H, Le^a, Le^b, I, K, k, Kp^a, P₁, D, C, E, c, e, C^w, M, N, Bg^a, S, s, Lu^a, Jk^a, Jk^b, Di^a, Fy^a, Fy^b, Xg^a, Do^a, Do^b, Yt^b, Co^b, Mg, Wr^a</p> <p>Requirement fulfilled <input checked="" type="checkbox"/> not fulfilled <input type="checkbox"/></p>
Avidity	tested
Microbial status:	This product was filtered through a filter with a pore size of 0.2 µm and subjected to a bacterial load test on Columbia agar, Tryptic Soy Agar and Sabouraud Dextrose Agar with an incubation time of 2 to 5 days, <1 CFU / ml.
Preservatives:	< 0.1 % Sodium-Azide
Storage:	+2 - +8°C
Origin:	Culture supernatant / monitored animal populations / material of human origin
Virus detection:	Materials of human origin have been tested and are negative for HIV 1 + 2, HCV, HBV (HbsAg). Adequate safety precautions are recommended when handling potentially infectious materials.
Note:	All cattle materials used come from veterinary controlled and certified animal herds. Products of bovine origin are assessed as having a low TSE risk.

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Anti-B	
Reference No.:	02110
LOT:	OBM178-5
Clone:	B-6F9
Titer:	1:512
Manufacture date:	2025-01
Expiry date:	2027-01
Purpose of use:	<i>In-vitro-diagnostics</i>
Colour:	Yellow / Clear
PH value:	7,2
Specificity:	<p>Untested antigens circled:</p> <p>A, B, H, Le^a, Le^b, I, K, k, Kp^a, P₁, D, C, E, c, e, C^w, M, N, Bg^a, S, s, Lu^a, Jk^a, Jk^b, Di^a, Fy^a, Fy^b, Xg^a, Do^a, Do^b, Yt^b, Co^b, Mg, Wr^a</p> <p>Requirements fulfilled <input checked="" type="checkbox"/> not fulfilled <input type="checkbox"/></p>
Avidity:	tested
Microbial status:	This product was filtered through a filter with a pore size of 0.2 µm and subjected to a bacterial load test on Columbia agar, Tryptic Soy Agar and Sabouraud Dextrose Agar with an incubation time of 2 to 5 days, <1 CFU / ml.
Preservatives:	< 0.1 % Sodium-Azide
Storage:	+2 - +8°C
Origin:	Culture supernatant / monitored animal populations / material of human origin
Virus detection:	Materials of human origin have been tested and are negative for HIV 1 + 2, HCV, HBV (HbsAg). Adequate safety precautions are recommended when handling potentially infectious materials.
Note:	All cattle materials used come from veterinary controlled and certified animal herds. Products of bovine origin are assessed as having a low TSE risk.

Neckargemünd, 25.02.2025

DNr.:102en 11.2021

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