Klaipėda,

2021 m. *06* mėn. *25* d.

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
Integrated Power Console (IPC™)



Models EC300 and 1898001

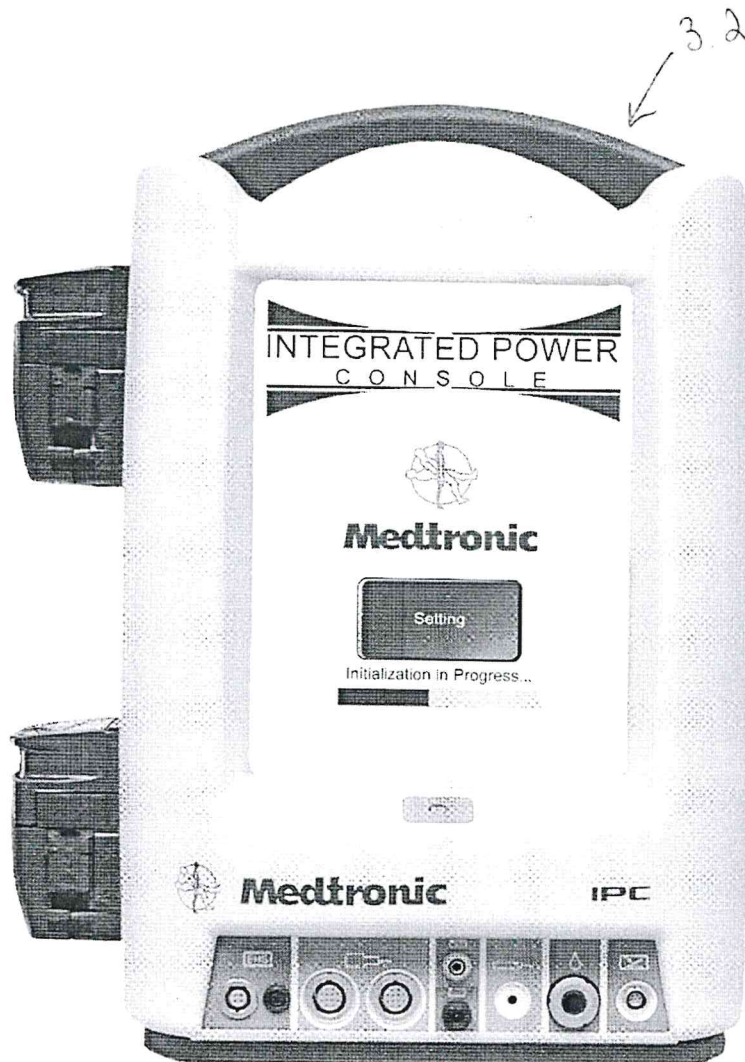Service Manual for Console and Attachments

Rx Only

- A small amount of irrigant is observed flowing at the tip of irrigation device(s).
 - Pump(s) turns off.
10. Confirm system operation.
- Verify:
- Pedal (Coolant) Starts handpiece and coolant flow (coolant pump continues to run for 1 minute after pedal is released).
 - Pedal (Irrigation) Starts and stops the handpiece and irrigation flow (At this step you should also verify that the characters on the SPEED display changed from white to yellow).
 - Pedal Buttons: Please refer to "Multifunction Foot Control Unit".
11. Depress the intraoperative button on the back of the console.
- Verify:
- Starts and stops the handpiece, irrigation and/or coolant flow.
12. Touch Screen
- Verify:
- Speed can be adjusted.
 - Mode can be changed.
 - In oscillate and cut modes check:
 - The "Blade Position" panel opens.
 - The clockwise and counterclockwise buttons move the position indicator and blade in the appropriate direction.

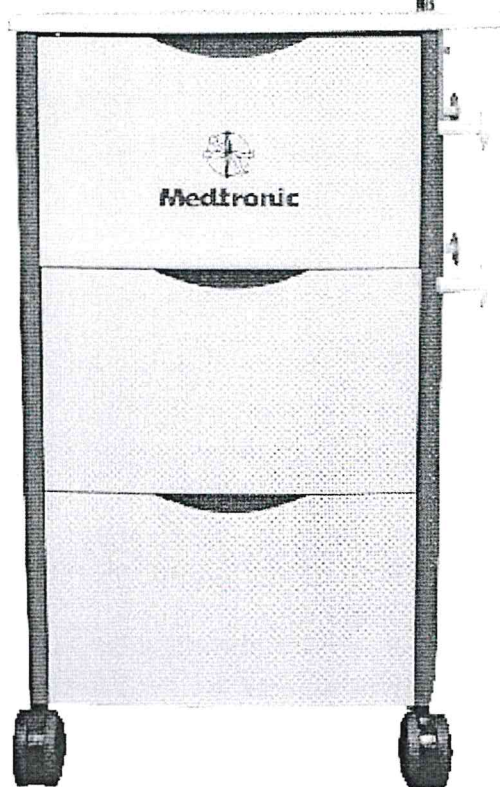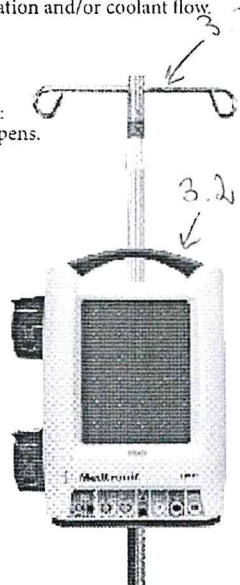- Depressing the 180° button moves the position indicator and blade 180°.
- Flow rate for irrigation is adjustable.

Power Down

1. Turn power switch Off.
2. Disconnect:
 - a. Accessories.
 - b. Suction, irrigation, and coolant tubing.
 - c. Power cord.
3. Discard disposables following health-care facility guidelines on contaminated materials.

NOTE: If any of these condition are different check your set-up, if still incorrect contact Customer Service.

Cleaning

IPC™, Foot Control Unit, and Endo-Scrub® 2 Footswitch

- Do not immerse or sterilize the units.
 - Do not use alcohol, other solvents, or abrasive cleaners.
1. Wipe down the IPC™, Foot Control Unit, and Endo-Scrub® 2 Footswitch with a cloth dampened with a neutral enzymatic detergent, pH 6.0-8.0 or phenol based disinfectant.

Non-Slip Pad ONLY

- 1a. Spray a neutral enzymatic detergent, pH 6.0 – 8.0, or a phenol based disinfectant, mixed to manufactures instructions, directly onto foot pad.
- 1b. Allow the solution to remain in contact with the surface for approximately 10 minutes.
- 1c. Wipe the solution or disinfectant off the foot pedal until visually clean.

2. Dry the units with a clean, non-abrasive cloth.

NOTE: If debris is found under the Foot Control Unit's boot, return for warranty service.

Console Specifications

Functional Standards for Electric Systems		
ANSI / AAMI: - ES 60601-1	Medical electrical equipment -- Part 1: General requirements for basic safety and essential performance	2005
IEC - 60601-1	Medical electrical equipment -- Part 1: General requirements for basic safety and essential performance	2005
EN - 60601-1	Medical electrical equipment -- Part 1: General requirements for basic safety and essential performance (IEC 60601-1:2005)	2006
IEC - 60601-1-4	Medical Electrical Equipment - Part 1: General Requirements for Safety, Part 4: Programmable Electrical Medical Systems	2000
EN - 60601-1-2	Medical Electrical Equipment - Part 1-2: General Requirements for Safety - Collateral Standard: Electromagnetic Compatibility - Requirements and Tests	2001/ A1: 2006
CSA - C22.2 No. 601.1	Medical Electrical Equipment - Part 1: General Requirements for Safety	2005

Physical Dimensions

Size: 277 mm W x 353 mm H x 267 mm D

Weight: 7.3 kg 3.2

Operational Environment

Temperature: +10°C to +33°C
 Humidity: 30% to 75% RH
 Barometric Pressure: 700 - 1060 hPa

Transport and Storage Environment

Temperature: -40°C to +70°C
 Humidity: 10% to 95% RH
 Barometric Pressure: 500 to 1060 hPa

Display / Touch Screen

Type: High contrast, digital, graphic Color, visible in complete darkness.
 Resolution: Display 21 cm diagonal, resolution 480 X 640 pixels

ACCESS /
CLOSURE

THE POWER TO DO MORE

THE INTEGRATED POWER CONSOLE (IPC®) PLATFORM

Inspired by listening to your needs, the IPC® platform delivers maximum torque and ease of use. It's the first multispecialty powered surgery console from Medtronic for cranial, spine, and ENT. Our innovative IPC system drives more handpieces with more power than ever before.

IntelliFlow™
irrigation
remote control

Legend EHS Stylus®
motor delivers
>50% more torque


Color coded
connectors
for ease of use

Helpful
graphics for
each handpiece

Easy-to-use
touchscreen

Connect up to
four handpieces
at once


Console Front

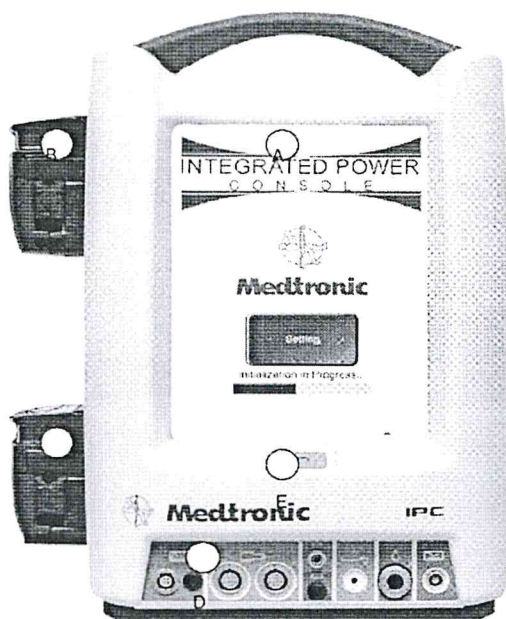- A. Touchscreen – User interface.
- B. Pump 1 – Coolant, lens cleaning, or irrigation.
- C. Pump 2 – Irrigation.
- D. Connector Panel – peripheral devices.
- E. Power Switch – System On/Off switch.

NOTE: 2 micromotors can operate independently of each other. 2

Connector Panel

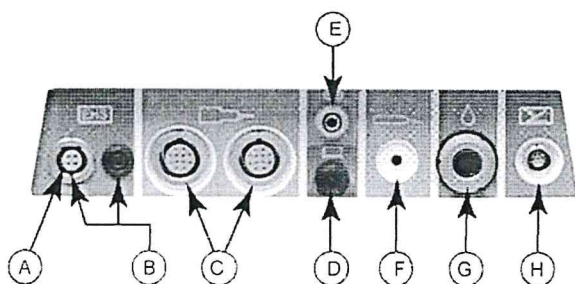Port #	Component	Quantity
A	Midas Rex® Legend EHS® motor	1
B	Midas Rex® Legend EHS® Stylus motor	1
C	StraightShot® M4 Microdebrider	2
	Midas Rex® Legend EHS® Stylus Touch™ motor	
	Midas Rex® SC1	
	StraightShot® Magnum® II and StraightShot® III Visao®	
D	Stimulus input from Patient Interface connection (NIM).	1
E	Stimulus output to STIM Bur Guard	1
F	Skeeter® Handpiece	1
G	Endo-Scrub® 2 Finger Switch	1
	Endo-Scrub® 2 Footpedal	
	IntelliFlow Irrigation Remote Control	
H	Foot Control Unit (FCU)	1

Connector Panel Cable Connection

Cable to console connection red/silver dot:


- Red or silver dot connections are multi pin and must be correctly aligned (oriented).

Cable to console connection without dot:

- Connectors without the red or silver dot are single pin and may be inserted without regard to orientation.


Connector Panel Cable Disconnection (multi pin) To Remove Midas Rex® :

To Remove Midas Rex® Legend EHS® Motor and Legend EHS Stylus® Motor, Cable from motor or console:

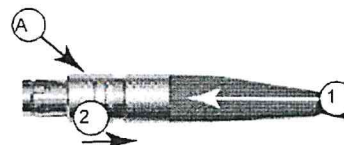
1. Push the cable towards the motor or console.
2. Then pull out by locking ring (A).

To Remove Midas Rex® Legend EHS Stylus® Cable from console:

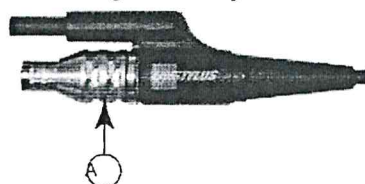Push the cable towards the console, then pull by locking ring (A)

To Remove cables (multi pin) with polymer insulating boots:

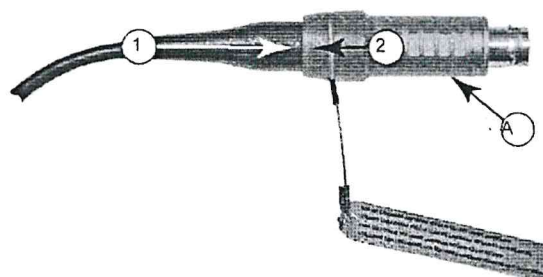
NOTE: Confirm handpieces contain polymer insulating boot (A). If handpiece contain missing or cracked polymer boots, contact Medtronic Customer Care for upgrade.

NOTE: If units with polymer insulating boots have debris under the insulator:

- Reclean according to Cleaning and Sterilization instructions.
- If debris was not removed return for warranty servicing. See warning W4.

3. Push the cable towards the console.
4. Then pull out by the polymer insulating boot (A).

To Remove cables (multi pin) with silicone insulating boots:

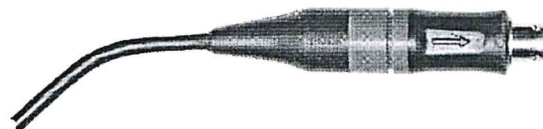Silicone insulated connectors do not have a locking device (ring) and may be removed by pulling straight out on the connector.

Cable Disconnection (single pin)

Single pin connectors do not have a locking device (ring) and may be removed by pulling straight out on the connector.

System Description

List used in Manual

This manual contains two basic types of list: Alphabetic list (A, B, C, etc.) and Numeric list (1, 2, 3). The alphabetic list contains general information about the part such as name or usages. Numeric list are instructions for completing a task.

Definitions

This section explains the essential meaning of a word or acronym as used in this manual. Also explains changes in words or phrases variations from one product generation to the next.

FCU -	Foot Control Unit
IPC -	Integrated Power Console
I.V. -	Intravenous
NIM -	Nerve Integrity Monitor - One or all of the following units: NIM-Response® 2.0, NIM-Neuro® 2.0, NIM-Response® 3.0, and NIM-Neuro® 3.0

Nomenclature	The act or process or an instance of naming
XPS -	Xomed Power System

When The System Arrives

Unpacking and Inspection

Check off the contents of the box against packing slip. If incomplete or damaged, notify Customer Care.

If container is damaged, or cushioning material shows stress, notify carrier and Customer Care. Keep shipping materials for carrier inspection.

After unpacking, save the cartons and packing material. If the instrument is to be shipped the shipping package will provide proper protection.

System Description

The IPC® System is a powered microdebrider, drill and saw system that will remove soft tissue, hard tissue, and bone during surgical procedures. The system consists of a power control console, footswitch, connection cables, and assorted handpieces to drive various burs, blades, drills, rasps, cannulae, and saws. It includes integrated irrigation pumps for irrigation of blades, burs and for motor coolant.

The Nerve Integrity Monitor (NIM®) is a separate device that stimulates and monitors the nerve. This system has connections that allow the NIM® to be connected with the Visao® handpiece and Stimulating Bur Guard enabling the NIM® to stimulate and monitor the nerve at the surgical site.

The system can be used to clear the end of a rigid rod endoscope in order to maintain good visualization of endoscopic procedures without having to remove the scope from the surgical site.

This device is intended for use by physicians trained in the procedures described.

Intended Use / Indications for use

The IPC® is indicated for the incision / cutting, removal, drilling, and sawing of soft and hard tissue and bone, and biomaterials in Neurosurgical (Cranial, Craniofacial, Orthopedic, Arthroscopic, Spinal, Stereotaxy, and General surgical procedures).

Contraindications

The IPC® system is contraindicated for arthroscopic microdissection in individuals with the following:

- Severe/progressive neurological deficits
- Cauda equine syndrome
- Active infection.

Arthroscopic microdissection is not indicated for individuals with sequestered disc fragments, discogenic pain, internal disc destruction, or lumbago.

Sales And Customer Care

Medtronic Powered Surgical Solutions

4620 North Beach Street Fort Worth, TX 76137 USA

www.medtronic.com

U.S. Help Line


(800) 468-9710

International Service


International customers should contact their Medtronic Neurology Technologies representative.

Console

Console Front

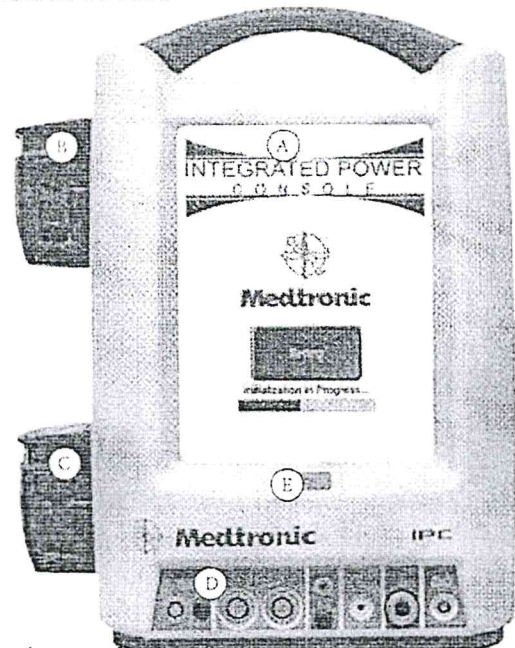- A. Touchscreen - User interface.
- B. Pump 1 - For coolant, lens cleaning, or irrigation.
- C. Pump 2 - For Irrigation.
- D. Connector Panel - peripheral devices.
- E. Power Switch - Turns system on or off.

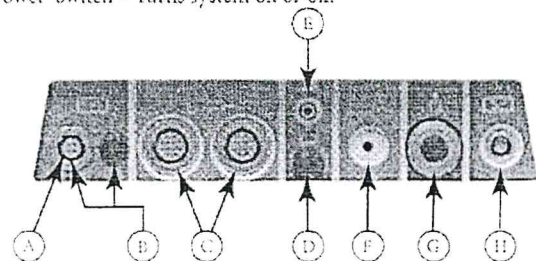Connector Panel

Port #	Component	Quantity
A	Midas Rex® Legend EHS® motor	1
B	Midas Rex® Legend EHS® Stylus motor	1
C	StraightShot® M4 Microdebrider	2
	Midas Rex® Legend EHS® Stylus Touch™ motor	1
	Midas Rex® SC1	
	Straightshot® Magnum® II and Straightshot® III	
	Visao®	
D	Stimulus input from Patient Interface connection (NIM).	1
E	Stimulus output to STIM Bur Guard	1
F	Skeeter® Handpiece	1
G	Endo-Scrub® 2 Finger Switch	1
	Endo-Scrub® 2 Footpedal	
	IntelliFlow Irrigation Remote Control	
H	Foot Control Unit (FCU)	1

Set-Up

General instructions for set-up, inspection, and use of the Integrated Power Console.

The IPC™ is designed to require no regular maintenance or service. The following recommended sustaining care may be performed annually or at more frequent intervals to extend the life of the IPC™.

NOTE: Use sterile water or saline for irrigation and cooling.


Inspect:

- Visually inspect:
 - the entire console for signs of cracks or other damage.
 - the front and back connector panel for damaged or loose connectors
 - air inlet and exhaust on the bottom and rear of console is clean and free of debris.
- Check that the handle is securely attached to console and slowly rotate the IPC™ to an inverted position while listening for loose material inside the enclosure.
- Ensure pole clamp is securely attached to console and knob spins freely.

Configure System:

- On Irrigation pole, mount IPC™ and irrigation/coolant bag(s).
Note: Irrigant and coolant bags should be placed above the console to ensure adequate flow.
- Position the IPC™ in a manner that does not obstruct the power inlet for the purpose of disconnecting the Mains voltage by the power cord. Plug unit into power source.
- Connect FCU.
- Connect the accessories to console (test all accessories available one at a time).
- Irrigation/Coolant Pumps
 - Connect tubing as needed (suction, cooling, irrigation).

Functional Test Splash Screen

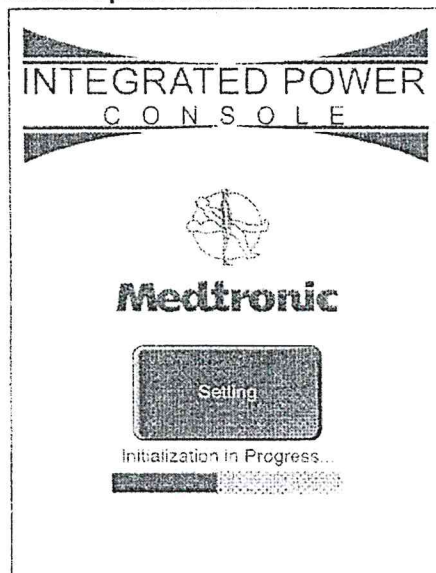
- Turn power switch ON and verify:
 - the power switch operates smoothly
 - switch lamp illuminates
 - fan on bottom of device is turning by listening for fan noise.
 - The Splash Screen is displayed while the system is starting up and executing its self tests
 - System passes self test
 - Default screen opens.

Main Screen

During the boot up/self-test operation the IPC™ will identify attachments such as handpieces (foot pedal etc). Actual screen displayed is dependent on attachment found.

The following generalized screen is meant to familiarize the technician with adjustments available.

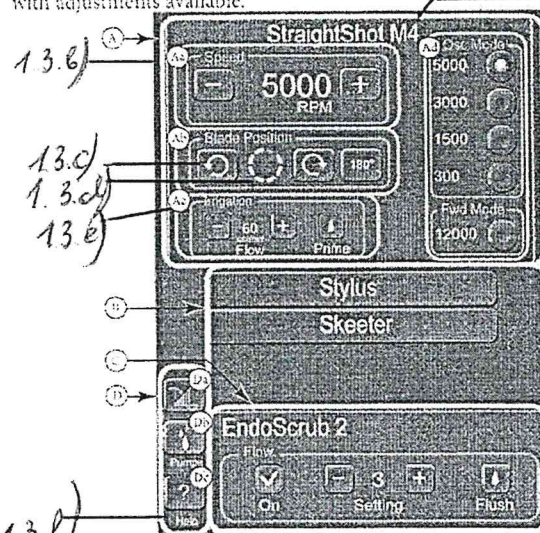
- Handpiece Panel - This panel shows the active handpiece.
 - Speed - Variable adjustment on most handpieces. Default value is handpiece specific.
 - Blade Position - Handpiece specific, allows the user to rotate the inner cutting tip of specially designed rotatable blades.
 - Irrigation - Adjust the flow rate of irrigation. Default value is handpiece specific.
 - Mode - Handpiece specific cutting mode selection.
- Augment Area - Shows supplemental information such as inactive handpiece(s), special function panel, pump panel etc.
- Special Function Panel - Shows Suction Irrigator or Endo-Scrub* 2 panel.
- Main Screen subsection
 - Foot Control Unit (FCU) Button - changes foot pedal from variable speed control to On/Off.
 - Pumps - Opens pump panel.
 - Help - Opens help screens.

Handpiece Touchscreen Default Value Table		RPMs		Mode		Pump 1		Pump 2
Handpiece Name	Touchscreen Display Name	RPM Range	Default RPMs	Fwd (Forward)	Rev (Reverse)	On (Oscillate)	Cut (1 revolution)	Cooling
Visio	Visio	200-80000	80000	X				X
Suction Irrigator Optional					X			X
Midas Rex SCI	SCI Handpiece	50-3000	5000		X			
		50-12000	12000	X				X
StraightShot M4	StraightShot M4	50-5000	5000		X			X
Microdebrider		50-12000	12000	X				
EndoScrub 2 Optional								X
StraightShot III	StraightShot M4	50-5000	5000		X			X
Microdebrider		50-12000	12000	X				
EndoScrub 2 Optional								X
StraightShot Magnum II	StraightShot M4	50-3000	5000		X			
Microdebrider		50-12000	12000	X				X
EndoScrub 2 Optional								X
Midas Rex Legend EHS* motor	EHS	200-75000	75000	X	X			X
Suction Irrigator Optional								X
Midas Rex Legend EHS Stylus motor	Stylus	200-75000	60000	X	X			X
Suction Irrigator Optional								X
Skeeter Handpiece	Skeeter	1000-16000	16000	X	X			
Suction Irrigator Optional								X


System figures

Figure 1-1. IPC™ Console Front

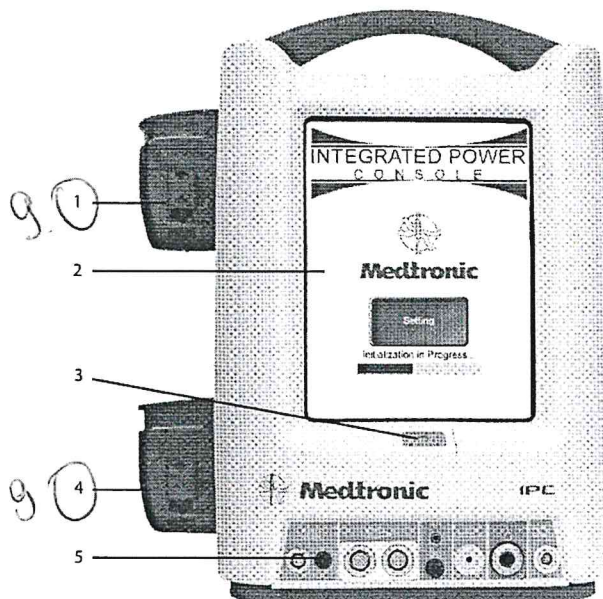
- 1 Pump 1: Coolant, irrigation
- 2 Touchscreen
- 3 Power on/off
- 4 Pump 2: Irrigation
- 5 Console connector panel for peripheral devices

Figure 1-2. IPC™ Console Back

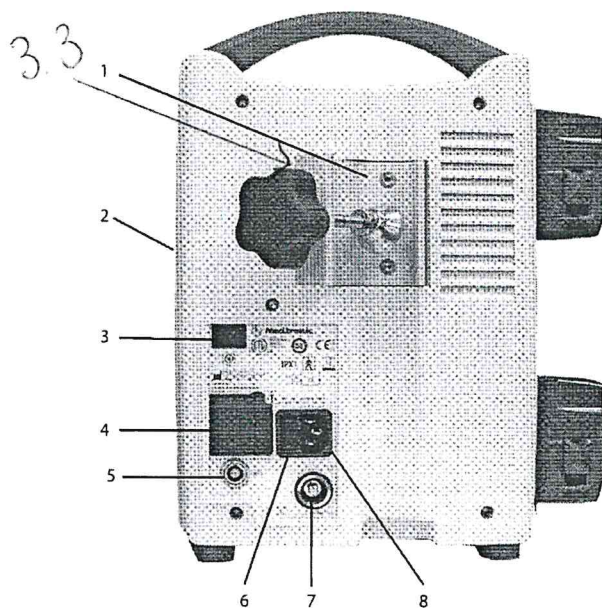
- 1 Pole clamp
- 2 Compact flash card port (Medtronic Use)
- 3 Manual start/stop
- 4 Auxiliary power outlet
- 5 Endo-Scrub 2 connector
- 6 Fuse access
- 7 Equipotential Ground Connector.
Apply potential equalization conductor.
- 8 Hospital grade power cord connector

Figure 1-3. IPC™ Console Connector Panel

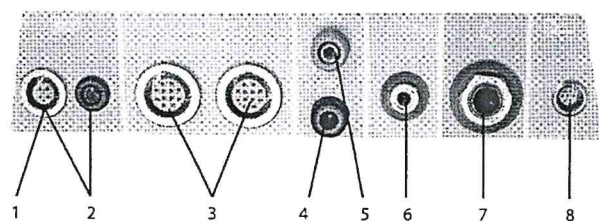- 1 Legend EHS™
- 2 Legend EHS™ Stylus, Stealth-Midas™
- 3 Midas Rex™ Stylus Touch™ motor, Microsaws, High-torque handpiece, and MR8 motors
- 4 Stimulus input from patient interface (NIM™ or NIM-Eclipse™)
- 5 Stimulus output to stim bur guard
- 6 Not applicable
- 7 IntelliFlow™ irrigation remote control
- 8 Multifunction foot pedal


Figure 1-4. Multifunction Foot Pedal & Y-Splitter

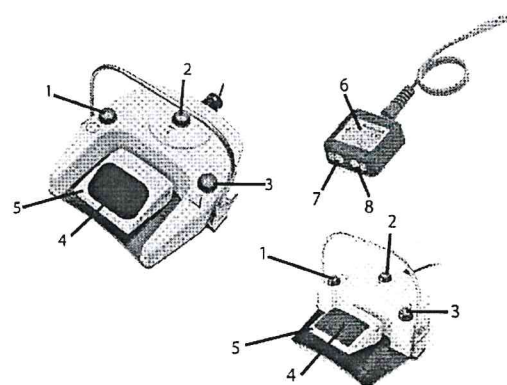
- 1 Mode button
- 2 Handpiece button
- 3 Control button
- 4 Slip-resistant foot pad
- 5 Foot pedal
- 6 Y-Splitter
- 7 Port 1
- 8 Port 2

Visao® Coolant Pump Set-Up


The Pump Cartridge snaps onto the lower section of pump # 1.

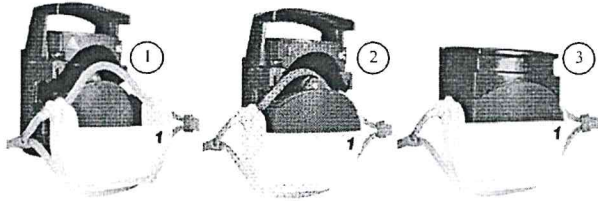Standard Pump Set-up

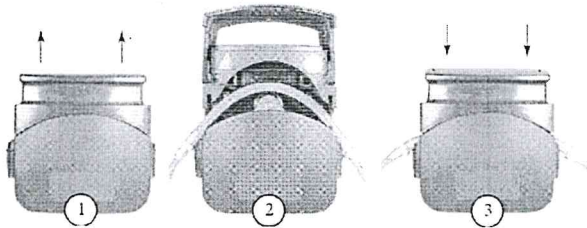
Tips on loading the pump

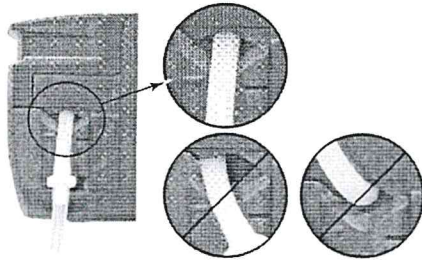Accessories/attachments

Multifunction Foot Control Unit (FCU)

Part No. 1898430 or EF200

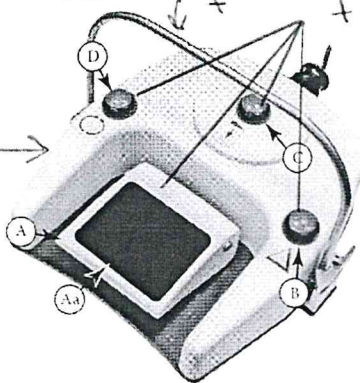Buttons and Pedal

NOTE: Each button must be depressed and held for a definable amount of time (100 mS by default).

Drills

- A. Foot Pedal - Start/Stop, Variable speed.
- Aa. Non-Slip Foot Pad.
- B. Right Button - Pedal function, (Start-Stop or Variable speed).
- C. Top Button - Active handpiece selection
- D. Left Button - Mode selection, (FWD/REV).

Microdebrider

- A. Pedal: Start/stop, variable speed (start/stop, or variable speed selectable via FCU button on Main Screen).
- B. Right Button: In OSC Mode this button will rotate inner blade 60°/180° (touchscreen selected). In FWD Mode this button will select Pedal function (Start/stop, or variable speed).
- C. Top Button: Active handpiece selection.
- D. Left Button: Mode/RPM selection -


SC1

- A. Pedal: Start/stop, variable speed (start/stop, or variable speed selectable via FCU button on main screen).
- B. Right Button: If mode is set to OSC this button will, rotate inner tube on blades 180°. If mode is set to FWD this button will, select pedal function (Start/stop, or variable speed).
- C. Top Button: Active handpiece selection.
- D. Left Button: Mode selection -FWD/OSC

NOTE: If any of these condition are different check your set-up, if still incorrect contact Customer Service.

Intelliflow Irrigation Remote Control

- A. Pause/On-Off:
 - Pause if used with handpiece irrigation (Flow rate will flash yellow).
 - On-Off/Pause if used with Suction Irrigator.
- B. Increase/Decrease:
 - Handpiece Irrigation - fine adjustment for irrigation rate.
 - Suction Irrigator - fine adjustment for irrigation rate.
- C. Increase/Decrease:
 - Handpiece Irrigation - coarse adjustment for irrigation rate.
 - Suction Irrigator - selects stainless steel (Fr) tubing size.

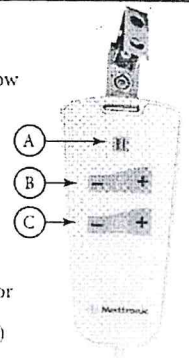
Endo-Scrub® 2

NOTE: Can be used only with a microdebrider.


The IPC™ System incorporates Endo-Scrub® 2 functionality by using irrigation pump number one (1) and controlling operation with the touch screen and an external footswitch or finger switch.

It is not to be used for infusion, for disinfection or sterilization of an endoscope, or for suction removal of blood and debris.

NOTE: Use the Endo-Scrub® 2 sheath only with an endoscope listed on the sheath product label, as malfunction or poor performance could result.

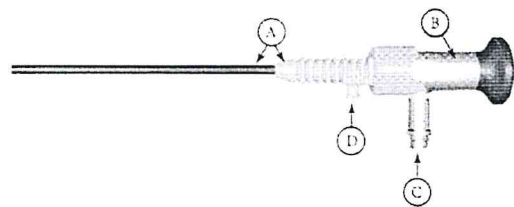- A. Endo-Scrub® 2 Sheath.
- B. Endoscope.
- C. Light source connection.
- D. Irrigation connection.

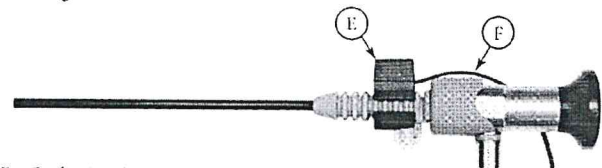- E. Endo-Scrub® 2 Finger switch.
- F. Finger switch cable.

YOUR
COMPLETE
SOLUTION.
ONEPARTNER.


In your daily practise you constantly strive for improved clinical outcome and optimized workflow.

Deformity Solutions
Degenerative Solutions
Fracture Solutions

ACCURACY. EFFICIENCY. ADDED CONFIDENCE.

Console Rear

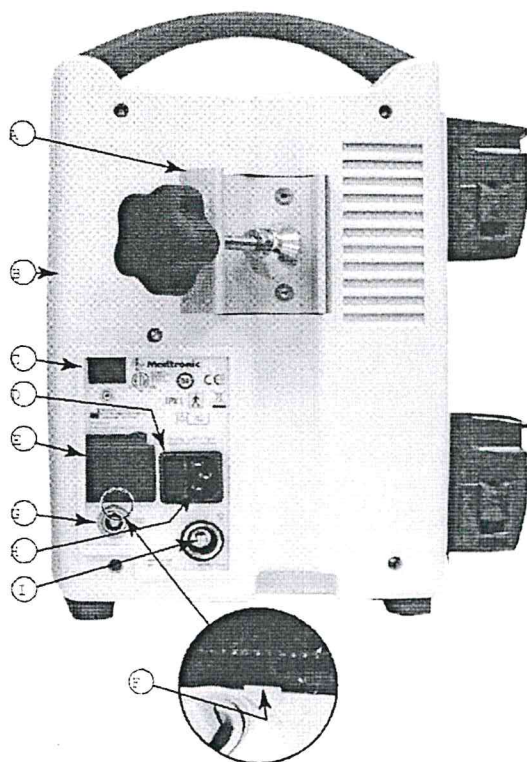
- E. Pole clamp.
- F. Compact flash card port (factory use only).
- G. Manual Start/stop button.
- H. Fuse Access – Replace only with 5 x 20 T. L. 5A, 250 V fuse.
- I. Auxiliary power outlet with protective cover:
 - For use at grid voltage < 120 VAC only.
 - HydroDebrider™, or Bone Mill consoles only. See warning W5.
- F. To remove cover, place small screwdriver in notch at bottom and pull/pry off.
- J. Endo-Scrub® 2 - power connector.
- K. Power cord connector: See appendix B for part numbers.
 - Hospital grade power cord connects here.
 - Means of disconnecting device from Mains voltage by the power cord.
- Equipotential:
 - Uniform potential.
 - Means for eliminating noise or interference with sensitive equipment by application of a POTENTIAL EQUALIZATION CONDUCTOR.

Power Cords

North America: USA, Barbados, Belize, Bolivia, Canada, Columbia, Ecuador, Venezuela Standard P/N EA600 or 1895820 6 meter P/N EA650 or 189721	United Kingdom, Ireland, Hong Kong, Malaysia, Singapore P/N EA606 or 1895821	Continental Europe: Austria, Belgium, Finland, France, Germany, Greece, Korea, Luxembourg, Netherlands, Norway, Portugal, Spain, Sweden P/N EA602 or 1895822
China P/N EA604	India, South Africa P/N EA607	Switzerland P/N EA601
Argentina P/N EA608	Israel P/N EA609	Denmark P/N EA610
Australia, New Zealand P/N EA605	Japan P/N EA603 or 1895823	Italy, Chile P/N EA611

Console Pump Designator

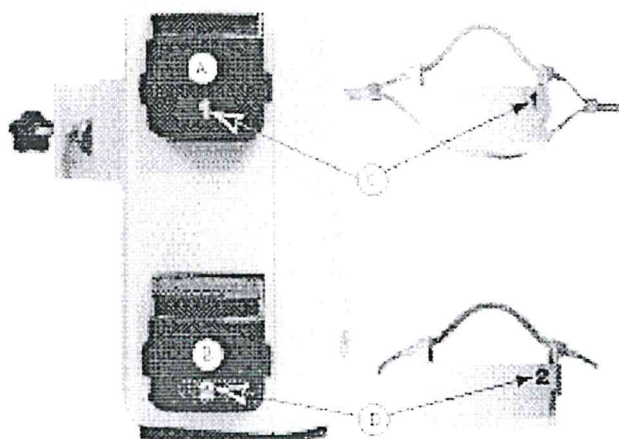
- A. Pump 1: Coolant, lens cleaning, or irrigation.
- B. Pump 2: Irrigation or lens cleaning.
- C. Pump 1 Designator – This designator number is used to coordinate the pump (by number) with the cartridge number and/or pump set-up screen number listed on the touch screen. When setting up the console these **numbers must match**.
- D. Pump 2 Designator.

NOTE: Not all Pump Cartridges have pump designator numbers. For these cartridges the operator should view the Pump Setup Screen prior to installing the cartridge. Irrigation speed 0.5-100 ml/min.


Irrigation/Coolant Pumps

Pump Cartridge Set-up

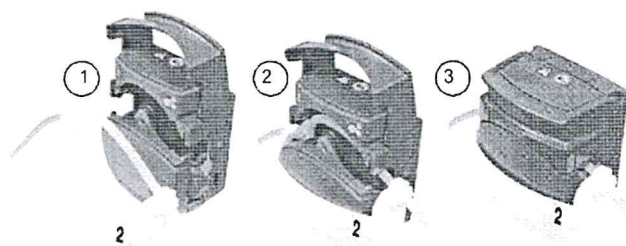Visao® Pump Cartridge


The Visao® Pump Cartridge has both a pump tube and a return tube.

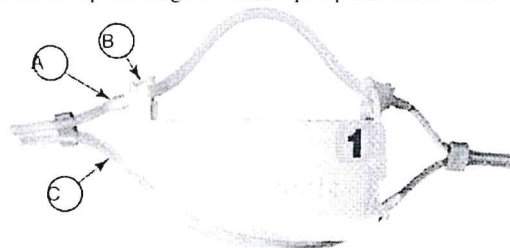- L. Pump tubing.
- M. Pump tubing is clipped into the Pump Cartridge.
- N. Return tube.

Legend EHS® Motor Cable


Connects the motor to the console.

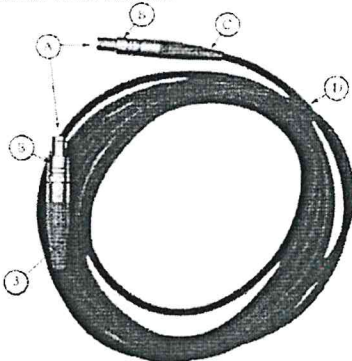- A. 4-pin connector.
- B. Locking sleeve.
- C. Green boot.
- D. Cable.

Motor Collet

Prior to installing an attachment, ensure that arrows on the motor collet are in proper alignment.

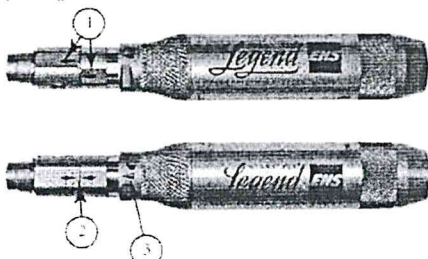1. Improperly aligned collets.
2. Properly aligned collets.
3. Motor side attachment alignment arrow.

If the arrows are not aligned, use the Motor Wrench to turn the rotational collet until its arrow is aligned with the arrow on the stationary collet.

Technical Specifications


Part No.	EM100-A
Speed	200-75000 RPM forward/reverse
Size	9.02 cm length x 2.03 cm diameter
Weight	180 g

Duty Cycle (To avoid overheating):

- For continuous use in operating room temperatures up to 40°C, the Legend EHS® Motor is rated for a cutting time of 3 minutes at 70,000 RPM.
- For normal operating room temperatures (typically 20°C) the Legend EHS® Motor is rated for a continuous cutting time of 10 minutes followed by 25 minutes of rest.
- The Legend EHS® Motor is rated for intermittent use of 20 seconds ON / 20 seconds OFF, indefinitely at 70,000 RPM.

Midas Rex® Legend EHS Stylus® Motor

A smaller compact high speed, high torque, reversible electric motor used to dissect bone and biomaterials at selectable speeds from 200 to 75,000 RPM. The Midas Rex® Legend EHS Stylus® Motor cable is integral with the Handpiece and is not removable from the motor.

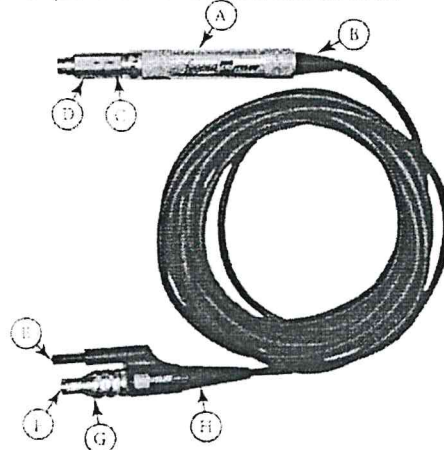- A. Midas Rex® Legend EHS Stylus® Motor.
- B. Cable.
- C. Rotational collet.
- D. Stationary collet.
- E. Ground connector.
- F. 4-pin connect.
- G. Locking sleeve.
- H. Black boot.

Technical Specifications

Legend EHS Stylus® Motor

Part No.	EM200
Speed	200-75000 RPM forward/reverse
Size	7.77 cm length x 1.65 cm diameter
Weight	90 g

Duty Cycle (To avoid overheating):

- For continuous use in operating room temperatures up to 40°C, the Legend Stylus® Motor is rated for 3 minutes at 60,000 RPM, followed by 25 minutes of rest.
- For normal operating room temperatures (typically 20°C) the Legend Stylus® Motor is rated for continuous cutting indefinitely at 60,000 RPM.

HIGH SPEED ELECTRIC MOTOR SPECS

		Legend Stylus	Legend EHS
Console		IPC/Legend EHS	IPC/Legend EHS
13.1	Max Speed (rpm)	75,000	75,000
13.2	Max Power (W)	160	146
Max Torque (Mn-m)		42	38
13.6	Noise (Db)	58	58
Weight (oz)		3.17	6.35
Length (in)		3.06	3.55
Diameter (in)		0.65	0.8
Foot/Finger		Both	Foot only
Special Features		Trans Nasal Skullbase burs	

LOW PROFILE ON DRILL IRRIGATION

Low-Profile Irrigation Tubing improves surgical visibility and is designed to connect easily and securely with small-bore attachment

- Thin fluid delivery tubing minimizes visual obstructions
- Quick clip design provides seamless line management within the surgical field
- Translucent clip design minimizes glare in surgical field
- Malleable tip to direct flow with changes in tool sizes

Newest upgrades of IPC® System have improved irrigation stop and start times and added continuous flow options at lower rate

Medtronic

ACCESS /
CLOSURE

TRANS-NASAL SKULLS BASE BURS

Trans-Nasal Skull Base Burs
allows you to perform endoscopic
approaches to anterior skull base.
These burs enhance visualization
with endoscope and can help
protect anatomy from a spinning
bur shaft.

TN30MFL TN30MCD TN40RFL


CLEARVIEW™ TOOLS

The Midas Rex ClearView™ line
of tools have been designed,
developed, verified and validated
specifically for spine surgery.


These burs allows better visibility,
have a protected shaft, various
shapes and lengths and integrated
irrigation.

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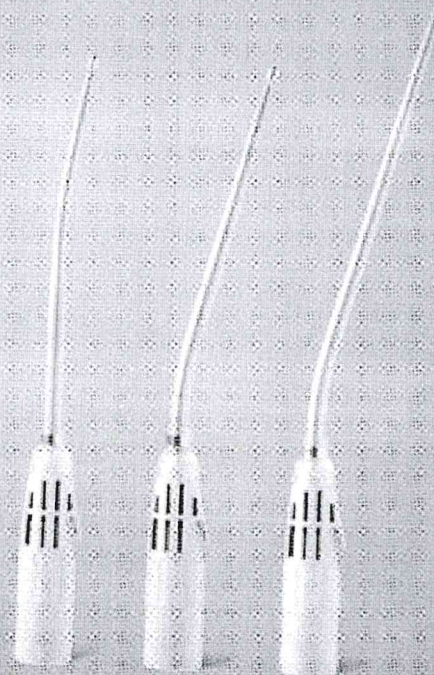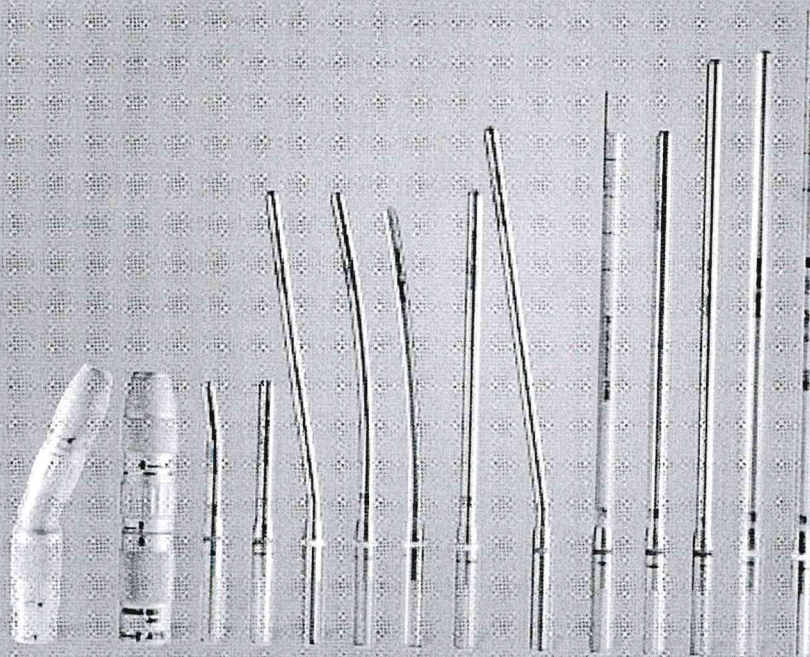
TELESCOPING TUBES

The Legend[®] system telescoping tubes offer a range of straight, angled, and curved tubes that improve operative site visualization and help protect surrounding soft tissues.

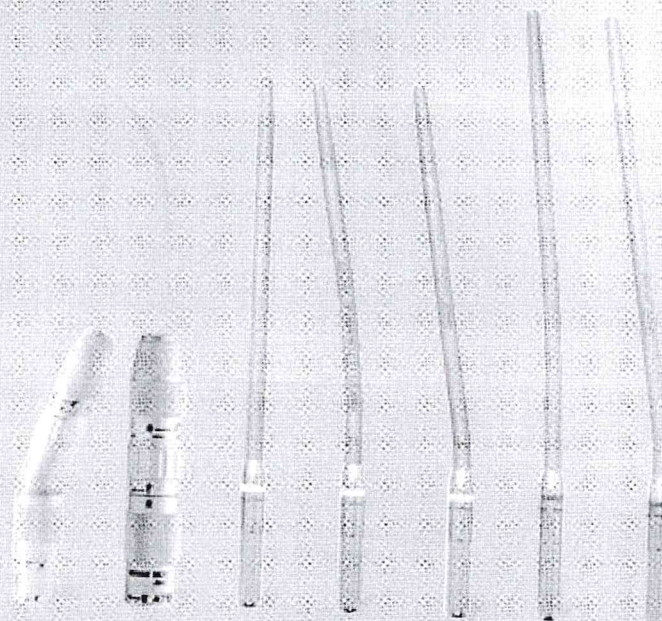
AT10A AT10 TT12AM5 TT12CM5 TT12AM15 TT12CM15 TT14AM15

TAPERED TELESCOPING TUBES

¹³ Variable exposure attachments allow you to vary dissecting tool exposure for optimal visibility. There is no need to touch the end of the dissecting tool, and they are adjustable up to 12 mm in 1 mm increments. Compatible with Midas Rex[®] Legend[®] tools and motors.

13

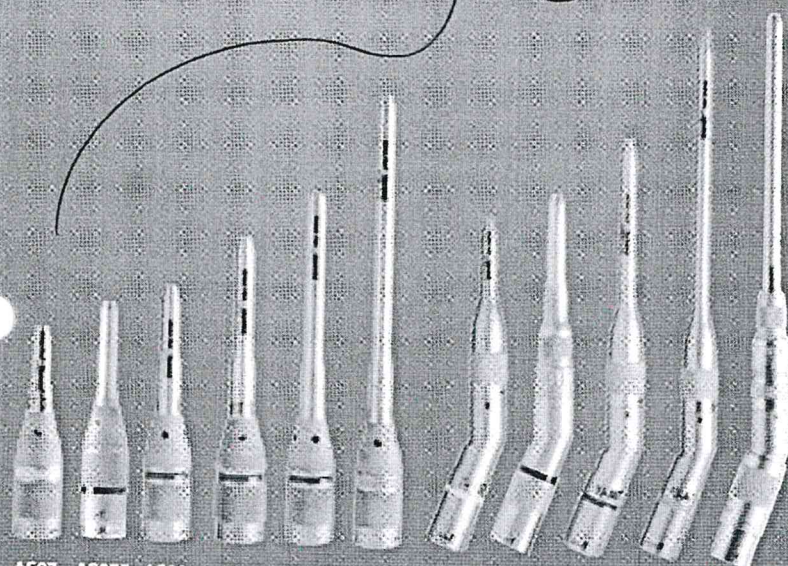
AS07 AS075 AS08 AS10 AS10S* AS12S AS15 AA07 AA075 AA10 AA10S* AA15 AVA15DK

SMALL BORE ATTACHMENTS

Accommodates Legend® 2.4 mm diameter tools

These attachments are used in various procedures for better visibility in small surgical sites. All Legend® attachments and tools utilize a patented 3-step locking process for increased ease of use and efficiency in the OR.

*The Bore has increased from 2.4 mm to 2.5 mm

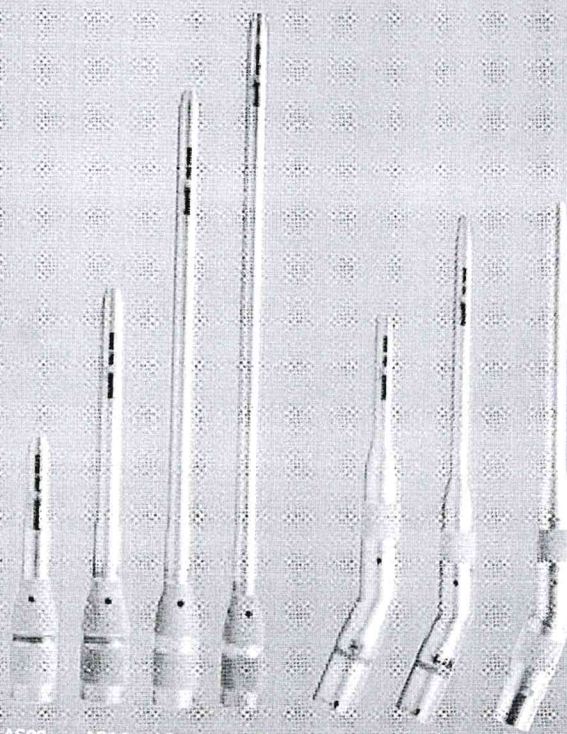AS09 AS09S* AS14 AS14S* AS21 AS26 AA09S AA14 AA14S* AA14DK

LARGE BORE ATTACHMENTS

Accommodates Legend® 3.2 mm diameter tools

These attachments contain the most complete line of tool shapes and sizes in the Legend line. The large bore attachments can accommodate a large number of clinical specialties. The Legend® System large bore straight and angled attachments are designed to be used with the larger 3.2 mm shaft tools.

*The Bore has increased from 3.2 mm to 3.5 mm

16

SURGICAL TECHNIQUE

THE EASYDRILL™ CRANIAL PERFORATORS

The EasyDrill Cranial Perforator device is a bone cutting and drilling instrument used in conjunction with a surgical motor and a Hudson Chuck - speed reducer attachment, to drill access holes through a patient's skull.

When properly used, the EasyDrill Cranial Perforator employs a clutch mechanism to automatically disengage once perforation is accomplished and as the drill ceases to find resistance to bone.

EQUIPMENT NEEDED

1. EasyDrill Cranial Perforator™

2. Attachment Driver (AD03)

3. Drill

EASYDRILL CRANIAL PERFORATOR FUNCTIONAL TEST


Perform the following function test before each trepanation. Use sterile disposable gloves to perform the test.

CAUTION: If the sterile protective cover is not available, use sterile gauze (a protective barrier) when handling the cutting end of the EasyDrill Cranial Perforator.


NOTE: Do not use the product and contact customer service if either of the two function tests below fail to present the indicated results.

RECOMMENDED SPEED RANGE:

We recommend a speed range of 800 to 1000 RPM for the drilling with EasyDrill Cranial Perforators.

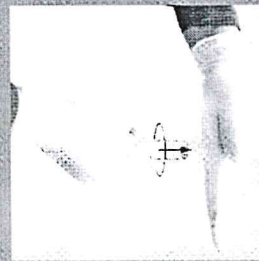
1. Firmly hold the EasyDrill Cranial Perforator's plastic drive shank (Hudson shank) and the protective cover, turning the cutting end of the EasyDrill Cranial Perforator device in a clockwise direction. The cutting end of the device should turn freely.

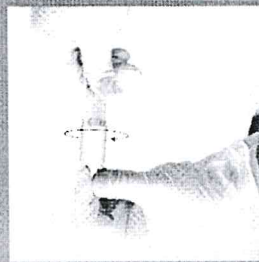
2. Next, grasp the EasyDrill Cranial Perforator by the plastic drive shank (Hudson shank) and press the protected cutting end of the device into the opposite hand; turn the drill clockwise until it locks to apply pressure.

CRANIAL BONE DRILLING

1. After the functional tests, attach the EasyDrill Cranial Perforator to the Hudson-type chuck of your motorized surgical motor. Verify that the EasyDrill Cranial Perforator is firmly fixed to the Hudson-type chuck before proceeding.
2. Remove the protection cover from the cutting end of the EasyDrill Cranial Perforator device before drilling.
3. Position and press the cutting end of the EasyDrill Cranial Perforator gently against the skull to be perforated.
4. Manually turn the EasyDrill Cranial Perforator to your desired drill orientation.
5. Position the EasyDrill Cranial Perforator perpendicular to the skull (perpendicular to the point of perforation). Supply power to the drilling system and begin trepanation.
6. Maintain a gentle and continuous pressure while drilling.
7. Stop the point of perforation with physiological saline while drilling to prevent bone necrosis.
8. As the drilling depth within the EasyDrill Cranial Perforator device ceases to find resistance to bone, it automatically disengages the cutting end of the device, which stops rotating. Note: At this stage, the plastic handle of the device may appear to continue to rotate.
9. Carefully remove the EasyDrill Cranial Perforator from the skull.
10. The entire skull's surface that is in the hole that must be carefully removed with appropriate surgical care. If the layer is removed to the brain, do not take appropriate care not to injure the brain.
11. Before performing another perforation, after the cutting end of the EasyDrill Cranial Perforator removes the red dust on the previous drilling. Perform testing procedures again.

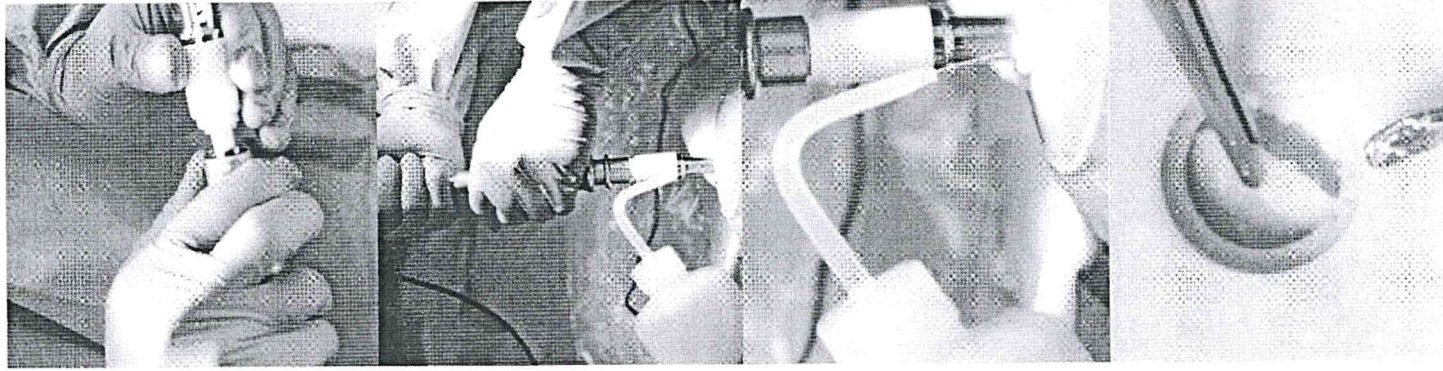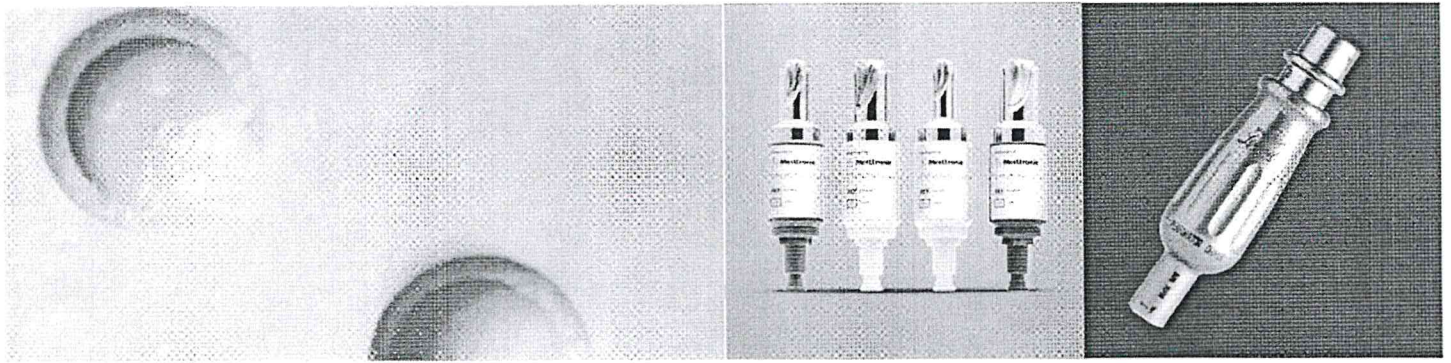SURGICAL TECHNIQUE

OTHER TIPS

The drill can cut or tear the dura mater when it unlocks. Check some conditions before performing the procedure, such as the existence of adhering dura mater, high intracranial pressure, or other anomalies adjacent to the drilling site.

- Do not use the drill in a hole or an area that has already been drilled, or the surrounding areas. The composition of reconstructed bone differs from that of healthy bone, increasing the risk of unwanted breakage. Leave a minimum spacing of 15 mm from one hole to another.

- It is essential to keep the drill perpendicular (90°) at the predetermined point of the skull to be drilled, as an excessive deviation from perpendicularity may cause the product to fail and lead to serious patient injury. The perforator can nick or tear the dura or brain if not held at a 90° angle to the skull.
- Select a drill bit suitable for the bone thickness according to the table below. This prevents the drill from tearing the bone or brain tissue (similar effect when the bit is not positioned at 90°).
- Take extra care with the pressure applied when drilling the skull of infants, children, and the elderly or pathological bone.

ORDERING INFORMATION

Item #	Colour	Outer Diam.	Inner Diam.	Shelf	Minimum Bone Thickness
DM0008FAA		9.0 mm	6.0 mm	6.0 mm	3.25 mm
DM0011FAA		11.0 mm	7.0 mm	3.0 mm	3.25 mm
DM0210FAA		14.0 mm	11.0 mm	1.5 mm	1.75 mm
DM0010FAA		14.0 mm	11.0 mm	3.0 mm	3.25 mm
Item #	Description				
*AD05	Midas Rex® Perforator Driver Attachment Manufactured by Medtronic				

Medtronic

Europe
Contact your local
representative
for more information
or to place an order.
www.Medtronic.com

United Kingdom/Ireland
Contact your local
representative
for more information
or to place an order.
www.Medtronic.com

Manufactured by
MICROMAR
Distributed by Medtronic

Read only. Refer to product instruction manual/
package insert for instructions, warnings,
precautions and contraindications.

For more information, contact your Medtronic
Neurosurgery sales representative or refer to
www.MedtronicNeurosurgery.com.

For a listing of indications, contraindications,
precautions, warnings, and potential adverse
events, please refer to the Instructions for Use.

www.medtronic.eu

MIDAS REX® Legend® System

Specialty Attachments

Product	Description
AR14	Right-Angle Attachment
AC16	Contra-Angle Attachment (ME)
ASMC	Metal Cutter Attachment
AD01	Perforator Driver Attachment — 800 rpm
AD03	Perforator Driver Attachment — 1000 rpm
AD02	Jacobs Attachment (5/32") with key
ASJL	J-Latch Attachment

AR14

AC16

ASMC

AD01

AD03

AD02

ASJL

SymmeTRI Dissecting Tools

Product	Description	Head Dia. (mm)
75BA30T	SymmeTRI 7.5 cm, 3 mm	3.0
75BA40T	SymmeTRI 7.5 cm, 4 mm	4.0
75BA50T	SymmeTRI 7.5 cm, 5 mm	5.0
75BA60T	SymmeTRI 7.5 cm, 6 mm	6.0
75BA75T	SymmeTRI 7.5 cm, 7.5 mm	7.5
9BA30T	SymmeTRI 9 cm, 3 mm	3.0
9BA40T	SymmeTRI 9 cm, 4 mm	4.0
9BA50T	SymmeTRI 9 cm, 5 mm	5.0
9BA60T	SymmeTRI 9 cm, 6 mm	6.0
10BA30T	SymmeTRI 10 cm, 3 mm	3.0
10BA40T	SymmeTRI 10 cm, 4 mm	4.0
12BA30T	SymmeTRI 12 cm, 3 mm	3.0
12BA40T	SymmeTRI 12 cm, 4 mm	4.0
14BA30T	SymmeTRI 14 cm, 3 mm	3.0
14BA40T	SymmeTRI 14 cm, 4 mm	4.0
14BA50T	SymmeTRI 14 cm, 5 mm	5.0
14BA60T	SymmeTRI 14 cm, 6 mm	6.0
14BA75T	SymmeTRI 14 cm, 7.5 mm	7.5
15BA30T	SymmeTRI 15 cm, 3 mm	3.0
15BA40T	SymmeTRI 15 cm, 4 mm	4.0

EASYDRILL™ CRANIAL PERFORATORS

The EasyDrill™ Cranial Perforators combine efficient bone cutting performance with the protection of AutoStop technology.

PROTECTION


Employs a clutch mechanism to automatically disengage once it ceases to find resistance to bone.

PERFORMANCE

Drill tip design helps prevent slipping, while sharp edges cut quickly through bone.

VARIETY


Four color-coded sizes to meet a wide range of cranium types.

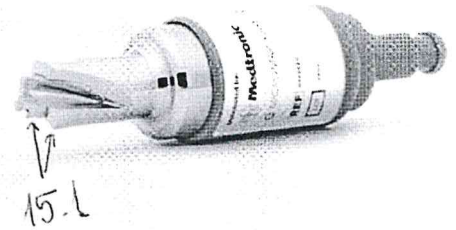ORDERING INFORMATION

Item #	Colour	Outer Diam.	Inner Diam.	Shelf	Minimum Bone Thickness
DM0008FAA	White	9.0 mm	6.0 mm	3.0 mm	3.25 mm
DM0011FAA	Blue	11.0 mm	7.0 mm	3.0 mm	3.25 mm
DM0210FAA	Green	14.0 mm	11.0 mm	1.5 mm	1.75 mm
DM0010FAA	Yellow	14.0 mm	11.0 mm	3.0 mm	3.25 mm

Item #	Description
AD03	Midas Rex® Perforator Driver Attachment Manufactured by Medtronic

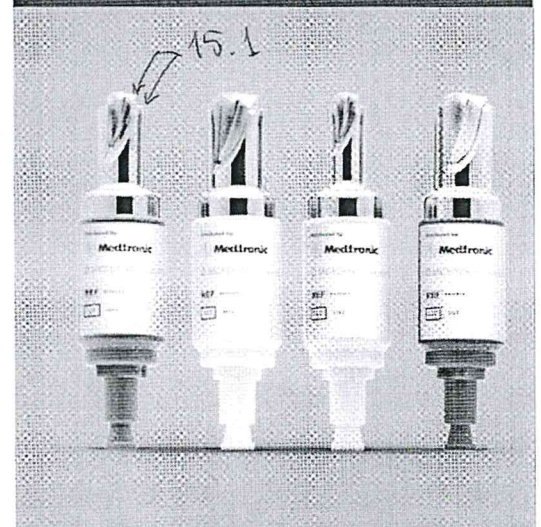Medtronic

Europe
15, rue de la Libération
92015 Neuilly-sur-Seine
France
Tél: +33 (0)1 47 00 10 00
Fax: +33 (0)1 47 00 10 01
www.medtronic.eu

United Kingdom/Ireland
The Medtronic Group
15, rue de la Libération
92015 Neuilly-sur-Seine
France
Tél: +33 (0)1 47 00 10 00
Fax: +33 (0)1 47 00 10 01
www.medtronic.eu

Manufactured by
MICROMAR
Distributed by Medtronic

15, rue de la Libération
92015 Neuilly-sur-Seine
France
Tél: +33 (0)1 47 00 10 00
Fax: +33 (0)1 47 00 10 01
www.medtronic.eu


Rx only. Refer to product instruction manual/
package insert for instructions, warnings,
precautions and contraindications.

For more information, contact your Medtronic
Neurosurgery sales representative or refer to
www.MedtronicNeurosurgery.com.

For a listing of indications, contraindications,
precautions, warnings, and potential adverse
events, please refer to the Instructions for Use.

www.medtronic.eu


20

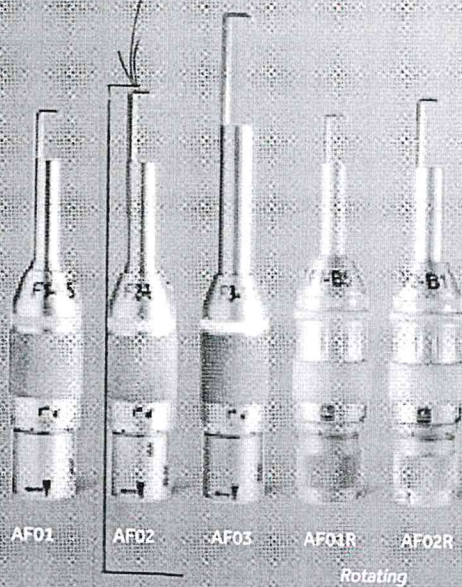FOOTED ATTACHMENTS


Also referred to as craniotome attachments, footed attachments are used in cranial and spinal procedures.

The footed attachments allow the surgeon to protect soft tissue when cutting through bone.

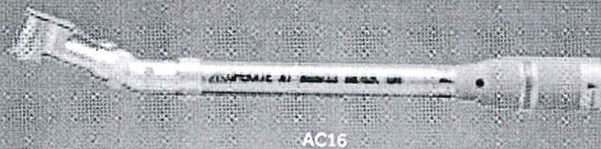
CONTRA ANGLED ATTACHMENT

The contra angled attachment is designed for limited access right-angled applications.

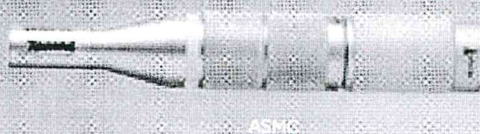
METAL CUTTING ATTACHMENT

This attachment can accommodate tungsten carbide or diamond tools for the trans-section of metal.

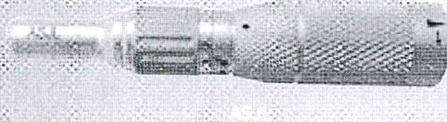
J-LATCH ATTACHMENT

ASJL universal attachment that accepts all manufacturer tools with a J shaped end; compatible with all Midas Rex motors.

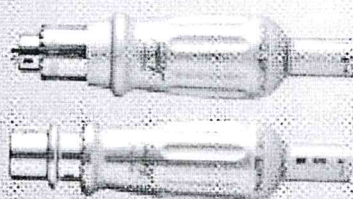PERFORATOR DRIVER AND JACOBS® CHUCK

The Perforator Driver is used to drive perforator bits with a Hudson shank, while the non-cannulated Jacobs Chuck drives pins and wires.

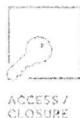F2/8TA23S

F3/9TA30

Legend[®] Tapered - Footed

TOOL TYPE	Head Diameter (mm)	Head Length (mm)	SHAFT LENGTH	
			8 cm small bore	9 cm large bore
	1.5	11.1	F1/8TA15	
	1.5	13.1	F1/8TA15S	
	2.3	15.9	F2/8TA23	
	2.3	16.4	F2/8TA23S	
	3.0	25.4		F3/9TA30
ATTACHMENT			 AS08 AVS08	 AS09, AS09S, AA09S AVS09, AVA09

Circle indicates other compatible attachments:

AF01 & AF01R

AF02 & AF02R

AF03

Legend[®] Tapered - Reverse Taper

TOOL TYPE	Head Diameter (mm)	Head Length (mm)	SHAFT LENGTH
			7 cm large bore
	6.4	6.1	14RT64
ATTACHMENT			 AS14, AS14S, AA14, AA14S, AVS14, AVA14

MIDAS REX® Legend® System

Straight Small Bore Attachments (2.4 mm)

Product	Description	Length (cm)
AS07	Straight Small Bore	7
AS075	Straight Small Bore	7.5
AS08	Straight Small Bore	8
AS10	Straight Small Bore	10
AS10S	Straight Small Bore	9.5
AS12S	Straight Small Bore	11.5
AS15	Straight Small Bore	15

AS07

Variable Exposure — 10 mm Adjustment

Product	Description	Length (cm)
AVS07	Straight Small Bore	7
AVS08	Straight Small Bore	8
AVS10	Straight Small Bore	10
AVS12	Straight Small Bore	12
AVS15	Straight Small Bore	15

AVS15

Straight Large Bore Attachments (3.2 mm)

Product	Description	Length (cm)
AS09	Straight Large Bore	9
AS09S	Straight Large Bore	8.5
AS14	Straight Large Bore	14
AS14P	Straight Large Bore — Precision Graft	14
AS14S	Straight Large Bore	13.5
AS21	Straight Large Bore	21
AS26	Straight Large Bore	26

AS09

Variable Exposure — 12 mm Adjustment

Product	Description	Length (cm)
AVS09	Straight Large Bore	9
AVS14	Straight Large Bore	14

AVS14

ACCESS /
CLOSURE

Legend® Ball - Fluted

21, 22, 25, 26, 27

TOOL TYPE	Head Diameter (mm)	SHAFT LENGTH							
		7 cm small bore	9 cm large bore	10 cm small bore	12 cm small bore	14 cm large bore	15 cm small bore	21 cm large bore	26 cm large bore
	0.5	7BA05							
	1.0	7BA10		10BA10					
	1.5	7BA15							
	2.0	7BA20		10BA20		14BA20			
	2.5	7BA25				14BA25			
	3.0	7BA30	9BA30	10BA30	12BA30	14BA30			
	4.0	7BA40	9BA40	10BA40	12BA40	14BA40	15BA40		26BA40
	5.0	7BA50	9BA50	10BA50	12BA50	14BA50	15BA50	21BA50	
	6.0	7BA60	9BA60	10BA60	12BA60	14BA60	15BA60	21BA60	26BA60
	7.5		9BA75			14BA75		21BA75	26BA75
	9.0		9BA90			14BA90		21BA90	26BA90
ATTACHMENT		AS07 AA07 AVS07 AVA07	AS09 AS09S AA09S AVS09 AVA09	AS10 AS10S AA10S AVS10 AVA10	AS12S AVS12 AVA12	AS14 AS14S AA14 AA14S AVS14 AVA14	AS15 AA15 AVS15 AVA15	AS21	AS26

ACCESS /
CLOSURE

Legend® Ball - Diamond

-23

23


TOOL TYPE	Head Di- ameter (mm)	SHAFT LENGTH						
		7 cm small bore	9 cm large bore	10 cm small bore	12 cm small bore	14 cm large bore	15 cm small bore	26 cm large bore
	1.0	7BA10D		10BA10D				
	1.5	7BA15D						
	2.0	7BA20D		10BA20D		14BA20D		
	2.5	7BA25D				14BA25D		
	3.0	7BA30D	9BA30D	10BA30D	12BA30D	14BA30D		
	4.0	7BA40D	9BA40D	10BA40D	12BA40D	14BA40D	15BA40D	26BA40D
	5.0	7BA50D	9BA50D	10BA50D	12BA50D	14BA50D	15BA50D	
	6.0	7BA60D	9BA60D	10BA60D	12BA60D	14BA60D	15BA60D	
	7.5		9BA75D			14BA75D		
ATTACHMENT		AS07 AA07 AVS07 AVA07	 AS09 AS09S AA09S AVS09 AVA09	 AS10 AS10S AA10S AVS10 AVA10	 AS12S AVS12 AVA12	 AS14 AS14S AA14 AA14S AVS14 AVA14	 AS15 AA15 AVS15 AVA15	AS26

MIDAS REX® Legend® System

Powered Surgical Products


Angled Small Bore Attachments (2.4 mm)

Product	Description	Length (cm)
AA07	Angled Small Bore	7
AA075	Angled Small Bore	7.5
AA10	Angled Small Bore	10
AA10S	Angled Small Bore	9.5
AA15	Angled Small Bore	15
AA07DK	Double Lock	7
AA10DK	Double Lock	10
AA10SDK	Double Lock	10
AA12DK	Double Lock	12
AA15DK	Double Lock	15

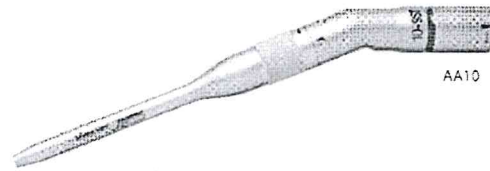Variable Exposure — 10 mm Adjustment

Product	Description	Length (cm)
AVA07	Angled Small Bore	7
AVA10	Angled Small Bore	10
AVA12	Angled Small Bore	12
AVA15	Angled Small Bore	15
AVA07DK	Double Lock	7
AVA10DK	Double Lock	10
AVA12DK	Double Lock	12
AVA15DK	Double Lock	15

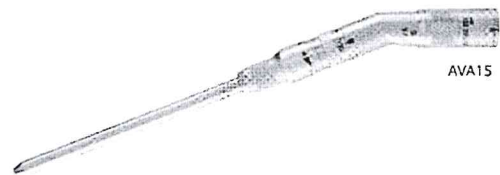Angled Large Bore Attachments (3.2 mm)

Product	Description	Length (cm)
AA09S	Angled Large Bore	8.5
AA14	Angled Large Bore	14
AA14S	Angled Large Bore	13.5
AA09DK	Double Lock	9
AA14DK	Double Lock	14
AA14SDK	Double Lock	14

Variable Exposure — 12 mm Adjustment

Product	Description	Length (cm)
AVA09	Angled Large Bore	9
AVA14	Angled Large Bore	14
AVA09DK	Double Lock	9
AVA14DK	Double Lock	14

1 AM-8 side cutting end SymmetRI designs


2 Shades and lengths for various applications


3 Midas Rex ClearView™ visibility and protected staff

4 Low-profile integrated irrigation

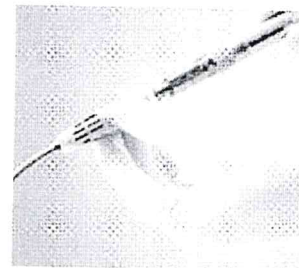5 Direct connect base for one-step assembly

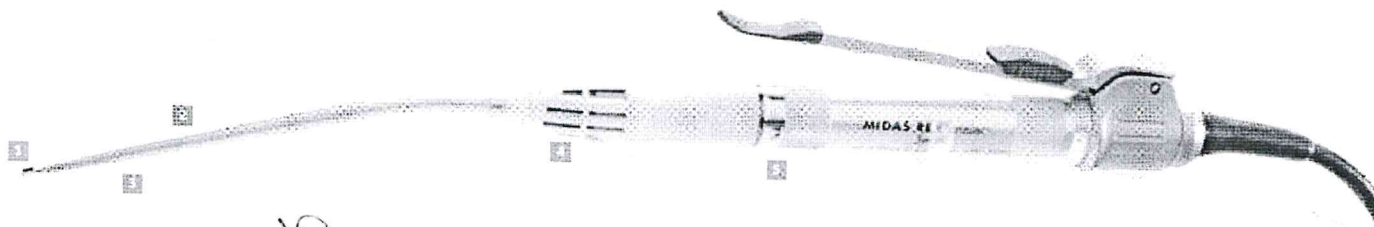28

ORDERING INFORMATION

12CM LENGTH, PROXIMAL BEND (SP12 PREFIX)

SP12MH25	2.5mm AM-8 Style Match Head
SP12MH25T	2.5mm SymmetRI Style Match Head
SP12MH25D	2.5mm Diamond Match Head
SP12MH30	3.0mm AM-8 Style Match Head
SP12MH30T	3.0mm SymmetRI Style Match Head
SP12MH30D	3.0mm Diamond Match Head
SP12BA30	3.0mm Fluted Ball
SP12BA40	4.0mm Fluted Ball

29

30

12CM LENGTH, DISTAL BEND (SD12 PREFIX)

SD12MH25	2.5mm AM-8 Style Match Head
SD12MH25T	2.5mm SymmetRI Style Match Head
SD12MH25D	2.5mm Diamond Match Head
SD12MH30	3.0mm AM-8 Style Match Head
SD12MH30T	3.0mm SymmetRI Style Match Head
SD12MH30D	3.0mm Diamond Match Head
SD12BA30	3.0mm Fluted Ball
SD12BA40	4.0mm Fluted Ball

14CM LENGTH, PROXIMAL BEND (SP14 PREFIX)

SP14MH25	2.5mm AM-8 Style Match Head
SP14MH25T	2.5mm SymmetRI Style Match Head
SP14MH25D	2.5mm Diamond Match Head
SP14MH30	3.0mm AM-8 Style Match Head
SP14MH30T	3.0mm SymmetRI Style Match Head
SP14MH30D	3.0mm Diamond Match Head
SP14BA30	3.0mm Fluted Ball
SP14BA40	4.0mm Fluted Ball

CLEARVIEW IRRIGATION

IRD370	Clear View Irrigation Tubing, Single
IRD375	Clear View Irrigation Tubing, 5 pack
IRD200	IntelliFlow™ Irrigation Remote Control, w5 pack

31

For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions of Use.

For more information please contact your Medtronic sales representative or refer to www.medtronicneurosurgery.com

Medtronic

Medtronic Limited
Building 9
Croxley Green Business Park
Watford
Hertfordshire WD18 8WW
UK
www.medtronic.co.uk
Tel: +44 (0)1923 212213
Fax: +44 (0)1923 241004

Medtronic Ireland Ltd
Swords Business Campus
Balheary Road
Swords
Co. Dublin
Ireland
Tel: +353 1 5111400


Europe
Medtronic International
Trading Sàrl
Route du Molliou 31
Case postale
CH-1151 Tolochenaz
www.medtronic.eu
Tel: +41 (0)21 802 70 00
Fax: +41 (0)21 802 79 00

www.medtronic.eu

27

32 PANA SPRAY Plus is a specially formulated, highly effective general handpiece and air motor lubricant. Frequent application of PANA SPRAY Plus will assure handpiece and air motor performance and durability. Apply PANA SPRAY Plus every time before autoclaving any handpiece or air motor.

PANA SPRAY Plus

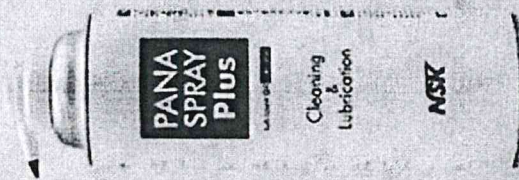Model PANA SPRAY Plus

Order Code Z182100

Model PANA SPRAY Plus Package*

Order Code Y900630

* Pack of 6

For High & Low speed handpieces, and air motors

Alcohol-contained lubricant

Contains 480 mL each



PANA SPRAY Plus is not intended for cleaning in the US

PANA SPRAY Plus is 510(k) cleared in the United States as a lubricant intended to be used during routine maintenance of dental and medical surgical handpieces after each patient use and prior to sterilization.

MIDAS REX[®] Legend[®] System

Accessories

Product	Description	Quantity
PC100	Pneumatic Control Unit FN2 - DISS (with 20 ft. hose)	1
PC101	Pneumatic Control Unit (Schrader)	1
PA100-A	Lubricant/Diffuser Pack	4
PA130	Legend Motor Wrench	1
PA200	Regulator (DISS conn)	1
PA300	1.2 mm Cleaning Brush (1/32")	1
PA305	1.2 mm Cleaning Brush	5
PA310	2.4 mm Cleaning Brush (3/32")	1
PA315	2.4 mm Cleaning Brush	5
PA320	3.2 mm Cleaning Brush (1/8")	1
PA325	3.2 mm Cleaning Brush	5
PA115	Tool Storage Rack	1
PC110	Adaptor FN2-DISS to M-Schr	1
PC120	Adaptor FN2-DISS to MA-DISS	1
PC130	WF4 Connector	1
PC140	Adaptor FN2-DISS to F-Schr	1
PC230	30 ft. Extension Hose, Black	1
PC231	30 ft. Extension Hose, Black (Schrader)	1
PC720	MR7 Triton [®] Adaptor	1
PC700	Midas Rex [®] MR7 Pneumatic Foot Control	1
PC710	Midas Rex MR7 Regulator Hose	1
PA700	Midas Rex MR7 Lubricant/Diffuse Cartridge	4-pk

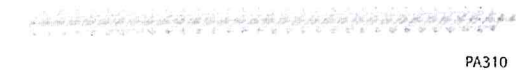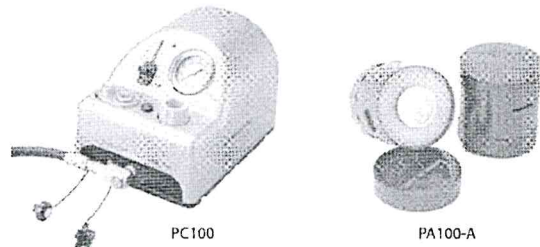
MIDAS REX® Legend® System

Powered Surgical Products

Instrument Case 3/4 DIN

Product	Description	Quantity
PA400		
<i>Case includes the following:</i>		
PA410	Instrument Case Base 3/4 DIN	1
PA415	Instrument Case Lid 3/4 DIN	1
PA420	Instrument Case Tray 3/4 DIN	1
PA421	Instrument Case Labels 3/4 DIN	1
PA425	Instrument Case Mat 3/4 DIN	1

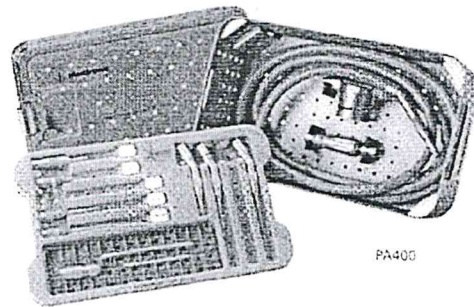
Instrument Case Kit Full DIN with Triton® Tray (not shown)

Product
PA500

Midas Rex® MR7 Instrument Case 3/4 DIN Case includes base, lid, tray and tray mat

Product	Description	Quantity
PA600	Midas Rex® MR7 Instrument Case 3/4 DIN	1
<i>Replacement components:</i>		
PA615	Midas Rex® MR7 Instrument Case Replacement (Lid)	1
PA620	Midas Rex® MR7 Instrument Case Replacement (Tray)	1
PA625	Midas Rex® MR7 Instrument Case Replacement (Tray Mat)	1

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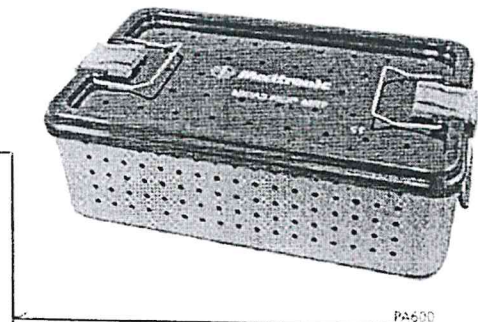ORDERING INFORMATION

Item #	Colour	Outer Diam.	Inner Diam.	Shelf	Minimum Bone Thickness
DM0008FAA		9.0 mm	6.0 mm	3.0 mm	3.25 mm
DM0011FAA		11.0 mm	7.0 mm	3.0 mm	3.25 mm
DM0210FAA		14.0 mm	11.0 mm	1.5 mm	3.75 mm
DM0010FAA		14.0 mm	11.0 mm	3.0 mm	3.25 mm
Item #	Description				
*AD003	Midas Rex™ Perforator Driver Attachment Manufactured by Medtronic				

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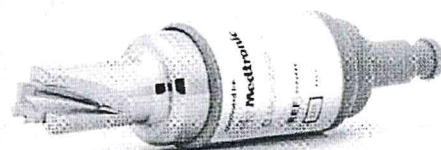
• AF01 & AF01R • AF02 & AF02R • AF03

Legend® Tapered - Reverse Taper

TOOL TYPE	Head Diameter (mm)	Head Length (mm)	SHAFT LENGTH
			7 cm large bore
	6.4	6.1	14RT64
ATTACHMENT			AS14, AS14S, AA14, AA14S, AVS14, AVA14

ACCESS /
CLOSURE

Legend® Ball - Fluted

4, 5, 7, 8

TOOL TYPE	Head Diameter (mm)	SHAFT LENGTH							
		7 cm small bore	9 cm large bore	10 cm small bore	12 cm small bore	14 cm large bore	15 cm small bore	21 cm large bore	26 cm large bore
	0.5	7BA05							
	1.0	7BA10		10BA10					
	1.5	7BA15							
	2.0	7BA20		10BA20		14BA20			
	2.5	7BA25				14BA25			
	3.0	7BA30	9BA30	10BA30	12BA30	14BA30			
	4.0	7BA40	9BA40	10BA40	12BA40	14BA40	15BA40		26BA40
	5.0	7BA50	9BA50	10BA50	12BA50	14BA50	15BA50	21BA50	
	6.0	7BA60	9BA60	10BA60	12BA60	14BA60	15BA60	21BA60	26BA60
	7.5		9BA75			14BA75		21BA75	26BA75
	9.0		9BA90			14BA90		21BA90	26BA90
ATTACHMENT		AS07 AA07 AVS07 AVA07	AS09 AS09S AA09S AVS09 AVA09	AS10 AS10S AA10 AA10S AVS10 AVA10	AS12S AVS12 AVA12	AS14 AS14S AA14 AA14S AVS14 AVA14	AS15 AA15 AVS15 AVA15	AS21	AS26

4, 5
7, 8, 9

ACCESS /
CLOSURE

Legend® Ball - Diamond

TOOL TYPE	Head Di- ameter (mm)	SHAFT LENGTH						
		7 cm small bore	9 cm large bore	10 cm small bore	12 cm small bore	14 cm large bore	15 cm small bore	26 cm large bore
	1.0	7BA10D		10BA10D				
	1.5	7BA15D						
	2.0	7BA20D		10BA20D		14BA20D		
	2.5	7BA25D				14BA25D		
	3.0	7BA30D	9BA30D	10BA30D	12BA30D	14BA30D		
	4.0	7BA40D	9BA40D	10BA40D	12BA40D	14BA40D	15BA40D	26BA40D
	5.0	7BA50D	9BA50D	10BA50D	12BA50D	14BA50D	15BA50D	
	6.0	7BA60D	9BA60D	10BA60D	12BA60D	14BA60D	15BA60D	
	7.5		9BA75D			14BA75D		
ATTACHMENT		AS07 AA07 AVS07 AVA07	AS09 AS09S AA09S AVS09 AVA09	AS10 AS10S AA10 AA10S AVS10 AVA10	AS12S AVS12 AVA12	AS14 AS14S AA14 AA14S AVS14 AVA14	AS15 AA15 AVS15 AVA15	AS26


1 AM-8 side cutting and SymmeTRI designs

2 Shapes and lengths for various applications

3 Midas Rex ClearView™ visibility and protected shaft

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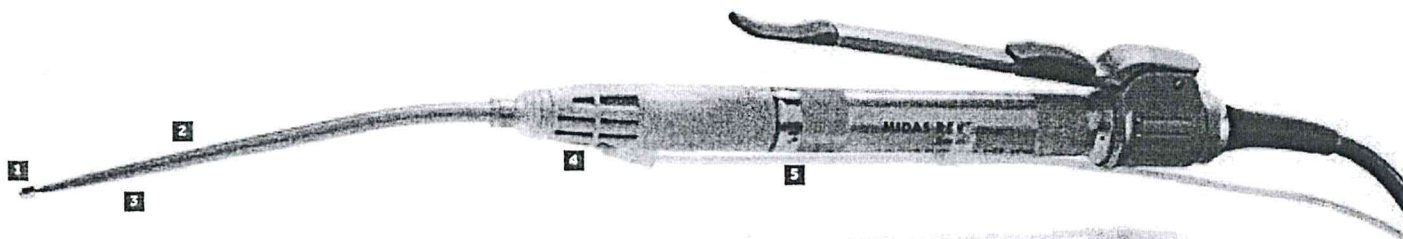ORDERING INFORMATION

12CM LENGTH, PROXIMAL BEND (SP12 PREFIX)

SP12MH25	2.5mm AM-8 Style Match Head
SP12MH25T	2.5mm SymmeTRI Style Match Head
SP12MH25D	2.5mm Diamond Match Head
SP12MH30	3.0mm AM-8 Style Match Head
SP12MH30T	3.0mm SymmeTRI Style Match Head
SP12MH30D	3.0mm Diamond Match Head
SP12BA30	3.0mm Fluted Ball
SP12BA40	4.0mm Fluted Ball

12CM LENGTH, DISTAL BEND (SD12 PREFIX)

SD12MH25	2.5mm AM-8 Style Match Head
SD12MH25T	2.5mm SymmeTRI Style Match Head
SD12MH25D	2.5mm Diamond Match Head
SD12MH30	3.0mm AM-8 Style Match Head
SD12MH30T	3.0mm SymmeTRI Style Match Head
SD12MH30D	3.0mm Diamond Match Head
SD12BA30	3.0mm Fluted Ball
SD12BA40	4.0mm Fluted Ball

14CM LENGTH, PROXIMAL BEND (SP14 PREFIX)

SP14MH25	2.5mm AM-8 Style Match Head
SP14MH25T	2.5mm SymmeTRI Style Match Head
SP14MH25D	2.5mm Diamond Match Head
SP14MH30	3.0mm AM-8 Style Match Head
SP14MH30T	3.0mm SymmeTRI Style Match Head
SP14MH30D	3.0mm Diamond Match Head
SP14BA30	3.0mm Fluted Ball
SP14BA40	4.0mm Fluted Ball

CLEARVIEW IRRIGATION

IRD370	ClearView Irrigation Tubing, Single
IRD375	ClearView Irrigation Tubing, 5 pack
IRD200	IntelliFlow™ Irrigation Remote Control, w5 pack

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Medtronic

Medtronic Limited
Building 9
Croxley Green Business Park
Watford
Hertfordshire WD18 8WW
UK
www.medtronic.co.uk
Tel: +44 (0)1923 212215
Fax: +44 (0)1923 241004

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