

# Pro Fusion

## Product Description

## Pro Fusion Product Description

v. 2.4

October 2021

All rights reserved.

Copyright © Tobii Pro AB (publ)

The information contained in this document is proprietary to Tobii Pro AB. Any reproduction in part or whole without prior written authorization by Tobii Pro AB is prohibited.

Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks.

While every precaution has been taken in the preparation of this document, the publisher and the author assume no responsibility for errors or omissions, or for damages resulting from the use of information contained in this document or from the use of programs and source code that may accompany it. In no event shall the publisher and the author be liable for any loss of profit or any other commercial damage caused or alleged to have been caused directly or indirectly by this document.

Content subject to change without notice.

Please visit [Tobii Pro Connect](#) for updated versions of this document.

# Table of Contents

<b>1 Overview of Tobii Pro Fusion</b>	<b>5</b>
1.1 Pro Fusion system requirements	5
1.2 Pro Fusion versions	5
1.3 Key features	6
1.4 What's in the Pro Fusion box	7
1.5 Tobii Pro Eye Tracker Manager	8
1.6 Tobii Pro accessories	8
1.7 Additional software options	9
<b>2 Technical Specifications</b>	<b>10</b>
2.1 Eye tracking specifications	10
2.2 Setup	11
2.3 Software compatibility	11
2.4 Hardware versions	11
2.5 Pro Fusion eye tracker specifications	11
<b>3 Compliance Information</b>	<b>13</b>
3.1 Pro Fusion compliance	13
3.1.1 FCC interference statement	13
3.1.2 Industry Canada statement	13
3.1.3 CE statement	13
3.1.4 UK statement	14
3.1.5 Safety compliance	14
3.1.6 EMC compliance	14
<b>4 Safety</b>	<b>15</b>
4.1 Mounting warning	15
4.2 Emergency warning	15
4.3 Epilepsy warning	15
4.4 Infrared warning	15
4.5 Magnetic field warning	16
4.6 Child safety	16
4.7 Electricity	16
4.8 Accessories	16
4.9 Third party	16
4.10 Power supply	16
<b>5 Product Care</b>	<b>17</b>
5.1 Operational environment	17
5.2 Transportation and storage	17
5.3 Cleaning	17
5.4 Disposal of the eye tracker	17
5.5 Limitation of liability	17
<b>Appendix A Support, Training, and Warranty</b>	<b>19</b>
A1 Customer Support	19
A1.1 Get help online	19
A2 Training and Education Services	19
A3 Warranty information	19

Appendix B Limitations and considerations .....	20
B1 Intended use .....	20
B2 Light conditions .....	20
B3 Eyelashes .....	20
B4 Droopy eyelids .....	20

# 1 Overview of Tobii Pro Fusion

Tobii Pro Fusion is the next generation of compact high-performance eye trackers from Tobii Pro. It provides a flexible solution that supports different research scenarios and study designs.

Pro Fusion has a slim design that integrates well with external monitors up to 24 inches and laptop screens. The eye tracker can also be mounted on a tripod and coupled with a scene camera so it can be used in study setups that use real-world stimuli such as physical objects or scenes. The sampling speed of up to 250 Hz allows you to capture data for a wide range of experimental paradigms (fixation, smooth pursuit, and saccade-based).

With high, selectable, sampling rates (up to 250 Hz), two eye tracking cameras and two pupil tracking modes (bright and dark pupil), Pro Fusion enables you to adapt your data collection setup to different research populations, scenarios, and data requirements.

Use Pro Fusion for studies on:

- Laptop and desktop computers
- PC monitors, 24" (16:9) or less
- Real-world setups such as physical objects and social interactions (using the Tripod Stand)
- Larger screens, projections, simulators (using the [Tripod Stand](#), available separately)

## 1.1 Pro Fusion system requirements

<b>Operating system:</b>	Windows 10 macOS 10.14 Mojave and 10.15 Catalina Linux Ubuntu 20.04
<b>CPU:</b>	1 GHz, 2 cores
<b>RAM:</b>	2 GB RAM memory* *Software for eye tracking research may require higher RAM
<b>Port:</b>	USB Type-A or USB Type-C

## 1.2 Pro Fusion versions

Pro Fusion is available in three different product versions: 60 Hz, 120 Hz, and 250 Hz.

The 60 Hz and 120 Hz versions can be upgraded to a faster version by contacting your sales representative. [Contact Tobii Pro](#).



All versions of the Pro Fusion eye tracker can also run in 30 Hz mode.

## 1.3 Key features

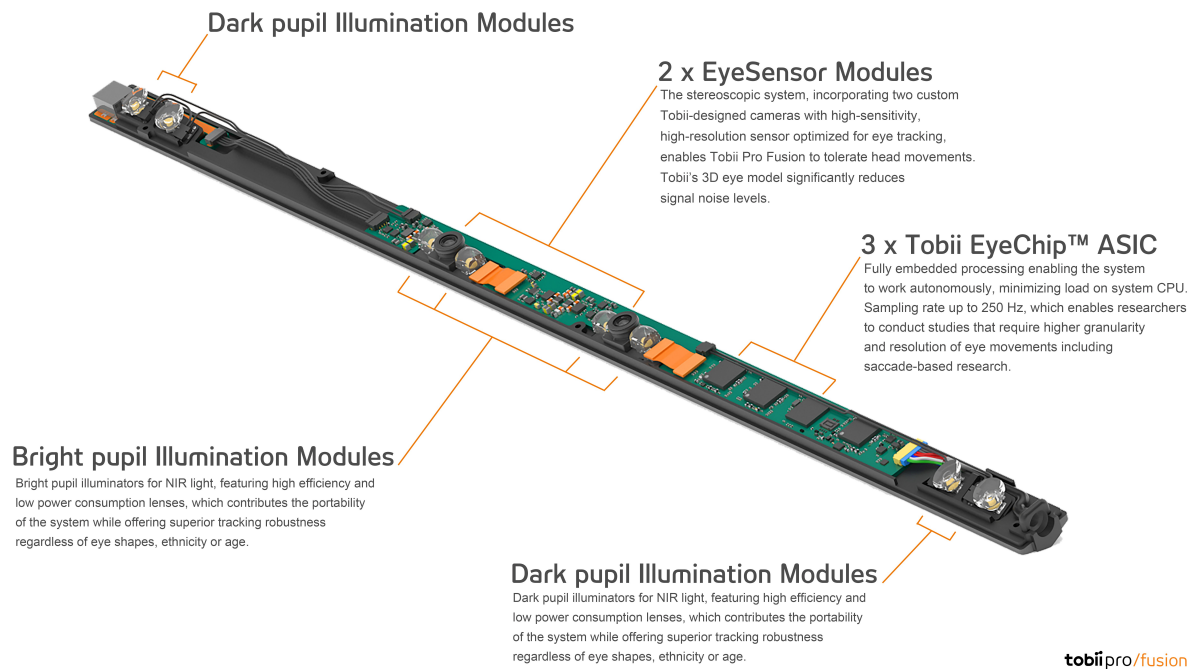


Figure 1. Pro Fusion's key features

## 1.4 What's in the Pro Fusion box



Figure 2. What's in the Pro Fusion box

#	Item	Additional information
1	Travel case	Stores, protects, and carries Pro Fusion and the included accessories
2	Pro Fusion eye tracker	Used for eye tracking

3	AC power adapter	Provides the additional power required when the USB Type-A to Type-C adapter is used
4	Safety and Compliance	Printed versions of important information
5	USB Type-A to Type-C adapter plug	Used when the computer does not have a USB Type-C port
6	6 International plug adapters	Used with the AC power adapter to match international outlets
7	Cleaning cloth	For cleaning the surface of the eye tracker
8	Four mounting plates with non-permanent adhesive	Attaches the Pro Fusion eye tracker to the screen bezel
9	Four cleaning wipes	Removes any dust and grease in the area around the screen before attaching a mounting plate

## 1.5 Tobii Pro Eye Tracker Manager

Tobii Pro Eye Tracker Manager is free software available on [Tobii Pro Connect](#).

Pro Eye Tracker Manager helps you manage your Tobii Pro eye tracker. It lets you see which eye tracker firmware and software versions you are running and also facilitates updates (where applicable). It provides:

- Drivers and Firmware installation or updates
- Display setup
- Positioning guide
- User calibration
- Gaze visualization<sup>1</sup>
- Eye images<sup>2</sup>



References in product documentation to Tobii Pro Eye Tracker Manager refer to the latest version available on Tobii Pro's website.

## 1.6 Tobii Pro accessories

Find out more about Tobii Pro products and accessories on [tobiipro.com](https://tobiipro.com).

<sup>1</sup>Gaze visualization in Pro Eye Tracker Manager lets you quickly evaluate the quality of the eye tracking data. The Gaze visualization overlaps data only on the Pro Eye Tracker Manager screen while Pro Eye Tracker Manager is open.

<sup>2</sup>Applicable when Pro Eye Tracker Manager is used with eye trackers that support eye images (Tobii Pro Spectrum and Tobii Pro Fusion).



## 1.7 Additional software options

Pro Lab <a href="https://tobiipro.com/product-listing/tobii-pro-lab/">tobiipro.com/product-listing/tobii-pro-lab/</a>	Pro Lab is a comprehensive research software platform for eye tracking designed to meet the highest demands on different research scenarios with exact timing accuracy. This software offers an efficient workflow, making it easy to design experiments, record data, analyze and visualize eye tracking data, and to sync this data with other biometric data streams.
Pro SDK <a href="https://tobiipro.com/product-listing/tobii-pro-sdk/">tobiipro.com/product-listing/tobii-pro-sdk/</a>	Pro SDK offers a broad set of tools that makes it simple to develop a variety of niche applications or scripts across multiple platforms, using a wide range of programming languages. This SDK gives the researcher access to the full set of relevant gaze data streams, such as 3D eye coordinates, raw data, pupil data, etc.
Third-party software	This term refers to any application built on Pro SDK.

## 2 Technical Specifications

### 2.1 Eye tracking specifications

Eye tracking technique	Video-based pupil- and corneal reflection eye tracking with dark and bright pupil illumination modes. <sup>1</sup> Two cameras capture stereo images of both eyes for robust accurate measurement of eye gaze and eye position in 3D space, as well as pupil diameter.
Sampling frequency	60, 120, & 250 Hz or 60 & 120 Hz, depending on the hardware version
Precision <sup>2</sup>	0.04° RMS at optimal conditions <sup>3</sup> 0.2° RMS in optimal conditions (raw signal)
Accuracy <sup>4</sup>	0.3° in optimal conditions
Binocular eye tracking	Yes
Total system latency <sup>5</sup>	Mean latency < 13 ms @ 250 Hz (SD < 1.2 ms)
Blink recovery time	1 frame (immediate)
Gaze recovery time	250 ms
Data sample output <sup>6</sup>	Timestamp Gaze origin Gaze point Pupil diameter
Eye image data stream	Eye image data stream frequency is approximately 4 Hz (one image per eye). Zoomed-in eye images available in tracking mode. Full-frame camera images are available in gaze recovery mode.
TTL input stream	Not Available
Tracker and client time synchronization	Integrated between the eye tracker time domain and the client computer time domain

---

<sup>1</sup>Dark pupil tracking is supported in all sample frequencies. Bright pupil tracking mode is supported at 60 and 120 Hz

<sup>2</sup>Tobii Pro uses an extensive test method to measure and report performance and quality of data. Please download the [Tobii Pro Fusion Eye Tracker Data Quality Test Report](#) for more detailed information.

<sup>3</sup>Applying Savitzky-Golay filtering (settings listed in the [Tobii Pro Fusion Eye Tracker Data Quality Test Report](#)).

<sup>4</sup>(see footnote 2)

<sup>5</sup>Visit [Tobii Pro Connect](#) for more information on how to optimize your Pro Fusion setup.

<sup>6</sup>For more information about output data and the supplementary data stream, read the [Tobii Pro SDK documentation](#).

## 2.2 Setup

Head movement tolerance*	Excellent. Dual-camera system with more images than a one camera system which gives a more accurate data calculation and the best level of precision and robustness for head movement.
Freedom of head movement * (at 65cm distance)	Width × height: 40 cm × 25 cm (15.7" × 9.84") (At least one eye tracked)
Freedom of head movement* (at 80 cm distance)	Width x height: 45 cm x 30 cm (17.7" x 11.8") (At least one eye tracked)
Operating distance (mounted on screen)	50-80 cm (19.69"-31.49") from the eye tracker
Tracker setup options	Tracker mounted on tripod, allows for even larger screens and physical objects to be tracked.
Optimal screen size	24" (16:9 aspect ratio)

\* Describes the region in space where the participant can move their head and still have at least one eye within the eye tracker's field of view (trackbox) at the specific distance.

## 2.3 Software compatibility

Software and framework compatibility	Tobii Pro Lab Tobii Pro Eye Tracker Manager Tobii Pro SDK Any application built on the Tobii Pro SDK
Operating system	Windows, macOS, Linux

## 2.4 Hardware versions

60 HZ
120 Hz
250 Hz

## 2.5 Pro Fusion eye tracker specifications

Dimensions	374 mm L × 18 mm H × 13.7 mm W (14.72" × 0.70" × 0.53")
Weight	168g (5.9 oz.)
Connectors	USB Type-C (USB Type-A to Type-C adapter provided) AC power
Data processing	3 Tobii EyeChip™ ASIC with fully embedded data processing
Eye tracking cameras	2 x Tobii EyeSensor™ Modules

Illuminators	Dark Pupil Illumination Modules Bright Pupil Illumination Modules
Power consumption	Typical power consumption: 4.3 W Max. rated power consumption: 9 W
Power options	Directly via USB Type-C connector, or with provided AC power adapter, when using a computer with USB 2.0 Type-A port

## 3 Compliance Information

### 3.1 Pro Fusion compliance



This Tobii Pro eye tracker is CE-marked, indicating compliance with applicable health and safety legal requirements set out in the European Union and EEA.



この装置は、現在設置されている場所で妨害波の測定がされた情報技術装置です。この場所以外で使用する場合は、その場所で、再び妨害波の測定が必要となります。



This Tobii Pro eye tracker is registered in the Australian National Equipment Registration System.

The RCM (Regulatory Compliance Mark) label signifies compliance with the Australian and New Zealand product safety requirement for electrical and electronic equipment according to EESS and radio communication and electromagnetic compatibility according to ACMA.



This Tobii Pro eye tracker is KC-marked which signifies compliance with Korea's product safety requirements for electrical and electronic equipment and components.

#### 3.1.1 FCC interference statement

Pro Fusion has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. Pro Fusion can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If Pro Fusion does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by increasing the separation between Pro Fusion and the receiver. Modifications not expressly approved by Tobii Pro AB could void the user's authority to operate the equipment under FCC rules.

#### 3.1.2 Industry Canada statement

This Class B digital apparatus complies with Canadian ICES-003.

Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.

#### 3.1.3 CE statement

Pro Fusion is CE-marked, certifying compliance with the legal requirements for health, safety and environmental protection requirements for the European market. Pro Fusion complies with the following directives:

- 2011/65/EU (RoHS 2) - Restriction of Hazardous Substances Directive
- 2014/30/EU (EMC) - Electromagnetic Compatibility Directive

- 2012/19/EU (WEEE) - Waste Electrical and Electronic Equipment Directive
- 2014/35/EU (LVD) - Low Voltage Directive

#### 3.1.4 UK statement

Pro Fusion is UKCA-marked, certifying compliance with the legal requirements for health, safety and environmental protection requirements for the European market. Pro Fusion complies with the following Statutory Instruments:

- The Electrical Equipment (Safety) Regulations 2016. UK Statutory Instruments 2016 No.1101.
- The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012. UK Statutory Instruments 2012 No.3032.
- The Electromagnetic Compatibility Regulations 2016. UK Statutory Instruments 2016 No.1091.

#### 3.1.5 Safety compliance

Pro Fusion complies with the following standards:

- IEC/EN 62471:2008, Photo Biological Safety of Lamps and Lamp Systems
- IEC 60950-1:2005 +A1 +A2, Safety of Information Technology Equipment
- IEC/EN 62368-1, including collateral standard IEC 62368-3 (European Standard including UK)
- CAN/CSA-C22.2 No. 62368-1/UL 62368-1 Edition No. 2, including collateral standard IEC 62368-3 (US and Canadian Standard)
- J62368-1 (2020) (Japanese Standard)
- GB 4943.1-2011 (Chinese Standard)
- AS/NZS 62368.1:2018 (New Zealand, Australia)

#### 3.1.6 EMC compliance

Pro Fusion complies with the following standards:

- EN 55032:2015, Electromagnetic compatibility of multimedia equipment. Emission requirements (European standard)
- EN 55035:2017 Electromagnetic compatibility of multimedia equipment - Immunity requirements
- FCC part 15 Class B, Regulations under which an intentional, unintentional, or incidental radiator that can be operated without an individual license (American Standard)
- ICES-003 Issue Class B, Interference-Causing Equipment Standard (Canadian Standard)
- CISPR32:2015, Electromagnetic compatibility of multimedia equipment - Emission requirements (International standard accepted in AS/NZS and Japan)
- GB17625.1-2012 and GB9254-2008 (Chinese national standard)
- KN32 (CISPR 32), KN35 (CISPR 35) (South Korean Standard)

# 4 Safety

## 4.1 Mounting warning



Tobii Pro eye trackers should be mounted according to Tobii Pro's instructions for approved mounts.

The mounting bracket is designed for adhesion to a wide variety of surfaces, such as monitors or laptop screens. It is also designed for clean removability. Failure to follow the recommended removal instructions can result in damage to both the monitor/laptop and the mounting bracket. Do not mount the eye tracker on monitors positioned above the head or face of a user, as it might fall.

## 4.2 Emergency warning



Tobii Pro eye trackers are designed to be used only for research purposes. Be aware that due to the low, but possible risk of failure or distraction, the eye tracker should not be relied upon or used in dangerous or critical situations.

## 4.3 Epilepsy warning



Some people with photosensitive epilepsy are susceptible to epileptic seizures or loss of consciousness when exposed to certain flashing lights or light patterns in everyday life. This may happen even if the person has no medical history of epilepsy or has never had any epileptic seizures.

A person with photosensitive epilepsy would also be likely to have problems with TV screens, some arcade games, and flickering fluorescent bulbs. Such people may have a seizure while watching certain images or patterns on a monitor, or even when exposed to the light sources of an eye tracker. It is estimated that about 3-5% of people with epilepsy have this type of photosensitive epilepsy. Many people with photosensitive epilepsy experience an "aura" or feel odd sensations before the seizure occurs. If you feel odd during use, move your eyes away from the eye tracker.

## 4.4 Infrared warning



When activated, the Tobii Pro eye tracker emits pulsed infrared (IR) light. Certain medical devices are susceptible to disturbance by IR light and/or radiation. Do not use the eye tracker in the vicinity of these kinds of susceptible medical devices, as their accuracy or proper functionality could be inhibited. Do not stare into the LED lights at a close distance (<15cm).

## 4.5 Magnetic field warning



This Tobii Pro eye tracker contains magnets. Magnetic fields may interfere with the function of cardiac pacemakers and implantable cardioverter defibrillators. As a general rule, maintain a minimum distance of 6 inches (15 centimeters) between any item with magnets and your heart device.

## 4.6 Child safety



A Tobii Pro eye tracker is an advanced computer system and electronic device. As such, it is composed of numerous separate, assembled parts. In the hands of a child, some of these parts have the possibility of being separated from the device, possibly resulting in a choking hazard or other danger to the child.

Young children should not have access to, nor use of the device or its accessories without parental or guardian supervision.

## 4.7 Electricity



Do not open the casing of the Tobii Pro eye tracker, this could expose you to potentially hazardous electrical voltage. The device contains no user-serviceable parts. Non-compliance will result in loss of warranty! Contact Tobii Pro Support if your eye tracker is not working properly.

## 4.8 Accessories



Only use accessories provided by or approved by Tobii Pro AB.

## 4.9 Third party



Any use of a Tobii Pro eye tracker outside the intended use and together with any third-party software or hardware that changes the intended use is a risk, and Tobii Pro AB can not take any responsibility in these situations.

## 4.10 Power supply



Only use the AC power adapter provided by Tobii Pro to power your Tobii Pro eye tracker. If the AC power adapter is damaged, it should be replaced only by Tobii Pro Service Personnel. If damaged, do not use the AC power adapter until it has been replaced. Tobii Pro and its agents are not liable for any damages or injuries to a person or property due to incorrect use of the provided accessories.



# 5 Product Care

## 5.1 Operational environment

The recommended range for temperature and humidity during operations is:

- **Temperature:** 10°C to 35°C (50°F to 95°F)
- **Humidity:** Max 20% to 95% (no condensation on the device)

## 5.2 Transportation and storage

Use the provided case when traveling with the eye tracker. For shipping and storage, use the original packing.

For transportation and storage, the recommended range for temperature and humidity for the device is as follows:

- **Temperature:** -40°C to 70°C (-40°F to 158°F)
- **Humidity:** 20% to 95% (no condensation on the device)



Tobii Pro eye trackers are not waterproof or water resistant. The eye tracker should not be kept in excessively humid, damp or wet conditions. Do not submerge the eye tracker in water or in any other liquid. Be careful not to spill liquids on the device.

## 5.3 Cleaning

- Before cleaning your Tobii Pro eye tracker, unplug its USB cable and electrical cable if there is one.
- Use the included cleaning cloth and clean the eye tracker's front gently.
- Avoid getting moisture in openings on the device.
- Do not use window cleaners, household cleaners, aerosol sprays, solvents, alcohol, ammonia, or abrasives.



Scratches on the front surface of the eye tracker may impair its performance.

## 5.4 Disposal of the eye tracker

Do not dispose of your Tobii Pro eye tracker in general household or office garbage receptacles. Follow your local regulations for the proper disposal of electrical and electronic equipment.

## 5.5 Limitation of liability

Except where prohibited by law, Tobii Pro AB is not liable for any loss or damage arising from the use or misuse of a Tobii Pro eye tracker or the mounting plate, whether direct, indirect, spe-

cial, incidental or consequential, regardless of the legal theory asserted, including warranty, contract, negligence, or strict liability.

# Appendix A      Support, Training, and Warranty

## A1      Customer Support

If you need help, please contact [Customer Support](#) at Tobii Pro. In order to receive assistance as quickly as possible, make sure you have access to your Tobii Pro device and, if possible, an Internet connection. You should also be able to supply the serial number of the device, which you will find on a sticker on the back or bottom of the device.

### A1.1      Get help online

Many questions can be answered by visiting Tobii Pro Connect. It contains the latest information about contacting Support, helpful articles and FAQs, links to downloads, and much more. Log in or register to see information about your account and to reach Customer Support on [Tobii Pro Connect](#).

## A2      Training and Education Services

If you are new to eye tracking, or want to extend your knowledge about eye tracking research, sign up for one of our learning programs and events on our website at [Training and Education Services](#).

## A3      Warranty information

Read more online about [Tobii Pro Care and Tobii Pro's eye tracker warranty](#).

# Appendix B      Limitations and considerations

## B1      Intended use



Tobii Pro eye trackers are intended to be used in research activities about human behavior including eye movements, in a dry and dust free indoor environment. The product should only be used as described in this User's Manual. Please read the User Manual and other supplied documentation thoroughly before using the product.

## B2      Light conditions



Tobii Pro recommends that eye tracking studies be performed in a controlled environment. Sunlight should be avoided since it contains high levels of infrared light which will interfere with the eye tracker system. Sunlight affects eye tracking performance severely and longer exposure can overheat the eye tracker. This eye tracker is not designed for exposure to (direct) sunlight. Eye tracking generally does not work in strong direct sunlight.

Shielding the eye tracker adequately from the sun may prevent sunlight from interfering with eye tracking.

## B3      Eyelashes



Long eyelashes can be obstructive when the participant's eyes are less open, especially if the participant is wearing mascara. In rare cases, eyelashes may completely block the view of the participant's pupils, making eye tracking impossible.

## B4      Droopy eyelids



Droopy eyelids or otherwise obstructive eyelids can block the view of the participant's pupils. In rare cases, such eyelids may completely block the view of the participant's pupils, making eye tracking impossible.



---

Copyright ©Tobii AB (publ). Not all products and services offered in each local market. Specifications are subject to change without prior notice. All trademarks are the property of their respective owners.

## Support for Your Tobii Pro Device

### Get Help Online

Visit Tobii Pro Connect for help with your Tobii Pro device. It contains the latest information about contacting Support, links to our Learning Center, and much more. Visit [connect.tobii.com](https://connect.tobii.com).

### Contact Your Solution Consultant or Reseller

For questions or problems with your product, contact your Tobii Pro sales representative or authorized reseller for assistance. They are most familiar with your personal setup and can best help you with tips and product training.

Visit [tobii.com/contact](https://tobii.com/contact).