EN
22

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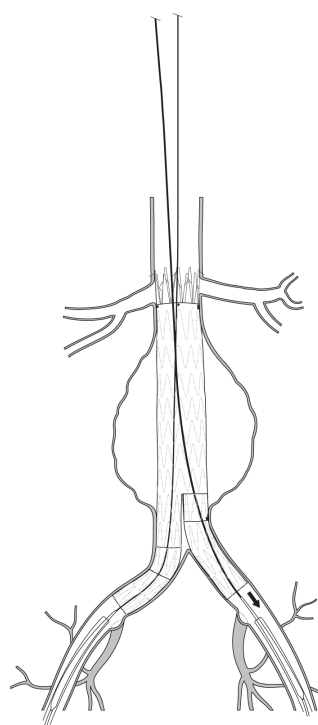
Fig. 45

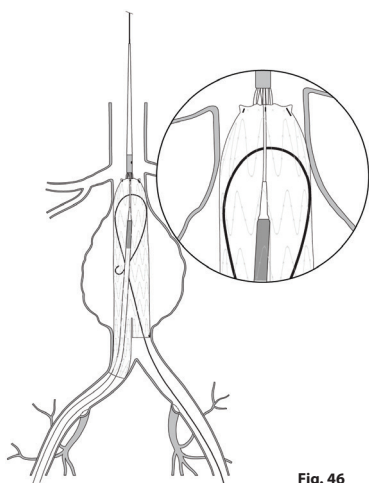
Fig. 46

ZENITH FLEX® AAA ENDOVASCULAR GRAFT WITH THE Z-TRAK™ INTRODUCTION SYSTEM

Read all instructions carefully. Failure to properly follow the instructions, warnings and precautions may lead to serious surgical consequences or injury to the patient.

CAUTION: Federal (U.S.A.) law restricts this device to sale by or on the order of a physician.

CAUTION: All contents of the outer pouch (including the introduction system and the endovascular grafts) are supplied sterile, for single use only.

For the Zenith product line there are four applicable Suggested Instructions for Use (IFU). This IFU describes the Suggested Instructions for Use for the Zenith Flex AAA Endovascular Graft (main body and iliac legs). For information on other Zenith components, please refer to the following Suggested Instructions for Use:

- Zenith AAA Endovascular Graft (Zenith AAA Endovascular Graft main body, iliac legs and occluder);
- Zenith AAA Endovascular Graft Ancillary Components (main body extension, iliac leg extension, converter and iliac plug);
- Zenith Renu® AAA Ancillary Graft (main body extension and converter configurations); and
- Coda® Balloon Catheter.

1 DEVICE DESCRIPTION

1.1 Aortic Main Body and Iliac Leg Components

The Zenith Flex AAA Endovascular Graft is a modular system consisting of three components, a bifurcated aortic main body and two iliac legs. **(Fig. 1)** The graft modules are constructed of full-thickness woven polyester fabric sewn to self-expanding stainless steel Cook-Z® stents with braided polyester and monofilament polypropylene suture. The modules are fully stented to provide stability and the expansile force necessary to open the lumen of the graft during deployment. Additionally, the Cook-Z stents provide the necessary attachment and seal of the graft to the vessel wall.

The bare suprarenal stent at the proximal end of the graft contains barbs that are placed at 3 mm increments for additional fixation of the device. To facilitate fluoroscopic visualization of the stent graft, gold radiopaque markers are positioned as follows: one on the lateral aspect of the most distal stent on the contralateral limb of the bifurcated section of the main body and four in a circumferential orientation within 2 mm of the most superior aspect of the graft material.

1.2 Main Body Delivery System

The Zenith Flex AAA Endovascular Graft main body is shipped preloaded onto the Z-Trak Introduction System. **(Fig. 2)** It has a sequential deployment method with built-in features to provide continuous control of the endovascular graft throughout the deployment procedure. The Z-Trak Introduction System is designed for precise positioning and allows readjustment of the final graft position before deployment of the barbed suprarenal stent.

The main body graft delivery system uses an 18, 20 or 22 French Z-Trak Introduction System. Dual trigger-wire release mechanisms lock the endovascular graft onto the delivery system until released by the physician. All systems are compatible with a .035 inch wire guide.

For added hemostasis, the Captor® Hemostatic Valve can be loosened or tightened for the introduction and/or removal of ancillary devices into and out of the sheath. The main body delivery systems feature a Flexor® introducer sheath which resists kinking and is hydrophilically coated. Both features are intended to enhance trackability in the iliac arteries and abdominal aorta.

1.3 Iliac Leg Delivery System

The Zenith Flex AAA Endovascular Graft iliac legs are shipped preloaded onto the Z-Trak Introduction System. **(Fig. 3)** The delivery system is designed for ease of use with minimal preparation. The iliac leg delivery system uses a 14 French or 16 French Z-Trak Introduction System. All systems are compatible with a .035 inch wire guide.

For added hemostasis, the Captor Hemostatic Valve can be loosened or tightened for the introduction and/or removal of ancillary devices into and out of the sheath.

1.4 Zenith AAA Endovascular Graft Ancillary Components

Additional ancillary endovascular components (main body extensions, iliac leg extensions, converters and iliac plugs) are available. **(Fig. 4)** Refer to the Zenith AAA Endovascular Graft Ancillary Components Instructions for Use for more information.

2 INDICATIONS FOR USE

The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is indicated for the endovascular treatment of patients with abdominal aortic or aortiliac aneurysms having morphology suitable for endovascular repair, including:

- Adequate iliac/femoral access compatible with the required introduction systems,
- Non-aneurysmal infrarenal aortic segment (neck) proximal to the aneurysm:
 - with a length of at least 15 mm,
 - with a diameter measured outer wall to outer wall of no greater than 32 mm and no less than 18 mm,
 - with an angle less than 60 degrees relative to the long axis of the aneurysm, and
 - with an angle less than 45 degrees relative to the axis of the suprarenal aorta.
- Iliac artery distal fixation site greater than 10 mm in length and 7.5-20 mm in diameter (measured outer wall to outer wall).

3 CONTRAINDICATIONS

The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is contraindicated in:

- Patients with known sensitivities or allergies to stainless steel, polyester, solder (tin, silver), polypropylene, or gold.
- Patients with a systemic infection who may be at increased risk of endovascular graft infection.

4 WARNINGS AND PRECAUTIONS

4.1 General

- Read all instructions carefully. Failure to properly follow the instructions, warnings and precautions may lead to serious consequences or injury to the patient.
- Always have a qualified surgery team available during implantation or reintervention procedures in the event that conversion to open surgical repair is necessary.

- The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction Systems should only be used by physicians and teams trained in vascular interventional techniques (catheter-based and surgical) and in the use of these devices. Specific training expectations are described in **Section 10.1, Physician Training**.
- Additional endovascular interventions or conversion to standard open surgical repair following initial endovascular repair should be considered for patients experiencing an enlarging aneurysm, unacceptable decrease in fixation length (vessel and component overlap) and/or endoleak. An increase in aneurysm size and/or persistent endoleak or migration may lead to aneurysm rupture.
- Patients experiencing reduced blood flow through the graft limb and/or leaks may be required to undergo secondary interventions or surgical procedures.

4.2 Patient Selection, Treatment and Follow-Up

- The Zenith Flex AAA Endovascular Graft is designed to treat aortic neck diameters no smaller than 18 mm and no larger than 32 mm. The Zenith Flex AAA Endovascular Graft is designed to treat proximal aortic necks (distal to the lowest renal artery) of at least 15 mm in length. Iliac artery distal fixation site greater than 10 mm in length and 7.5 - 20 mm in diameter (measured outer wall to outer wall) is required. These sizing measurements are critical to the performance of the endovascular repair.
- Key anatomical elements that may affect successful exclusion of the aneurysm include severe proximal neck angulation (>60 degrees for infrarenal neck to axis of AAA or >45 degrees for suprarenal neck relative to the immediate infrarenal neck); short proximal aortic neck (<15 mm); an inverted funnel shape (greater than 10% increase in diameter over 15 mm of proximal aortic neck length); and circumferential thrombus and/or calcification at the arterial implantation sites, specifically the proximal aortic neck and distal iliac artery interface. In the presence of anatomical limitations, a longer neck may be required to obtain adequate sealing and fixation. Irregular calcification and/or plaque may compromise the attachment and sealing at the fixation sites. Necks exhibiting these key anatomical elements may be more conducive to graft migration or endoleak.
- Adequate iliac or femoral access is required to introduce the device into the vasculature. Access vessel diameter (measured inner wall to inner wall) and morphology (minimal tortuosity, occlusive disease and/or calcification) should be compatible with vascular access techniques and delivery systems of a 16 French to 22 French vascular introducer sheath. Vessels that are significantly calcified, occlusive, tortuous or thrombus-lined may preclude placement of the endovascular graft and/or may increase the risk of embolization. A vascular conduit technique may be necessary to achieve success in some patients.
- The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is not recommended in patients who cannot tolerate contrast agents necessary for intraoperative and postoperative follow-up imaging. All patients should be monitored closely and checked periodically for a change in the condition of their disease and the integrity of the endoprosthesis.
- The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is not recommended in patients exceeding weight and/or size limits which compromise or prevent the necessary imaging requirements.
- Inability to maintain patency of at least one internal iliac artery or occlusion of an indispensable inferior mesenteric artery may increase the risk of pelvic/bowel ischemia.
- Multiple large, patent lumbar arteries, mural thrombus and a patent inferior mesenteric artery may all predispose a patient to Type II endoleaks. Patients with uncorrectable coagulopathy may also have an increased risk of Type II endoleak or bleeding complications.
- The safety and effectiveness of the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System has not been evaluated in the following patient populations:
 - traumatic aortic injury
 - leaking, pending rupture or ruptured aneurysms
 - mycotic aneurysms
 - pseudoaneurysms resulting from previous graft placement
 - revision of previously placed endovascular grafts
 - uncorrectable coagulopathy
 - indispensable mesenteric artery
 - genetic connective tissue disease (e.g., Marfans or Ehlers-Danlos Syndromes)
 - concomitant thoracic aortic or thoracoabdominal aneurysms
 - active systemic infections
 - pregnant or nursing females
 - morbidly obese patients
 - less than 18 years of age
 - patients with less than 15 mm in length or greater than 60 degrees angulation of the proximal aortic neck relative to the long axis of the aneurysm.
- Successful patient selection requires specific imaging and accurate measurements; please see **Section 4.3 Pre-Procedure Measurement Techniques and Imaging**.
- All lengths and diameters of the devices necessary to complete the procedure should be available to the physician, especially when preoperative case planning measurements (treatment diameters/lengths) are not certain. This approach allows for greater intraoperative flexibility to achieve optimal procedural outcomes.

4.3 Pre-Procedure Measurement Techniques and Imaging

- Lack of non-contrast CT imaging may result in failure to appreciate iliac or aortic calcification, which may preclude access or reliable device fixation and seal.
- Pre-procedure imaging reconstruction thicknesses >3 mm may result in sub-optimal device sizing, or in failure to appreciate focal stenoses from CT.
- Clinical experience indicates that contrast-enhanced spiral computed tomographic angiography (CTA) with 3-D reconstruction is the strongly recommended imaging modality to accurately assess patient anatomy prior to treatment with the Zenith Flex AAA Endovascular Graft. If contrast-enhanced spiral CTA with 3-D reconstruction is not available, the patient should be referred to a facility with these capabilities.
- Clinicians recommend positioning the x-ray C-arm during procedural angiography such that the origins of the renal arteries, and particularly the lowest patent renal artery, are well demonstrated prior to deployment of the proximal edge of the graft material (sealing stent) of the main body. Additionally, angiography should demonstrate the iliac artery bifurcations such that the distal common iliacs are well defined relative to the origin of the internal iliac arteries bilaterally, prior to deployment of the iliac leg components.

Diameters

Utilizing CT, diameter measurements should be determined from the outer wall to outer wall vessel diameter (not lumen measurement) to help with proper device sizing and device selection. The contrast-enhanced spiral CT scan must start 1 cm superior to the celiac axis and continue through the femoral heads at an axial thickness slice of 3 mm or less.

Lengths

Utilizing CT, length measurements should be determined to accurately assess infrarenal proximal neck length as well as planning main body sizes and leg components for the Zenith Flex AAA Endovascular Graft. These reconstructions should be performed in sagittal, coronal , and 3-D.

- **The long-term performance of endovascular grafts has not yet been established. All patients should be advised that endovascular treatment requires life-long, regular follow-up to assess their health and the performance of their endovascular graft.** Patients with specific clinical findings (e.g., endoleaks, enlarging aneurysm or changes in the structure or position of the endovascular graft) should receive enhanced follow-up. Specific follow-up guidelines are described in **Section 12, IMAGING GUIDELINES AND POSTOPERATIVE FOLLOW-UP.**
- The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is not recommended in patients unable to undergo, or who will not be compliant with, the necessary preoperative and postoperative imaging and implantation studies as described in **Section 12, IMAGING GUIDELINES AND POSTOPERATIVE FOLLOW-UP.**
- After endovascular graft placement, patients should be regularly monitored for perigraft flow, aneurysm growth or changes in the structure or position of the endovascular graft. At a minimum, annual imaging is required, including: 1) abdominal radiographs to examine device integrity (separation between components, stent fracture or barb separation) and 2) contrast and non-contrast CT to examine aneurysm changes, perigraft flow, patency, tortuosity and progressive disease. If renal complications or other factors preclude the use of image contrast media, abdominal radiographs and duplex ultrasound may provide similar information.

4.4 Device Selection

- Strict adherence to the Zenith Flex AAA Endovascular Graft IFU sizing guide is strongly recommended when selecting the appropriate device size (**Tables 10.5.1 through 10.5.2**). Appropriate device oversizing has been incorporated into the IFU sizing guide. Sizing outside of this range can result in endoleak, fracture, migration, device infolding or compression.

4.5 Implant Procedure

(Refer to **Section 11, DIRECTIONS FOR USE**)

- Appropriate procedural imaging is required to successfully position the Zenith Flex AAA Endovascular Graft and assure accurate apposition to the aortic wall.
- Do not bend or kink the delivery system. Doing so may cause damage to the delivery system and the Zenith Flex AAA Endovascular Graft.
- To avoid any twist in the endovascular graft, during any rotation of the delivery system, be careful to rotate all of the components of the system together (from outer sheath to inner cannula).
- Do not continue advancing any portion of the delivery system if resistance is felt during advancement of the wire guide or delivery system. Stop and assess the cause of resistance; vessel, catheter or graft damage may occur. Exercise particular care in areas of stenosis, intravascular thrombosis or in calcified or tortuous vessels.
- Inadvertent partial deployment or migration of the endoprosthesis may require surgical removal.
- Unless medically indicated, do not deploy the Zenith Flex AAA Endovascular Graft in a location that will occlude arteries necessary to supply blood flow to organs or extremities. Do not cover significant renal or mesenteric arteries (exception is the inferior mesenteric artery) with the endoprosthesis. Vessel occlusion may occur. During the clinical study, this device was not studied in patients with two occluded internal iliac arteries.
- Do not attempt to re-sheath the graft after partial or complete deployment.
- Repositioning the stent graft distally after partial deployment of the covered proximal stent may result in damage to the stent graft and/or vessel injury.
- Inaccurate placement and/or incomplete sealing of the Zenith Flex AAA Endovascular Graft within the vessel may result in increased risk of endoleak, migration or inadvertent occlusion of the renal or internal iliac arteries. Renal artery patency must be maintained to prevent/reduce the risk of renal failure and subsequent complications.
- Inadequate fixation of the Zenith Flex AAA Endovascular Graft may result in increased risk of migration of the stent graft. Incorrect deployment or migration of the endoprosthesis may require surgical intervention.
- Systemic anticoagulation should be used during the implantation procedure based on hospital and physician preferred protocol. If heparin is contraindicated, an alternative anticoagulant should be considered.
- To activate the hydrophilic coating on the outside of the Flexor introducer sheath, the surface must be wiped with sterile gauze pads soaked in saline solution. Always keep the sheath hydrated for optimal performance.
- Minimize handling of the constrained endoprosthesis during preparation and insertion to decrease the risk of endoprosthesis contamination and infection.
- Maintain wire guide position during delivery system insertion.
- Fluoroscopy should be used during introduction and deployment to confirm proper operation of the delivery system components, proper placement of the graft, and desired procedural outcome.
- The use of the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction Systems requires administration of intravascular contrast. Patients with pre-existing renal insufficiency may have an increased risk of renal failure postoperatively. Care should be taken to limit the amount of contrast media used during the procedure and to observe preventative methods of treatment to decrease renal compromise (e.g., adequate hydration).
- As the sheath and/or wire guide is withdrawn, anatomy and graft position may change. Constantly monitor graft position and perform angiography to check position as necessary.
- The Zenith Flex AAA Endovascular Graft incorporates a suprarenal stent with fixation barbs. Exercise extreme caution when manipulating interventional and angiographic devices in the region of the suprarenal stent.
- Use caution during manipulation of catheters, wires and sheaths within an aneurysm. Significant disturbances may dislodge fragments of thrombus, which can cause distal embolization or rupture of the aneurysm.
- Avoid damaging the graft or disturbing graft positioning after placement in the event reinstrumentation (secondary intervention) of the graft is necessary.
- Before deployment of the suprarenal stent, verify that the position of the access wire guide extends just distal to the aortic arch.
- Verify that the predetermined contralateral iliac leg is selected for insertion on the contralateral side of the patient before implantation.

4.6 Molding Balloon Use

- Do not inflate the balloon in vessel outside of the graft, as doing so may cause damage to the vessel. Use the balloon in accordance with its labeling.
- Use care in inflating the balloon within the graft in the presence of calcification, as excessive inflation may cause damage to the vessel.
- Confirm complete deflation of the balloon prior to repositioning.
- For added hemostasis, the Captor Hemostatic Valve can be loosened or tightened to accommodate the insertion and subsequent withdrawal of a molding balloon.

4.7 MRI Information

Non-clinical testing has demonstrated that the Zenith AAA Endovascular Graft is MR Conditional. A patient with this endovascular graft can be scanned safely immediately after placement under the following conditions:

Static Magnetic Field

- Static magnetic field of 3.0 Tesla or less.
- Highest spatial magnetic gradient field of 720 Gauss/cm

Non-clinical evaluation was conducted in a 3 Tesla MR system (General Electric Excite) with a maximum spatial magnetic gradient field of 720 Gauss/cm as measured with a gaussmeter in the position of the static magnetic field pertinent to the patient (i.e., outside of scanner covering, accessible to a patient or individual).

MRI-Related Heating

1.5 Tesla Systems:

- Static magnetic field of 1.5 Tesla
- Maximum whole-body-averaged specific absorption rate (SAR) of 2 W/kg for 15 minutes of scanning (i.e., per scanning sequence)

In non-clinical testing, the Zenith AAA Endovascular Graft produced a temperature rise of less than or equal to 1.4 °C at a MR system reported whole-body-averaged specific absorption rate (SAR) of 2.8 W/kg, for 15 minutes of MR scanning in a 1.5 Tesla Magnetom, Siemens Medical Magnetom, Numaris/4 Software, Version Syngo MR 2002B DHHS MR Scanner. The maximum whole-body-averaged specific absorption rate (SAR) was 2.8 W/kg, which corresponds to a calorimetry measured value of 1.5 W/kg.

3.0 Tesla Systems:

- Static magnetic field of 3.0 Tesla
- Maximum whole-body-averaged specific absorption rate (SAR) of 2 W/kg for 15 minutes of scanning (i.e., per scanning sequence)

In non-clinical testing, the Zenith AAA Endovascular Graft produced a temperature rise of less than or equal to 1.9 °C at a MR system reported whole-body-averaged specific absorption rate (SAR) of 3.0 W/kg, for 15 minutes of MR scanning in a 3.0 Tesla Excite, GE Electric Healthcare, G3.0-052B Software, MR Scanner. The maximum whole-body-averaged specific absorption rate (SAR) was 3.0 W/kg, which corresponds to a calorimetry measured value of 2.8 W/kg.

Image Artifact

The image artifact extends throughout the anatomical region containing the device, obscuring the view of immediately adjacent anatomical structures within approximately 20 cm of the device, as well as the entire device and its lumen, when scanned in nonclinical testing using the sequence: Fast spin echo, in a 3.0 Tesla, Excite, GE Electric Healthcare, with G3.0-052B Software, MR system with body radiofrequency coil.

For all scanners, the image artifact dissipates as the distance from the device to the area of interest increases. MR scans of the head and neck and lower extremities may be obtained without image artifact. Image artifact may be present in scans of the abdominal region and upper extremities, depending on distance from the device to the area of interest.

Clinical information is available for seventeen patients who received MRI scans after stent-graft implantation. There have been no reported adverse events or device problems in any of these patients as a result of having received an MRI. Additionally, there have been well over 100,000 Zenith AAA Endovascular Grafts implanted worldwide, in which there have been no reported adverse events or device problems as a result of MRI.

Cook recommends that the patient register the MR conditions disclosed in this IFU with the MedicAlert Foundation. The MedicAlert Foundation can be contacted in the following manners.

Mail:	MedicAlert Foundation International 2323 Colorado Avenue Turlock, CA 95382
Phone:	888-633-4298 (toll free) 209-668-3333 from outside the US
Fax:	209-669-2450
Web:	www.medicalert.org

5 ADVERSE EVENTS

5.1 Observed Adverse Events

A U.S. multicenter, prospective study of a previous version of the device (Zenith AAA Endovascular Graft) conducted at 15 centers which included 352 endovascular patients (200 standard risk, 100 high risk and 52 roll-in) and 80 control patients provides the basis of the observed adverse event rates in **Table 5.1.1**. Patients were enrolled in the standard risk arm if they were physiologically capable of withstanding an open or endovascular repair and had anatomy suitable for treatment with the Zenith AAA Endovascular Graft with the H&L-B One-Shot™ Introduction System (a previous version of the introduction system). Patients with suitable anatomy, but at higher risk of morbidity or mortality with open repair were enrolled into the high risk arm. Initial patients treated in the study were enrolled in the roll-in arm. The control group included patients whose vascular anatomy may not have been suitable for endovascular AAA repair.

Table 5.1.1 Death and Rupture from Clinical Study

Death and Rupture	Zenith Standard Risk ¹		Surgical Standard Risk		P value	Zenith High Risk		Zenith Roll-in	
All death (0-30 days) ²	0.5%	(1/199)	2.5%	(2/80)	.20	2.0%	(2/100)	1.9%	(1/52)
(31-365 days) ³	3.0%	(6/198)	1.3%	(1/78)	.68	7.1%	(7/98)	9.8%	(5/51)
AAA-related	0.0%	(0/198)	1.3%	(1/78)	.29	3.1%	(3/98)	0.0%	(0/51)
Non-AAA-related	3.0%	(6/198)	0.0%	(0/78)	.19	4.1%	(4/98)	9.8%	(5/51)
(0-365 days) ^{2,3,4}	3.5%	(7/199)	3.8%	(3/80)	>.99	9.0%	(9/100)	11.5%	(6/52)
AAA-related	0.5%	(1/199)	3.8%	(3/80)	.07	5.0%	(5/100)	1.9%	(1/52)
Non-AAA-related	3.0%	(6/199)	0.0%	(0/80)	.19	4.0%	(4/100)	9.6%	(5/52)
Rupture									
(0-30 days)	0.0%	(0/199)	n/a	n/a	n/a	0.0%	(0/100)	0.0%	(0/52)
(31-365 days)	0.0%	(0/198)	n/a	n/a	n/a	1.0%	(1/98)	0.0%	(0/51)
(0-365 days)	0.0%	(0/199)	n/a	n/a	n/a	1.0%	(1/100)	0.0%	(0/52)

¹Denominator of 199 because one standard risk patient did not receive a device.

²All deaths (0-30 days) were considered AAA and procedure related.

³Of the deaths (31-365 days), four were considered AAA related: 1 surgical (septic shock from ischemic colitis) and 3 high risk (pancreatitis with renal failure and sepsis, hemorrhage from upper abdominal aneurysm [not treated AAA] and multiple system failure).

⁴Of the deaths (0-365 days), ten were considered AAA related: 1 standard risk (cardiac failure), 3 surgical (massive hemorrhage, mesenteric ischemia and septic shock from ischemic colitis), 5 high risk (respiratory failure, cardiac failure with pulmonary embolism, pancreatitis with renal failure and sepsis, hemorrhage from upper abdominal aneurysm [not treated AAA] and multiple system failure) and 1 roll-in (suspected cardiac failure).

Table 5.1.2 Adverse Events¹ in Clinical Study

	Zenith Standard Risk		Surgical Standard Risk		P value	Zenith High Risk		Zenith Roll-in	
Freedom from Morbidity									
(0-30 days)	80%	(160/200)	58%	(46/80)	<.001	68%	(68/100)	73%	(38/52)
Cardiovascular ²	3.0%	(6/200)	11%	(9/80)	.02	14%	(14/100)	1.9%	(1/52)
Pulmonary ³	1.0%	(2/200)	15%	(12/80)	<.001	2.0%	(2/100)	0.0%	(0/52)
Renal ^{4, 9}	2.5%	(5/200)	10%	(8/80)	.01	6.0%	(6/100)	5.8%	(3/52)
Bowel ⁵	1.0%	(2/200)	3.8%	(3/80)	.14	1.0%	(1/100)	1.9%	(1/52)
Wound ⁶	4.5%	(9/200)	7.5%	(6/80)	.38	2.0%	(2/100)	3.8%	(2/52)
Neurologic ⁷	0.0%	(0/200)	2.5%	(2/80)	.08	0.0%	(0/100)	0.0%	(0/52)
Vascular ^{8,11}	11%	(21/200)	31%	(25/80)	<.001	20%	(20/100)	19%	(10/52)
Freedom from Morbidity									
(31-365 days)	91%	(181/198)	86%	(67/78)	.25	79%	(77/98)	86%	(44/51)
Cardiovascular ²	2.5%	(5/198)	3.8%	(3/78)	.69	5.1%	(5/98)	2.0%	(1/51)
Pulmonary ³	0.5%	(1/198)	1.3%	(1/78)	.49	4.1%	(4/98)	0.0%	(0/51)
Renal ^{4,10}	0.5%	(1/198)	0.0%	(0/78)	>.99	3.1%	(3/98)	0.0%	(0/51)
Bowel ⁵	0.5%	(1/198)	1.3%	(1/78)	.49	0.0%	(0/98)	0.0%	(0/51)
Wound ⁶	2.0%	(4/198)	5.1%	(4/78)	.23	3.1%	(3/98)	2.0%	(1/51)
Neurologic ⁷	1.0%	(2/198)	0.0%	(0/78)	>.99	1.0%	(1/98)	3.9%	(2/51)
Vascular ⁸	3.0%	(6/198)	3.8%	(3/78)	.72	8.2%	(8/98)	5.9%	(3/51)
Freedom from Morbidity									
(0-365 days)	76%	(151/200)	49%	(39/80)	<.001	55%	(55/100)	62%	(32/52)
Cardiovascular ²	5.0%	(10/200)	14%	(11/80)	.02	19%	(19/100)	3.8%	(2/52)
Pulmonary ³	1.5%	(3/200)	16%	(13/80)	<.001	6.0%	(6/100)	0.0%	(0/52)
Renal ^{4,9,10}	2.5%	(5/200)	10%	(8/80)	.01	9.0%	(9/100)	5.8%	(3/52)
Bowel ⁵	1.5%	(3/200)	3.8%	(3/80)	.36	1.0%	(1/100)	1.9%	(1/52)
Wound ⁶	5.5%	(11/200)	13%	(10/80)	.08	5.0%	(5/100)	5.8%	(3/52)
Neurologic ⁷	1.0%	(2/200)	2.5%	(2/80)	.32	1.0%	(1/100)	3.8%	(2/52)
Vascular ^{8,11}	12%	(24/200)	33%	(26/80)	<.001	25%	(25/100)	23%	(12/52)

¹From the morbidity index.

²Cardiovascular included: Q-wave and non-Q-wave myocardial infarctions, congestive heart failure, arrhythmias requiring new medication or treatment, cardiac ischemia requiring intervention, inotropic support, medically intractable hypertension.

³Pulmonary included: reintubation or ventilation >24 hours, pneumonia requiring antibiotics, supplemental oxygen at discharge.

⁴Renal included: dialysis in patients with normal preoperative renal function, creatinine rise >30% from baseline on two or more follow-up tests.

⁵Bowel included: bowel obstruction, bowel ischemia, aorto-enteric fistula, paralytic ileus >4 days.

⁶Wound included: infection requiring antibiotic treatment, hernia, lymph fistula, dehiscence, necrosis requiring debridement.

⁷Neurologic included: stroke, TIA, spinal cord ischemia/paralysis.

⁸Vascular included: limb thrombosis, distal embolization resulting in tissue loss or requiring intervention, transfusion post-procedure (resulting from pseudoaneurysm, vascular injury, aneurysm leak or other procedure-related causes), pseudoaneurysm, vascular injury (such as inadvertent occlusion, dissection or other procedure related causes), aneurysm leak or rupture, increase in aneurysm size by more than 0.5 cm relative to the smallest of any prior measurement.

⁹Investigators reported one additional high risk patient to have occlusion of an accessory renal artery and one additional high risk patient to have chronic renal insufficiency as “other” adverse events.

¹⁰Investigators reported one additional roll-in patient to have renal insufficiency and one additional surgical patient to have renal insufficiency as “other” adverse events.

¹¹Investigators reported one additional roll-in patient to have experienced intraoperative aortic plaque rupture resulting in renal artery occlusion as an “other” adverse event.

5.2 Potential Adverse Events

Adverse events that may occur and/or require intervention include, but are not limited to:

- Amputation
- Anesthetic complications and subsequent attendant problems (e.g., aspiration)
- Aneurysm enlargement
- Aneurysm rupture and death
- Aortic damage, including perforation, dissection, bleeding, rupture and death
- Arterial or venous thrombosis and/or pseudoaneurysm
- Arteriovenous fistula
- Bleeding, hematoma or coagulopathy
- Bowel complications (e.g., ileus, transient ischemia, infarction, necrosis)
- Cardiac complications and subsequent attendant problems (e.g., arrhythmia, myocardial infarction, congestive heart failure, hypotension, hypertension)
- Claudication (e.g., buttock, lower limb)
- Death
- Edema
- Embolization (micro and macro) with transient or permanent ischemia or infarction
- Endoleak
- Endoprosthesi: improper component placement; incomplete component deployment; component migration; suture break; occlusion; infection; stent fracture; graft material wear; dilatation; erosion; puncture; perigraft flow; barb separation and corrosion
- Fever and localized inflammation
- Genitourinary complications and subsequent attendant problems (e.g., ischemia, erosion, fistula, incontinence, hematuria, infection)
- Graft or native vessel occlusion
- Hepatic failure

- Impotence
- Infection of the aneurysm, device or access site, including abscess formation, transient fever and pain
- Lymphatic complications and subsequent attendant problems (e.g., lymph fistula)
- Neurologic local or systemic complications and subsequent attendant problems (e.g., stroke, transient ischemic attack, paraplegia, paraparesis, paralysis)
- Pulmonary/respiratory complications and subsequent attendant problems (e.g., pneumonia, respiratory failure, prolonged intubation)
- Renal complications and subsequent attendant problems (e.g., artery occlusion, contrast toxicity, insufficiency, failure)
- Surgical conversion to open repair
- Vascular access site complications, including infection, pain, hematoma, pseudoaneurysm, arteriovenous fistula
- Vessel damage
- Wound complications and subsequent attendant problems (e.g., dehiscence, infection)
- Vascular spasm or vascular trauma (e.g., iliofemoral vessel dissection, bleeding, rupture, death)

Device-Related Adverse Event Reporting

Any adverse event (clinical incident) involving the Zenith Flex AAA Endovascular Graft should be reported to COOK immediately. For customers in the United States, to report an incident, call the Customer Relations Department at 1-800-457-4500 (24 hour) or 1-812-339-2235. For customers outside the United States, please call your distributor.

6 SUMMARY OF CLINICAL STUDIES

6.1 Objectives

The primary objective of the clinical study was to evaluate the safety and effectiveness of a previous version of the device (Zenith AAA Endovascular Graft) as an alternative to open surgical repair in the primary treatment of infrarenal abdominal aortic aneurysms. Study hypotheses examined whether standard risk patients experienced less 30-day morbidity, non-inferior 30-day survival rates, non-inferior 12-month survival rates, non-inferior 12-month treatment success and improved clinical utility measures compared to surgical control patients. Safety was determined by evaluating whether the Zenith AAA Endovascular Graft subjects would have reduced 30-day morbidity, non-inferior 30-day and 12-month survival and non-inferior 12-month treatment success compared to the subjects treated with open surgical treatment. Effectiveness was based on exclusion of the aneurysm including the absence of any endoleak, the absence of aneurysm enlargement (≥5 mm) and the absence of major device adverse events evaluated through one year follow-up. Secondary objectives included an assessment of clinical benefit and quality-of-life measures.

6.2 Study Design

The U.S. clinical study was a multicenter, non-randomized study comparing standard medical risk patients who received an endovascular graft to an open surgical control. There were two additional study arms for high medical risk and roll-in treatment groups. Fifteen centers enrolled 200 standard risk, 80 surgical control, 100 high risk and 52 roll-in patients. The control group included patients whose vascular anatomy may not have been suitable for endovascular AAA repair. Follow-up evaluations were scheduled for pre-discharge, 1 month, 6 months, 12 months and 24 months. Patient follow-up and accountability at 1 month, 12 months and 24 months are presented in **Table 6.2.1**, as these were the primary data analysis time points. Imaging

data provided in this summary are based on findings from an independent centralized image analysis laboratory (Core Lab), which reviewed CT scans and abdominal X-rays to assess aneurysm diameter changes, device and relative component migration, device integrity (wire and graft) and the presence and type of endoleaks. Clinical events were adjudicated by an independent clinical events committee, and safety was monitored by a data safety monitoring committee.

Surgical and Zenith standard risk patients met identical pathophysiologic risk criteria. The endovascular groups excluded circumferential thrombus in the proximal neck, proximal neck less than 15 mm in length, outer wall to outer wall proximal neck diameter less than 18 mm or greater than 28 mm, severe proximal neck angulation, outer wall to outer wall iliac artery diameter less than 7.5 mm or greater than 20 mm at distal fixation site or iliac artery distal fixation site less than 10 mm in length.

Patients were considered at higher risk for surgical repair if they had age greater than 80, baseline creatinine >2.0 mg/dl, home oxygen therapy, FEV₁ <1 liter, ejection fraction <25%, disabling COPD, New York Heart Classification 3 or 4, hostile abdomen, dialysis, MI within last 6 months, medically intractable hypertension, previous stroke with residual deficit, cultural objection to receipt of blood or blood products, previous renal bypass surgery or inflammatory aneurysm.

Before enrolling patients into the pivotal trial, centers without Zenith AAA Endovascular Graft experience were required to treat initial patients under the supervision of a proctor. These roll-in patients were a combination of standard and high risk patients and were followed according to the same schedule as the patients in the pivotal trial.

Table 6.2.1 Patient Follow-Up and Accountability¹

Treatment Interval	Zenith Standard Risk			Surgical Standard Risk		
	1 mo.	12 mo.	24 mo.	1 mo.	12 mo.	24 mo.
No device	1	1²	1²	0	0	n/a
Conversion to open repair	0	2²	2²	n/a	n/a	n/a
Expired	1	7	17	2	3	n/a
Withdrawn/lost to follow-up	0	0	2	0	4	n/a
Available	198	190	178	78	73	n/a
Site CT imaging	191	168	110	69	58	n/a
Core lab CT imaging	190	165	99	69	59	n/a
Site KUB imaging	179	153	108	n/a	n/a	n/a
Core lab KUB imaging	178	149	93	n/a	n/a	n/a
Site evaluated for endoleak	187	163	107	n/a	n/a	n/a
Core lab evaluated for endoleak	161	148	92	n/a	n/a	n/a
Site evaluated for aneurysm enlargement	n/a	149	104	n/a	n/a	n/a
Core lab evaluated for aneurysm enlargement	n/a	151	94	n/a	n/a	n/a

¹Data analysis sample size varies for each of the time points above and in the following tables. This variability is due to patient availability for follow-up, as well as quantity and quality of images available from specific time points for evaluation. For example, the number and quality of images available for evaluation of endoleak at 12 months is different than the number and quality of images available at 24 months due to variation in the number of image exams performed, the number of images provided from the clinical site to the Core Lab and/or the number of images with acceptable evaluation quality. Totals at time points are not cumulative, unless otherwise noted.

²Totals at time points are cumulative.

6.3 Patient Demographics

Tables 6.3.1 and **6.3.2** compare the subject characteristics and initial aneurysm diameter of the Zenith AAA Endovascular Graft and open surgical population, respectively.

Table 6.3.1 Comparison of Subject Characteristics

Item	Zenith Standard Risk		Surgical Standard Risk		P value	Zenith High Risk		Zenith Roll-in	
Age (years)	71 ± 7		69 ± 7		.03	77 ± 7		74 ± 8	
Gender male	94%	(187/200)	89%	(71/80)	.22	92%	(92/100)	90%	(47/52)
Current medical conditions									
Peripheral vascular disease	16%	(31/195)	25%	(19/76)	.12	24%	(23/96)	9.6%	(5/52)
Hypertension	64%	(127/200)	83%	(65/78)	.001	68%	(67/99)	67%	(35/52)
Renal failure	0.0%	(0/197)	0.0%	(0/79)	>.99	5.2%	(5/97)	1.9%	(1/52)
COPD	20%	(39/199)	18%	(14/78)	.87	34%	(33/98)	22%	(11/51)
Thromboembolic event	4.5%	(9/199)	7.7%	(6/78)	.38	7.1%	(7/99)	1.9%	(1/52)
Liver disease	2.1%	(4/192)	5.1%	(4/79)	.24	1.0%	(1/99)	1.9%	(1/52)
Diabetes mellitus	12%	(24/199)	15%	(12/79)	.55	17%	(17/99)	14%	(7/51)
Insulin-dependent	17%	(4/24)	8.3%	(1/12)	.65	24%	(4/17)	43%	(3/7)
Previous medical conditions									
MI	39%	(74/192)	29%	(23/80)	.13	35%	(34/98)	35%	(18/52)
Congestive heart failure	5.0%	(10/199)	12%	(9/78)	.07	16%	(16/100)	10%	(5/50)
Angina	49%	(98/198)	39%	(31/79)	.14	45%	(44/98)	44%	(23/52)
Arrhythmia	20%	(40/197)	22%	(17/78)	.87	28%	(27/98)	24%	(12/51)
Cerebrovascular disease	9.5%	(19/199)	16%	(13/79)	.14	20%	(20/99)	9.8%	(5/51)
Systemic infection	1.0%	(2/196)	0.0%	(0/78)	>.99	3.1%	(3/97)	0.0%	(0/49)
Cancer	22%	(43/200)	19%	(15/80)	.74	31%	(31/99)	29%	(15/51)
Family history of aneurysmal disease	16%	(24/150)	27%	(17/63)	.09	14%	(11/77)	26%	(10/38)
Previous surgery at site	10%	(20/200)	15%	(12/79)	.22	10%	(10/99)	14%	(7/51)
Previous radiation at site	0.5%	(1/197)	0.0%	(0/79)	>.99	2.0%	(2/100)	2.0%	(1/51)
Excessive alcohol use	3.6%	(7/193)	10%	(8/77)	.04	3.1%	(3/96)	4.0%	(2/50)
Tobacco use									
Never smoked	10%	(20/193)	5.0%	(4/80)		14%	(13/96)	24%	(12/49)
Past smoker	69%	(133/193)	60%	(48/80)	.03	69%	(66/96)	57%	(28/49)
Still smokes	21%	(40/193)	35%	(28/80)		18%	(17/96)	18%	(9/49)

Due to inclusion criteria, high risk patients were older ($P<.001$), had more renal failure ($P=.004$), COPD ($P=.01$), congestive heart failure ($P=.004$) and cerebrovascular disease ($P=.02$) than standard risk patients.

Table 6.3.2 Aneurysm Diameter Distribution

Diameter Range	Zenith Standard Risk		Surgical Standard Risk		Zenith High Risk		Zenith Roll-in	
<30 mm	0.0%	(0/199)	0.0%	(0/78)	0.0%	(0/100)	0.0%	(0/52)
30-39 mm	0.5%	(1/199)	0.0%	(0/78)	1.0%	(1/100)	0.0%	(0/52)
40-49 mm	23%	(45/199)	7.7%	(6/78)	15%	(15/100)	13%	(7/52)
50-59 mm	48%	(95/199)	33%	(26/78)	47%	(47/100)	40%	(21/52)
60-69 mm	24%	(47/199)	29%	(23/78)	27%	(27/100)	42%	(22/52)
70-79 mm	3.0%	(6/199)	21%	(16/78)	5.0%	(5/100)	1.9%	(1/52)
80-89 mm	2.5%	(5/199)	6.4%	(5/78)	1.0%	(1/100)	0.0%	(0/52)
≥90 mm	0.0%	(0/199)	2.6%	(2/78)	1.0%	(1/100)	0.0%	(0/52)

Aneurysm diameter distribution was not assessed in three high-risk and one roll-in patient.

6.4 Results

Data gathered in **Tables 6.4.1** through **6.4.3** were collected by the clinical study sites and Core Lab. Where available, 24-month data are provided. Control patients were not followed beyond 12 months and some data have not yet been adjudicated beyond 12 months. Therefore, some results are presented to 12 months while other results are presented to 24 months in this section. **Table 6.4.1** describes the devices implanted in clinical study patients. **Table 6.4.2** describes the primary results from the clinical study. **Figs. 6.4.1** and **6.4.2** are the Kaplan-Meier plots of all-cause and AAA-related survival to 24 months, respectively. An independent clinical events committee adjudicated all deaths for possible relationship to aneurysm repair. All early deaths (0-30 days) were considered AAA-related. Deaths after 30 days were considered AAA-related if AAA disease or device involvement was confirmed. **Table 6.4.3** presents Success Measures and **Fig. 6.4.3** is the Kaplan-Meier plot of Freedom from Morbidity.

Table 6.4.1 Devices Implanted

Item	Zenith Standard Risk		Zenith High Risk		Zenith Roll-in	
Main body and legs	99.5%	(199/200)*	100%	(100/100)†	100%	(52/52)
Main body extension²	1.5%	(3/199)	1.0%	(1/100)	5.8%	(3/52)
Ipsilateral iliac leg extension³,⁵	9.5%	(19/199)	11%	(11/100)	0.0%	(0/52)
Contralateral iliac leg extension⁴,⁵	11%	(21/199)	11%	(11/100)	13.5%	(7/52)
Converter	0.5%	(1/199)**	0.0%	(0/100)	0.0%	(0/52)
Occluder	0.0%	(0/199)	0.0%	(0/100)	0.0%	(0/52)

*One standard risk patient did not receive a device due to tortuosity and calcification of the access vessel.

**Converter was used without occluder.

†One device was custom.

²Two standard risk and one high risk patient received main body extensions post-procedure; one standard risk patient received two main body extensions.

³Two standard risk and one high risk patient received ipsilateral leg extensions post-procedure.

⁴Four standard risk, two high risk and one roll-in patient received contralateral leg extensions post-procedure.

⁵Three standard risk and three high risk patients received both ipsilateral and contralateral extensions during the procedure; one standard risk received both ipsilateral and contralateral extensions post-procedure.

Table 6.4.2 Primary Results

Item	Zenith Standard Risk¹		Surgical Standard Risk		P value	Zenith High Risk		Zenith Roll-in	
All death (0-30 days)²	0.5%	(1/199)	2.5%	(2/80)	.20	2.0%	(2/100)	1.9%	(1/52)
All death									
(31-365 days)²,³	3.0%	(6/199)	1.3%	(1/80)	.68	7.0%	(7/100)	9.6%	(5/52)
AAA-related	0.0%	(0/199)	1.3%	(1/80)	.29	3.0%	(3/100)	0.0%	(0/52)
Non-AAA-related	3.0%	(6/199)	0.0%	(0/80)	.19	4.0%	(4/100)	9.6%	(5/52)
All death									
(0-365 days)²,³,⁴	3.5%	(7/199)	3.8%	(3/80)	>.99	9.0%	(9/100)	11.5%	(6/52)
AAA-related	0.5%	(1/199)	3.8%	(3/80)	.07	5.0%	(5/100)	1.9%	(1/52)
Non-AAA-related	3.0%	(6/199)	0.0%	(0/80)	.19	4.0%	(4/100)	9.6%	(5/52)
Rupture									
(0-30 days)	0.0%	(0/199)	n/a	n/a	n/a	0.0%	(0/100)	0.0%	(0/52)
(31-365 days)	0.0%	(0/199)	n/a	n/a	n/a	1.0%	(1/100)	0.0%	(0/52)
(0-365 days)	0.0%	(0/199)	n/a	n/a	n/a	1.0%	(1/100)	0.0%	(0/52)
Conversion									
(0-30 days)	0.0%	(0/199)	n/a	n/a	n/a	0.0%	(0/100)	0.0%	(0/52)
(31-365 days)⁵	1.0%	(2/199)	n/a	n/a	n/a	1.0%	(1/100)	0.0%	(0/52)
(0-365 days)⁵	1.0%	(2/199)	n/a	n/a	n/a	1.0%	(1/100)	0.0%	(0/52)
Adverse events									
(0-30 days)⁶	20%	(40/200)	43%	(34/80)	<.001	32%	(32/100)	27%	(14/52)
(31-365 days)⁶	8.6%	(17/199)	14%	(11/78)	.25	21%	(21/98)	14%	(7/51)
(0-365 days)⁶	25%	(49/200)	51%	(41/80)	<.001	45%	(45/100)	38%	(20/52)

¹Denominator of 199 because one standard risk patient did not receive a device.

²All deaths (0-30 days) were considered AAA and procedure related.

³Of the deaths (31-365 days), four were considered AAA related: 1 surgical (septic shock from ischemic colitis) and 3 high risk (pancreatitis with renal failure and sepsis, hemorrhage from upper abdominal aneurysm [not treated AAA] and multiple system failure).

⁴Of the deaths (0-365 days), ten were considered AAA related: 1 standard risk (cardiac failure), 3 surgical (massive hemorrhage, mesenteric ischemia and septic shock from ischemic colitis), 5 high risk (respiratory failure, cardiac failure with pulmonary embolism, pancreatitis with renal failure and sepsis, hemorrhage from upper abdominal aneurysm [not treated AAA] and multiple system failure) and 1 roll-in (suspected cardiac failure).

⁵Standard risk patients underwent conversions due to a persistent, proximal Type I endoleak and a new suprarenal aortic aneurysm. Three surgical patients had massive hemorrhages, of which 2 required re-operation and one died.

⁶Adverse events included in the morbidity index.

Zenith AAA Endovascular Graft patients exhibited no significant differences between males and females for survival and freedom from major adverse events. (Error bars in **Figs. 6.4.1**, **6.4.2** and **6.4.3** represent 95% confidence limits.) **Fig. 6.4.1** presents all-cause survival to 24 months. The accompanying table presents the Kaplan-Meier analysis at 1, 6, 12 and 24 months.

Fig. 6.4.1 Survival at 24 Months


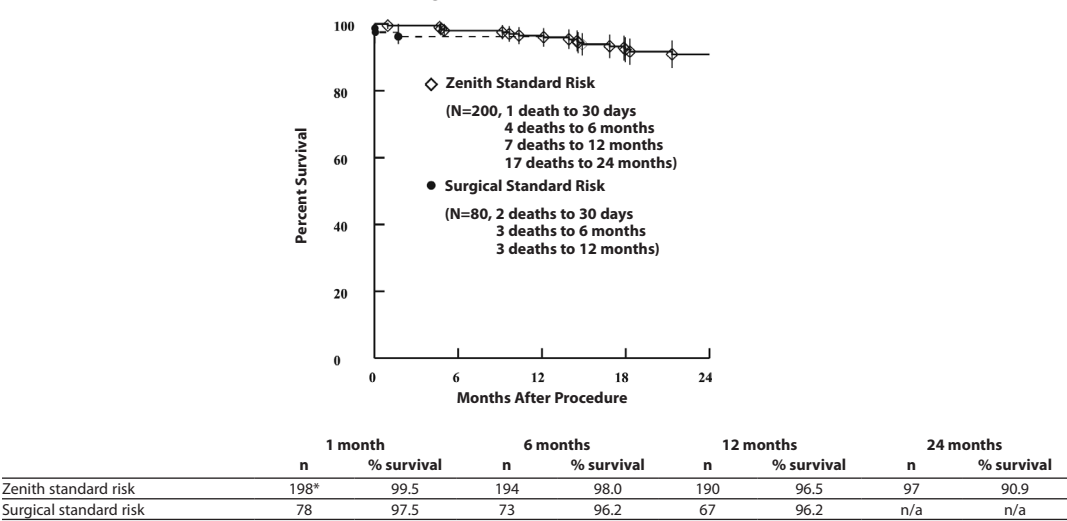
Fig. 6.4.2 presents AAA-related survival (determined by Clinical Events Committee) to 24 months. The accompanying table presents the Kaplan-Meier analysis at 1, 6, 12 and 24 months.

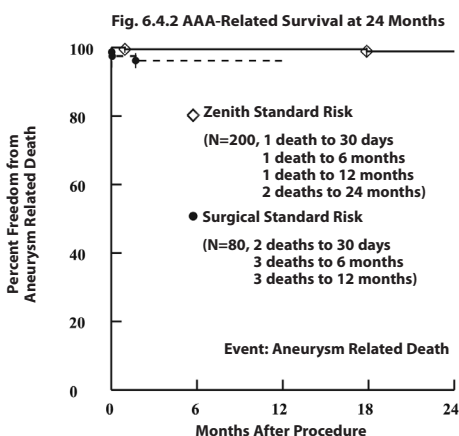	1 month		6 months		12 months		24 months	
	n	% survival	n	% survival	n	% survival	n	% survival
Zenith standard risk	198*	99.5	194	99.5	190	99.5	97	99.0
Surgical standard risk	78	97.5	73	96.2	67	96.2	n/a	n/a

n = Patients alive and available for follow-up at the end of the interval

P = .04

*One patient expired before 1 month and one patient did not receive a device.

Table 6.4.3 Success Measures


Item	Zenith Standard Risk		Surgical Standard Risk		Zenith High Risk		Zenith Roll-in	
Technical success ¹	99.5%	(199/200)	98.8%	(79/80)	100%	(100/100)	100%	(52/52)
Procedural success at 30 days ²	95.1%	(155/163)	88%	(60/68)	86%	(70/81)	91%	(30/33)
Treatment success at 12 months ³	89%	(122/137)	85%	(52/61)	70%	(44/63)	87%	(26/30)

¹Patent graft following deployment.

²Technical success with no major complications, patent graft and no Type I or Type III endoleaks at 30 days.

³Procedural success extended to 12 months with no aneurysm enlargement (>5 mm).

Fig. 6.4.3 presents freedom from morbidity (events in the morbidity index) to 12 months. The accompanying table presents the Kaplan-Meier analysis at 1, 6 and 12 months.

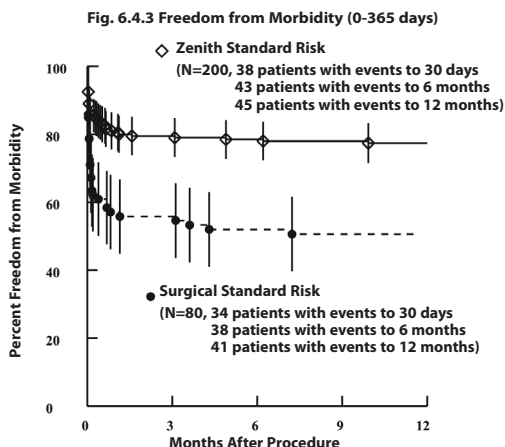	1 month		6 months		12 months	
	n	%	n	%	n	%
Zenith standard risk	162	81.0	154	78.5	133	77.4
Surgical standard risk	45	57.1	39	52.0	33	47.8

n = Patients alive and free of morbidity at the end of the interval

P = <.001

Tables 6.4.4 through **6.4.7** describe results of the Zenith AAA Endovascular Graft subjects as reported by the Core Lab. Device performance factors analyzed by the Core Lab include device integrity (**Table 6.4.4**), device patency (**Table 6.4.5**), migration (**Table 6.4.6**) and limb separation (**Table 6.4.7**).

Table 6.4.4 Abdominal Radiographic Findings – Device Integrity

Item	Zenith Standard Risk		Zenith High Risk		Zenith Roll-in	
Stent Fractures ¹						
Pre-discharge	0.0%	(0/172)	0.0%	(0/81)	0.0%	(0/39)
30 day	0.0%	(0/172)	0.0%	(0/83)	0.0%	(0/43)
6 month	0.0%	(0/166)	0.0%	(0/78)	0.0%	(0/35)
12 month	0.0%	(0/148)	0.0%	(0/60)	0.0%	(0/28)
24 month	0.0%	(0/93)	0.0%	(0/42)	0.0%	(0/19)
Barb Separation ²						
Pre-discharge	0.0%	(0/176)	0.0%	(0/86)	0.0%	(0/39)
30 day	0.0%	(0/178)	0.0%	(0/86)	0.0%	(0/43)
6 month	1.2%	(2/167)	2.5%	(2/80)	0.0%	(0/35)
12 month	2.0%	(3/149)	1.7%	(1/60)	0.0%	(0/28)
24 month	1.1%	(1/93)	0.0%	(0/42)	0.0%	(0/19)
Graft material rupture						
Pre-discharge	0.0%	(0/176)	0.0%	(0/86)	0.0%	(0/39)
30 day	0.0%	(0/178)	0.0%	(0/86)	0.0%	(0/43)
6 month	0.0%	(0/167)	0.0%	(0/80)	0.0%	(0/35)
12 month	0.0%	(0/149)	0.0%	(0/60)	0.0%	(0/28)
24 month	0.0%	(0/93)	0.0%	(0/42)	0.0%	(0/19)

¹Stent fracture percentages are for main body. There were also no right iliac leg, left iliac leg, occluder, converter, left iliac extension, right iliac extension or main body extension fractures observed by the Core Lab.

²Patients with separation of 1 or 2 barbs (of 10 or 12 total); no adverse clinical sequelae.

Table 6.4.5 CT Findings – Graft Patency							
Item	Zenith Standard Risk			Zenith High Risk		Zenith Roll-in	
Graft patency							
30 day	100%	(185/185)		99%	(85/86)	100%	(47/47)
6 month	99%	(183/184)		100%	(74/74)	100%	(39/39)
12 month	99%	(153/155)		100%	(62/62)	100%	(30/30)
24 month	100%	(96/96)		100%	(33/33)	100%	(25/25)

Table 6.4.6 CT Findings – Graft (Main Body) Migration							
Item	Zenith Standard Risk			Zenith High Risk		Zenith Roll-in	
Graft migration (>5 mm) at 12 months							
with clinical sequelae ¹ or intervention	0.0%	(0/162)		0.0%	(0/71)	0.0%	(0/34)
without clinical sequelae ¹ or intervention	2.5%	(4/162)		2.8%	(2/71)	0.0%	(0/34)
Graft migration (>10 mm)	0.0%	(0/162)		0.0%	(0/71)	0.0%	(0/34)

¹Migration with clinical sequelae would include endoleak, conversion, rupture or AAA-related death.

Table 6.4.7 Abdominal Radiograph Findings – Limb Separation							
Item	Zenith Standard Risk			Zenith High Risk		Zenith Roll-in	
Limb separation							
Pre-discharge	0.0%	(0/176)		0.0%	(0/86)	0.0%	(0/39)
30 day	0.0%	(0/178)		0.0%	(0/86)	0.0%	(0/43)
6 month	0.0%	(0/167)		0.0%	(0/80)	0.0%	(0/35)
12 month	0.0%	(0/149)		0.0%	(0/60)	0.0%	(0/28)
24 month	0.0%	(0/93)		0.0%	(0/42)	0.0%	(0/19)

6.5 Endoleak Management

During the clinical study Type I endoleaks were treated during the initial procedure by use of additional balloon seating or if unsuccessful, additional prostheses. Type II endoleaks were observed for a period of one to six months to determine if they would spontaneously thrombose, or in the absence of enlarging aneurysms, they were treated with endovascular techniques at the discretion of the practicing physician. If the aneurysm enlarged, treatment by embolization or ligation was considered and, in some cases, performed. Type III endoleaks caused by graft defects, inadequate seal or disconnection of the modular components were treated with additional ballooning or prostheses. As reported by the angiographic core lab, there were no Type IV endoleaks during the U.S. clinical study. The graft material used to manufacture the Zenith AAA Endovascular Graft is of standard thickness and is the same material used in open surgical procedures. **Table 6.5.1** presents the incidence of endoleaks by evaluation interval, as identified by the Core Lab for the standard risk, high risk and roll-in patients, respectively.

Table 6.5.1 Endoleaks (All Types, New and Persistent)							
Item	Zenith Standard Risk			Zenith High Risk		Zenith Roll-in	
Endoleaks							
Pre-discharge	15%	(23/153)		14%	(11/78)	12%	(3/26)
30 day ¹	9.9%	(16/161)		12%	(9/75)	6.3%	(2/32)
6 month ¹	8.7%	(15/172)		11%	(8/70)	8.6%	(3/35)
12 month ¹	7.4%	(11/148)		8.8%	(5/57)	3.4%	(1/29)

¹Includes both persistent endoleaks and new observations.

Tables 6.5.2 – 6.5.4 present the incidence of first occurrence of an endoleak according to evaluation interval, as identified by the Core Lab at or before the 30-day, 6-month and 12-month exams for the standard risk, high risk and roll-in patients, respectively. The number of patients who are leak-free thereafter is also given.

Table 6.5.2 First Occurrence of Endoleak ¹ for Standard Risk Patients									
Item	One-Month Exam N=179			Six-Month Exam N=172			Twelve-Month Exam ³ N=148		
	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²
Endoleaks	17	31	17	2.3	4	3	3.4	5	2
Proximal Type I	2.8	5	4	0.0	0	0	0.0	0	0
Distal									
Type I	1.7	3	1	0.0	0	0	0.7	1	1
Type II	9.5	17	9	2.3	4	3	1.4	2	1
Type III	1.1	2	2	0.0	0	0	0.7	1	0
Type IV	0.0	0	0	0.0	0	0	0.0	0	0
Multiple	1.1	2	1	0.0	0	0	0.0	0	0
Unknown	1.1	2	0	0.0	0	0	0.7	1	0

¹Identified by Core Lab.

²Subsequent endoleaks may have been of different type than original.

³Only 2 patients had new endoleaks after 12 months; follow-up after 24 months not available.

Table 6.5.3 First Occurrence of Endoleak ¹ for High Risk Patients									
Item	One-Month Exam N=88			Six-Month Exam N=70			Twelve-Month Exam ³ N=57		
	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²
Endoleaks	18	16	6	2.9	2	0	3.5	2	1
Proximal Type I	2.3	2	1	0.0	0	0	0.0	0	0
Distal									
Type I	1.1	1	1	0.0	0	0	0.0	0	0
Type II	9.1	8	3	1.4	1	0	1.8	1	0
Type III	0.0	0	0	1.4	1	0	1.8	1	1
Type IV	0.0	0	0	0.0	0	0	0.0	0	0
Multiple	4.5	4	0	0.0	0	0	0.0	0	0
Unknown	1.1	1	1	0.0	0	0	0.0	0	0

¹Identified by Core Lab.

²Subsequent endoleaks may have been of different type than original.

³No endoleaks after 12 months; follow-up after 24 months not available.

Table 6.5.4 First Occurrence of Endoleak ¹ for Roll-In Patients									
Item	One-Month Exam N=36			Six-Month Exam N=35			Twelve-Month Exam ³ N=29		
	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²	%	Endoleak ¹	Leak-free Thereafter ²
Endoleaks	11	4	2	2.9	1	0	0.0	0	0
Proximal Type I	0.0	0	0	0.0	0	0	0.0	0	0
Distal									
Type I	0.0	0	0	0.0	0	0	0.0	0	0
Type II	5.6	2	1	2.9	1	0	0.0	0	0
Type III	2.8	1	0	0.0	0	0	0.0	0	0
Type IV	0.0	0	0	0.0	0	0	0.0	0	0
Multiple	0.0	0	0	0.0	0	0	0.0	0	0
Unknown	2.8	1	1	0.0	0	0	0.0	0	0

¹Identified by Core Lab.

²Subsequent endoleaks may have been of different type than original.

³No endoleaks after 12 months; follow-up after 24 months not available.

6.6 Aneurysm Change

Tables 6.6.1 - 6.6.3 present the change in aneurysm diameter for the endovascular patients, as identified by the Core Lab. **Table 6.6.1** presents maximum aneurysm diameter change by interval. **Tables 6.6.2** and **6.6.3** present aneurysm change and endoleak at 12 and 24 months, respectively.

Table 6.6.1 Change in Maximum Aneurysm Diameter by Interval*						
Item	Zenith Standard Risk		Zenith High Risk		Zenith Roll-in	
From pre-discharge 30-day						
Decrease >5 mm	1.7%	(3/180)	4.8%	(4/84)	0.0%	(0/40)
Unchanged	97%	(174/180)	94%	(79/84)	97.5%	(39/40)
Increase >5 mm	1.7%	(3/180)	1.2%	(1/84)	2.5%	(1/40)
Change from pre-discharge 6-month						
Decrease >5 mm	36%	(63/173)	41%	(30/73)	49%	(18/37)
Unchanged	62%	(108/173)	59%	(43/73)	51%	(19/37)
Increase >5 mm	1.2%	(2/173)	0.0%	(0/73)	0.0%	(0/37)
Change from pre-discharge 12-month						
Decrease >5 mm	68%	(102/151)	63%	(39/62)	67%	(20/30)
Unchanged	31%	(47/151)	35%	(22/62)	33%	(10/30)
Increase >5 mm	1.3%	(2/151)	1.6%	(1/62)	0.0%	(0/30)
Change from pre-discharge 24-month						
Decrease >5 mm	78%	(74/94)	75%	(27/36)	71%	(17/24)
Unchanged	19%	(18/94)	25%	(9/36)	25%	(6/24)
Increase >5 mm	2.1%	(2/94)	0.0%	(0/36)	4.2%	(1/24)

*Only includes subjects with interpretable films and measurements of aneurysm change from 1 to 24 months.

Table 6.6.2 Change in Aneurysm Size and Endoleak at 12 Months									
	Zenith Standard Risk N=140			Zenith High Risk N=53			Zenith Roll-in N=27		
	Endoleak			Endoleak			Endoleak		
Item	N	n	%	N	n	%	N	n	%
Aneurysm size change from pre-discharge to 12 months									
Decrease >5 mm	96	3	3.0	35	3	9.0	20	0	0.0
Unchanged	42	8	19	17	2	12	7	1	14
Increase >5 mm	2	0	0.0	1	0	0.0	0	0	0.0

Table 6.6.3 Change in Aneurysm Size and Endoleak at 24 Months									
Item	Zenith Standard Risk N=90			Zenith High Risk N=31			Zenith Roll-in N=21		
	Endoleak			Endoleak			Endoleak		
	N	n	%	N	n	%	N	n	%
Aneurysm size change from pre-discharge to 24 months									
Decrease >5 mm	71	3	4.0	24	1	4.0	16	0	0.0
Unchanged	18	1	6.0	7	3	43	5	0	0.0
Increase >5 mm	1	1	100	0	0	0.0	0	0	0.0

6.7 AAA-Related Secondary Interventions

AAA-related secondary interventions within the first year were performed in 11% of the Zenith standard risk, 13% Zenith high risk and 5.8% Zenith roll-in subjects as shown in **Table 6.7.1**. Greater than 50% of the secondary interventions involved catheterization to treat an endoleak. AAA-related secondary interventions within the second year were performed in 4.2% of the Zenith standard risk, 2.2% Zenith high risk and 2.3% Zenith roll-in subjects as shown in **Table 6.7.2**.

Table 6.7.1 Secondary Interventions (to 12 Months)						
Intervention	Zenith Standard Risk N=199		Zenith High Risk N=100		Zenith Roll-in N=52	
	n	%	n	%	n	%
Conversion to open repair	2	1.0	1	1.0	0	0.0
Subjects with ≥1 intervention	21	11	13	13	3	5.8
Treat an endoleak						
Embolization	5	2.5	3	3.0	1	2.0
Ancillary component	6	3.0	4	4.0	1	2.0
Stent	1	0.5	0	0.0	0	0.0
Angioplasty	0	0.0	1	1.0	0	0.0
Treat an aneurysm increase						
Embolization	0	0.0	0	0.0	0	0.0
Treat a limb occlusion	1	0.5	3	3.0	1	2.0
Treat a limb stenosis	1	0.5	0	0.0	0	0.0
Treat a renal artery	5	2.5	0	0.0	0	0.0
Treat infra-inguinal ischemia	1	0.5	1	1.0	0	0.0
Treat multiple events	1	0.5	1	1.0	0	0.0

Table 6.7.2 Secondary Interventions (>12 to 24 Months)						
Intervention	Zenith Standard Risk N=190		Zenith High Risk N=90		Zenith Roll-in N=44	
	n	%	n	%	n	%
Conversion to open repair	1	0.5	1	1.1	0	0.0
Subjects with ≥1 intervention	8	4.2	2	2.2	1	2.3
Treat an endoleak						
Embolization	3 ¹	1.6	0	0.0	1 ²	2.3
Ancillary component	1	0.5	0	0.0	0	0.0
Stent	0	0.0	1	1.1	0	0.0
Angioplasty	0	0.0	0	0.0	0	0.0
Treat an aneurysm increase						
Embolization	1	0.5	0	0.0	1 ²	2.3
Treat a limb occlusion	1	0.5	0	0.0	0	0.0
Treat a limb stenosis	0	0.0	0	0.0	0	0.0
Treat a renal artery	0	0.0	1	1.1	0	0.0
Treat infra-inguinal ischemia	1	0.5	0	0.0	0	0.0
Treat multiple events	1	0.5	0	0.0	0	0.0

¹ Patient also received ancillary component.

² Patient underwent intervention to treat aneurysm increase and endoleak.

6.8 Secondary Outcome Measures

As described in **Table 6.8.1**, treatment of AAA with the Zenith AAA Endovascular Graft compared to the surgical control group demonstrated significant benefits in recovery and quality of life measures.

Table 6.8.1 Secondary Outcomes by Treatment Group

Item	Zenith Standard Risk	Surgical Standard Risk	P value	Zenith High Risk	Zenith Roll-in
Anesthesia time (min)	221.6 ± 67.3	304.5 ± 102.7	<.001	218.9 ± 69.6	213.9 ± 57.7
Procedure time (min)	153.2 ± 56.3	238.7 ± 92.2	<.001	153.5 ± 58.6	155.9 ± 43.2
Blood bank products received	5.0% (10/200)	84% (67/80)	<.001	12% (12/100)	3.8% (2/52)
Blood loss (cc)	299 ± 324	1676 ± 1676	<.001	356 ± 514	265 ± 226
Days in ICU	0.4 ± 0.9	3.4 ± 4.6	<.001	0.5 ± 1.2	0.5 ± 0.9
Days to discharge	2.6 ± 1.7	8.8 ± 5.6	<.001	3.0 ± 2.8	2.7 ± 1.5
Days to oral fluids	0.5 ± 0.8	3.9 ± 2.5	<.001	0.5 ± 0.6	0.7 ± 0.5
Days to normal diet	1.3 ± 1.2	6.6 ± 4.9	<.001	1.3 ± 0.8	1.1 ± 0.7
Days to normal bowel function	2.6 ± 1.4	4.2 ± 2.1	<.001	2.6 ± 1.5	2.0 ± 1.2
Days to ambulation	1.2 ± 0.7	3.5 ± 3.4	<.001	1.2 ± 0.7	1.2 ± 0.6
Hours of intubation	1.9 ± 2.2	11.7 ± 13.6	<.001	1.2 ± 1.7	2.6 ± 4.6
Maximum temperature	101.1 ± 1.3	100.7 ± 1.2	.06	100.8 ± 1.1	101.0 ± 1.2

7 PATIENT SELECTION AND TREATMENT

(See **Section 4, WARNINGS AND PRECAUTIONS**)

7.1 Individualization of Treatment

Cook recommends that the Zenith Flex AAA Endovascular Graft component diameters be selected as described in **Tables 10.5.1** and **10.5.2**. The length of the Zenith Flex AAA Endovascular Graft should extend from the lowest renal artery to just above the internal iliac (hypogastric) artery bifurcation. All lengths and diameters of the devices necessary to complete the procedure should be available to the physician, especially when preoperative case planning measurements (treatment diameters/lengths) are not certain. This approach allows for greater intraoperative flexibility to achieve optimal procedural outcomes. The risks and benefits previously described in **Section 6, SUMMARY OF CLINICAL STUDIES** should be carefully considered for each patient before use of the Zenith Flex AAA Endovascular Graft. Additional considerations for patient selection include but are not limited to:

- Patient's age and life expectancy
- Co-morbidities (e.g., cardiac, pulmonary or renal insufficiency prior to surgery, morbid obesity)
- Patient's suitability for open surgical repair
- Patient's anatomical suitability for endovascular repair
- The risk of aneurysm rupture compared to the risk of treatment with the Zenith Flex AAA Endovascular Graft
- Ability to tolerate general, regional or local anesthesia
- Iliofemoral access vessel size and morphology (minimal thrombus, calcium and/or tortuosity) should be compatible with vascular access techniques and accessories of the delivery profile of a 16 French to 22 French vascular introducer sheath
- Non-aneurysmal infrarenal aortic segment (neck) proximal to the aneurysm:
 - with a length of at least 15 mm,
 - with a diameter measured outer wall to outer wall of no greater than 32 mm and no less than 18 mm,
 - with an angle less than 60 degrees relative to the long axis of the aneurysm, and
 - with an angle less than 45 degrees relative to the axis of the suprarenal aorta
- Iliac artery distal fixation site greater than 10 mm in length and 7.5 to 20 mm in diameter (measured outer wall to outer wall)
- Freedom from significant femoral/iliac artery occlusive disease that would impede flow through the endovascular graft

The final treatment decision is at the discretion of the physician and patient.

8 PATIENT COUNSELING INFORMATION

The physician and patient (and/or family members) should review the risks and benefits when discussing this endovascular device and procedure, including:

- Risks and differences between endovascular repair and surgical repair
- Potential advantages of traditional open surgical repair
- Potential advantages of endovascular repair
- The possibility that subsequent interventional or open surgical repair of the aneurysm may be required after initial endovascular repair

In addition to the risks and benefits of an endovascular repair, the physician should assess the patient's commitment and compliance to postoperative follow-up as necessary to ensure continuing safe and effective results. Listed below are additional topics to discuss with the patient as to expectations after an endovascular repair:

- **The long-term performance of endovascular grafts has not yet been established. All patients should be advised that endovascular treatment requires lifelong, regular follow-up to assess their health and the performance of their endovascular graft.** Patients with specific clinical findings (e.g., endoleaks, enlarging aneurysms or changes in the structure or position of the endovascular graft) should receive enhanced follow-up. Specific follow-up guidelines are described in **Section 12, IMAGING GUIDELINES AND POSTOPERATIVE FOLLOW-UP**.
- Patients should be counseled on the importance of adhering to the follow-up schedule, both during the first year and at yearly intervals thereafter. Patients should be told that regular and consistent follow-up is a critical part of ensuring the ongoing safety and effectiveness of endovascular treatment of AAAs. At a minimum, annual imaging and adherence to routine postoperative follow-up requirements are required and should be considered a lifelong commitment to the patient's health and well-being.
- The patient should be told that successful aneurysm repair does not arrest the disease process. It is still possible to have associated degeneration of vessels.
- Physicians must advise each patient that it is important to seek prompt medical attention if they experience signs of limb occlusion, aneurysm enlargement or rupture. Signs of graft limb occlusion include pain in the hip(s) or leg(s) during walking or at rest or discoloration or coolness of the leg. Aneurysm rupture may be asymptomatic, but usually presents as: pain; numbness; weakness in the legs; any back, chest, abdominal or groin pain; dizziness; fainting; rapid heartbeat or sudden weakness.
- Due to the imaging required for successful placement and follow-up of endovascular devices, the risks of radiation exposure to developing tissue should be discussed with women who are or suspect they are pregnant. Men who undergo endovascular or open surgical repair may experience impotence.

Physicians should refer patients to the *Patient Guide* regarding risks occurring during or after implantation of the device. Procedure-related risks include cardiac, pulmonary, neurologic, bowel and bleeding complications. Device-related risks include occlusion, endoleak, aneurysm enlargement, fracture, potential for reintervention and open surgical conversion, rupture and death (see **Section 5.1, Observed Adverse Events** and **Section 5.2, Potential Adverse Events**). The physician should complete the *Patient I.D. Card* and give it to the patient so that he/she can carry it with him/her at all times. The patient should refer to the card anytime he/she visits additional health practitioners, particularly for any additional diagnostic procedures (e.g., MRI).

9 HOW SUPPLIED

- The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is sterilized by ethylene oxide gas and pre-loaded in peel-open packages.
- The device is intended for single use only. Do not re-sterilize the device.
- The product is sterile if the package is unopened or undamaged. Inspect the device and packaging to verify that no damage has occurred as a result of shipping. Do not use this device if damage has occurred or if the sterilization barrier has been damaged or broken. If damage has occurred, do not use the product and return to COOK.
- Prior to use, verify correct devices (quantity and size) have been supplied for the patient by matching the device to the order prescribed by the physician for that particular patient.
- The main body device is loaded onto an 18, 20 or 22 French Flexor introducer sheath. The sheath's surface is treated with a hydrophilic coating that, when hydrated, enhances trackability. To activate the hydrophilic coating, the surface must be wiped with a sterile gauze pad soaked in saline solution.
- Do not use after the "USE BY" (expiration) date printed on the label.
- Store in a cool, dry place.

10 CLINICAL USE INFORMATION

10.1 Physician Training

CAUTION: Always have a qualified surgery team available during implantation or re-intervention procedures in the event that conversion to open surgical repair is necessary.

CAUTION: The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System should only be used by physicians and teams trained in vascular interventional techniques and in the use of this device. The recommended skill/knowledge requirements for physicians using the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System are outlined below:

Patient Selection:

- Knowledge of the natural history of abdominal aortic aneurysms (AAA) and co-morbidities associated with AAA repair.
- Knowledge of radiographic image interpretation, device selection, planning and sizing.

A multidisciplinary team that has combined procedural experience with:

- Femoral cutdown, arteriotomy and repair
- Percutaneous access and closure techniques
- Non-selective and selective wire guide and catheter techniques
- Fluoroscopic and angiographic image interpretation
- Embolization
- Angioplasty
- Endovascular stent placement
- Snare techniques
- Appropriate use of radiographic contrast material
- Techniques to minimize radiation exposure
- Expertise in necessary patient follow-up modalities

10.2 Inspection Prior to Use

Inspect the device and packaging to verify that no damage has occurred as a result of shipping. Do not use this device if damage has occurred or if the sterilization barrier has been damaged or broken. If damage has occurred, do not use the product and return to COOK. Prior to use, verify correct devices (quantity and size) have been supplied for the patient by matching the device to the order prescribed by the physician for that particular patient.

10.3 Materials Required

(Not included in 3-piece modular system)

- Zenith AAA Endovascular Graft Ancillary Kit
- Fluoroscope with digital angiography capabilities (C-arm or fixed unit)
- Contrast media
- Syringe
- Heparinized saline solution
- Sterile gauze pads

10.4 Materials Recommended

(Not included in 3-piece modular system)

The following products are recommended for implantation of any component in the Zenith product line. For information on the use of these products, refer to the individual product's *Suggested Instructions for Use*.

- .035 inch (0.89 mm) extra stiff wire guide, 260 cm; for example:
 - Cook Underquist Extra Stiff Wire Guides (LES)
- .035 inch (0.89 mm) standard wire guide; for example:
 - Cook .035 inch wire guides
 - Cook Nimble™ Wire Guides
- Molding Balloons; for example:
 - Cook Coda Balloon Catheter
- Introducer sets; for example:
 - Cook Check-Flo® Introducer Sets
 - Cook Extra Large Check-Flo Introducer Sets
 - Cook Flexor Balkin Up & Over® Contralateral Introducers
- Sizing catheter; for example:
 - Cook Auros® Centimeter Sizing Catheters
- Angiographic radiopaque tip catheters; for example:
 - Cook Beacon® Tip Angiographic Catheters
 - Cook Beacon Tip Royal Flush Catheters
- Entry needles; for example:
 - Cook single wall entry needles
- Endovascular Dilators; for example:
 - Cook Endovascular Dilator Sets

10.5 Device Sizing Guidelines

The choice of diameter should be determined from the outer wall to outer wall vessel diameter and **not** the lumen diameter. Undersizing or oversizing may result in incomplete sealing or compromised flow.

Table 10.5.1 Main Body Graft Sizing Guide*

Intended Aortic Vessel Diameter ^{1,2} (mm)	Main Body Diameter ³ (mm)	Overall Length to Contralateral Limb/ Overall Length to Ipsilateral Limb (mm)	Introducer Sheath (Fr)
18-19	22	82/112, 96/126, 111/141, 125/155, 140/170	18
20-21	24	82/112, 96/126, 111/141, 125/155, 140/170	18
22	26	82/112, 96/126, 111/141, 125/155, 140/170	18
23-24	28	82/112, 96/126, 111/141, 125/155, 140/170	20
25-26	30	82/112, 96/126, 111/141, 125/155, 140/170	20
27-28	32	82/112, 96/126, 111/141, 125/155, 140/170	20
29-32	36	95/125, 113/143, 131/161, 149/179	22

¹Maximum diameter along the proximal fixation site.

²Round measured aortic diameter to nearest mm.

³Additional considerations may affect choice of diameter.

*All dimensions are nominal.

Table 10.5.2 Iliac Leg Graft Sizing Guide*

Intended Iliac Vessel Diameter ^{1,2} (mm)	Iliac Leg Diameter ³ (mm)	Iliac Leg Working Length ⁴ (mm)	Introducer Sheath (Fr)
<8	8	37, 54, 71, 88, 105, 122	14
8-9	10	37, 54, 71, 88, 105, 122	14
10-11	12	39, 56, 73, 90, 107, 124	16
12-13	14	39, 56, 73, 90, 107, 124	16
14-15	16	39, 56, 73, 90	16
16-17	18	39, 56, 73, 90	16
18	20	39, 56, 73, 90	16
19	22	39, 56, 73, 90	16
20	24	39, 56, 73, 90	16

¹Maximum diameter along the distal fixation site.

²Round measured iliac diameter to nearest mm.

³Additional considerations may affect choice of diameter.

⁴Overall leg length = working length + 22 mm docking stent.

*All dimensions are nominal.

11 DIRECTIONS FOR USE

Anatomical Requirements

- Iliofemoral access vessel size and morphology (minimal thrombus, calcium and/or tortuosity) should be compatible with vascular access techniques and accessories. Arterial conduit techniques may be required.
- Proximal aortic neck lengths should be a minimum of 15 mm with a diameter measured outer wall to outer wall of 18 – 32 mm.
- Iliac artery distal fixation site should be greater than 10 mm in length and 7.5 – 20 mm in diameter (measured outer wall to outer wall).

Prior to use of the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System, review this *Suggested Instructions for Use* booklet. The following instructions embody a basic guideline for device placement. Variations in the following procedures may be necessary. These instructions are intended to help guide the physician and do not take the place of physician judgment.

General Use Information

Standard techniques for placement of arterial access sheaths, guiding catheters, angiographic catheters and wire guides should be employed during use of the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System. The Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is compatible with .035 inch diameter wire guides.

Endovascular stent grafting is a surgical procedure, and blood loss from various causes may occur, infrequently requiring intervention (including transfusion) to prevent adverse outcomes. It is important to monitor blood loss from the hemostatic valve throughout the procedure, but is specifically relevant during and after manipulation of the gray positioner. After the gray positioner has been removed, if blood loss is excessive, consider placing an uninflated molding balloon or an introduction system dilator within the valve, restricting flow.

Pre-Implant Determinants

Verify from pre-implant planning that the correct device has been selected. Determinants include:

1. Femoral artery selection for introduction of the main body system (i.e., define respective contralateral and ipsilateral iliac arteries).
2. Angulation of aortic neck, aneurysm and iliac arteries.
3. Quality of the aortic neck.
4. Diameters of infrarenal aortic neck and distal iliac arteries.
5. Distance from renal arteries to the aortic bifurcation.
6. Length from the aortic bifurcation to the internal iliac arteries/attachment site(s).
7. Aneurysm(s) extending into the iliac arteries may require special consideration in selecting a suitable graft/artery interface site.
8. Consider the degree of vascular calcification.

Patient Preparation

1. Refer to institutional protocols relating to anesthesia, anticoagulation and monitoring of vital signs.
2. Position patient on imaging table allowing fluoroscopic visualization from the aortic arch to the femoral bifurcations.
3. Expose both common femoral arteries using standard surgical technique.
4. Establish adequate proximal and distal vascular control of both femoral vessels.

11.1 Bifurcated System (Figs. 2 and 3)

11.1.1 Bifurcated Main Body Preparation/Flush

1. Remove gray-hubbed shipping stylet (from the inner cannula) and dilator tip protector (from the dilator tip). Remove Peel-Away™ sheath from back of hemostatic valve. (**Fig. 5**) Elevate distal tip of system and flush through the stopcock on the hemostatic valve until fluid emerges from the sideport near the tip of the introducer sheath. (**Fig. 6**) Continue to inject a full 20 cc of flushing solution through the device. Discontinue injection and close stopcock on connecting tube.

NOTE: Graft flushing solution of heparinized saline is often used.

2. Attach syringe with heparinized saline to the hub on the inner cannula. Flush until fluid exits the dilator tip. (**Fig. 7**)

NOTE: When flushing system, elevate distal end of system to facilitate removal of air.

3. Soak sterile gauze pads in saline solution and use to wipe Flexor introducer sheath to activate the hydrophilic coating. Hydrate both sheath and dilator liberally.

11.1.2 Contralateral Iliac Leg Preparation/Flush

1. Remove gray-hubbed shipping stylet (from the inner cannula) and dilator tip protector (from the dilator tip). Remove Peel-Away sheath from back of the hemostatic valve. (**Fig. 8**) Elevate distal tip of system and flush through the stopcock on the hemostatic valve until fluid emerges from the sideport near the tip of the introducer sheath. (**Fig. 6**) Continue to inject a full 20 cc of flushing solution through the device. Discontinue injection and close stopcock on connecting tube.

NOTE: Graft flushing solution of heparinized saline is often used.

2. Attach syringe with heparinized saline to the hub on the inner cannula. Flush until fluid exits the dilator tip. (**Fig. 7**)

NOTE: When flushing system, elevate distal end of system to facilitate removal of air.

11.1.3 Ipsilateral Iliac Leg Preparation/Flush

Follow the previous Section **11.1.2, Contralateral Iliac Leg Preparation/Flush** instructions to ensure proper flushing of the ipsilateral iliac leg graft.

11.1.4 Vascular Access and Angiography

1. Puncture the selected common femoral arteries using standard technique with an 18 or 19 UT gage arterial needle. Upon vessel entry, insert:
 - Wire guides — standard .035 inch diameter, 145 cm long, J tip or Bentson Wire Guide
 - Appropriate size sheaths (e.g., 6 or 8 French)
 - Flush catheter (often radiopaque sizing catheters – e.g., Centimeter Sizing Catheter or straight flush catheter)
2. Perform angiography to identify level(s) of renals, aortic bifurcation and iliac bifurcations.

NOTE: If fluoroscope angulation is used with an angulated neck, it may be necessary to perform angiograms using various projections.

11.1.5 Main Body Placement

1. Ensure the delivery system has been flushed with heparinized saline and that all air is removed from the system.
2. Give systemic heparin and check flushing solutions. Flush after each catheter and/or wire guide exchange.

NOTE: Monitor the patient’s coagulation status throughout the procedure.

3. On ipsilateral side, replace J wire with stiff wire guide (LES) .035 inch, 260 cm long, and advance through catheter and up to the thoracic aorta. Remove flush catheter and sheath. Maintain wire guide position.
4. Before insertion, position main body delivery system on patient’s abdomen under fluoroscopy to determine the orientation of the contralateral limb radiopaque marker. The sidearm of the hemostatic valve may serve as an external reference to the contralateral limb radiopaque marker.
5. Insert main body delivery system over the wire, into the femoral artery with attention to sidearm reference.

CAUTION: Maintain wire guide position during delivery system insertion.

CAUTION: To avoid any twist in the endovascular graft, during any rotation of the delivery system, be careful to rotate all of the components of the system together (from outer sheath to inner cannula).

6. Advance delivery system until the four gold radiopaque markers (which are positioned 2 mm from the most proximal segment of the graft material) (**Fig. 9, Illustration 1**) are just inferior to the most inferior renal orifice.
7. Verify position of wire guide in the thoracic aorta. Ensure the graft system is oriented such that the contralateral limb is positioned above and anterior to the origin of the contralateral iliac. If the contralateral limb radiopaque marker is not properly aligned, rotate the entire system until it is correctly positioned halfway between a lateral and an anterior position on the contralateral side.
 - A marker formation of a ✓ indicates an anterior position of the short (contralateral) limb. (**Fig. 9, Illustration 4**)
 - A marker formation of a ✗ indicates a posterior position of the short (contralateral) limb. (**Fig. 9, Illustration 5**)
 - A marker formation of a | indicates a lateral position of the short (contralateral) limb. (**Fig. 9, Illustration 6**)

NOTE: Radiopaque cannula on each stent between the renal arteries and the contralateral limb align with the contralateral limb radiopaque marker.

8. Repeat the angiogram to verify the four gold radiopaque markers are 2 mm or more below the most inferior renal orifice.
9. Ensure the Captor Hemostatic Valve on the Flexor introducer sheath is turned to the open position. **(Fig. 10)**
10. Use the gripper to stabilize the gray positioner (the shaft of the delivery system) while withdrawing the sheath. Deploy the first two covered stents by withdrawing the sheath while monitoring device location. **(Fig. 18)**
11. Without moving the table, decrease magnification to check position of the contralateral limb radiopaque marker and location of renal arteries. Proceed with deployment until the contralateral limb is fully deployed. **(Fig. 11)** Stop withdrawing sheath.

NOTE: Verify contralateral limb is at least 5 mm above the aortic bifurcation and in desired location for cannulation.

11.1.6 Contralateral Iliac Wire Guide Placement

1. Manipulate catheter and wire guide through open end of contralateral limb into body of the graft. Advance the wire guide until it curves inside the body of the graft. AP and oblique fluoroscopic views can aid in verification of device cannulation.
2. After cannulation, advance the angiographic catheter over the wire into the body of the endovascular graft. Remove wire and perform angiography to confirm position. Advance wire guide until it curves inside the body of the graft. **(Fig. 12)** Remove angiographic catheter.

11.1.7 Main Body Proximal (Top) Deployment

1. Perform angiography through an angiographic catheter to verify position of the endovascular graft with respect to the renal arteries. If necessary, carefully reposition the covered portion of the endovascular graft with respect to the renal arteries. (Repositioning can only take place over a small range of distance at this stage.)

NOTE: Ensure patency of renal arteries by confirming that the proximal graft markers are 2 mm or more below the lowest patent renal artery.

CAUTION: During proximal trigger-wire removal, top cap advancement, and subsequent suprarenal stent deployment, verify that the position of the main body wire guide extends just distal to the aortic arch and that support of the system is maximized.

2. Remove the safety lock from the top stent trigger-wire release mechanism. Under fluoroscopy, withdraw and remove the trigger-wire by sliding the top stent trigger-wire release mechanism off the handle and then remove via its slot over the inner cannula. **(Fig. 13)**

If resistance is felt or system bowing is noticed, the trigger-wire is under tension. Excessive force may cause the graft position to be altered. If excessive resistance or delivery system movement is noted, stop and assess the situation. If unable to remove the top stent trigger-wire release mechanism from the top cap, perform the following steps under fluoroscopy:

- a. Remove tension on the trigger-wire by loosening the pin vise and slightly **pulling** the inner cannula to move the top cap down over the suprarenal stent. Avoid compressing the Zenith Flex main body.
- b. Retighten the pin vise.
- c. Remove the top stent trigger-wire release mechanism.
- d. Continue with Step 3 in **Section 11.1.7, Main Body Proximal (Top) Deployment**.

NOTE: If still unable to remove the top stent trigger-wire release mechanism from the top cap, see **Section 14.1, Trigger-Wire Release Troubleshooting**.

3. Loosen the pin vise. **(Fig. 14)** Control the position of the graft by stabilizing the gray positioner of the introducer.
4. Deploy the suprarenal stent by advancing the top cap inner cannula 1 to 2 mm at a time while controlling the position of the main body until the top stent is fully deployed. **(Figs. 15 and 16)** Advance the top cap cannula an additional 1 to 2 cm and then retighten the pin vise to avoid contact with the deployed suprarenal stent.

NOTE: If unable to fully deploy suprarenal stent by advancing the top cap inner cannula, see **Section 14.2, Suprarenal Stent Deployment Troubleshooting**.

NOTE: Once the barbed suprarenal stent has been deployed, further attempts to reposition the graft are not recommended.

WARNING: The Zenith Flex AAA Endovascular Graft incorporates a suprarenal stent with fixation barbs. Exercise extreme caution when manipulating interventional devices in the region of the suprarenal stent.

5. Advance the contralateral wire guide into the thoracic aorta.

11.1.8 Contralateral Iliac Leg Placement and Deployment

CAUTION: Verify that the predetermined contralateral iliac leg is selected for insertion on the contralateral side of the patient before implantation.

1. Position the image intensifier to show both the contralateral internal iliac artery and contralateral common iliac artery.
2. Prior to the introduction of the contralateral iliac leg delivery system, inject contrast through the contralateral femoral sheath to locate the contralateral internal iliac artery.
3. Introduce the contralateral iliac leg delivery system into the artery. Advance slowly until the iliac leg graft overlaps at least one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the contralateral limb of the main body. **(Fig. 17)** If there is any tendency for the main body graft to move during this maneuver, hold it in position by stabilizing the gray positioner on the ipsilateral side.

NOTE: If difficulty is encountered advancing the iliac leg delivery system, exchange to a more supportive wire guide. In tortuous vessels the anatomy may alter significantly with the introduction of the rigid wires and sheath systems.

4. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure both internal iliac patency and a minimum overlap of one full iliac leg stent (i.e., proximal stent of iliac leg graft, maximum overlap of 1.5 stents) within the main body endovascular graft.

NOTE: Ensure Captor Hemostatic Valve on the iliac leg introducer sheath is turned to the open position. **(Fig. 10)**

5. To deploy, hold the iliac leg graft in position with the gripper on the gray positioner while withdrawing the sheath. **(Figs. 18 and 19)** Ensure one stent overlap is maintained.
6. Stop withdrawing the sheath as soon as the distal end of the iliac leg graft is released.
7. Under fluoroscopy and after verification of iliac leg graft position, loosen pin vise and retract inner cannula to dock tapered dilator to gray positioner. Tighten pin vise. Maintain sheath position while withdrawing gray positioner with secured inner cannula.
8. Close the Captor Hemostatic Valve on the iliac leg introducer sheath by turning it in a clockwise direction until it stops. **(Fig. 26)**
9. Re-check the position of the wire guide.

11.1.9 Main Body Distal (Bottom) Deployment

1. Return to the ipsilateral side.
2. Fully deploy the ipsilateral limb of the main body by withdrawing the sheath until the most distal stent has expanded. **(Figs. 18 and 20)** Stop withdrawing sheath.

NOTE: The distal stent is still secured by the trigger-wire.

3. Remove the safety lock from the ipsilateral limb trigger-wire release mechanism. Withdraw and remove the trigger-wire by sliding the ipsilateral limb trigger-wire release mechanism off the handle and then remove via its slot over the inner cannula. **(Fig. 21)**

11.1.10 Docking of Top Cap

1. Loosen the pin vise. **(Fig. 22)**
2. Secure sheath and inner cannula to avoid any movement of these components.
3. Advance the gray positioner over the inner cannula until it docks with the top cap. **(Figs. 23, 24 and 25)**

NOTE: If resistance occurs, slightly rotate gray positioner and continue to gently advance.

4. Retighten the pin vise and withdraw the entire top cap and gray positioner through the graft and through the sheath by pulling on the inner cannula. Leave the sheath and wire guide in place.

NOTE: Maintain position of sheath and wire guide.

5. Close the Captor Hemostatic Valve on the main body introducer sheath by turning it in a clockwise direction until it stops. **(Fig. 26)**

11.1.11 Ipsilateral Iliac Leg Placement and Deployment

NOTE: Ensure the Captor Hemostatic Valve on the main body introducer sheath is turned to the open position. **(Fig. 27)**

1. Utilize the main body graft wire and sheath assembly to introduce the ipsilateral iliac leg graft. Advance dilator and sheath assembly into the main body sheath.

NOTE: In tortuous vessels, the position of the internal iliac arteries may alter significantly with the introduction of the rigid wires and sheath systems.

2. Advance slowly until the ipsilateral iliac leg graft overlaps a minimum of one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the ipsilateral limb of the main body. **(Fig. 28)**

NOTE: If an overlap of greater than three iliac leg stents is required (greater than two iliac leg stents for 37, 39, 54, and 56 mm leg lengths), it may be necessary to consider use of a leg extension in the bifurcation area of the opposite side.

3. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure internal iliac patency.

NOTE: Ensure the Captor Hemostatic Valve on the iliac leg introducer sheath is turned to the open position. **(Fig. 10)**

4. To deploy, stabilize the iliac leg graft with the gripper on the gray positioner while withdrawing the iliac leg sheath. **(Figs. 29 and 30)** If necessary, withdraw the main body sheath.
5. Under fluoroscopy and after verification of iliac leg graft position, loosen pin vise, and retract inner cannula to dock tapered dilator to gray positioner. Tighten pin vise. Maintain main body sheath position while withdrawing iliac leg sheath and gray positioner with secured inner cannula.
6. Close the Captor Hemostatic Valve on the main body introducer sheath by turning it in a clockwise direction until it stops.
7. Re-check the position of the wire guides. Leave sheath and wire guides in place.

11.1.12 Molding Balloon Insertion

1. Prepare molding balloon as follows:
 - Flush wire lumen with heparinized saline.
 - Remove all air from balloon.
2. In preparation for the insertion of the molding balloon, open the Captor Hemostatic Valve by turning counter-clockwise.
3. Advance the molding balloon over the wire guide and through the Captor Hemostatic Valve of the main body introduction system to the level of the renal arteries. Maintain proper sheath position.
4. Tighten the Captor Hemostatic Valve around the molding balloon with gentle pressure by turning it clockwise.

CAUTION: Do not inflate balloon in vessel outside of graft.

5. Expand the molding balloon with diluted contrast media (as directed by the manufacturer) in the area of the most proximal covered stent and the infrarenal neck, starting proximally and working in the distal direction. **(Fig. 31)**

CAUTION: Confirm complete deflation of balloon prior to repositioning.

CAUTION: Captor Hemostatic Valve must be open prior to repositioning of molding balloon.

6. Withdraw the molding balloon to the ipsilateral limb overlap and expand.

CAUTION: Captor Hemostatic Valve must be open prior to repositioning of molding balloon.

7. Withdraw the molding balloon to the ipsilateral distal fixation site and expand.

CAUTION: Do not inflate balloon in vessel outside of graft.

CAUTION: Captor Hemostatic Valve must be open prior to repositioning of molding balloon.

8. Deflate and remove molding balloon. Transfer the molding balloon onto the contralateral wire guide and into the contralateral iliac leg introduction system. Advance molding balloon to the contralateral limb overlap and expand.

CAUTION: Confirm complete deflation of balloon prior to repositioning.

CAUTION: Captor Hemostatic Valve must be open prior to repositioning of molding balloon.

9. Withdraw the molding balloon to the contralateral iliac leg/vessel distal fixation and expand. **(Fig. 31)**

CAUTION: Do not inflate balloon in vessel outside of graft.

10. Remove molding balloon and replace it with an angiographic catheter to perform completion angiograms.
11. Remove or replace all stiff wire guides to allow iliac arteries to resume their natural position.

Final Angiogram

1. Position angiographic catheter just above the level of the renal arteries. Perform angiography to verify that the renal arteries are patent and that there are no endoleaks. Verify patency of internal iliac arteries.
2. Confirm there are no endoleaks or kinks and verify position of proximal gold radiopaque markers. Remove the sheaths, wires and catheters.

NOTE: If endoleaks or other problems are observed, refer to the *Suggested Instructions for Use* for the Zenith AAA Endovascular Graft Ancillary Components.

3. Repair vessels and close in standard surgical fashion.

12 IMAGING GUIDELINES AND POSTOPERATIVE FOLLOW-UP

12.1 General

- **The long-term performance of endovascular grafts has not yet been established. All patients should be advised that endovascular treatment requires lifelong, regular follow-up to assess their health and performance of their endovascular graft.** Patients with specific clinical findings (e.g., endoleaks, enlarging aneurysms or changes in the structure or position of the endovascular graft) should receive additional follow-up. Patients should be counseled on the importance of adhering to the follow-up schedule, both during the first year and at yearly intervals thereafter. Patients should be told that regular and consistent follow-up is a critical part of ensuring the ongoing safety and effectiveness of endovascular treatment of AAAs.

- Physicians should evaluate patients on an individual basis and prescribe their follow-up relative to the needs and circumstances of each individual patient. The recommended imaging schedule is presented in **Table 12.1.1**. This schedule was used in the pivotal trial, and is recommended even in the absence of clinical symptoms (e.g., pain, numbness, weakness). Patients with specific clinical findings (e.g., endoleaks, enlarging aneurysms or changes in the structure or position of the stent graft) should receive follow-up at more frequent intervals.
- Annual imaging follow-up should include abdominal radiographs and both contrast and non-contrast CT examinations. If renal complications or other factors preclude the use of image contrast media, abdominal radiographs, non-contrast CT and duplex ultrasound may be used.
- The combination of contrast and non-contrast CT imaging provides information on aneurysm diameter change, endoleak, patency, tortuosity, progressive disease, fixation length and other morphological changes.
- The abdominal radiographs provide information on device integrity (separation between components, stent fracture and barb separation).
- Duplex ultrasound imaging may provide information on aneurysm diameter change, endoleak, patency, tortuosity and progressive disease. In this circumstance, a non-contrast CT should be performed to use in conjunction with the ultrasound. Ultrasound may be a less reliable and sensitive diagnostic method compared to CT.

Table 12.1.1 lists the minimum imaging follow-up for patients with the Zenith Flex AAA Endovascular Graft. Patients requiring enhanced follow-up should have interim evaluations.

Table 12.1.1 Recommended Imaging Schedule for Endograft Patients

	Angiogram	CT (Contrast and non-contrast)	Abdominal Radiographs
Pre-procedure	X ¹	X ¹	
Procedural	X		
Pre-discharge (within 7 days)		X ^{2,3,4}	X
1 month		X ^{2,3,4}	X
3 month		X ^{2,4,5}	
6 month		X ^{2,4}	X
12 month (annually thereafter)		X ^{2,4}	X

¹Imaging should be performed within 6 months before the procedure.

²Duplex ultrasound may be used for those patients experiencing renal failure or who are otherwise unable to undergo contrast enhanced CT scan. With ultrasound, non-contrast CT is still recommended.

³Either pre-discharge or 1 month CT recommended.

⁴If Type I or III endoleak, prompt intervention and additional follow-up post-intervention recommended. See **Section 12.6, Additional Surveillance and Treatment**.

⁵Recommended if endoleak reported at pre-discharge or 1 month.

12.2 Contrast and Non-Contrast CT Recommendations

- Film sets should include all sequential images at lowest possible slice thickness (≤3 mm). DO NOT perform large slice thickness (>3 mm) and/or omit consecutive CT images/film sets, as it prevents precise anatomical and device comparisons over time.
- All images should include a scale for each film/image. Images should be arranged no smaller than 20:1 images on 14 inch x 17 inch sheets if film is used.
- Both non-contrast and contrast runs are required, with matching or corresponding table positions.
- Pre-contrast and contrast run slice thickness and interval must match.
- DO NOT change patient orientation or re-landmark patient between non-contrast and contrast runs.

Non-contrast and contrast enhanced baseline and follow-up imaging are important for optimal patient surveillance. It is important to follow acceptable imaging protocols during the CT exam. **Table 12.2.1** lists examples of acceptable imaging protocols.

Table 12.2.1 Acceptable Imaging Protocols

	Non-Contrast	Contrast
IV contrast	No	Yes
Acceptable machines	Spiral CT or high performance MDCT capable of >40 seconds	Spiral CT or high performance MDCT capable of >40 seconds
Injection volume	n/a	Per institutional protocol
Injection rate	n/a	>2.5 cc/sec
Injection mode	n/a	Power
Bolus timing	n/a	Test bolus: SmartPrep, C.A.R.E. or equivalent
Coverage - start	Diaphragm	1 cm superior to celiac axis
Coverage - finish	Proximal femur	Profunda femoris origin
Collimation	<3 mm	<3 mm
Reconstruction	2.5 mm throughout – soft algorithm	2.5 mm throughout – soft algorithm
Axial DFOV	32 cm	32 cm
Post-injection runs	None	None

12.3 Abdominal Radiographs

The following views are required:

- Four films: supine-frontal (AP), cross-table lateral, 30 degree LPO and 30 degree RPO views centered on umbilicus.
- Record the table-to-film distance and use the same distance at each subsequent examination.

Ensure entire device is captured on each single image format lengthwise.

If there is any concern about the device integrity (e.g., kinking, stent breaks, barb separation, relative component migration), it is recommended to use magnified views. The attending physician should evaluate films for device integrity (entire device length including components) using 2-4X magnification visual aid.

12.4 Ultrasound

Ultrasound imaging may be performed in place of contrast CT when patient factors preclude the use of image contrast media. Ultrasound may be paired with non-contrast CT. A complete aortic duplex is to be videotaped for maximum aneurysm diameter, endoleaks, stent patency and stenosis. Included on the videotape should be the following information as outlined below:

- Transverse and longitudinal imaging should be obtained from the level of the proximal aorta demonstrating mesenteric and renal arteries to the iliac bifurcations to determine if endoleaks are present utilizing color flow and color power angiography (if accessible).
- Spectral analysis confirmation should be performed for any suspected endoleaks.
- Transverse and longitudinal imaging of the maximum aneurysm should be obtained.

12.5 MRI Information

Non-clinical testing has demonstrated that the Zenith AAA Endovascular Graft is MR Conditional. A patient with this endovascular graft can be scanned safely immediately after placement under the following conditions:

Static Magnetic Field

- Static magnetic field of 3.0 Tesla or less.
- Highest spatial magnetic gradient field of 720 Gauss/cm

Non-clinical evaluation was conducted in a 3 Tesla MR system (General Electric Excite) with a maximum spatial magnetic gradient field of 720 Gauss/cm as measured with a gaussmeter in the position of the static magnetic

field pertinent to the patient (i.e., outside of scanner covering, accessible to a patient or individual).

MRI-Related Heating

1.5 Tesla Systems:

- Static magnetic field of 1.5 Tesla
- Maximum whole-body-averaged specific absorption rate (SAR) of 2 W/kg for 15 minutes of scanning (i.e., per scanning sequence)

In non-clinical testing, the Zenith AAA Endovascular Graft produced a temperature rise of less than or equal to 1.4 °C at a MR system reported whole-body-averaged specific absorption rate (SAR) of 2.8 W/kg, for 15 minutes of MR scanning in a 1.5 Tesla Magnetom, Siemens Medical Magnetom, Numaris/4 Software, Version Syngo MR 2002B DHHS MR Scanner. The maximum whole-body-averaged specific absorption rate (SAR) was 2.8 W/kg, which corresponds to a calorimetry measured value of 1.5 W/kg.

3.0 Tesla Systems:

- Static magnetic field of 3.0 Tesla
- Maximum whole-body-averaged specific absorption rate (SAR) of 2 W/kg for 15 minutes of scanning (i.e., per scanning sequence)

In non-clinical testing, the Zenith AAA Endovascular Graft produced a temperature rise of less than or equal to 1.9 °C at a MR system reported whole-body-averaged specific absorption rate (SAR) of 3.0 W/kg, for 15 minutes of MR scanning in a 3.0 Tesla Excite, GE Electric Healthcare, G3.0-052B Software, MR Scanner. The maximum whole-body-averaged specific absorption rate (SAR) was 3.0 W/kg, which corresponds to a calorimetry measured value of 2.8 W/kg.

Image Artifact

The image artifact extends throughout the anatomical region containing the device, obscuring the view of immediately adjacent anatomical structures within approximately 20 cm of the device, as well as the entire device and its lumen, when scanned in non-clinical testing using the sequence: Fast spin echo, in a 3.0 Tesla, Excite, GE Electric Healthcare, with G3.0-052B Software, MR system with body radiofrequency coil.

For all scanners, the image artifact dissipates as the distance from the device to the area of interest increases. MR scans of the head and neck and lower extremities may be obtained without image artifact. Image artifact may be present in scans of the abdominal region and upper extremities, depending on distance from the device to the area of interest.

Clinical information is available for seventeen patients who received MRI scans after stent-graft implantation. There have been no reported adverse events or device problems in any of these patients as a result of having received an MRI. Additionally, there have been well over 100,000 Zenith AAA Endovascular Grafts implanted worldwide, in which there have been no reported adverse events or device problems as a result of MRI.

Cook recommends that the patient register the MR conditions disclosed in this IFU with the MedicaAlert Foundation. The MedicaAlert Foundation can be contacted in the following manners.

Mail: MedicaAlert Foundation International
2323 Colorado Avenue
Turlock, CA 95382
Phone: 888-633-4298 (toll free)
209-668-3333 from outside the US
Fax: 209-669-2450
Web: www.medicalert.org

12.6 Additional Surveillance and Treatment

Additional surveillance and possible treatment is recommended for:

- Aneurysms with Type I endoleak
- Aneurysms with Type III endoleak
- Aneurysm enlargement, ≥5 mm of maximum diameter (regardless of endoleak status)
- Migration
- Inadequate seal length

Consideration for reintervention or conversion to open repair should include the attending physician's assessment of an individual patient's co-morbidities, life expectancy and the patient's personal choices. Patients should be counseled that subsequent reinterventions including catheter based and open surgical conversion are possible following endograft placement.

13 PATIENT TRACKING INFORMATION

In addition to these *Instructions for Use*, the Zenith Flex AAA Endovascular Graft with the Z-Trak Introduction System is packaged with a *Device Tracking Form* which the hospital staff is required to complete and forward to COOK for the purposes of tracking all patients who receive the Zenith Flex AAA Endovascular Graft (as required by U.S. Federal Regulation).

14 TROUBLESHOOTING

NOTE: Technical assistance from a Cook product specialist may be obtained by contacting your local Cook representative.

14.1 Trigger Wire Release Troubleshooting

CAUTION: The following steps should be performed only if unable to remove the proximal trigger-wire as described in Section 11.1.7.

14.1.1 Main Body Proximal (Top) Deployment

1. If the top stent trigger-wire release mechanism cannot be removed from the handle, cut the wire adjacent to the release mechanism (**Fig. 32**) and remove the release mechanism from the handle.
2. Stabilize the gray positioner while withdrawing the sheath to fully deploy the ipsilateral limb.
3. Remove the safety lock from the ipsilateral limb trigger-wire release mechanism.
4. Withdraw and remove the trigger-wire by sliding the ipsilateral limb trigger-wire release mechanism off the handle, and then remove via its slot over the device inner cannula.
5. Using locking forceps, clamp and secure the cut end of the top cap trigger-wire. (**Fig. 33**)
6. Loosen the pin vise and, while maintaining inner cannula and trigger-wire position, advance the gray positioner and sheath into the graft until the tip of the gray positioner is approximately 2 cm from the gold markers. (**Fig. 34**) The advanced gray positioner provides added support to the inner cannula.

NOTE: Maintain gentle tension on the trigger-wire to remove any slack in the wire as the gray positioner and sheath are being advanced.

NOTE: Ensure that the tip of the gray positioner is *not* advanced into the top cap.

7. Lock the pin vise. Confirm that the trigger-wire is secured to the forceps.
8. Stabilize the gray positioner and slowly advance the sheath until the sheath tip is 2 mm from the gold markers. (**Fig. 35**)

NOTE: Take care not to advance the graft itself during sheath advancement.

9. Stabilize the sheath and retract the gray positioner with inner cannula to draw the top cap over the suprarenal stent. (**Fig. 36**)
10. Verify position of the gold markers inferior to the renal arteries.
11. Remove the trigger-wire.
12. Withdraw the sheath until the tapered tip of the gray positioner is exposed.
13. Advance a molding balloon through the contralateral limb of the main body and position it just superior to the graft bifurcation.
14. Inflate the balloon to the full diameter of the graft. (**Fig. 37**)
15. Loosen the pin vise.
16. Stabilize the gray positioner and balloon catheter, and advance the inner cannula to deploy the top cap.
17. Tighten the pin vise.
18. Deflate the balloon and advance the contralateral wire guide into the thoracic aorta.

NOTE: Due to early ipsilateral leg deployment and trigger-wire release, it is suggested to leave the molding balloon in or just above the contralateral limb for graft stabilization during placement of the ipsilateral limb.

14.1.2 Docking of Top Cap

1. Loosen the pin vise. (**Fig. 22**)
2. Secure sheath and inner cannula to avoid any movement of these components.
3. Advance the gray positioner over the inner cannula until it docks with the top cap. (**Figs. 23, 24 and 38**)

NOTE: If resistance occurs, slightly rotate gray positioner and continue to gently advance.

4. Retighten the pin vise and withdraw the entire top cap and gray positioner through the graft and through the sheath by pulling on the inner cannula. (**Fig. 39**) Leave the sheath and wire guide in place.
5. Close the Captor Hemostatic Valve on the Flexor introducer sheath by turning it in a clockwise direction until it stops. (**Fig. 26**)

14.1.3 Ipsilateral Iliac Leg Placement and Deployment

NOTE: Ensure that the Captor Hemostatic Valve on the introducer sheath is turned to the open position. (**Fig. 27**)

1. Utilize the main body graft wire and sheath assembly to introduce the ipsilateral leg graft.

NOTE: Due to top cap release modification, the main body sheath assembly must be withdrawn to a point 1-2 cm inside the proximal ipsilateral limb. Advance the dilator sheath assembly of the ipsilateral limb into the main body sheath.

NOTE: The molding balloon may be inflated in the contralateral limb of the main body graft if additional graft stabilization is necessary.

CAUTION: Do not inflate the balloon outside of the graft.

NOTE: In tortuous vessels, the position of the internal iliac arteries may alter significantly with the introduction of the rigid wires and sheath systems.

2. Advance slowly until the ipsilateral iliac leg graft overlaps a minimum of one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the ipsilateral limb of the main body. (**Fig. 40**)

NOTE: If an overlap of greater than three iliac leg stents is required (greater than two iliac leg stents for 37, 39, 54 and 56 mm leg lengths), it may be necessary to consider use of a leg extension in the bifurcation area of the opposite side.

3. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure internal iliac patency.
4. To deploy, stabilize the iliac leg graft with the gray positioner while withdrawing the iliac leg sheath and main body sheath together. (**Figs. 29 and 41**) If necessary, withdraw the main body sheath.
5. Under fluoroscopy and after verification of iliac leg graft position, loosen the pin vise, and retract the inner cannula to dock the tapered dilator to the gray positioner. Tighten pin vise. Maintain sheath position while withdrawing gray positioner with secured inner cannula. (**Fig. 42**)
6. Close the Captor Hemostatic Valve on the Flexor introducer sheath by turning it in a clockwise direction until it stops.
7. Re-check the position of the wire guides. Leave sheath and wire guides in place.
8. Remove deflated balloon from contralateral side.

14.1.4 Contralateral Iliac Leg Placement and Deployment

CAUTION: Verify the predetermined contralateral iliac leg is selected for insertion on the contralateral side of the patient before implantation.

1. Position the image intensifier to show both the contralateral internal iliac artery and contralateral common iliac artery.
2. Prior to introduction of contralateral iliac leg delivery system, inject contrast through the contralateral femoral sheath to locate the contralateral internal iliac artery.
3. Introduce the contralateral iliac leg delivery system into the artery. Advance slowly until the iliac leg graft overlaps at least one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the contralateral limb of the main body. (**Fig. 43**)

NOTE: If difficulty is encountered advancing the iliac leg delivery system, exchange to a more supportive wire guide. In tortuous vessels, the anatomy may alter significantly with the introduction of the rigid wires and sheath systems.

4. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure both internal iliac patency and a minimum overlap of one full iliac leg stent (i.e., proximal stent of iliac leg graft, maximum overlap of 1.5 stents) within the main body endovascular graft.
5. To deploy, hold the iliac leg graft in position with the gray positioner while withdrawing the sheath. (**Figs. 18 and 44**) Ensure one stent overlap is maintained.
6. Stop withdrawing the sheath as soon as the distal end of the iliac leg graft is released.
7. Under fluoroscopy and after verification of the iliac leg graft position, loosen the pin vise, and retract the inner cannula to dock the tapered dilator to the gray positioner. Tighten pin vise. Maintain sheath position while withdrawing gray positioner with secured inner cannula. (**Fig. 45**)
8. Re-check the position of the wire guides.
9. Perform molding balloon insertion and final angiogram as described in Section 11.1.12, **Molding Balloon Insertion**.

14.2 Suprarenal Stent Deployment Troubleshooting

CAUTION: The following steps should be performed only if unable to deploy the suprarenal stent as described in Section 11.1.7.

14.2.1 Main Body Proximal (Top) Deployment

1. If the suprarenal stent cannot be fully deployed by advancing the top cap inner cannula, advance a molding balloon through the contralateral limb of the main body and position it just superior to the bifurcation of the stent graft.
2. To add support to the inner cannula, inflate the balloon to the full diameter of the graft.
3. Loosen the pin vise. (**Fig. 14**)
4. Stabilize the gray positioner and balloon catheter, and advance the inner cannula to deploy the suprarenal stent.

If the suprarenal stent is fully deployed:

5. Tighten the pin vise; deflate and withdraw the balloon.
6. Advance the contralateral wire guide into the thoracic aorta and return to the appropriate step in the Instructions for Use to complete the procedure.

If still unable to fully deploy the suprarenal stent:

5. Tighten the pin vise and deflate the balloon. While maintaining balloon position, stabilize the gray positioner and withdraw the sheath to fully deploy the ipsilateral limb.
6. Remove the safety lock from the ipsilateral limb trigger-wire release mechanism.
7. Withdraw and remove the trigger-wire by sliding the ipsilateral limb trigger-wire release mechanism off the handle, and remove via its slot over the device inner cannula.
8. Loosen the pin vise (**Fig. 22**) and, while maintaining inner cannula position, advance the gray positioner and sheath into the graft until the tip of the gray positioner is approximately 2 cm inferior to the proximal gold markers. (**Fig. 46**) The advanced gray positioner provides added support to the inner cannula.
NOTE: Take care not to advance the graft during gray positioner and sheath advancement.
NOTE: Ensure that the tip of the gray positioner is not advanced into the top cap.
9. Lock the pin vise.
10. Verify position of the gold markers inferior to the renal arteries.
11. Reposition molding balloon so that it is seated against the bifurcation.
12. Inflate the balloon to the full diameter of the graft. (**Fig. 37**)
13. Loosen the pin vise.
14. Stabilize the gray positioner and balloon catheter, and advance the inner cannula to deploy the suprarenal stent.
15. Tighten the pin vise.

16. Deflate the balloon and advance the contralateral wire guide into the thoracic aorta.
NOTE: Due to early ipsilateral leg deployment and trigger-wire release, it is suggested to leave the molding balloon in or just above the contralateral limb for use in graft stabilization during placement of the ipsilateral limb.

14.2.2 Docking of Top Cap

1. Loosen the pin vise. (**Fig. 22**)
2. Secure sheath and inner cannula to avoid any movement of these components.
3. Advance the gray positioner over the inner cannula until it docks with the top cap. (**Figs. 23, 24 and 38**)
NOTE: If resistance occurs, slightly rotate gray positioner and continue to gently advance.
4. Retighten the pin vise and withdraw the entire top cap and gray positioner through the graft and through the sheath by pulling on the inner cannula. (**Fig. 39**) Leave the sheath and wire guide in place.
5. Close the Captor Hemostatic Valve on the Flexor introducer sheath by turning it in a clockwise direction until it stops. (**Fig. 26**)

14.2.3 Ipsilateral Iliac Leg Placement and Deployment

NOTE: Ensure that the Captor Hemostatic Valve on the introducer sheath is turned to the open position. (**Fig. 27**)

1. Position the image intensifier to show both the ipsilateral internal iliac artery and ipsilateral common iliac artery.
2. After the introduction of the ipsilateral iliac leg delivery system, inject contrast through the main body sheath to locate the ipsilateral internal iliac artery.
3. Utilize the main body graft wire and sheath assembly to introduce the ipsilateral leg graft.
NOTE: Due to top cap release modification, the main body sheath assembly must be withdrawn to a point 1-2 cm inside the proximal ipsilateral limb. Advance the dilator sheath assembly of the ipsilateral limb into the main body sheath.
NOTE: The molding balloon may be inflated in the contralateral limb of the main body graft if additional graft stabilization is necessary.
CAUTION: Do not inflate the balloon outside of the graft.
NOTE: In tortuous vessels, the position of the internal iliac arteries may alter significantly with the introduction of the rigid wires and sheath systems.
4. Advance the ipsilateral iliac leg delivery system slowly until the ipsilateral iliac leg graft overlaps a minimum of one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the ipsilateral limb of the main body. (**Fig. 40**)
NOTE: If an overlap of greater than three iliac leg stents is required (greater than two iliac leg stents for 37, 39, 54 and 56 mm leg lengths), it may be necessary to consider use of a leg extension into the bifurcation area of the opposite side.
5. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure internal iliac patency.
6. To deploy, stabilize the iliac leg graft with the gray positioner while withdrawing the iliac leg sheath and main body sheath together. (**Figs. 29 and 41**)
7. Under fluoroscopy and after verification of iliac leg graft position, loosen the pin vise, and retract the inner cannula to dock the tapered dilator to the gray positioner. Tighten pin vise. Maintain sheath position while withdrawing gray positioner with secured inner cannula. (**Fig. 42**)
8. Close the Captor Hemostatic Valve on the Flexor introducer sheath by turning it in a clockwise direction until it stops.
9. Re-check the position of the wire guides. Leave sheath and wire guides in place.
10. Remove deflated balloon from contralateral side.

14.2.4 Contralateral Iliac Leg Placement and Deployment

CAUTION: Verify the predetermined contralateral iliac leg is selected for insertion on the contralateral side of the patient before implantation.

1. Position the image intensifier to show both the contralateral internal iliac artery and contralateral common iliac artery.
2. Prior to introduction of contralateral iliac leg delivery system, inject contrast through the contralateral femoral sheath to locate the contralateral internal iliac artery.
3. Introduce the contralateral iliac leg delivery system into the artery. Advance slowly until the iliac leg graft overlaps at least one full iliac leg stent (i.e., proximal stent of iliac leg graft) inside the contralateral limb of the main body. (**Fig. 43**)
NOTE: If using a Spiral-Z Iliac Leg, a radiopaque marker is positioned 30 mm from the proximal end of the iliac leg graft to identify the maximum amount of overlap. To position this graft, advance slowly until the iliac leg graft overlaps at least one stent and not past the radiopaque marker band.
NOTE: If difficulty is encountered advancing the iliac leg delivery system, exchange to a more supportive wire guide. In tortuous vessels, the anatomy may alter significantly with the introduction of the rigid wires and sheath systems.
4. Confirm position of distal end of the iliac leg graft. Reposition the iliac leg graft if necessary to ensure both internal iliac patency and a minimum overlap of one full iliac leg stent (i.e., proximal stent of iliac leg graft, maximum overlap of 1.5 stents) within the main body endovascular graft.
5. To deploy, hold the iliac leg graft in position with the gray positioner while withdrawing the sheath. (**Figs. 18 and 44**) Ensure one stent overlap is maintained.
6. Stop withdrawing the sheath as soon as the distal end of the iliac leg graft is released.
7. Under fluoroscopy and after verification of the iliac leg graft position, loosen the pin vise, retract the inner cannula to dock the tapered dilator to the gray positioner. Tighten pin vise. Maintain sheath position while withdrawing gray positioner with secured inner cannula. (**Fig. 45**)
8. Re-check the position of the wire guides.
9. Perform molding balloon insertion and final angiogram as described in the **Molding Balloon Insertion** section of the product's Instructions for Use.

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© COOK 2019
Bloomington, IN 47404 U.S.A.
Cook Medical Europe Ltd