

Produkto pavadinimas	UMDNS paaiškinimas	Artikelio numeros	UMDN S kodas	MDD klasė
		01801, 01802, 01803, 01804, 01805, 01806, 01807, 01808, 01809, 01810, 01811, 01812, 01813, 01814, 01815, 01816, 01817, 01818, 01819, 01820, 01821, 01822, 01823, 01824, 01825, 01826, 01827, 01828, 01829, 01830, 01831, 01832, 01833, 01834, 01835, 01836, 01837, 01838, 01839, 01840, 01841, 01842, 01843, 01844, 01845, 01846, 01847, 01848, 01850, 01851, 01852, 01853, 01854, 01855, 01856, 01857, 01858, 01859, 01860, 01861, 01862, 01863, 01865, 01866, 01867, 01868, 01869, 01870, 01871, 01872, 01873, 01874, 01875, 01876, 01877, 01878, 01879, 01880, 01881, 01884, 01885, 01886, 01887, 01888, 01889, 01890, 01891, 01892, 01893, 30000, 30001, 30002, 30003, 30003, 30004, 30005, 30016, 30017, 30006, 30007, 30008, 30009, 30010, 30011, 30012, 30013, 30014, 30015, 30016, 30017		
Tip protectors Instrumentenschutztaschen Apsaugos	-	01469, 01470, 01471, 01476, 01477, 01478, 01479, 01492, 01496, 01499	16-222	1
Lubrinol	-	02004	-	1
Container Konteineriai	Sterilization container Konteineriai sterilizacijai	28060, 28061, 28062, 28063, 28064, 28066, 28067, 28068, 28069, 28070, 28071, 28072, 28073, 28074, 28075, 28076, 28077, 28078, 28079, 28080, 28081, 28106, 28107, 28108, 28109, 28110, 28111, 28112, 28113, 28114, 28115, 28166, 28167, 28213, 28228, 28311, 28351, 28421	13-730	1
Cleaning brushes Reinigungsbürsten Šepetėliai instrumentams	Brushes Šepetėliai	09049, 09050, 09051, 09052, 09053, 09054, 09055, 09056, 09058, 09059, 09064, 09066, 09068, 09083, 09084, 09095, 09096, 09098, 09303, 09306, 09307, 09308, 09309, 09310, 09311, 09312, 09313, 09314, 09315, 09316, 09317, 09318, 09319, 09320, 09322, 09324, 09331, 09332, 09333, 09334, 09335, 09336, 09338, 09339, 09340, 09341, 09342, 09343, 09344, 09345, 09346, 09349, 09365, 09366, 09380, 09388, 09458, 09459, 09460, 09478, 09458, 09490, 09491, 09492, 09499, 09745, 09746, 09747, 26098, 26099, 26322, 26323, 26324, 26325, 26326, 26327, 26328, 26329, 26334, 26335, 26336, 26337, 26342, 26441, 26480, 26485, 26491, 26492, 26475, 26476, 26477, 26478, 45906, 45908, 45910, 45912, 45914, 96280.2, 96280.3, 26280.4, 96280.6, 96280.8, 96280.9, 96281.4, 96281.7	10-500	1
	Brush, endoscope cleaning Šepetelis endoskopų valymui	26022, 26020, 26021, 26482, 26025, 26028, 26481, 26027, 26090, 26305, 26026, 09356, 09353, 09352	17-143	1
Tape-Ex	-	02111	-	1
Textile Filter-Tekstiliniai filtrai	Filter, Steam	01070, 01088, 01097, 01108, 01109,	15-035	1

Produkto pavadinimas	UMDNS paaiškinimas	Artikelio numeros	UMDN S kodas	MDD klasė
	Filtras, garai	01112, 01113, 01171, 01173, 01174, 01175, 01179, 01600, 01601, 01610, 01656, 01654, 01659, 01661, 01677, 01667, 01669, 01682, 01683, 01692, 01691, 01694, 02305, 02306, 02314, 01067, 01606, 01616, 01617, 01618, 01619, 01620, 01629, 01631, 01636, 01642, 01665, 02315, 01060, 01061, 01062, 01063, 01064, 01065, 01066, 01071, 01072, 01073, 01074, 01075, 01076, 01077, 01078, 01079, 01080, 01081, 01082, 01083, 01084, 01085, 01086, 01087, 01089, 01090, 01091, 01092, 01093, 01094, 01095, 01096, 01098, 01099, 01100, 01101, 01102, 01103, 01104, 01105, 01106, 01107, 01111, 01114, 01115, 01170, 01172, 01176, 01177, 01178, 01602, 01603, 01604, 01605, 01607, 01608, 01609, 01611, 01612, 01655, 01068, 01653, 01651, 01652, 01657, 01658, 01660, 01177, 01662, 01664, 01106, 01099, 01083, 01078, 01666, 01676, 01674, 01673, 01672, 01668, 01671, 01675, 01684, 01685, 01686, 01687, 01688, 01690, 01693, 01695, 02300, 01696, 02303, 02302, 02307, 01096, 01678, 01663, 01063, 01681, 01680, 02301, 02304, 01670, 02312, 02313, 02311, 02309, 02310, 02308, 01069, 01613, 01614, 01615, 01621, 01622, 01623, 01624, 01625, 01626, 01627, 01628, 01630, 01632, 01633, 01634, 01635, 01637, 01638, 01639, 01640, 01641, 01643, 01644, 01645, 01646, 01647, 01648, 01649, 01650, 01679, 01689, 02316, 02317, 02318, 71170		

Patvirtinta:

Parašas:

Vardas:

Pareigos:



83 -

Instrumentų plovimo šepetėlis (pagal katalogą)

Paskirtis: skirtas profesionaliam naudojimui sveikatos priežiūros įstaigose, medicinos instrumentų valymui.

Naudojimas: instrumentams valyti naudojami tik švarus, nesuteršti šepetėliai. Instrumentų valymui parenkami šepetėliai atitinkantys instrumento sudėtingumą ir konstrukcijos ypatumus.

Priežiūra: po panaudojimo šepetėlį būtina plauti tekančiu vandeniu, kad būtų pašalinti teršalai, esantys tarp šerelių ir ant šepetėlio pagrindo. Išplovus šepetėlį rekomenduojama dezinfekuoti pagal įstaigoje numatytą tvarką. Esant reikalui, šepetėlį galima sterilizuoti (žiūrėti etiketę) garuose (134°C) pagal įstaigoje numatytą tvarką.

Šalinimas: šepetėliai, kurių korpusas nepataisomai išsikreipė, šereliai prarado lankstumą, iškrito, išsilydė yra laikomi netinkamais naudoti, nes gali neigiamai paveikti valomą instrumentą. Netinkami naudoti šepetėliai šalinami kaip nepavojingos medicininės atliekos.

Klasifikavimas: šepetėliai klasifikuojami kaip I klasės medicinos priemonė pagal MDD 93/42/EEC.



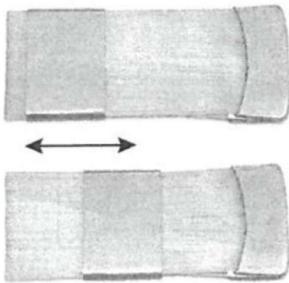
Gamintojas/Tiekėjas: Key Surgical GmbH, Zum Windpark 1, 23738 Lensahn, Vokietija

Platintojas: L.R. Tamulio firma "Meditalika", Radvilų dvaro g.4, LT-48320 Kaunas, Lietuva
tel.: +370 37 22 22 23 e-paštas: info@meditalika.lt www.meditalika.lt

Brush with sliding cover and nylon bristles

REF

pack. unit



26029

length: 55 mm, width: 23 mm

3 pieces

hardness of bristles: adjustable

Cleaning brushes with nylon bristles for MIS instruments

heat resistant: up to 134 °C

washable in washer-disinfector

	brush diameter	brush length	total length	
✓ 09339	2.0 mm	100 mm	300 mm	5 pieces
09342	2.0 mm	60 mm	300 mm	5 pieces
09308	2.5 mm	50 mm	300 mm	5 pieces
09340	2.5 mm	100 mm	300 mm	5 pieces
09343	2.5 mm	80 mm	400 mm	5 pieces
09050	2.5 mm	30 mm	500 mm	5 pieces
09313	2.5 mm	50 mm	610 mm	5 pieces
09341	3.0 mm	100 mm	300 mm	5 pieces
✓ 09344	3.0 mm	100 mm	450 mm	5 pieces
09335	3.5 mm	100 mm	250 mm	5 pieces
09346	3.5 mm	35 mm	500 mm	5 pieces
09052	4.0 mm	100 mm	300 mm	5 pieces
09345	4.0 mm	100 mm	500 mm	5 pieces
09336 ✓	5.0 mm	100 mm	250 mm	5 pieces
09306	5.0 mm	100 mm	300 mm	5 pieces
09318	5.0 mm	80 mm	405 mm	5 pieces
09324	5.0 mm	35 mm	500 mm	5 pieces
✓ 09312	5.0 mm	100 mm	610 mm	5 pieces
09317	5.0 mm	100 mm	650 mm	5 pieces
09365	5.0 mm	100 mm	850 mm	5 pieces
09053	6.0 mm	100 mm	300 mm	5 pieces
09307	7.0 mm	100 mm	300 mm	5 pieces
09051	7.0 mm	80 mm	350 mm	5 pieces
09068	7.0 mm	100 mm	500 mm	5 pieces
09311	7.0 mm	100 mm	610 mm	5 pieces
09316	7.0 mm	100 mm	650 mm	5 pieces
09054	8.0 mm	100 mm	300 mm	5 pieces
09055	10.0 mm	100 mm	300 mm	5 pieces
✓ 09310	10.0 mm	100 mm	610 mm	5 pieces
09315	10.0 mm	100 mm	650 mm	5 pieces
✓ 09309	12.0 mm	100 mm	610 mm	5 pieces
09314	12.0 mm	100 mm	650 mm	5 pieces
09056	15.0 mm	100 mm	300 mm	5 pieces
09066	15.0 mm	100 mm	800 mm	5 pieces
09332	20.0 mm	100 mm	500 mm	5 pieces
09059	30.0 mm	80 mm	300 mm	5 pieces
09334	40.0 mm	60 mm	500 mm	5 pieces

Cleaning brush with synthetic bristles for instruments, double sided

REF		pack. unit
09064	bristle length: 15 and 12 - 15 mm brush length: 10 and 30 mm total length: 155 mm	12 pieces



Cleaning brush, with extra strong synthetic bristles for stubborn soils

09320	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 12 mm brush length: 30 and 45 mm total length: 165 mm	10 pieces
-------	---	-----------



Cleaning brush, with extra strong bristles for stubborn soils

	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 15 mm brush length: 75 mm	
09460	total length: 215 mm	2 pieces
09485	handle of brushes: angled, total length: 235 mm	2 pieces



Cleaning brush, with soft synthetic bristles

09478	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 15 mm bristle material: plastic brush length: 75 mm total length: 215 mm	1 piece
-------	---	---------



Cleaning brush with synthetic bristles for instruments, double sided

REF		pack. unit
09064	bristle length: 15 and 12 - 15 mm brush length: 10 and 30 mm total length: 155 mm	12 pieces



Cleaning brush, with extra strong synthetic bristles for stubborn soils

09320	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 12 mm brush length: 30 and 45 mm total length: 165 mm	10 pieces
-------	---	-----------



Cleaning brush, with extra strong bristles for stubborn soils

	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 15 mm brush length: 75 mm	
09460	total length: 215 mm	2 pieces
09485	handle of brushes: angled, total length: 235 mm	2 pieces



Cleaning brush, with soft synthetic bristles

09478	blue heat resistant: up to 134 °C washable in washer-disinfector bristle length: 15 mm bristle material: plastic brush length: 75 mm total length: 215 mm	1 piece
-------	---	---------



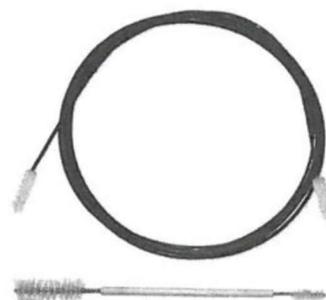
Disposable cleaning brushes for instruments, double sided, single packed

REF		pack. unit
26022	brush diameter: 3.0 mm for flexible scopes and instruments with a channel size of: 2.0 - 2.8 mm brush length: 20 mm total length: 230 cm blue	50 pieces
26020	brush diameter: 5.0 mm for flexible scopes and instruments with a channel size of: 2.8 - 4.2 mm brush length: 20 mm total length: 230 cm red	50 pieces
26021	brush diameter: 6.0 mm for flexible scopes and instruments with a channel size of: 2.8 - 4.5 mm brush length: 20 mm total length: 230 cm turquoise	50 pieces
26482	brush diameter: 6.0 mm for flexible scopes and instruments with a channel size of: 2.8 - 4.5 mm brush length: 20 mm total length: 360 cm green	50 pieces



Disposable set of cleaning brushes for endoscopes, non-sterilizable

09356	Set consisting of: <ul style="list-style-type: none"> • double-sided endoscope brush, blue, total length: 257 cm, brush diameter: 6.5 mm, brush length: 20 mm • double-sided valve cleaning brush, yellow, total length: 150 mm, brush diameter: 4.5 mm, brush length: 18 mm / brush diameter: 10 mm, brush length: 30 mm 	1 set
-------	---	-------



Disposable instrument cleaning brush, double-sided, with pipe cleaner and nylon brush

26305	brush diameter: 3.0 mm For flexible scopes and instruments with a channel size of: 2.0 - 2.8 mm brush length: 20 mm (nylon) / 300 mm (pipe cleaner) total length: 230 cm blue	50 pieces
-------	--	-----------



Disposable instrument cleaning brush, double-sided, with pipe cleaner and nylon brush

pack. unit

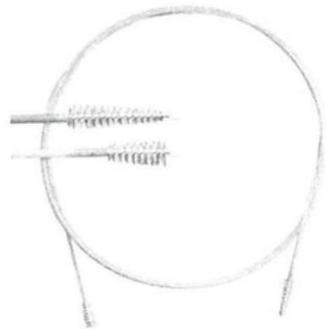
REF



26026

brush diameter: 5.0 mm
 for flexible scopes and instruments
 with a channel size of: 2.8 - 4.0 mm
 brush length: 20 mm (nylon) / 310 mm (pipe
 cleaner)
 total length: 265 cm
 yellow

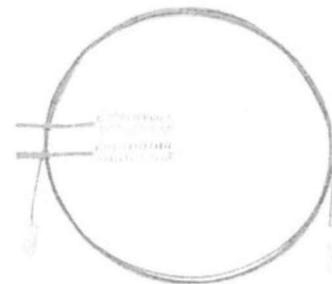
50 pieces

Disposable cleaning brushes for instruments, double sided, tapered brush heads, single packed

26025

brush diameter: 4 - 6 mm / 3 - 5 mm (conical)
 for flexible scopes and instruments
 with a channel size of: 2.8 - 4.5 mm
 brush length: 20 mm
 total length: 230 cm
 yellow

50 pieces

Disposable cleaning brushes for instruments, double sided, with one valve and one channel cleaning brush

26028

brush diameter:
 channel brush 5.5 mm / valve brush 10.0 mm
 for flexible scopes and instruments
 with a channel size of: 2.8 - 4.5 mm
 brush diameter: 5.5 mm / 10.0 mm
 total length: 230 cm
 turquoise

50 pieces

Disposable valve cleaning brush, double-sided

26027

brush diameter: 5.2 mm / 11.0 mm
 total length: 156 mm
 red
 for single use

100 pieces



LUBRINOL

lubrikantas medicinos prietaisams

Dėmesio

ši naudojimo instrukcija skirta profesionaliems naudotojams. Profesionalūs naudotojai privalo taikyti etiketėje ir saugos duomenų lape nurodytas darbų saugos ir sveikatos bei atliekų tvarkymo priemones.

Savybės

Lubrikantas LUBRINOL yra pagamintas sintetinės alyvos pagrindu.

Skirtas nuplautiems medicinos prietaisams (jų jungtims, funkcinėms dalims) sutepti prieš naudojimą arba sterilizacijos procesus.

Priemonė pasižymi optimaliomis lubrikanto savybėmis, tolygiai pasiskirsto ant paviršiaus. Lubrinol yra fiziologiškai saugi priemonė.

Naudojimas

Lubrinol paskirstomas ant medicinos prietaisų paviršiaus purškiant. Prieš purškiant Lubrinol indelį būtina supurtyti.

Lubrinol purškiamas ant medicinos prietaisų paviršiaus per saugų 20-30 cm atstumą.

Lubrinol purškiamas tik ant išvalytų medicinos prietaisų paviršių.

Saugos ir rizikos frazės

nurodytos ant pakuotės

Laikymas

Lubrinol rekomenduojama laikyti kambario temperatūroje ir tik gamintojo pakuotėje.

Pakuotė

400 ml talpos indeliai

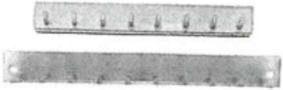
Klasifikavimas: klasifikuojamas kaip I klasės medicinos priemonė pagal MDD 93/42/EEC.



Gamintojas/Tiekėjas: Key Surgical GmbH, Zum Windpark 1, 23738 Lensahn, Vokietija

Platintojas: L.R. Tamulio firma "Meditalika", Radvilų dvaro g.4, LT-48320 Kaunas, Lietuva
tel.: +370 37 22 22 23 e-paštas: info@meditalika.lt www.meditalika.lt

Racks for brushes to hang up brushes

REF		pack. unit	
	for wall mounting, incl. screws and dowels		
06011	white length: 260 mm racks for brushes: 8-fold material: plastic	1 piece	
06012	length: 300 mm racks for brushes: 10-fold material: stainless steel	1 piece	
06021	length: 600 mm racks for brushes: 20-fold material: stainless steel	1 piece	

Leonardo-Wet-Set

	For transport storage of EndoWrist® instruments, keeps the distal end moist and facilitates the reprocessing. Protects the instrument tip during transport.		
50055	Set consisting of: • 3 tubes • 3 caps	1 set	
50058	Set consisting of: • 20 tubes • 20 caps	1 set	

Lubrinol

02004	Lubrinol maintenance oil for surgical instruments has a high creeping effect, excellent lubricating qualities for joints and is physiologically safe. The cleaned instruments have to be sprayed at a distance of 20 - 30 cm and then be sterilized. Lubrinol is suitable for steam sterilization at 134 °C. content: 400 ml	1 aerosol		
-------	---	-----------	---	---

Lubrifiuid

02008	lubricant for maintenance of all Bien-Air handpieces, contra-angles and air motors content: 500 ml	1 aerosol	
-------	---	-----------	---