




ETHICON SECURESTRAP™

5mm Absorbable strap fixation device. For fixation of prosthetic material to soft tissues in various laparoscopic and open surgical procedures, such as Hernia Repair.

Unique design delivers secure and consistent fixation for as long as you need it also in various angles then absorbs - leaving less implanted foreign material.

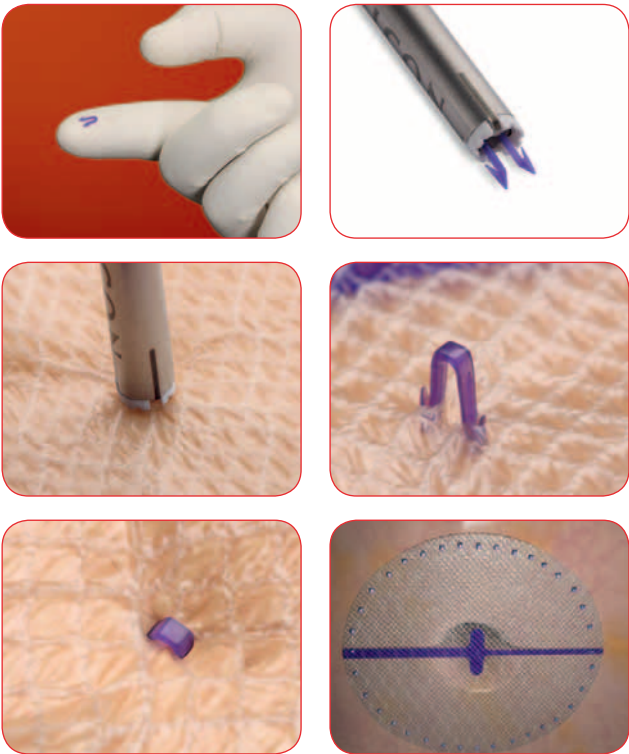

Consistent Strap Deployment

When the trigger is depressed a spring inside the handle is compressed. At full stroke, the spring is released, deploying the strap with the same consistency - reducing user variability.



Absorbable Dyed Straps

Synthetic absorbable straps dyed purple for clear visibility in the abdominal cavity. Minimal absorption in the first two weeks post implantation delivers secure mesh fixation during the tissue integration period. Degradation by hydrolysis significant after the critical tissue integration period has occurred and the polymer absorption is essentially complete at approx. 12 months.

Ref: Prediction of In Vivo Total Absorption Time for ORION Implant Using In Vitro/In Vivo Correlation Degradation Model (AST-2009-0332)



Dimensions	Packaging	Product Code
5mm singel use device with 25 absorbable STRAPS	6 device per box	STRAP25
5mm singel use device with 25 absorbable STRAPS	1 device per box	STRAP25R
5mm singel use device with 12 absorbable STRAPS	6 device per box	STRAP12