


## PrimeLine System



Reusable filter protects against external influences in transport and storage. Lid made of impact-resistant, thermostable plastic.

1.2.5.

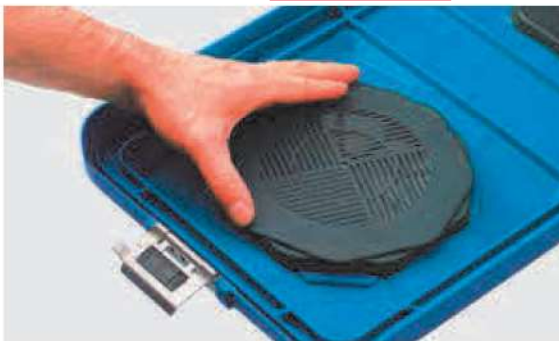

Detachable inner lamella plate for protecting the reusable filter against mechanical damage during the washing process.



The core of the permanent germ barrier is a reusable filter made of PTFE, validated for 5,000 sterilization cycles acc. to EN 868 (2009), Part 8, Annex H.



**Thanks to their material,  
PrimeLine lids are:**

- impact-resistant
- resistant to deformation
- colorfast

This means additional protection and easier organization for the sterile supplies cycle.

**Important processing advice  
for PrimeLine:**

- Do not use solvents such as acetone or benzine
- Final rinse with fully demineralised water
- Make certain that there are no residues from the cleaning process left on the lid
- Do not use rinsing agents in the final rinse.