

SAFE IV THERAPY WITH NEEDLELESS CONNECTORS

SPYDER RANGE

MD class: IIa, CE₀₄₅₉. Please read the instruction for use

Spyder with non-return valves



Non-return valves are one-way valves:

- » it is possible to use the device by gravity, with an IV pump and syringe pump, but withdrawal is impossible
- » they prevent mixing and fluid backflow between the different ways

Warning! The non-return valves are open systems. The tube must be clamped downstream before any disconnection/ connection in order to prevent air entry, especially on the central line⁽³⁾.


REFERENCE	DESCRIPTION	INNER BOX QTY.
KR370	Non-return valve luer lock male luer lock female	200

šakų galuose- besisukančios luer jungtys prijungimui prie kitų sistemų

Dvišakis konektorius sterilus, be adatos


Spyder 2-way extension tubing with male luer lock mobile ring connector

REFERENCE	TUBE MATERIAL	TUBE DIAMETER (MM)	NUMBER OF VALVES	RESIDUAL VOLUME (ML)	INNER BOX QTY.	OUTER BOX QTY.
PY2101KR	PVC	1 x 2.5	2 NRV*	0.53	25	400
PY2301KR		2.5 x 4.1		1.45		

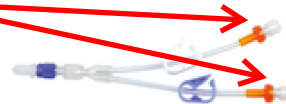


Spyder 3-way extension tubing with male luer lock mobile ring connector

REFERENCE	TUBE MATERIAL	TUBE DIAMETER (MM)	NUMBER OF VALVES	RESIDUAL VOLUME (ML)	INNER BOX QTY.	OUTER BOX QTY.
PY31015KR	PVC	1 x 2.5	3 NRV*	0.88	25	400
PY3101KR				0.76		
PY3301KR		2.5 x 4.1		2.26		

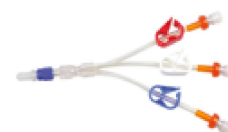

Spyder 4-way extension tubing with male luer lock mobile ring connector

REFERENCE	TUBE MATERIAL	TUBE DIAMETER (MM)	NUMBER OF VALVES	RESIDUAL VOLUME (ML)	INNER BOX QTY.	OUTER BOX QTY.
PY2101KR	PVC	1 x 2.5	4 NRV*	0.96	20	320

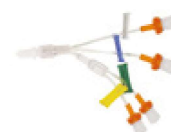

* NRV = Non-return valve
** BDV = Bi-directional valve