



AMPLATZER™ TORQVUE™ AND AMPLATZER™ TREVISIO™ INTRAVASCULAR DELIVERY SYSTEMS



Information contained herein for **DISTRIBUTION outside of the U.S. ONLY.**
Always check the regulatory status of the device in your region.



REINFORCED
POLYMER DELIVERY
SYSTEMS PROVIDING
**VERSATILE, CONTROLLED
DEVICE DELIVERY**

CONTROL

Amplatzer™ TorqVue™ and Amplatzer™ Trevisio™ Intravascular Delivery Systems are designed to provide a controlled pathway to the occlusion site.

Reinforced Polymer Delivery Sheath

- Engineered to provide added kink resistance through an internal stainless steel braid¹

Shaft Material Transitions From a Stiffer Proximal Section to a Softer Distal Section

- Built to ease pushability and maneuverability in tortuous anatomy

Soft, Radiopaque Tip

- Designed to minimize the risk of vessel damage and improve visualization of position

CONFIDENCE

Sheath performance is critical. Amplatzer TorqVue and Amplatzer Trevisio Intravascular Delivery Systems are designed to ease positioning and deployment of Amplatzer™ Occluders.

Radiopaque Distal Tip

- Facilitates precise positioning and accurate device placement

Low-friction PTFE Sheath Lining

- Allows for smooth, controlled delivery

Semi-clear Loader


- Assists in the confirmation of device collapse, flushing and air removal

Delivery Cables and Wires in Varying Gages and Flexibility

- Are specifically designed to aid in the delivery of each Amplatzer Occluder

1. Test(s) performed by and data on file at Abbott.

Information contained herein for **DISTRIBUTION outside of the U.S. ONLY**. Always check the regulatory status of the device in your region.

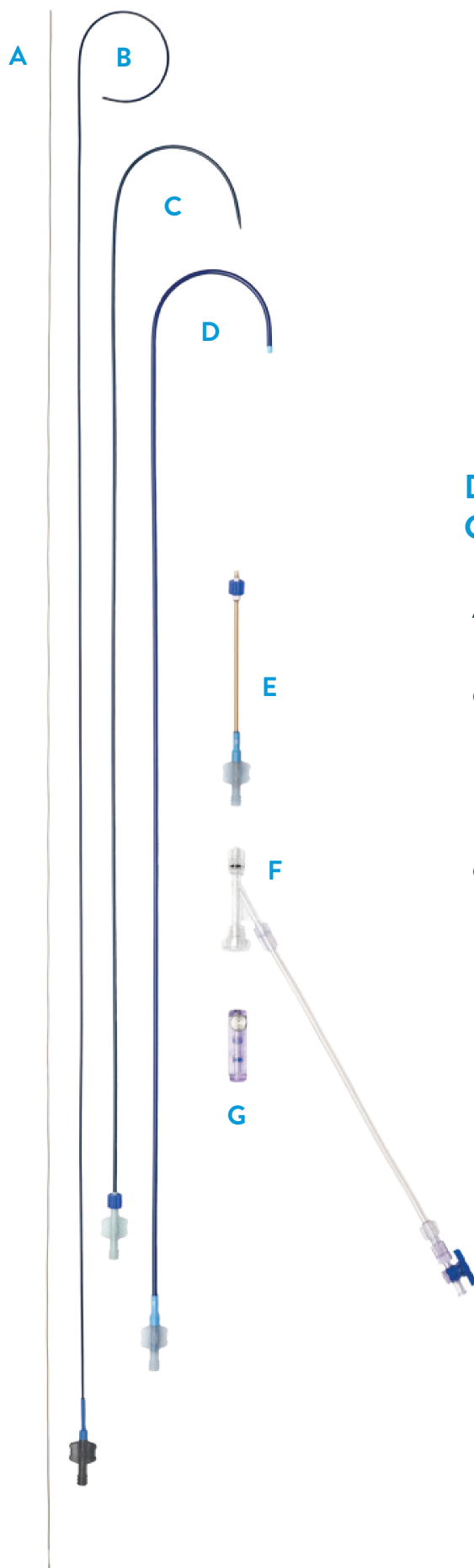

DELIVERY SYSTEM COMPONENTS*

- A. Delivery Cable
- B. Pusher Catheter
- C. Dilator
- D. Delivery Sheath
- E. Loader
- F. Hemostasis Valve
- G. Plastic Vise

*Components vary by delivery system. For a list of components that make up each Amplatzer™ TorqVue™ Amplatzer™ Trevisio™ Intravascular Delivery System, please see the ordering information on the following pages.



AMPLATZER™ TORQVUE™ 45° AND 180° DELIVERY SYSTEMS

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-ITV06F45/60	6	45°	60	Amplatzer™ Septal Occluder and Muscular VSD Occluder
9-ITV07F45/60	7	45°	60	Amplatzer™ Septal Occluder and Muscular VSD Occluder
9-ITV07F45/80	7	45°	80	Amplatzer™ Septal Occluder and Muscular VSD Occluder
9-ITV08F45/60	8	45°	60	Amplatzer™ Septal Occluder, Cribiform Occluder, PFO Occluder and Muscular VSD Occluder
9-ITV08F45/80	8	45°	80	Amplatzer™ Septal Occluder, Cribiform Occluder, PFO Occluder and Muscular VSD Occluder
9-ITV09F45/80	9	45°	80	Amplatzer™ Septal Occluder, Cribiform Occluder, PFO Occluder, Muscular VSD Occluder and P.I. Muscular VSD Occluder
9-ITV10F45/80	10	45°	80	Amplatzer™ Septal Occluder and P.I. Muscular VSD Occluder
9-ITV12F45/80	12	45°	80	Amplatzer™ Septal Occluder
9-ITV13F45/80a	13	45°	80	Amplatzer™ Septal Occluder
9-ITV05F180/60	5	180°	60	Amplatzer™ Duct Occluder and Muscular VSD Occluder
9-ITV06F180/60	6	180°	60	Amplatzer™ Duct Occluder and Muscular VSD Occluder
9-ITV06F180/80	6	180°	80	Amplatzer™ Duct Occluder and Muscular VSD Occluder
9-ITV07F180/80	7	180°	80	Amplatzer™ Duct Occluder and Muscular VSD Occluder
9-ITV08F180/80	8	180°	80	Amplatzer™ Muscular VSD Occluder
9-ITV09F180/80	9	180°	80	Amplatzer™ Muscular VSD Occluder and P.I. Muscular VSD Occluder

Note: Amplatzer™ TorqVue™ 45° and 180° Delivery Systems include a delivery cable, delivery sheath, dilator, hemostasis valve, loader and plastic vise.



AMPLATZER™ TREVISIO™ 45° DELIVERY SYSTEM

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-ATV06F45/60	6	45°	60	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV07F45/60	7	45°	60	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV07F45/80	7	45°	80	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV08F45/60	8	45°	60	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV08F45/80	8	45°	80	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV09F45/80	9	45°	80	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV10F45/80	10	45°	80	Amplatzer™ Septal Occluder, PFO Occluder, Muscular VSD Occluder
9-ATV12F45/80	12	45°	80	Amplatzer™ Septal Occluder
9-ATV13F45/80	13	45°	80	Amplatzer™ Septal Occluder

ORDERING INFORMATION



AMPLATZER™ TORQVUE™ EXCHANGE SYSTEMS

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-EITV09F45/80	9	45°	80	Amplatzer™ Septal Occluder, Cribriform Occluder, PFO Occluder and Muscular VSD Occluder
9-EITV12F45/80	12	45°	80	Amplatzer™ Septal Occluder, Cribriform Occluder, PFO Occluder and Muscular VSD Occluder
9-EITV06F180/80	6	180°	80	Amplatzer™ Duct Occluder and Muscular VSD Occluder
9-EITV08F180/80	8	180°	80	Amplatzer™ Duct Occluder and Muscular VSD Occluder

Note: Amplatzer™ TorqVue™ Exchange Systems include a delivery sheath, dilator, exchange wire, hemostasis valve, loader and plastic vise.



AMPLATZER™ TORQVUE™ LP DELIVERY SYSTEM^A

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	DELIVERY CABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-TVLP4F90/060	4	90°	60	160	Amplatzer™ Duct Occluder II
9-TVLP4F90/080	4	90°	80	195	Amplatzer™ Duct Occluder II
9-TVLP5F90/060	5	90°	60	160	Amplatzer™ Duct Occluder II
9-TVLP5F90/080	5	90°	80	195	Amplatzer™ Duct Occluder II

Note: Amplatzer™ TorqVue™ LP Delivery System includes a delivery cable, delivery catheter, hemostasis valve, loader and plastic vise.



AMPLATZER™ TORQVUE™ LP CATHETER^A

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	DELIVERY CABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-TVLPC4F90/080	4	90°	80	190	Amplatzer Piccolo™ Occluder

Note: Amplatzer™ TorqVue™ Exchange LP Catheter includes a delivery catheter, hemostasis valve, loader and self-sealing hemostasis valve.



AMPLATZER™ TORQVUE™ DELIVERY SYSTEM WITH PUSHER CATHETER

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
9-ITVP07F180/80	7	180°	80	Amplatzer™ Membranous VSD Occluder
9-ITVP08F180/80	8	180°	80	Amplatzer™ Membranous VSD Occluder
9-ITVP09F180/80	9	180°	80	Amplatzer™ Membranous VSD Occluder

Note: Amplatzer™ TorqVue™ Delivery System with Pusher Catheter includes a delivery cable, delivery sheath, dilator, hemostasis valve, loader, plastic vise and pusher catheter.



AMPLATZER™ TORQVUE™ 2 DELIVERY SHEATH

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED USE
9-TV2-05F120	5	Straight	120	Intended to provide a pathway through which devices are introduced within the chambers and coronary vasculature of the heart or in the peripheral vasculature.
9-TV2-06F120	6	Straight	120	
9-TV2-07F120	7	Straight	120	

Note: The Amplatzer™ TorqVue™ 2 Delivery Sheath includes a delivery sheath and dilator.



AMPLATZER™ AMULET™ DELIVERY SHEATH

MODEL/REORDER NUMBER	SHEATH SIZE (F)	TIP ANGLE	USEABLE LENGTH (CM)	RECOMMENDED FOR USE WITH
DS-TV45X45-12F-080	12	45° x 45°	80	Amplatzer™ Left Atrial Appendage Occluder
DS-TV45X45-14F-080	14	45° x 45°	80	Amplatzer™ Left Atrial Appendage Occluder

Note: The Amplatzer™ Amulet™ Delivery Sheath includes a delivery sheath and dilator. The 12 and 14 F delivery sheaths include a flush adapter to facilitate attachment of additional device components.

CAUTION: This product is intended for use by or under the direction of a physician. Prior to use, reference the Instructions for Use, inside the product carton (when available) or at eifu.abbottvascular.com or at medical.abbott/manuals for more detailed information on Indications, Contraindications, Warnings, Precautions and Adverse Events.

Information contained herein for **DISTRIBUTION outside of the U.S. ONLY**. Always check the regulatory status of the device in your region.

Illustrations are artist's representations only and should not be considered as engineering drawings or photographs.

Abbott

3200 Lakeside Dr., Santa Clara, CA 95054, USA Tel: 1 800 227 9902
www.cardiovascular.abbott

™ Indicates a trademark of the Abbott group of companies.

‡ Indicates a third party trademark, which is property of its respective owner.

© 2020 Abbott. All Rights Reserved.

MAT-2002247 v2.0 | Item approved for Global OUS use only.

