


STRUCTURAL HEART THERAPY

# AMPLATZER™ SEPTAL OCCLUDERS



Information contained herein for **DISTRIBUTION outside of the U.S. ONLY**.  
Always check the regulatory status of the device in your region.

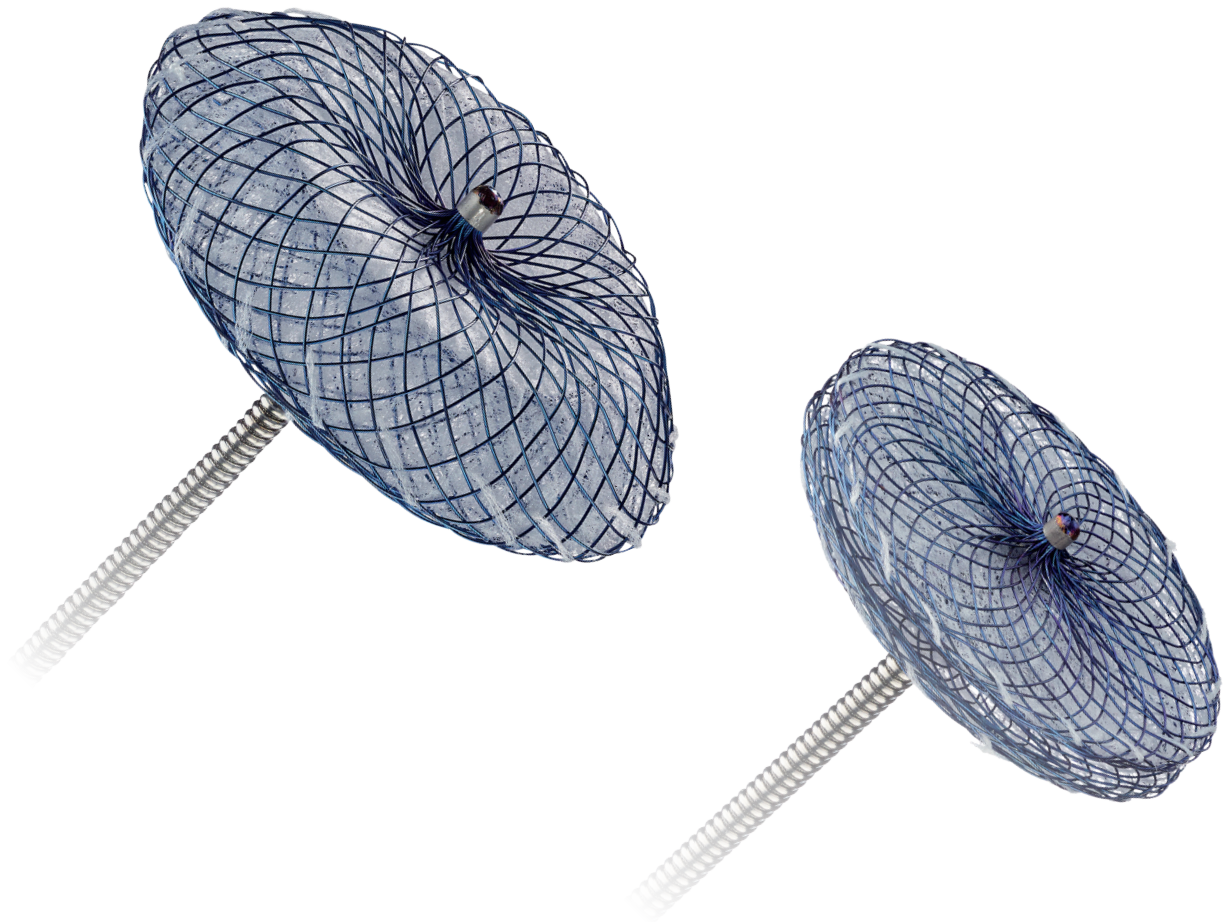

THE PROVEN STANDARD

# FOR TRANSCATHETER ATRIAL SEPTAL DEFECT CLOSURE<sup>1, 2</sup>

Amplatzer™ Septal Occluders are the standard of care for minimally invasive atrial septal defect (ASD) closure. These double-disc occluders are comprised of Nitinol mesh with polyester fabric. They are designed to securely appose the septal wall on each side of the defect and create a platform for tissue in-growth after implantation.

Information contained herein for **DISTRIBUTION outside of the U.S. ONLY**.  
Always check the regulatory status of the device in your region.

## WHEN SIMPLICITY MATTERS

The primary treatment option for closure of secundum atrial septal defects is with transcatheter devices.<sup>2</sup> This is when trust and simplicity matter. The Amplatzer Septal Occluder is the proven standard of care in transcatheter ASD closure.<sup>1,2</sup>

### SAFETY IN NUMBERS

**98.5% Closure Rate:** With no significant residual shunt (less than 2 mm) at 1-year follow-up<sup>1</sup>

**Low Major and Minor Complication Rates:**

Postprocedure and long-term based on 5-year follow-up<sup>1</sup>


### SPECIFICALLY DESIGNED FOR ASD CLOSURE

**Wide Waist:** Centers device and fills the ASD<sup>3</sup>

**Shape-memory Nitinol Mesh:** Designed to securely appose both sides of the septal wall<sup>3</sup>

**Polyester Material:** Promotes occlusion and tissue in-growth<sup>3</sup>

**Precise Placement:** Device can be easily recaptured and redeployed<sup>3</sup>

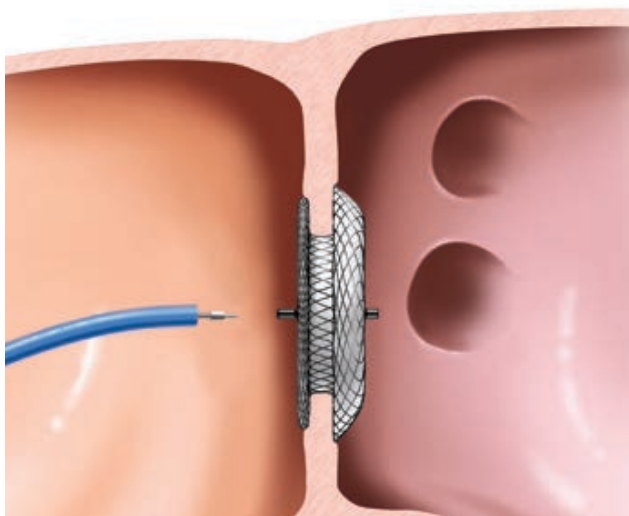

The waist of the Amplatzer™ Septal Occluder fills the defect for optimal occlusion.

## WHEN FLEXIBILITY MATTERS

Not all septal defects are the same. They come in a variety of shapes, sizes and often there can be multiple communications between the left and right atria. These types of atrial septal defects require a special device designed with the flexibility to meet multiple needs.

The addition of the Amplatzer Multifenestrated Septal Occluder - “Cribriform” to the Amplatzer family of occluders enables transcatheter closure for the majority of atrial septal defects.


### SPECIFICALLY DESIGNED FOR MULTI-FENESTRATED ASD CLOSURE

**Narrow Waist:** Allows for placement through a central defect<sup>3</sup>

**Matched Disc Diameters:** Maximizes coverage of multiple fenestrations<sup>3</sup>

**Polyester Material:** Promotes occlusion and tissue in-growth<sup>3</sup>

**Precise Placement:** Device can be easily recaptured and redeployed<sup>3</sup>

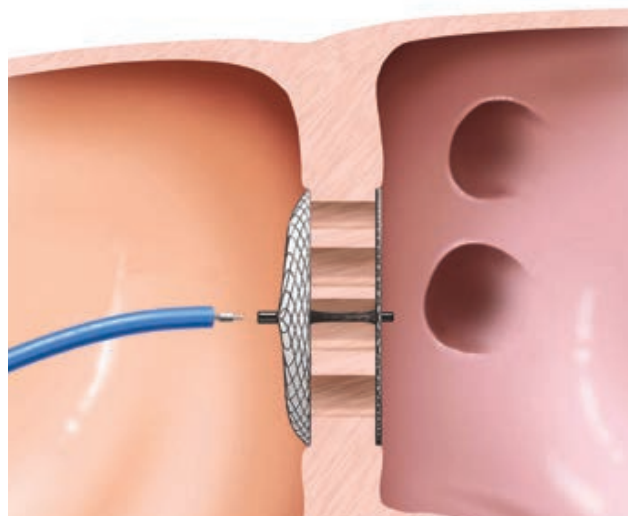

The Amplatzer™ Cribriform Multi-Fenestrated Septal Occluder enables occlusion of the defect by covering the fenestrations with a single device.

Stocking both **Amplatzer™ Septal and Cribriform Occluders** provides you the ability to choose the appropriate device for a wide size-range of defects.

# ORDERING INFORMATION

## AMPLATZER™ SEPTAL OCCLUDER

Model / Reorder Number	Waist Diameter (mm)	Waist Width (mm)	Right Atrial Disc Diameter (mm)	Left Atrial Disc Diameter (mm)
9-ASD-004	4	3	12	16
9-ASD-005	5	3	13	17
9-ASD-006	6	3	14	18
9-ASD-007	7	3	15	19
9-ASD-008	8	3	16	20
9-ASD-009	9	3	17	21
9-ASD-010	10	3	18	22
9-ASD-011	11	4	21	25
9-ASD-012	12	4	22	26
9-ASD-013	13	4	23	27
9-ASD-014	14	4	24	28
9-ASD-015	15	4	25	29
9-ASD-016	16	4	26	30
9-ASD-017	17	4	27	31
9-ASD-018	18	4	28	32
9-ASD-019	19	4	29	33
9-ASD-020	20	4	30	34
9-ASD-022	22	4	32	36
9-ASD-024	24	4	34	38
9-ASD-026	26	4	36	40
9-ASD-028	28	4	38	42
9-ASD-030	30	4	40	44
9-ASD-032	32	4	42	46
9-ASD-034	34	4	44	50
9-ASD-036	36	4	46	52
9-ASD-038	38	4	48	54
9-ASD-040	40	4	50	56

## AMPLATZER™ CRIBRIFORM MULTI-FENESTRATED SEPTAL OCCLUDER

Model / Reorder Number	Right & Left Atrium Disc Diameter (mm)	Waist Width (mm)
9-ASD-MF-018	18	3
9-ASD-MF-025	25	3
9-ASD-MF-030	30	3
9-ASD-MF-035	35	3
9-ASD-MF-040	40	3

### REFERENCE:

1. Amplatzer Septal Occluder Instructions for Use.
2. Kashour TS, Latroche B, Elhoury ME, et al. Successful Percutaneous Closure of a Secundum Atrial Septal Defect through Femoral Approach in a Patient with Interrupted Inferior Vena Cava. Congenital Heart Disease. 2010;5(6):620-623.
3. Test(s) performed by and data on file at Abbott.

**CAUTION:** This product is intended for use by or under the direction of a physician. Prior to use, reference the Instructions for Use, inside the product carton (when available) or at [eifu.abbottvascular.com](http://eifu.abbottvascular.com) or at [medical.abbott/manuals](http://medical.abbott/manuals) for more detailed information on Indications, Contraindications, Warnings, Precautions and Adverse Events.

Information contained herein for **DISTRIBUTION outside of the U.S. ONLY**. Always check the regulatory status of the device in your region.

Illustrations are artist's representations only and should not be considered as engineering drawings or photographs. Photo(s) on file at Abbott.

### Abbott

3200 Lakeside Dr, Santa Clara, CA 95054, USA  
[www.cardiovascular.abbott](http://www.cardiovascular.abbott)

™ Indicates a trademark of the Abbott group of companies.

© 2020 Abbott. All Rights Reserved. MAT-2004327 v1.0 | Item approved for Global OUS use only.

