



Sofia® PLUS
Soft torqueable catheter *Optimized For Intracranial Access*

Sofia® *Distal Access Catheter*

Soft torqueable catheter *Optimized For Intracranial Access*

EXTREMELY TRACKABLE
ASPIRATION AND

LARGE LUMEN
INTRACRANIAL CATHETERS

Distaliai matmenys – OD ne mažiau 0,068 / 1,7;
pasirinktinai 0,0815 / 2,1mm;
Proksimaliniai matmenys – OD ne mažiau 0,068 / 1,7;
pasirinktinai 0,0825 / 2,1mm;
Vidinis diametras ID 0,055" ir 0,070";
Distalinio galiuko forma – tiesi;
Distalinio galiuko ilgis 17 cm ir 19 cm;

6F

Proksimalaus galo ilgis: 96 cm, 98 cm kai kateterio darbinis ilgis 115 cm; 108 cm, 106 cm kai kateterio darbinis ilgis 125 cm; 112 cm kai kateterio darbinis ilgis 131 cm.
Darbinis kateterio ilgis: 115 cm; 125 cm; 131cm;

8. Intermediate kateteriai

Large inner diameter for improved aspiration power and enhanced softness for superior distal navigation


Sofia[®] PLUS
Soft torqueable catheter Optimized For Intracranial Access

Sofia[®] Distal Access Catheter
Soft torqueable catheter Optimized For Intracranial Access


One unit per box, includes shaping mandrel and introducer sheath

PRODUCT CODE	CATHETER SIZE	PROXIMAL OD	DISTAL OD	ID	WORKING LENGTH	DISTAL LENGTH	PROXIMAL LENGTH	DISTAL TIP SHAPE
DA6115ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	115cm	19cm	96cm	Straight
DA6125ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	125cm	19cm	106cm	Straight
DA6131ST	6F	0.0825" / 2.1mm	0.0815" / 2.1mm	0.070"	131cm	19cm	112cm	Straight

SOFIA[®] 6F 115cm

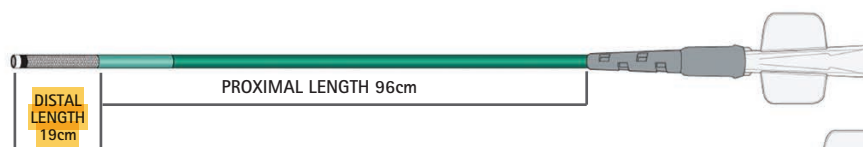

SOFIA[®] 6F 125cm

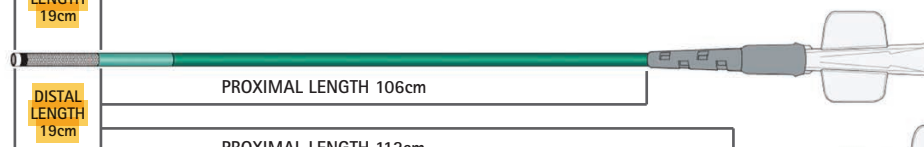

SOFIA[®] 6F 131cm



Design Features

FEATURES	BENEFITS
0.070" Lumen	Wide Inner Lumen for Increased Aspiration Ability
Exceptionally Soft Distal Tip	Allows Smooth Bypass of the Ophthalmic Artery
Hybrid Braid and Coil Design	Provides Superior 1:1 Push Response
Steam Shapeable Tip and Torqueable Shaft	Ability to Steer Distal Tip Past Difficult Bifurcations
Enhanced Kink Resistance	Maintains Distal and Proximal Lumen Integrity

INDICATIONS FOR USE: The SOFIA[®] Catheter is indicated for general intravascular use, including the neuro and peripheral vasculature. The SOFIA[®] Catheter can be used to facilitate introduction of diagnostic or therapeutic agents. The SOFIA[®] Catheter is not intended for use in coronary arteries. Moreover, the SOFIA[®] Catheter is intended for use in removal/aspiration of emboli and thrombi from selected blood vessels in the arterial system, including the peripheral and neuro vasculatures.

MicroVention, Inc., Worldwide Headquarters, 1311 Valencia Avenue, Tustin, CA 92780 USA
Phone: 714.247.8000 Customer Service: 800.990.8368

MicroVention[®]
TERUMO

microvention.com