

72 pirkinys dalis

CoCr stent designed in Japan

CoCr Coronary Stent System  
**Kaname**™

 **TERUMO**®



## Innovative Stent Design

**Cobalt chromium technology**  
& the **unique stent design**  
offer an excellent deliverability & crossability  
to treat even the most difficult lesions

### Enhanced stent flexibility



#### OPEN CELL STENT DESIGN AND THIN STRUTS CONTRIBUTE TO OUTSTANDING FLEXIBILITY

- ▶ **Uniform and homogeneous cell geometry** for enhanced stent flexibility and vessel conformability, to treat distal lesions and tortuous anatomy.
- ▶ **80  $\mu$ m thin struts** and a unique stent design provide enhanced flexibility, whilst maintaining radial strength and reducing stent recoil to provide long term outcomes.



## Excellent lesion access




### SUPERIOR TECHNOLOGY FOR EXCELLENT LESION ACCESS

- ▶ **Low entry profile 0.43 mm** ensures excellent lesion access to treat the most challenging cases and tight lesions.
- ▶ Innovative **enfold crimping technology** improves stent crossing profile and enhances lesion access in challenging anatomies.




- ▶ **The Ryujin™ shaft** provides improved pushability and allows easier access to complex lesions.
- ▶ **Hydrophilic coating** facilitates smooth passage to treat tortuous vessels and tight lesions.

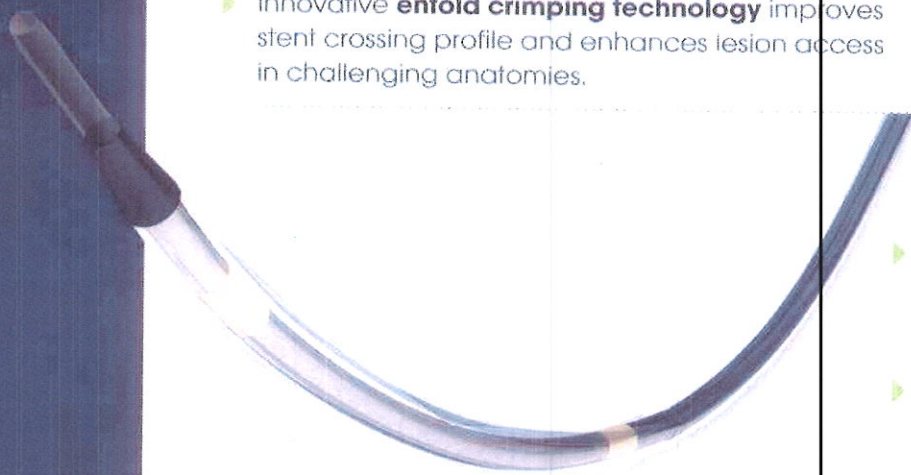

## Terumo's expertise for superior product performance



- ▶ **Open cell area** provides wide side-branch access to treat bifurcation lesions.
- ▶ **CoCr technology** offers outstanding radiopacity allowing a better stent visibility.
- ▶ **Low compliant balloon** for optimal stent expansion to improve long term outcomes.

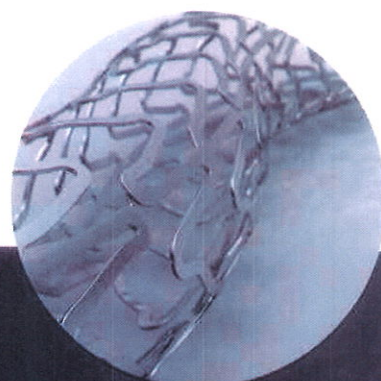

## Specifications

**Stent material:** Cobalt Chromium L605

**Guidewire diameter:** 0.014"/0.36 mm

**Balloon material:** Nylon elastomer


**Nominal pressure:** 9 atm

**Rated burst pressure:** 16 atm – 2.5 to 3.0 mm

14 atm – 3.5 to 4.0 mm

**Average stent crossing profile:** 0.041" for 3.0 mm diameter

**Usable length:** 145 cm

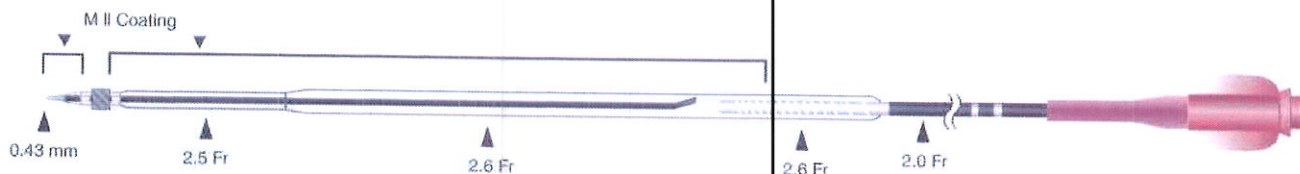

## Ordering information

Diameter (mm)	2.50	2.75	3.00	3.50	4.00
Length (mm)					
9	CS-RD2509KW	CS-RD2709KW	CS-RD3009KW	CS-RD3509KW	CS-RD4009KW
12	CS-RD2512KW	CS-RD2712KW	CS-RD3012KW	CS-RD3512KW	CS-RD4012KW
15	CS-RD2515KW	CS-RD2715KW	CS-RD3015KW	CS-RD3515KW	CS-RD4015KW
18	CS-RD2518KW	CS-RD2718KW	CS-RD3018KW	CS-RD3518KW	CS-RD4018KW
24	CS-RD2524KW	CS-RD2724KW	CS-RD3024KW	CS-RD3524KW	CS-RD4024KW
28	CS-RD2528KW	CS-RD2728KW	CS-RD3028KW	CS-RD3528KW	CS-RD4028KW



**TERUMO CORPORATION**  
441, 2-chome  
Hakogai  
Shinjuku-ku  
Tokyo 161-8522  
Japan  
Tel: +81 (3) 3374 8111  
Fax: +81 (3) 3374 8100

**European Headquarters and Manufacturing**

**TERUMO EUROPE N.V.**  
interleuven 40  
3001 Leuven  
Belgium  
Tel: +32 (0) 30 12 11 11  
Fax: +32 (0) 30 12 11 10

**European Business Unit**

**TERUMO EUROPE  
INTERVENTIONAL SYSTEMS**  
interleuven 40  
Avenue Jean-Pierre Tardieu 1  
CS 10721  
22007 Saint-Quentin-Yvelines Cedex  
France  
Tel: +33 (1) 30 30 30 30  
Fax: +33 (1) 30 30 30 30

**European Sales Offices**

**TERUMO EUROPE N.V.  
BENELUX SALES DIVISION**  
interleuven 40  
3001 Leuven  
Belgium  
Tel: +32 (0) 30 25 25 20  
Fax: +32 (0) 30 25 25 10  
The Netherlands  
Tel: 0800 002 05 00  
Fax: 0800 002 05 01

**LABORATOIRES TERUMO  
FRANCE S.A.**  
Pav. Avenir - Parc Biotech  
1 rue Malraux-Bouquet  
92015 Nanterre Cedex  
France  
Tel: +33 (1) 30 30 30 30  
Fax: +33 (1) 30 30 30 30

**TERUMO DEUTSCHLAND GMBH**  
Wegscheiderstrasse 46  
D-42699 Solingen  
Germany  
Tel: +49 (0) 21 24 19 10 10  
Fax: +49 (0) 21 24 19 10 11

**TERUMO EUROPE ESPANA SL**  
Avenida Juan Carlos I, 15 - 1<sup>a</sup> Planta  
Edificio Torre La Garena  
28904 Alcala de Henares (Madrid)  
Spain  
Tel: +34 (91) 611 11 11  
Fax: +34 (91) 611 11 11

**TERUMO ITALIA SRL**  
Via Salaria Mares 140/141  
I-00132 Roma  
Italy  
Tel: +39 (0) 6 49 11 430  
Fax: +39 (0) 6 49 11 430

**TERUMO DEUTSCHLAND GMBH**  
Friedrichstrasse 47  
D-50769 Eschborn  
Germany  
Tel: +49 (0) 61 93 93 00  
Fax: +49 (0) 61 93 93 00

**TERUMO SWEDEN AB**  
Myraveden 11  
S-205 21 Malmö  
Sweden  
Tel: +46 (0) 40 28 28 28  
Fax: +46 (0) 40 28 28 28

**TERUMO DENMARK**

**Filial of TERUMO SWEDEN AB**  
Tel: +45 7020 00 00  
Fax: +45 7020 00 00

**TERUMO UK LTD**  
Canvey, The Canvey  
Island  
Surrey  
TW20 9AW  
United Kingdom  
Tel: +44 (0) 184 4 70 500  
Fax: +44 (0) 184 4 70 501

**Terumo Representative  
Office in Russia**

**TERUMO EUROPE N.V.**  
Moscow Representative Office  
4, Krasnaya St.  
125109 Moscow  
Russia  
Tel: +7 (495) 244 1000  
Fax: +7 (495) 244 1000

© Registered trademarks  
distributed by Terumo Europe N.V.

[www.terumo-europe.com](http://www.terumo-europe.com)

Illustrations are subject to change  
without notice.