

92 posicija

NaviCross®

0.035" / 0.89 mm Support Catheter

Your first choice peripheral support catheter

Stainless Steel Double Braided Shaft


- ▶ Provides excellent steerability and efficient push transition in complex lesions
- ▶ Enhances torqueability and prevents kinking

Minimal Crossing Profile and Tapered Tip

- ▶ Guarantees a seamless transition between guide wire and catheter facilitating successful lesion access and crossing
- ▶ Angulated tip can be used for selecting the true lumen and navigating through bifurcated vessels

Unique Three Radiopaque Shaft-Markers


- ▶ Feasibility of accurate intraluminal measurement (e.g. assessment of treatment dimensions/positions of balloons and stents), because embedded shaft markers are positioned 1 mm from distal tip
 - 40mm and 60mm from precedent marker band

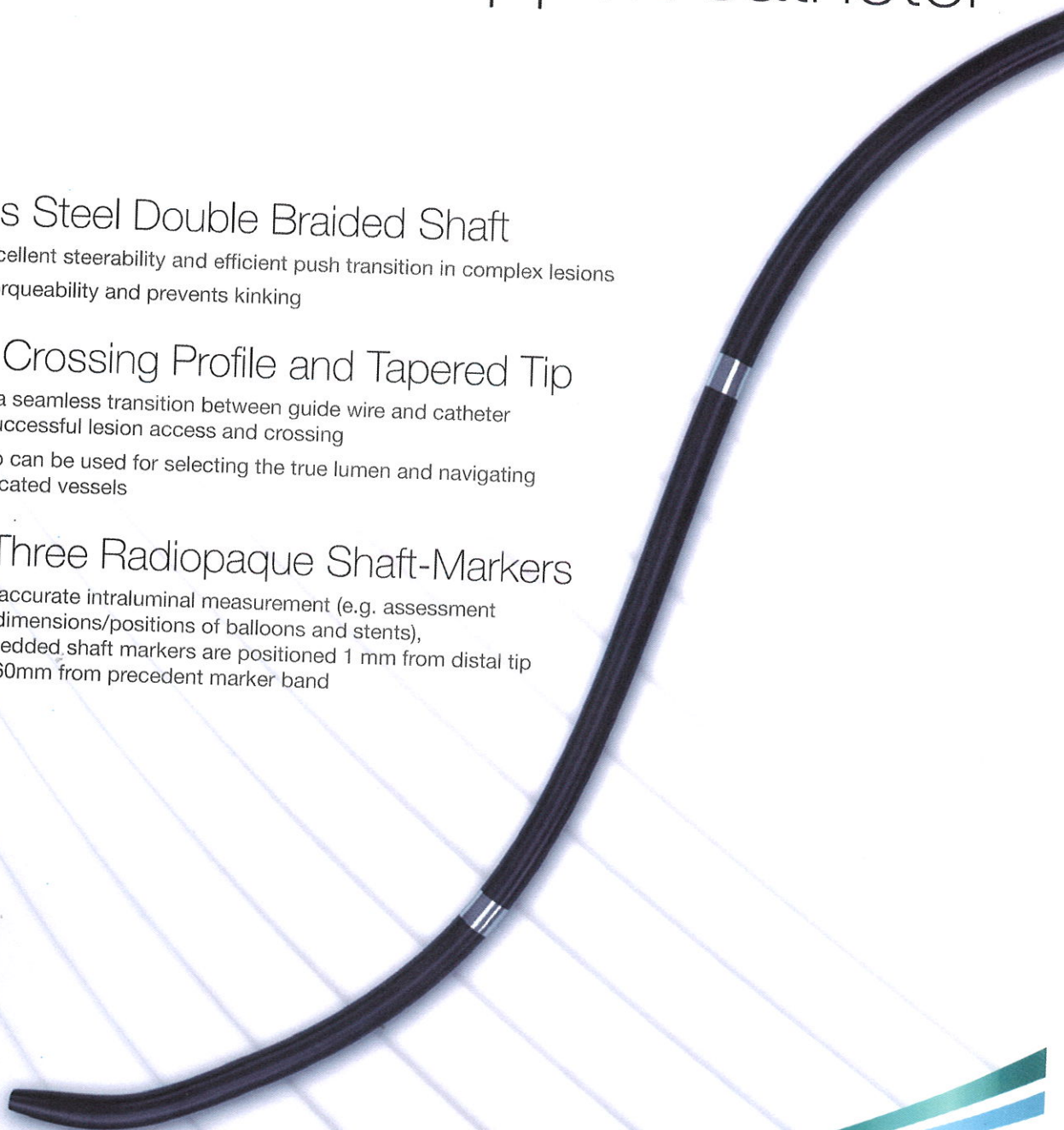

NaviCross®

Specifications

- Double braided stainless steel shaft
- 65, 90, 135, 150 cm lengths
- 0.035" / 0.89 mm wire compatibility
- 1 embedded and 2 swaged radiopaque markers
- Straight and 30° angled tip
- 4 Fr sheath compatibility

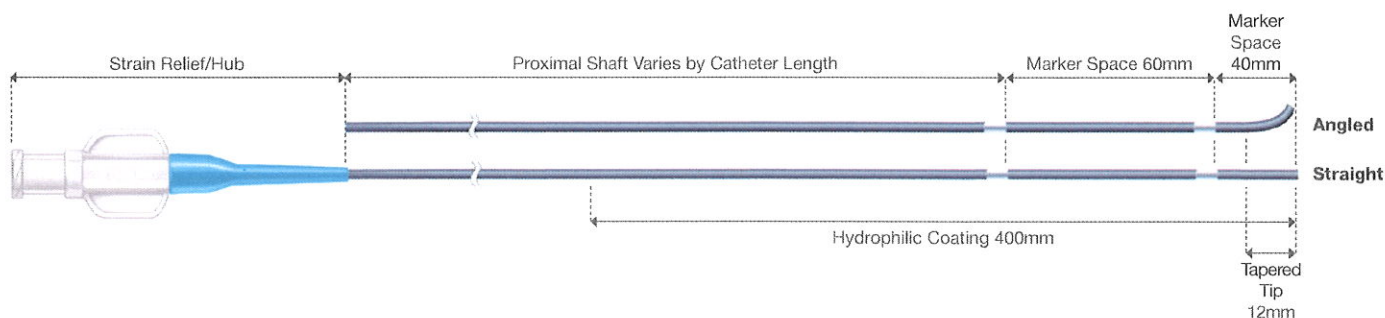

Ordering information

	0.035" / 0.89 mm Wire Compatibility	65 cm	90 cm	135 cm	150 cm
Tip shape	Straight	WS*NS350G3HM	WS*NS35093HM	WS*NS350N3HM	WS*NS35153HM
	30° Angle	WS*NA350G3HM	WS*NA35093HM	WS*NA350N3HM	WS*NA35153HM



Terumo Corporation
44-1, 2-chome, Hatagaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 38 12 11
Fax: +32 16 40 02 49

**Terumo Europe
Interventional Systems
Liaison Office**
Immeuble Krebs
8, rue Jean-Pierre Timbaud
CS 10720
78067 Saint Quentin
en Yvelines Cedex,
France
Tel.: +33 139 30 95 75
Fax: +33 139 30 95 70

European Offices

**Terumo Europe N.V.
Benelux Sales Division**
Interleuvenlaan 40
3001 Leuven,
Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex,
France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Italia S.r.l.
Via Simone Martini, 143 / 145
00142 Roma RM, Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo UK Ltd
Tamesis, The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid), Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda, Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

**Terumo Denmark
Filial of Terumo Sweden AB**
13Fr., 10 Testovskaya Street
Moscow 123317, Russia
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

**Terumo Europe N.V.
Moscow Representative Office**
BC Northern Tower
13Fr., 10 Testovskaya Street
Moscow 123317, Russia
Tel.: +7 495 988 4740
Fax: +7 495 988 4739

**Addresses subject to change,
please consult our Website.**

© Registered Trademark
Published by Terumo Europe N.V.