


*99 piri uo doli's*

# CLIMBER<sup>TM</sup>

PTCA Guiding Catheter

Route to success

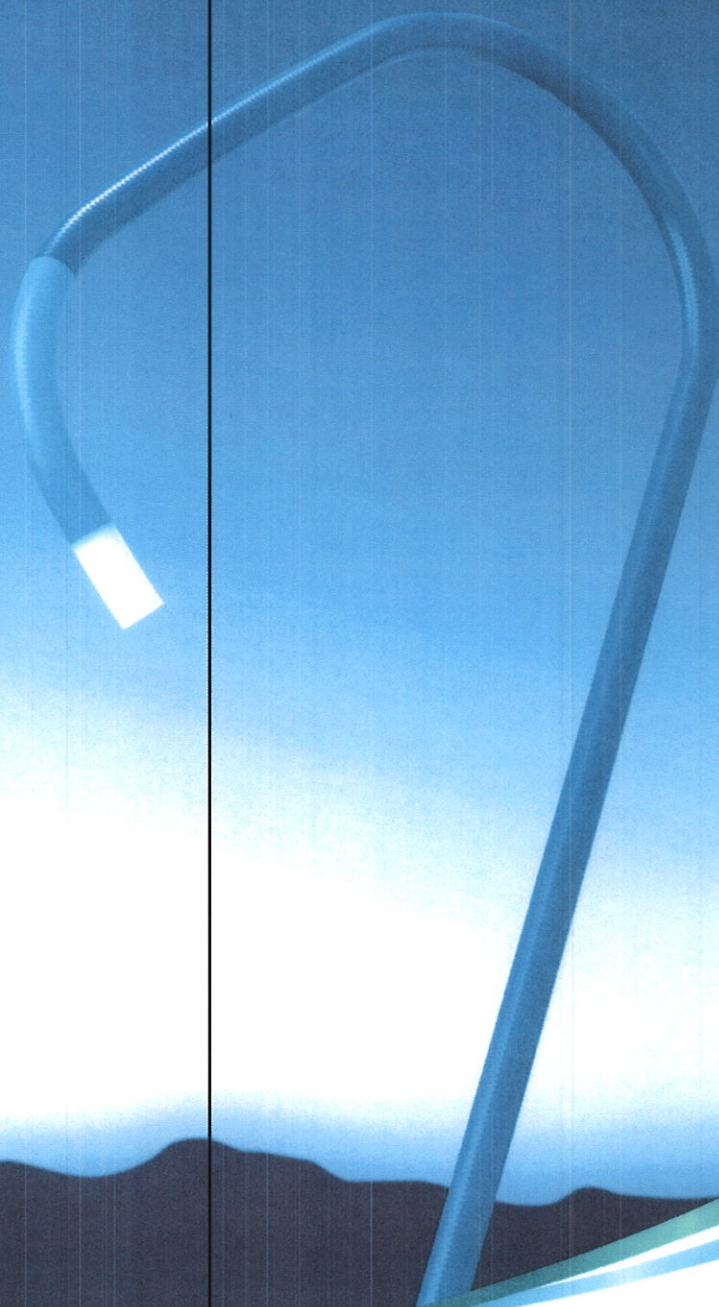

 **TERUMO®**




# New aspects in your daily PCI procedure

## **Hydrophilic coated new guiding catheter**

- **Frictionless insertion** even in tortuous or calcified anatomy
- **Reduce risk of spasm** in transradial procedure
- **Provide tactile 1:1 torque control**
- Distal 7 cm and proximal 25 cm uncoated parts assure a **safe and comfortable procedure** while maintaining the benefit of **hydrophilic coating**



**Hydrophilic coating**

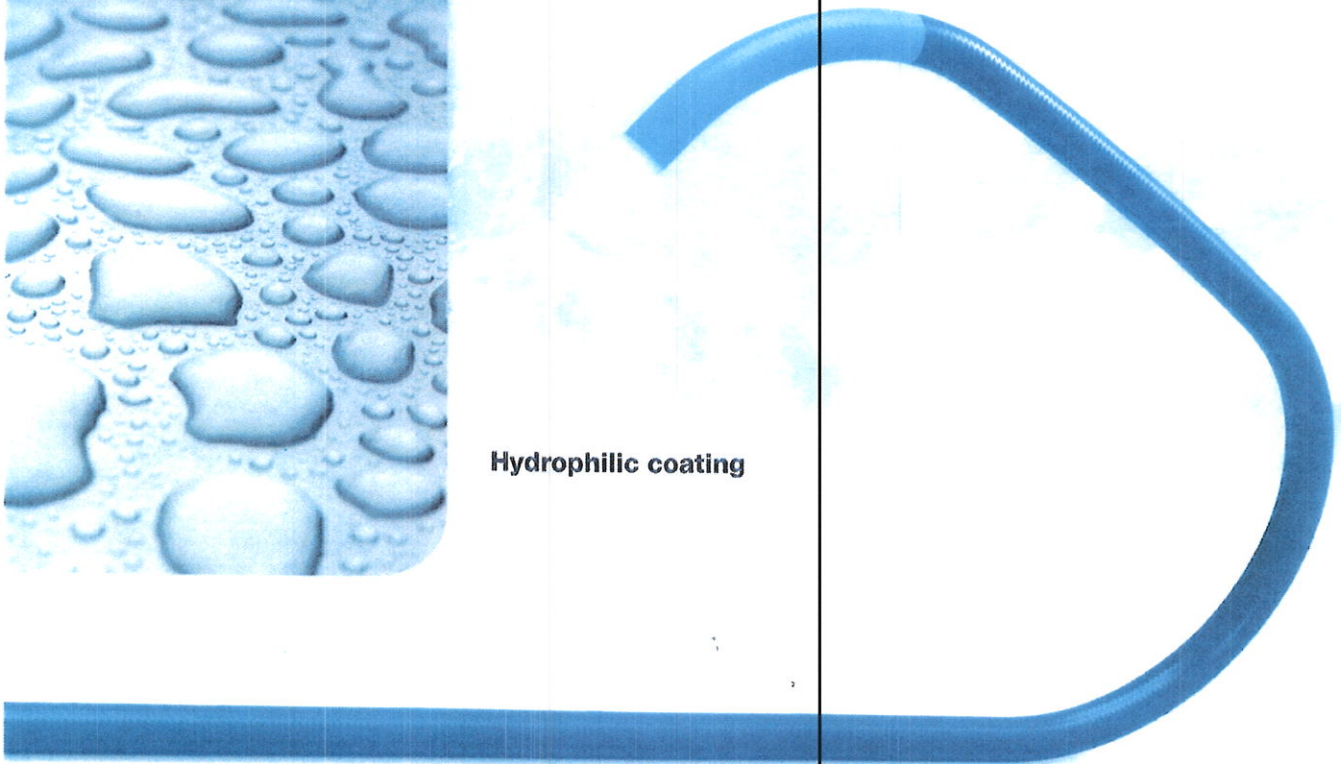

## Ideal combination of largest inner lumen and passive backup support shaft



5Fr I.D.  
0.058"/1.47 mm



6Fr I.D.  
0.071"/1.80 mm



7Fr I.D.  
0.082"/2.08 mm



8Fr I.D.  
0.091"/2.31 mm

- **Extra large 6 Fr inner lumen** allows to keep its system throughout procedure which contribute saving cost and contrast media
- **Flat braiding equipped inside shaft material** provides reliable backup support, pushability and excellent kink resistance
- **Polyamide outer shaft material** ensures excellent shape retention and minimal temperature softening of the catheter even in complex case treatment

## Broad choices of curves maximizes your success

				5Fr / 1.76 mm	6Fr / 2.08 mm	7Fr / 2.37 mm	7Fr / 2.37 mm with Sidehole	8Fr / 2.68 mm	8Fr / 2.68 mm with Sidehole
Standard curves	Amplatz Left	AL	0.75	GC-P5AL070N	GC-P6AL070N	GC-P7AL070N	GC-P7AL070A	GC-P8AL070N	GC-P8AL070A
		AL	1	GC-P5AL100N	GC-P6AL100N	GC-P7AL100N	GC-P7AL100A	GC-P8AL100N	GC-P8AL100A
		AL	1.5	GC-P5AL150N	GC-P6AL150N	GC-P7AL150N	GC-P7AL150A	GC-P8AL150N	GC-P8AL150A
		AL	2	GC-P5AL200N	GC-P6AL200N	GC-P7AL200N	GC-P7AL200A	GC-P8AL200N	GC-P8AL200A
	Amplatz Right	AR	1	GC-P5AR100N	GC-P6AR100N	GC-P7AR100N	GC-P7AR100A	GC-P8AR100N	GC-P8AR100A
		AR	2	GC-P5AR200N	GC-P6AR200N	GC-P7AR200N	GC-P7AR200A	GC-P8AR200N	GC-P8AR200A
		JL	3	GC-P5JL300N	GC-P6JL300N	GC-P7JL300N	GC-P7JL300A	GC-P8JL300N	GC-P8JL300A
		JL	3.5	GC-P5JL350N	GC-P6JL350N	GC-P7JL350N	GC-P7JL350A	GC-P8JL350N	GC-P8JL350A
	Judkins Left	JL	4	GC-P5JL400N	GC-P6JL400N	GC-P7JL400N	GC-P7JL400A	GC-P8JL400N	GC-P8JL400A
		JL	4.5	GC-P5JL450N	GC-P6JL450N	GC-P7JL450N	GC-P7JL450A	GC-P8JL450N	GC-P8JL450A
		JL	5	GC-P5JL500N	GC-P6JL500N	GC-P7JL500N	GC-P7JL500A	GC-P8JL500N	GC-P8JL500A
		JL	6	GC-P5JL600N	GC-P6JL600N	GC-P7JL600N	GC-P7JL600A	GC-P8JL600N	GC-P8JL600A
	Judkins Right	JR	3	GC-P5JR300N	GC-P6JR300N	GC-P7JR300N	GC-P7JR300A	GC-P8JR300N	GC-P8JR300A
		JR	3.5	GC-P5JR350N	GC-P6JR350N	GC-P7JR350N	GC-P7JR350A	GC-P8JR350N	GC-P8JR350A
		JR	4	GC-P5JR400N	GC-P6JR400N	GC-P7JR400N	GC-P7JR400A	GC-P8JR400N	GC-P8JR400A
		JR	4.5	GC-P5JR450N	GC-P6JR450N	GC-P7JR450N	GC-P7JR450A	GC-P8JR450N	GC-P8JR450A
		JR	5	GC-P5JR500N	GC-P6JR500N	GC-P7JR500N	GC-P7JR500A	GC-P8JR500N	GC-P8JR500A
		JR	6	GC-P5JR600N	GC-P6JR600N	GC-P7JR600N	GC-P7JR600A	GC-P8JR600N	GC-P8JR600A



# CLIMBER<sup>TM</sup>

			5Fr / 1.76 mm	6Fr / 2.08 mm	7Fr / 2.37 mm	7Fr / 2.37 mm with Sidehole	8Fr / 2.68 mm	8Fr / 2.68 mm with Sidehole
Extra Backup curves	ExtraBackUp for left coronary	EBU 3	GC-P5EB300N	GC-P6EB300N	GC-P7EB300N	GC-P7EB300A	GC-P8EB300N	GC-P8EB300A
		EBU 3.25	GC-P5EB320N	GC-P6EB320N	GC-P7EB320N	GC-P7EB320A	GC-P8EB320N	GC-P8EB320A
		EBU 3.5	GC-P5EB350N	GC-P6EB350N	GC-P7EB350N	GC-P7EB350A	GC-P8EB350N	GC-P8EB350A
		EBU 3.75	GC-P5EB370N	GC-P6EB370N	GC-P7EB370N	GC-P7EB370A	GC-P8EB370N	GC-P8EB370A
		EBU 4	GC-P5EB400N	GC-P6EB400N	GC-P7EB400N	GC-P7EB400A	GC-P8EB400N	GC-P8EB400A
		EBU 4.25	GC-P5EB420N	GC-P6EB420N	GC-P7EB420N	GC-P7EB420A	GC-P8EB420N	GC-P8EB420A
		EBU 4.5	GC-P5EB450N	GC-P6EB450N	GC-P7EB450N	GC-P7EB450A	GC-P8EB450N	GC-P8EB450A
		EBU 4.75	GC-P5EB470N	GC-P6EB470N	GC-P7EB470N	GC-P7EB470A	GC-P8EB470N	GC-P8EB470A
		EBU 5	GC-P5EB500N	GC-P6EB500N	GC-P7EB500N	GC-P7EB500A	GC-P8EB500N	GC-P8EB500A
		XBRCA 3	GC-P5XR300N	GC-P6XR300N	GC-P7XR300N	GC-P7XR300A	GC-P8XR300N	
Extra BackUp for Right coronary	XBRCA 3.5	GC-P5XR350N	GC-P6XR350N	GC-P7XR350N	GC-P7XR350A	GC-P8XR350N		
	XBRCA 3.75	GC-P5XR370N	GC-P6XR370N	GC-P7XR370N	GC-P7XR370A	GC-P8XR370N		
	XBRCA 4	GC-P5XR400N	GC-P6XR400N	GC-P7XR400N	GC-P7XR400A	GC-P8XR400N		
	XBRCA 4.25	GC-P5XR420N	GC-P6XR420N	GC-P7XR420N	GC-P7XR420A	GC-P8XR420N		
	XBRCA 4.5	GC-P5XR450N	GC-P6XR450N	GC-P7XR450N	GC-P7XR450A	GC-P8XR450N		
Other curves	Multipurpose	MPA -	GC-P5MP000N	GC-P6MP000N	GC-P7MP000N	GC-P7MP000A		
		HockeyStick HS -	GC-P5HS000N	GC-P6HS000N	GC-P7HS000N	GC-P7HS000A		
	IMA	IM -	GC-P5IM000N	GC-P6IM000N	GC-P7IM000N	GC-P7IM000A		
		Bypass	LCB -	GC-P5BP010N	GC-P6BP010N	GC-P7BP010N	GC-P7BP010A	
	RCB -		GC-P5BP020N	GC-P6BP020N	GC-P7BP020N	GC-P7BP020A		
Curves for Radial access	Tiger	TIG 3	GC-P5TR300N	GC-P6TR300N				
		TIG 3.5	GC-P5TR350N	GC-P6TR350N				
		TIG 4	GC-P5TR400N	GC-P6TR400N				
		TIG 4.5	GC-P5TR450N	GC-P6TR450N				
		TIG 5	GC-P5TR500N	GC-P6TR500N				

## Specifications

Outer coating: Hydrophilic

Catheter length: 100 cm

Non coated at Distal 7 cm and Proximal 25 cm



**Terumo Corporation**  
44-1, 2-chome, Hatagaya  
Shibuya-ku, Tokyo 151-0072,  
Japan  
Tel.: +81 3 3374 8111  
Fax: +81 3 3374 8196

### European Headquarters

**Terumo Europe N.V.**  
Interleuvenlaan 40  
3001 Leuven, Belgium  
Tel.: +32 16 38 12 11  
Fax: +32 16 40 02 49

**Terumo Europe  
Interventional Systems  
Liaison Office**  
Immeuble Krebs  
8, rue Jean-Pierre Timbaud  
CS 10720  
78067 Saint Quentin  
en Yvelines Cedex,  
France  
Tel.: +33 139 30 95 75  
Fax: +33 139 30 95 70

### European Offices

**Terumo Europe N.V.  
Benelux Sales Division**  
Interleuvenlaan 40  
3001 Leuven, Belgium  
Tel.: +32 16 39 25 80  
Fax: +32 16 39 25 99  
The Netherlands:  
Tel.: 0800 0220396  
Fax: 0800 0220414

**Laboratoires Terumo France S.A.**  
Parc Ariane - Bât. Uranus  
1, rue Hélène Boucher  
78284 Guyancourt Cedex, France  
Tel.: +33 130 96 13 00  
Fax: +33 130 43 60 85

**Terumo Italia S.r.l.**  
Via Simone Martini, 143/145  
00142 Roma RM, Italy  
Tel.: +39 0651 96 14 20  
Fax: +39 0650 30 407

**Terumo UK Ltd**  
Tamesis, The Causeway  
Egham, Surrey  
TW20 9AW, United Kingdom  
Tel.: +44 1784 476 200  
Fax: +44 1784 476 201

**Terumo Deutschland GmbH**  
Ludwig-Erhard-Strasse 6  
65760 Eschborn, Germany  
Tel.: +49 6196 80 230  
Fax: +49 6196 80 23 200

**Terumo Deutschland GmbH  
Zweigniederlassung Österreich**  
Liebmannstrasse F10 301  
2345 Brunn am Gebirge, Österreich  
Tel.: +43 2236 378 020  
Fax: +43 2236 378 020 100

**Terumo Deutschland GmbH  
Zweigniederlassung Spreitenbach**  
Bodenackerstrasse 3  
8957 Spreitenbach, Switzerland  
Tel.: +41 56 419 10 10  
Fax: +41 56 419 10 11

**Terumo Europe España SL**  
Avda. Juan Carlos I, Nº 13 - 7ª Planta  
Edificio Torre La Garena  
28806 Alcalá de Henares (Madrid), Spain  
Tel.: +34 9021 01 298  
Fax: +34 9021 01 358

**Terumo Denmark  
Filial of Terumo Sweden AB**  
Tel.: +45 7020 93 80  
Fax: +45 7020 94 80

**Terumo Sweden AB**  
Nya Varvet, Byggnad 90  
Sven Källfells Gata 18  
426 71 Västra Frölunda, Sweden  
Tel.: +46 3174 85 880  
Fax: +46 3174 85 890

**Terumo Russia LLC**  
BC Northern Tower  
13 Floor, 10 Testovskaya Street  
Moscow 123317, Russia  
Tel.: +7 495 988 4740  
Fax: +7 495 988 4739

**Terumo Corporation Dubai Branch  
(Middle East and Africa)**  
Al Masraf Tower, 22nd Floor  
P.O. Box 20291 Dubai  
U.A.E.  
Tel.: +971 4 221 2220  
Fax: +971 4 221 3330

*Addresses subject to change,  
please consult our Website.*

© Registered Trademark  
Published by Terumo Europe N.V.