


Single Use Polypectomy Snares

SnareMaster[™]
SnareMaster Soft[™]

More snares, more variety, more choices

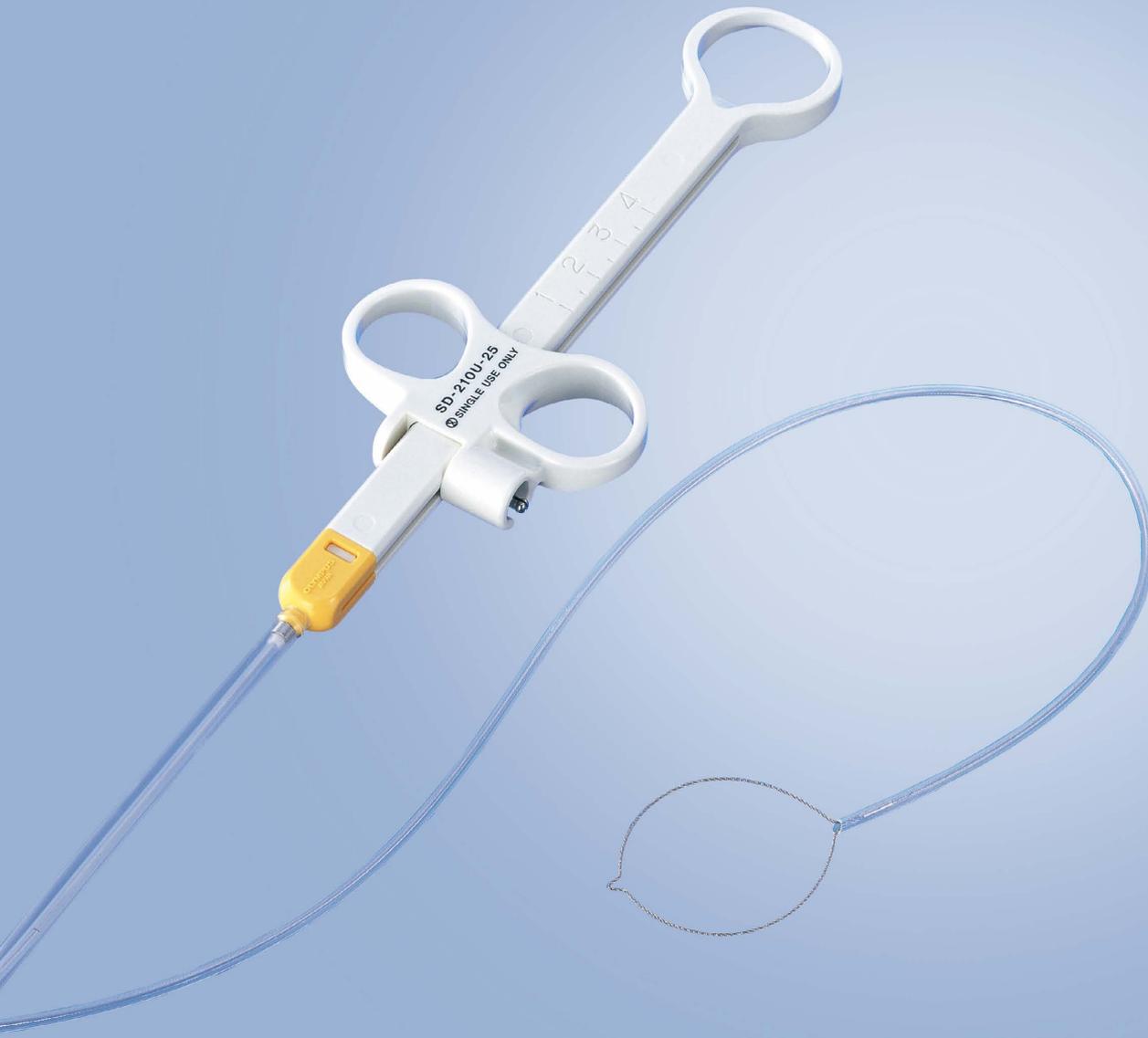

A Wider Choice of Loop Designs and Sizes for Varying Clinical Situations

SnareMaster™
SnareMaster Soft™

Oval Snares

SD-210U-10/15/25

Optimal stiffness and shape to precisely capture polyps

SnareMaster™

Firm snare

Resists warping when pressed against the mucosa, supporting accurate placement and a firm grip of the lesion

Wide-opening loop


Wide loop easily captures a larger resectable area

Thick braided wire

Optimal coagulation to help reduce bleeding

Diverse loop sizes

Suitable for various procedures and polyp sizes

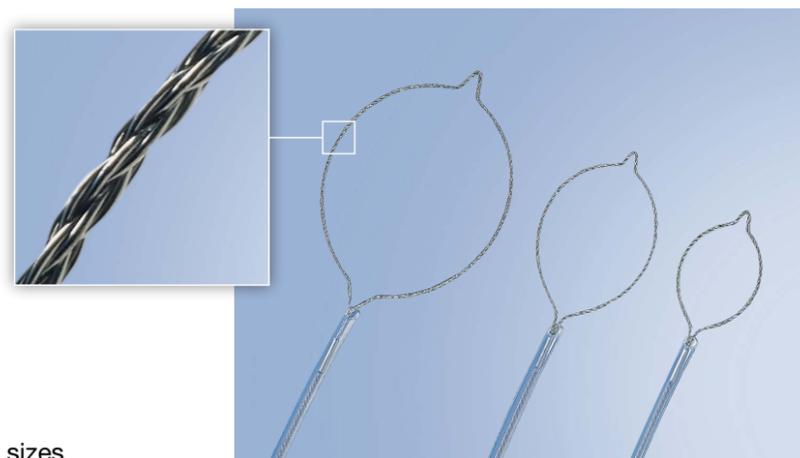

Soft Wire Snares

SD-240U-1U-10/15/25

Soft and flexible for smooth manipulation

SnareMaster Soft™

Thin, pliable 0.4mm wire


Loop opens and closes with less effort for easy manipulation and is designed for sharper, quicker cutting

Smaller distal end (triangle tip)

Reduces burning of non-target areas

Diverse loop sizes

Suitable for various procedures and polyp sizes

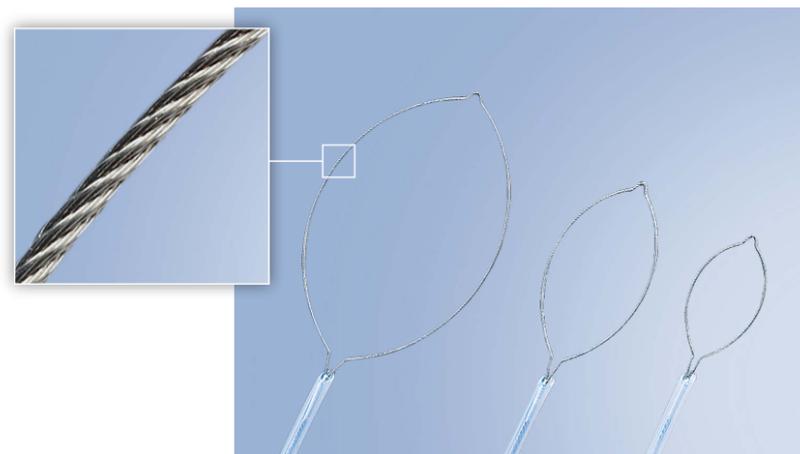

Snares designed for special purposes

Spiral Snare

SD-230U-20

Useful for EMR in lower GI tract

SnareMaster™

Twisted wire

Spiral wire designed to grasp flat mucosal lesions with less slippage and improved coagulation

Firm snare


Resists warping when pressed against the mucosa, supporting accurate placement and a firm grip of the lesion

Wide-opening loop

Wide loop easily captures the a larger resectable area

Thick wire

Opens and closes smoothly for slow, controlled cutting, optimal coagulation and less bleeding

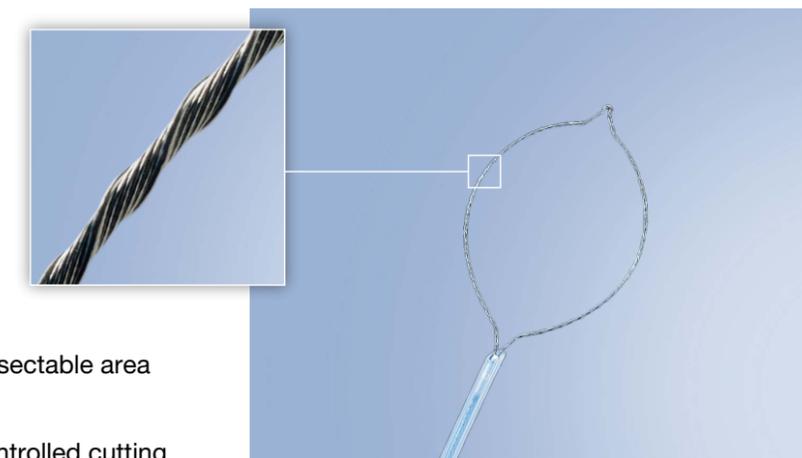

Crescent Snare

SD-221U-25

Specially designed for EMRC

SnareMaster™


Thin, pliable wire

Thin, flexible loop wire fits easily on rim within distal EMRC cap



Highly adjustable size

Special crescent design enables loop size adjustment to grasp lesions of any size

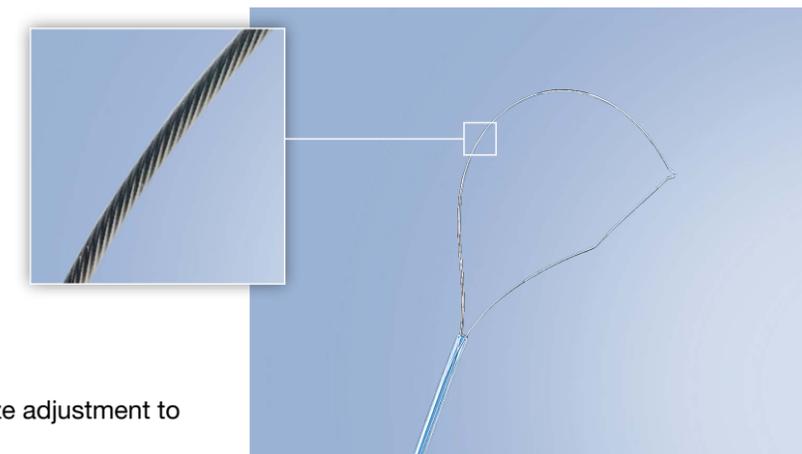

DISPOSABLE ELECTROSURGICAL SNARE

Specifications

Model	Loop design	Working length	Loop diameter	Minimum channel diameter	Diameter of snare loop wire
SD-210U-10	Oval	2300 mm	10 mm	2.8 mm	0.47 mm
SD-210U-15			15 mm		
SD-210U-25			25 mm		
SD-221L-25	Crescent	1650 mm	25 mm	2.0 mm	0.3 mm
SD-221U-25		2300 mm			
SD-230U-20	Oval Spiral	2300 mm	20 mm	2.8 mm	0.48 mm
SD-240U-10	Oval Soft	2300 mm	10 mm	2.8 mm	0.4 mm
SD-240U-15			15 mm		
SD-240U-25			25 mm		

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.



OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of
sales and distribution locations visit:
www.olympus.com