

Safeflow

Premium Safety



IV Accessories

The valve that makes infusion safety simple

B | BRAUN
SHARING EXPERTISE

When infusion safety really matters

Safe, simple and cost-effective needle-free infusion system

Safety first


Safeflow needle-free infusion system provides safe options for the preparation and administration of IV therapies. Safeflow's valve promotes clinical safety first and optimizes clinical infusion flow rates.

Versatile and innovative system

Safeflow needle-free infusion system provides needle-free products for all your IV therapy needs hospital wide, including IV administration sets and extension lines.

Safeflow products feature easy access for Luer-Slip and Luer-Lock connections providing clinicians with needle-free access for faster, easier drug delivery, blood sampling and infusions.


Perfect for your clinical success!

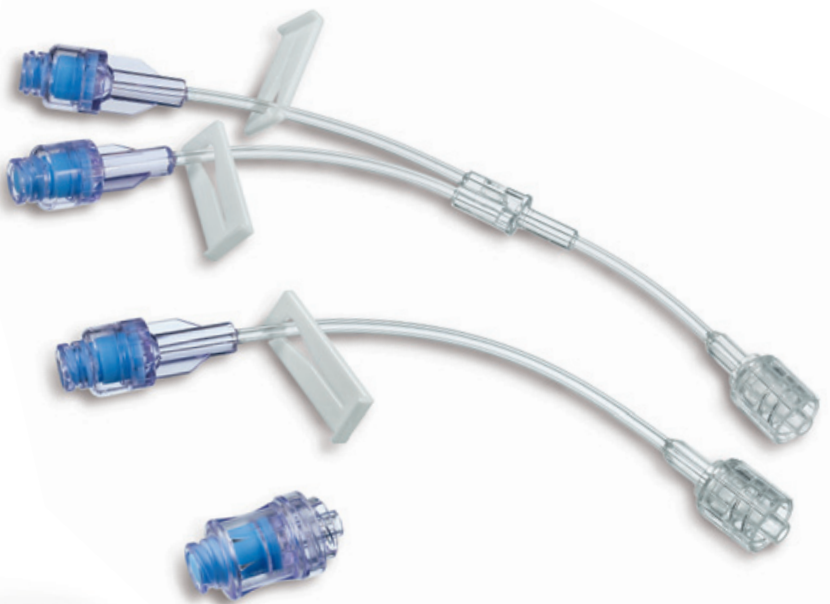

Extension lines

Safeflow is always a closed system. The smallbore extension line allows manipulation away from the injection side, which ensures treatment of patient by infusion, transfusion, aspiration of blood

samples and intermittent injection. Tubings with smaller internal diameter are minimizing the contained volume and can be used for pediatric/neonate patients.

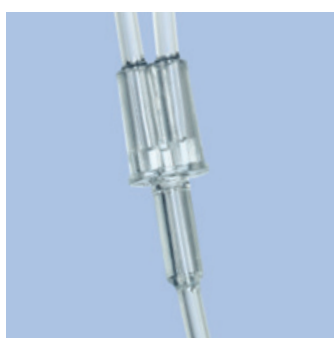

PVC/DEHP-free Neutrapur® tubing


Particularly patient-friendly.
No plasticizer leachable.
Ecological disposal.



Safeflow valve

Offers convenient and safe
needle-free access for injection,
aspiration or parallel
infusions as a closed system.

Ploščia
dezinfekuojava
membrana

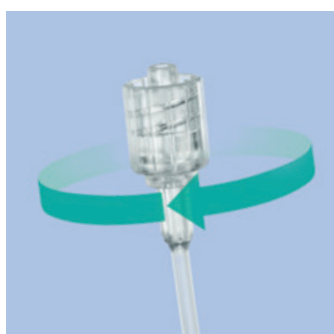

Rotating Luer-Lock connector

Permits safer and more convenient
connection without
twisting the tube.




Back check valve

Prevents backflow and therefore it offers higher safety
during parallel pressure and
gravity infusions.

User benefits

- Swabbable
- Compatible with Luer-Slip and Luer-Lock
- Allows multiple access for direct aspiration or fluid injection
- Suitable for all IV Systems

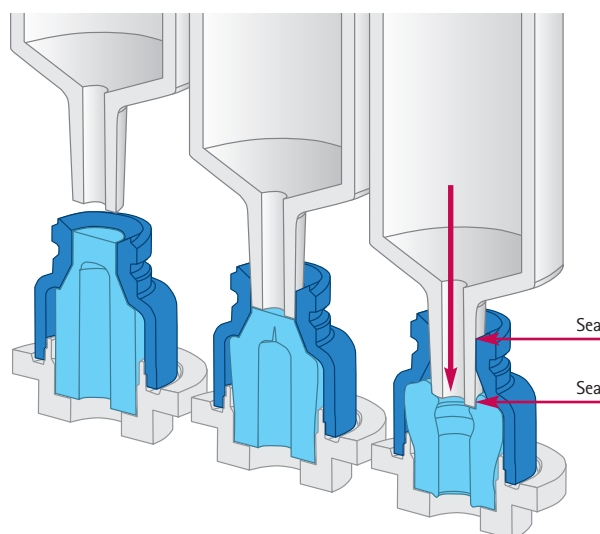







Optimal flow rates

The Safeflow Luer activated valve features an innovative design that provides a straight fluid pathway resulting in uninterrupted fluid flow providing flow rates of 360 ml/min.

Making infusion safety simple

Vamzdelio ilgis

IV Accessories

Product Type	Description		Latex free	DEHP free	PVC free		Sales unit/pcs.	Code No. (REF)
	Safeflow valve		•	•	•		50	409100H
	Extension Sets	Tubing inner Ø (mm)	Latex free	DEHP free tubing	PVC free tubing	Tubing length (cm)	Sales unit/pcs.	Code No. (REF)
	Smallbore Extension Set with Safeflow valve	1.0	•	•	•	10	50	4097154
	Smallbore Y-Extension Set with two Safeflow valves	1.0	•	•	•	12	50	4097145
	Smallbore Y-Extension Set with two Safeflow valves and back check valve	1.0	•	•	•	12	50	4097148
	IV Administration Sets							
	Intrafix® SafeSet with Safeflow	3.0	•			210	100	4063004
	Intrafix® Primline Comfort with Safeflow	3.0	•			180	100	4062158
	Infusomat® Space Line with Safeflow	3.0	•	•	•	230	100	8700110SP
	Infusomat® Space Line, Dosifix®	3.0	•	•	•	270	40	8250245SP



Luer-Lock



Luer-Slip

Product Data of Safeflow

Flow rate	360 ml/min
Latex-free	✓
Luer-Slip and Luer-Lock	✓
Complete system	✓
Priming volume	0.09 ml
Number of accesses	up to 200 accesses
Time in use	24 hours or as per hospital protocol

Tėkmės greitis

B | BRAUN
SHARING EXPERTISE

B. Braun Melsungen AG

34209 Melsungen
Germany
Tel +49(0) 56 61 71-0
www.bbraun.com