

# RENO Plasma Sterilizer


EXPERT IN INFECTION CONTROL

Low Temperature  
Plasma Sterilizer

## RENO Non Woven sheet

RENOSEM offers Non Woven sheet to wrap device containers or medical instruments for H<sub>2</sub>O<sub>2</sub> sterilization process.

Non-woven sheet package is useful for safe handling and storage of medical devices after sterilization process. The sheets of sterilization wrap are square or rectangular fabric produced using a three-layer SMS (spunbond-meltblown-spunbond) process. The wrap fabric is composed of polypropylene with addition of less than 2 % by weight of phthalocyanine blue pigment, less than 1 % by weight of titanium dioxide pigment, and less than 0,009 % by weight of a potassium phosphate anti-static treatment. Products contains no cellulose. The white sheet has the same material composition but contains no blue pigment. The wrap allows a sterilized package to be opened aseptically.

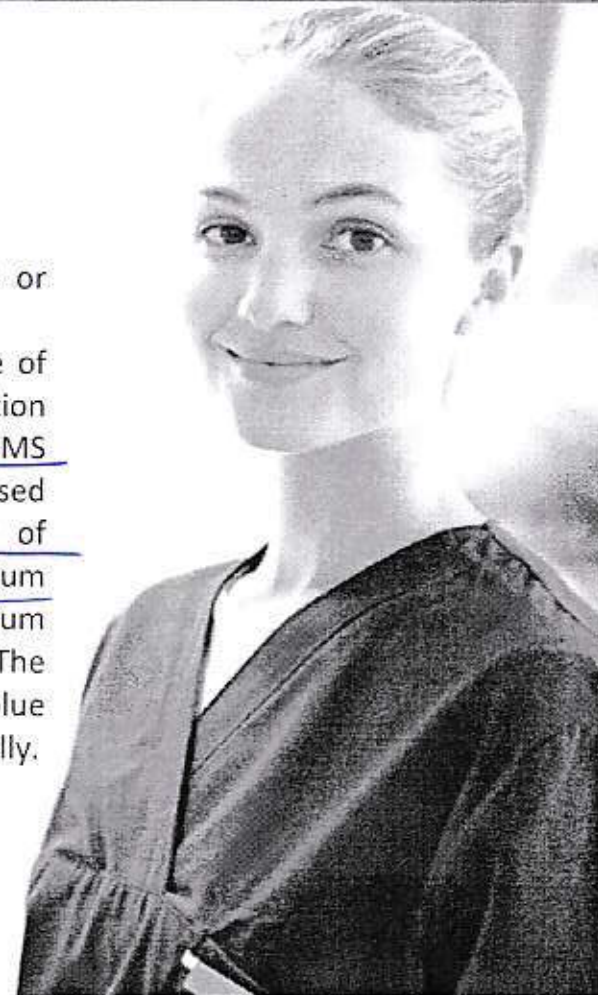

2022-10-13  
KOPIJA TIKRA

Vadybininkas  
Andrius Markūnas



[www.renosem.com](http://www.renosem.com)

# RENO Plasma Sterilizer


EXPERT IN INFECTION CONTROL

Low Temperature  
Plasma Sterilizer

## SPECIFICATION

Material

Composed of polypropylene

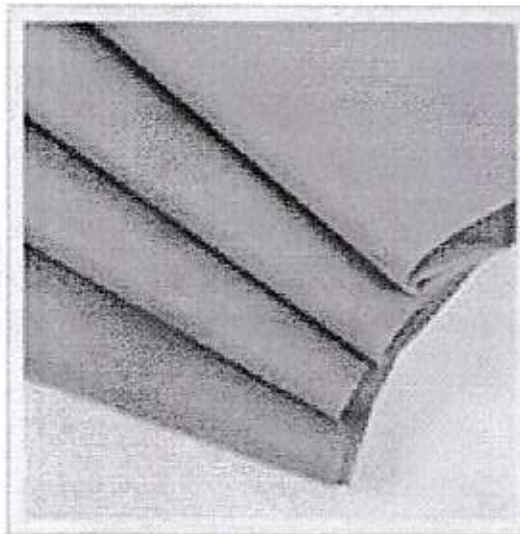

Parts No.

Product Description

Packing

Non-Woven Sheet (38\*38 cm)

38 cm x 38 cm

25/pack

Non-Woven Sheet (45\*45 cm)

45 cm x 45 cm

25/pack

Non-Woven Sheet (60\*60 cm)

60 cm x 60 cm

25/pack

Non-Woven Sheet (91\*91 cm)

91 cm x 91 cm

25/pack

Non-Woven Sheet (1140\*1140 cm)

114 cm x 114 cm

25/pack



2022-10-13  
KOPIJA TIKRA

Vadybininkas  
Andrius Markūnas




[www.renosem.com](http://www.renosem.com)

# RENO Plasma Sterilizer

EXPERT IN INFECTION CONTROL

## HOW TO PACK TYVEK POUCHES

1. Prepare the suitable tray for medical devices and place it in the tray.  
And put a CI Strip inside of the tray.

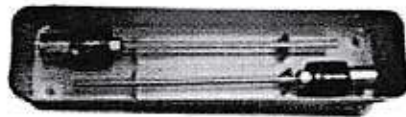


2. Put the tray including medical devices on the Non-woven sheet. And pack a tray as below.



2022-10-13  
KOPĖJA TIKRA

Vadybininkas  
Andrius Markūnas

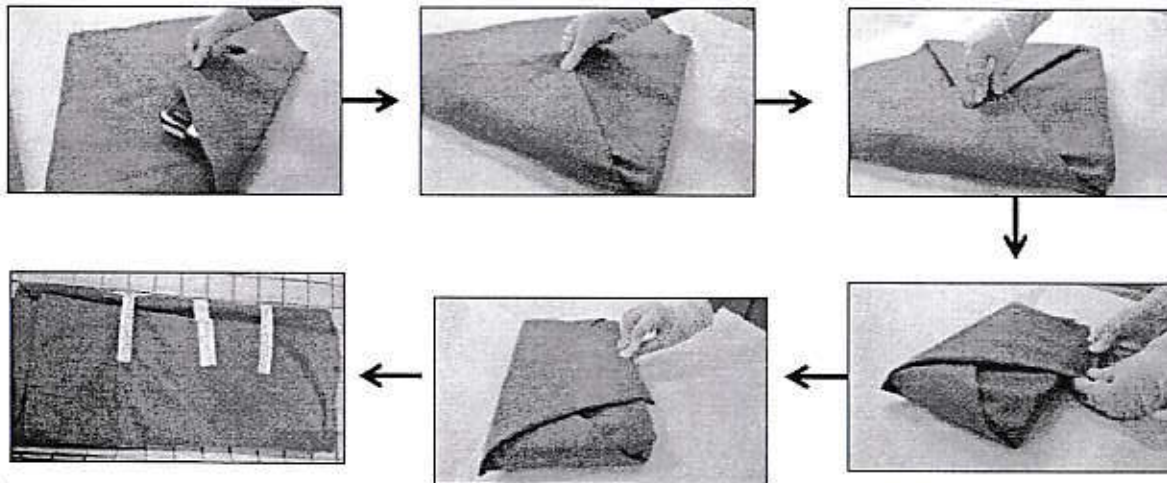

[www.renosem.com](http://www.renosem.com)