


Wire Selection Chart

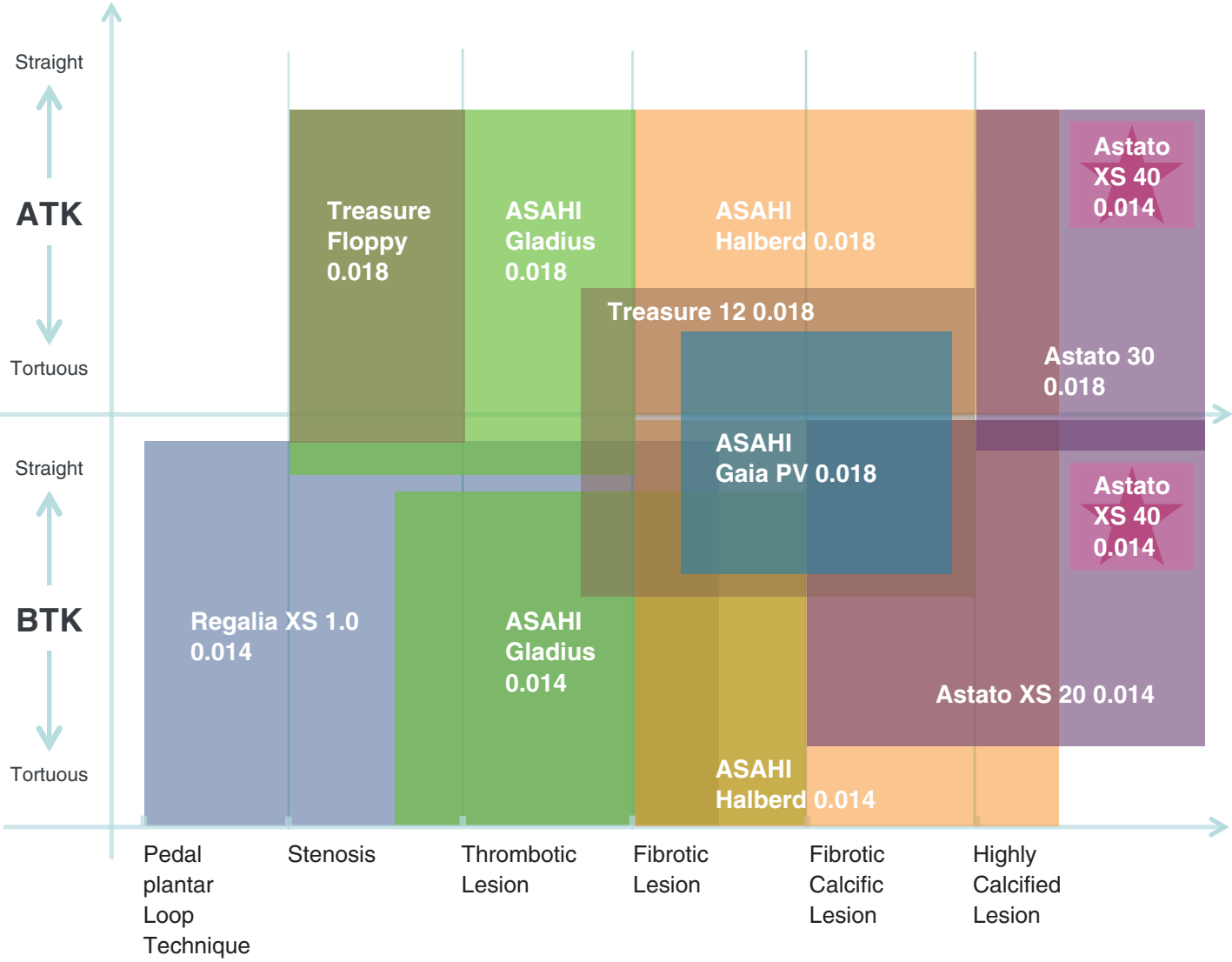

Ordering Information

... viel ilgiai - 200 cm ir 300 cm...

Product	Diameter	Catalog No.	Length	SION TECC	Tip Shape
ASAHI Gladius	0.36mm (0.014inch)	PPW14R100S	200cm	SION TECC	Straight
		PPW14R200S	235cm		
		PPW14R300S	300cm		
	0.36mm (0.014inch)	PPW14R100P	200cm		Pre-shaped
		PPW14R200P	235cm		
		PPW14R300P	300cm		
ASAHI Halberd	0.46mm (0.018inch)	PPW18R100S	200cm	SION TECC	Straight
		PPW18R200S	235cm		
		PPW18R300S	300cm		
	0.36mm (0.014inch)	PHW14R101S	200cm		Straight
		PHW14R201S	235cm		
		PHW14R301S	300cm		
	0.36mm (0.014inch)	PHW14R101P	200cm		Pre-shaped
		PHW14R201P	235cm		
		PHW14R301P	300cm		
	0.46mm (0.018inch)	PHW18R101S	200cm		Straight
		PHW18R201S	235cm		
		PHW18R301S	300cm		
ASAHI Gaia PV	0.46mm (0.018inch)	PHW18R101P	200cm	SION TECC	Pre-shaped
		PHW18R201P	235cm		
		PHW18R301P	300cm		
	0.46mm (0.018inch)	PHW18R102S	200cm		Straight
		PHW18R202S	235cm		
		PHW18R302S	300cm		
Regalia XS 1.0	0.36mm (0.014inch)	PHW18R102P	200cm	SION TECC	Pre-shaped
		PHW18R202P	235cm		
Astato XS 20	0.36mm (0.014inch)	PHW18R302P	300cm	SION TECC	Pre-shaped
		PAGP140000	180cm		
Astato XS 40	0.36mm (0.014inch)	PAGP140300	300cm	SION TECC	Pre-shaped
		PAGH143092	180cm		
Astato XS 30	0.46mm (0.018inch)	PAGH143392	300cm	SION TECC	Pre-shaped
		PAGHW143094	200cm		
Treasure Floppy	0.46mm (0.018inch)	PAGHW143394	300cm	SION TECC	Pre-shaped
		PAGH18M071	180cm		
Treasure 12	0.46mm (0.018inch)	PAGH18M371	300cm	SION TECC	Pre-shaped
		PAGH18M072	190cm		
Astato XS 20	0.36mm (0.014inch)	PAGH18M372	300cm	SION TECC	Pre-shaped
		PAGH18M070	180cm		
Astato XS 40	0.36mm (0.014inch)	PAGH18M370	300cm	SION TECC	Pre-shaped

Indications for Use
This Product is intended to facilitate the placement and exchange of diagnostic and therapeutic devices during intravascular procedures.
This device is intended for peripheral vascular use only.
Not all shapes and sizes are available in all markets. Please check availability with your sales representative.
Coated with SLIP-COAT® coating.



© 2016 ASAHI INTECC CO., LTD.
"ASAHI", "Astato", "ASAHI Gaia", "ASAHI Gladius" and "Treasure" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
Contact phone number: +81-561-86-9101
<http://www.asahi-intecc.com>

Distributed by:

ASAHI

Peripheral Guide Wire

New Guide Wires

ASAHI **Gladius** ASAHI **Halberd**
ASAHI **Gaia PV** ASAHI **Astato XS 40**




ASAHI's unique guide wire technology with Composite Core
SION TECC
Leading to the NEXT


ASAHI

Peripheral Guide Wire

Exceptional guide wires for precise interventions.

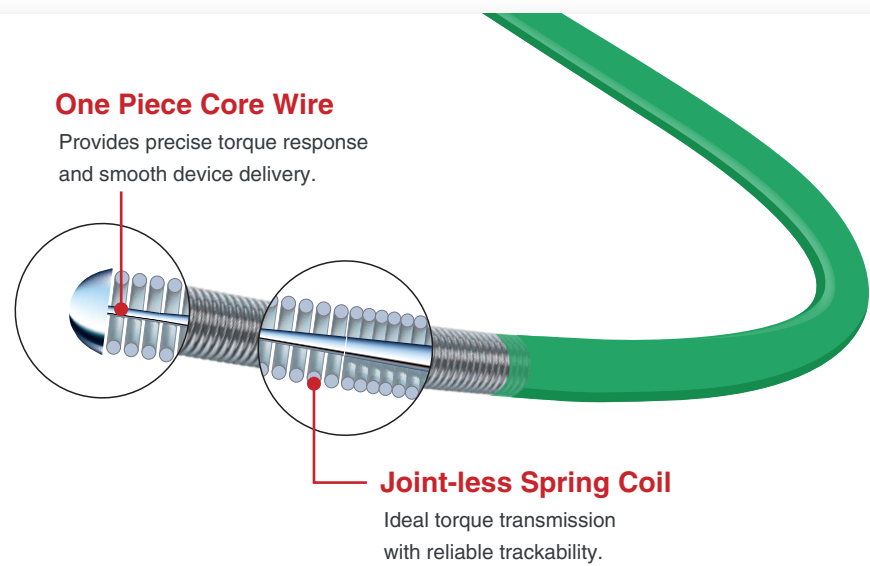

Regalia
Treasure
AstatO

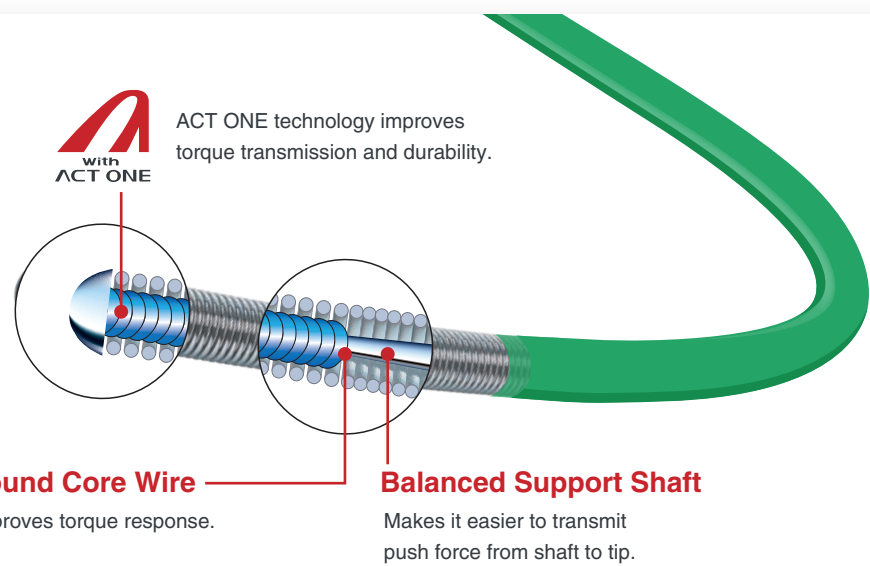

ASAHI's unique guide wire technology with Composite Core
SION TECC
Leading to the NEXT

ASAHI Gladius
ASAHI Halberd
ASAHI Gaia PV

Frontline

Regalia XS 1.0

0.36mm (0.014inch) Sub-total/Tortuous




Support ● ● ○ ○
Torque ● ○ ○ ○
Penetration ● ○ ○ ○

Balanced lubricity and trackability for BTK procedures.



Treasure Floppy

0.46mm (0.018inch) Workhorse

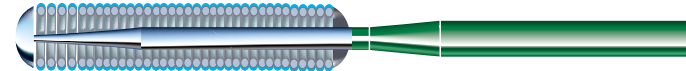

Support ● ● ○ ○
Torque ● ● ○ ○
Penetration ● ● ○ ○

Precise handling with a flexible tip for tortuous anatomy.



Specialty


ASAHI Halberd

0.36mm (0.014inch) Complex Lesion




Support ● ● ● ○
Torque ● ● ● ○
Penetration ● ● ● ●

Composite core improves torque for navigation, and the micro-cone tip improves penetration.

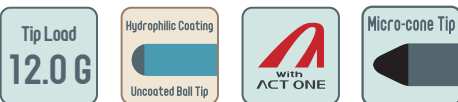

Treasure 12

0.46mm (0.018inch) Complex Lesion

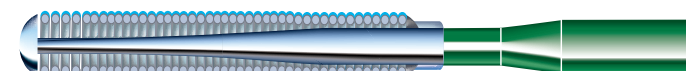


Support ● ● ○ ○
Torque ● ● ● ○
Penetration ● ● ○ ○

Consistent response and tactile feedback for complex lesion navigation.



AstatO XS 20

0.36mm (0.014inch) Highly Calcified Lesion

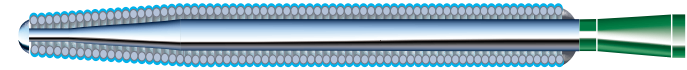


Support ● ● ○ ○
Torque ● ● ● ○
Penetration ● ● ● ●

Tapered tip provides excellent penetration.



AstatO 30

0.46mm (0.018inch) Highly Calcified Lesion

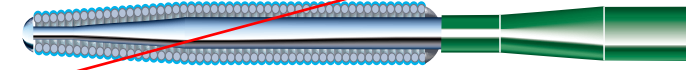

Support ● ● ● ○
Torque ● ● ● ●
Penetration ● ● ● ●

Tapered tip provides excellent penetration.



ASAHI Gladius

0.36mm (0.014inch) Workhorse



Support ● ● ○ ○
Torque ● ● ● ○
Penetration ● ● ○ ○

Composite core improves tip durability and torque to treat lesions in multiple vessels.



ASAHI Gladius

0.46mm (0.018inch) Workhorse



Support ● ● ● ●
Torque ● ● ● ○
Penetration ● ● ○ ○

Composite core improves torque response with increased support for device delivery.




ASAHI Halberd

0.46mm (0.018inch) Complex Lesion




Support ● ● ● ○
Torque ● ● ● ○
Penetration ● ● ● ○

Composite core provides precise torque for navigation with more support and micro-cone tip for proximal cap penetration.

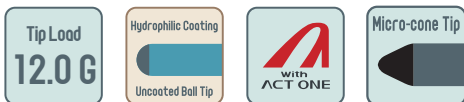

ASAHI Gaia PV

0.46mm (0.018inch) Complex Lesion

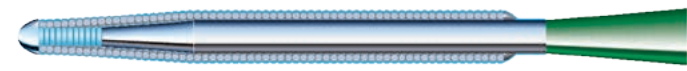

Support ● ● ● ○
Torque ● ● ● ●
Penetration ● ● ● ○

Composite core with a lower tip load improves navigation and torque. The tapered coil and micro-cone tip makes it easier to advance.



I tin ms visiškoms okliuzijoms smailiu galiuku: galiuko kietumas nuo 9,0 g iki 20,0 g; galiukas 17 - 20 cm rentgenokontrastiškas;