



Description



- Cold snare for small polyps resection,
- Thin and stiff braided wire,
- Rotable or not.

Primary packaging : 1 individual pouch
Secondary packaging : box of 10 / cardboard of 50
Minimum Delivery Quantity (MDQ) in units : 10

References	Description	Snare		Catheter		Minimum diameter of operating channel (mm)
		Width (mm)	Length (mm)	Length (cm)	Diameter (mm)	
CJ-CS-23-230-D10	Cold snare	10	20	230	2.3	≥ 2.8
CJ-CSR-23-230-D10	Rotative cold snare	10	20	230	2.3	≥ 2.8
CJ-CR-23-230-D15*	Cold snare	15	20	230	2.3	≥ 2.8
CJ-CSR-23-230-D15	Rotative cold snare	15	20	230	2.3	≥ 2.8

*Available on demand

HYBRID HOT & COLD SNARES



Braided wire



Monofilament wire

BRAIDED SNARES



Description

- Polypectomy snare for hot & cold resection*,
- Braided wire

Primary packaging : 1 individual pouch
 Secondary packaging : box of 10 / cardboard of 50
 Minimum Delivery Quantity in units : 10

References	Snare			Catheter		Minimum diameter of operating channel (mm)
	Shape	Width (mm)	Length (mm)	Length (cm)	Diameter (mm)	
CJ-ADR-23-230-HCO06 ⁽¹⁾	Oval	6	10	230	2.3	≥ 2.8
CJ-ADR-23-230-HCO10 ⁽¹⁾		10	15	230	2.3	≥ 2.8
CJ-ADR-23-230-HCO15 ⁽¹⁾		15	20	230	2.3	≥ 2.8

MONOFILAMENT SNARES



Description

- Polypectomy snare for hot & cold resection*,
- Monofilament wire

Primary packaging : 1 individual pouch
 Secondary packaging : box of 10 / cardboard of 50
 Minimum Delivery Quantity in units : 10

References	Snare			Catheter		Minimum diameter of operating channel (mm)
	Shape	Width (mm)	Length (mm)	Length (cm)	Diameter (mm)	
CJ-ADR-23-230-MFO10	Oval	6	10	230	2.3	≥ 2.8
CJ-ADR-23-230-MFO15		10	15	230	2.3	≥ 2.8

*Knowledge of the technical principles and clinical indications of polyp resections with or without electro-coagulation current is required prior to use of these products.

Prospective, randomized comparison of 2 methods of cold snare polypectomy for small colorectal polyps, Graham &co, Gastrointest Endosc 2015;82:686-92 - Gastrointest Endosc 2015; 82, No. 4 : 2015

HOT POLYPECTOMY SNARES



Désignation

- Polypectomy snare for hot resection,
- Braided wire,
- 3 available shapes : Oval, Hexagonal ou Crescent,
- Rotatable oval snare.



Primary packaging : 1 individual snare
 Secondary packaging : box of 10 / cardboard of 50
 Minimum Delivery Quantity in units : 10

References	Snare			Catheter		Minimum diameter of operating channel (mm)	Rotatable
	Shape	Width (mm)	Length (mm)	Length (cm)	Diameter (mm)		
CJ-ADR-23-230-O10	Oval	10	15	230	2.3	≥ 2.8	NO
CJ-ADR-23-230-O15		15	20	230	2.3	≥ 2.8	
CJ-ADR-23-230-O20		20	25	230	2.3	≥ 2.8	
CJ-ADR-23-230-O25		25	35	230	2.3	≥ 2.8	
CJ-ADR-23-230-O30		30	45	230	2.3	≥ 2.8	
CJ-ADR-23-230-O35		35	50	230	2.3	≥ 2.8	
CJ-ADR-23-230-O40	Oval big size	40	60	230	2.3	≥ 2.8	NO
CJ-ADR-23-230-O45		45	65	230	2.3	≥ 2.8	
CJ-ADR-23-230-O50		50	70	230	2.3	≥ 2.8	
CJ-ADR2-23-230-O10	Oval	10	15	230	2.3	≥ 2.8	YES
CJ-ADR2-23-230-O15		15	20	230	2.3	≥ 2.8	
CJ-ADR2-23-230-O20		20	25	230	2.3	≥ 2.8	
CJ-ADR2-23-230-O25		25	35	230	2.3	≥ 2.8	
CJ-ADR2-23-230-O30		30	45	230	2.3	≥ 2.8	
CJ-ADR2-23-300-O30		30	45	300	2.3	≥ 2.8	
CJ-ADR-23-230-H15	Hexagonal	15	20	230	2.3	≥ 2.8	NO
CJ-ADR-23-230-H25		25	35	230	2.3	≥ 2.8	
CJ-ADR-23-230-H30		30	45	230	2.3	≥ 2.8	
CJ-ADR-23-230-H35		35	50	230	2.3	≥ 2.8	
CJ-ADR-23-230-C15	Crescent	15	/	230	2.3	≥ 2.8	NO
CJ-ADR-23-230-C25		25	/	230	2.3	≥ 2.8	

HOT MUCOSECTOMY SNARES



SpiralSnare



Two in one

SPIRALSNARE

Description

- Polypectomy snare for hot resection of smooth lesion
- Additional spiral wire around the snare to prevent slipping on the mucosa
- Rotatable



Primary packaging : 1 individual pouch
Secondary packaging : box of 10 / cardboard of 50
Minimum Delivery Quantity in units : 10

References	Snare			Catheter		Minimum diameter of operating channel (mm)
	Shape	Width (mm)	Length (mm)	Length (cm)	Diameter (mm)	
CJ-ADR2-23-230-NG20	Oval	20	25	230	2.3	≥ 2.8
CJ-ADR2-23-230-NG25		25	35	230	2.3	≥ 2.8
CJ-ADR2-23-230-NG30		30	48	230	2.3	≥ 2.8

MONOFILAMENT SNARE

Description


- Polypectomy snare for hot resection,
 - Stiff monofilament snare,
 - Memory shape spring at the distal end.
- Two In One :
- Double size adjustable according to the polyp.



Primary packaging : 1 individual pouch
Secondary packaging : box of 10 / cardboard of 50
Minimum Delivery Quantity in units : 10

References	Snare			Catheter		Minimum diameter of operating channel (mm)
	Shape	Width (mm)	Length (mm)	Length (cm)	Diameter (mm)	
CJ-DADR-23-230-O25B	Double size	/	15 25	230	2.3	≥ 2.8
CJ-ADR-23-230-MFO20	Ovale	20	25	230	2.3	≥ 2.8

PEDIATRIC POLYPECTOMY SNARES

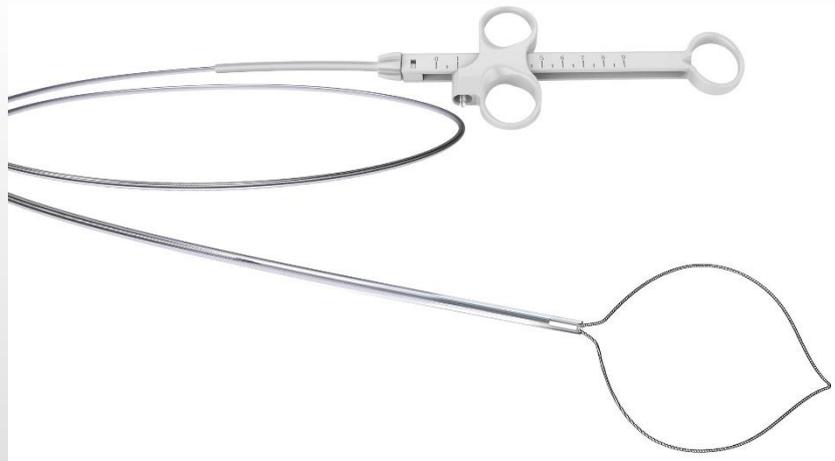

Description



- Polypectomy snare for hot resection,
- Braided wire.

Primary packaging : 1 individual pouch
Secondary packaging: box of 10 / cardboard of 50
Minimum Delivery Quantity in units : 10

References	Snare		Catheter		Minimum diameter of operating channel (mm)
	Shape	Width (mm)	Length (cm)	Diameter (mm)	
CJ-ADR-18-160-O10*	Oval	10	160	1.8	≥ 2.0 mm
CJ-ADR-18-160-O20*	Oval	20	160	1.8	≥ 2.0 mm

*On special order

Indications

All these products are used combined with endoscope for cutting polyps or other redundant tissues in digestive tract. Hybrid and hot snares require the use of high frequency current.

Technical informations

Materials :

Hot and Hybrid Hot & Cold Snare

Handle : ABS
Catheter : PTFE
Inner wire : Stainless steel
Snare wire : Stainless steel

Cold snare

Handle : ABS
Catheter : HDPE
Inner wire : Stainless steel
Snare wire : Stainless steel

Contains :

Latex : Yes No
Product of animal origin : Yes No
PVC : Yes No
Phtalate et DEHP : Yes No

Bisphenol A : Yes No
Nanoparticles : Yes No
Lead : Yes No
Cobalt : Yes No



Sterility :

Sterile : Yes No
Sterilization mode : Ethylene oxide gas
Number of use : Single use



Storage :

Storage conditions : Product should be stored in a cool, dry, clean, well ventilated, no corrosive gas environment.
Use by date : 2 years

Other data :

GMDN code : 35622 (Hot & cold polypectomy snares) ; 62615 (Cold snare)
Barcode : No GS1-128 standard HIBC standard
With the informations : Reference Lot number Use by date

Use :

Correct and safe use : Cf. Instruction for use
If yes, format Electronical Paper
User advices : Instructions for use, precautions and contraindications – cf Instruction for use



Manufacturer :

Manufacturer : Jiangsu Vedkang Medical Science & Technology Co., Ltd.
⁽¹⁾Changzhou JiuHong Medical Instrument
Country : People's Republic of China

CE MARK :

Notified body : TUV Rheinland /Certifying body n° : 0197
Regulation : European directive 93/42/CEE
Rule : Classe IIb medical device, according to the rule 9 (Hot, Hot & cold, pediatric snares)
Classe IIa medical device, according to the rule 5 (Cold snares)

IMPORTER

VYTIL

161 avenue Gallieni – 93170 Bagnole, France

Phone : 01 49 88 01 01

www.vytil.com / info@vytil.com

Order :

E-mail : sales@vytil.com

Vigilance :

E-mail : quality@vytil.com

