



6 pirkimo dalies PTVA balioniniai kateteriai vingiuotoms kraujagyslėms ir distalioms stenozėms atitikimai

CORONARY

DESIGNED FOR CHALLENGING ENVIRONMENTS

TREK™ & MINI TREK™

Coronary Dilatation Catheters

Available in RX and OTW

MINI TREK™ IS CTO INDICATED



Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

©2021 Abbott. All rights reserved. MAT-2111306 v1.0

TRANSITIONLESS DESIGN FOR CHALLENGING ANATOMY

TREK™ Coronary Dilatation Catheter's unique design provides smooth transitions from hub to tip.

- Transitionless tip
- Flexible distal shaft[^]
- Multi-layer CrossFlex^{2*} balloon and Slim Seal^{**} technology

MINI TREK™ INDICATED FOR CHRONIC TOTAL OCCLUSIONS

MINI TREK™ Coronary Dilatation Catheter's ultra low profile design enables lesion access.

- Small chassis designed distinctively for MINI TREK™ Coronary Dilatation Catheter
- Available in sizes as small as 1.20 mm




[^]Tests performed by and data on file at Abbott.

*CrossFlex² technology available on 2.25 mm – 5.0 mm balloon sizes.

**Slim Seal technology available on 2.0 mm – 3.75 mm balloon sizes.

Data on File at Abbott.

Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

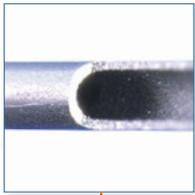

HYPOTUBE DESIGN



Designed for excellent vessel tracking and lesion crossing[^]. Skive design provides smooth transition between hypotube and distal shaft.

TRANSITIONLESS TIP
Smooth, streamlined tip for excellent push efficiency[^].



SLIM SEAL TECHNOLOGY**
Provides reduced wall thickness for flexibility and low tip crossing profiles[^].

FLEXIBLE TUNGSTEN MARKERS†
Radiopaque flexible markers for confident placement.

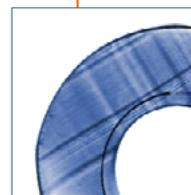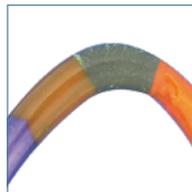

MULTI-LAYER CROSSFLEX^{2†} BALLOON
Thin, multi-layer balloon technology provides low crossing profiles[^] for access to complex lesions and a tri-fold design^{††} for excellent re-wrap.

BROAD SIZE MATRIX INCLUDING 1.20 mm DIAMETERS.

Available in RX and OTW.

[^]Tests performed by and data on file at Abbott.

^{**}Slim Seal technology available on 2.0 mm – 3.75 mm balloon sizes.

[†]Single marker on all 6 mm lengths and 1.20 – 1.50 mm diameters (all lengths). Dual markers on all other sizes.

^{*}CrossFlex2 technology available on 2.25 mm – 5.0 mm balloon sizes.

^{††}Tri-fold design for 2.00 mm - 4.00 mm diameters. Bi-fold design for 1.20 mm & 1.50 mm diameters.

Data on File at Abbott.

Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

TECHNICAL INFORMATION

| TREK™ & MINI TREK™ CORONARY DILATATION CATHETERS | |
|--|---|
| NOMINAL | 8 atm |
| RBP | 14 atm |
| TIP ENTRY PROFILE [^] | 0.017 in (3.00 mm Ø) |
| TIP CROSSING PROFILE | 0.021 in average [^] (3.00 mm Ø) 0.023 in maximum |
| FOLDED BALLOON CROSSING PROFILE | 0.032 in average [^] (3.00 mm Ø) 0.036 in maximum (3.00 mm Ø) |
| REFOLDED BALLOON CROSSING PROFILE | 0.056 in maximum (3.00 mm Ø) |
| RX NOTCH DIAMETER (FOR RX PLATFORM) | 2.5 French average [^] (3.00 mm Ø) 2.7 French maximum (all Ø) |
| TIP LENGTH | <3mm |
| BALLOON MATERIAL | Pebax [‡] , semi-compliant |
| MARKER MATERIAL | Flexible Tungsten/Pebax [‡] |
| SHAFT LENGTH | 145 cm |
| BALLOON COATING | Hydrophilic |
| MINIMUM GUIDE CATHETER | 5 French |
| MAXIMUM GUIDE WIRE COMPATIBILITY | 0.014 in |
| LATEX FREE | Yes |



[^]Tests performed by and data on file at Abbott.
Data on File at Abbott.

Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

ORDERING INFORMATION

BROAD SIZE MATRIX WITH 77 RX SIZES AND 75 OTW SIZES

| | BALLOON DIAMETER (mm) | BALLOON LENGTH (mm) | | | | | | | |
|--------------|-----------------------|---------------------|-------------|-------------|-------------|-------------|-------------|------------|---------------|
| | | 6 | 8 | 12 | 15 | 20 | 25 | 30 | |
| MINITREK™ RX | 1.20 | 1012268-06 | 1012268-08 | 1012268-12 | 1012268-15 | 1012268-20 | - | - | INDICATED CTO |
| | 1.50 | 1012269-06 | 1012269-08 | 1012269-12 | 1012269-15 | 1012269-20 | - | - | |
| | 2.00 | 1012270-06 | 1012270-08 | 1012270-12 | 1012270-15 | 1012270-20 | 1012270-25 | 1012270-30 | |
| TREK™ RX | 2.25 | 1012271-06 | 1012271-08 | 1012271-12 | 1012271-15 | 1012271-20 | 1012271-25 | 1012271-30 | |
| | 2.50 | 1012272-06 | 1012272-08 | 1012272-12 | 1012272-15 | 1012272-20 | 1012272-25 | 1012272-30 | |
| | 2.75 | 1012273-06 | 1012273-08 | 1012273-12 | 1012273-15 | 1012273-20 | 1012273-25 | 1012273-30 | |
| | 3.00 | 1012274-06 | 1012274-08 | 1012274-12 | 1012274-15 | 1012274-20 | 1012274-25 | 1012274-30 | |
| | 3.25 | 1012275-06 | 1012275-08 | 1012275-12 | 1012275-15 | 1012275-20 | 1012275-25 | 1012275-30 | |
| | 3.50 | 1012276-06 | 1012276-08 | 1012276-12 | 1012276-15 | 1012276-20 | 1012276-25 | 1012276-30 | |
| | 3.75 | 1012277-06 | 1012277-08 | 1012277-12 | 1012277-15 | 1012277-20 | 1012277-25 | 1012277-30 | |
| | 4.00 | 1012278-06 | 1012278-08 | 1012278-12 | 1012278-15 | 1012278-20 | 1012278-25 | 1012278-30 | |
| | 4.50 | - | - | 1012279-12 | 1012279-15 | - | - | - | |
| | 5.00 | - | - | 1012280-12 | 1012280-15 | - | - | - | |
| | MINITREK™ II OTW | 1.20 | 1012401-06A | 1012401-08A | 1012401-12A | 1012401-15A | 1012401-20A | - | - |
| 1.50 | | 1012402-06A | 1012402-08A | 1012402-12A | 1012402-15A | 1012402-20A | - | - | |
| 2.00 | | 1012403-06A | 1012403-08A | 1012403-12A | 1012403-15A | 1012403-20A | - | - | |
| TREK™ OTW | 2.25 | 1012404-06 | 1012404-08 | 1012404-12 | 1012404-15 | 1012404-20 | 1012404-25 | 1012404-30 | |
| | 2.50 | 1012405-06 | 1012405-08 | 1012405-12 | 1012405-15 | 1012405-20 | 1012405-25 | 1012405-30 | |
| | 2.75 | 1012406-06 | 1012406-08 | 1012406-12 | 1012406-15 | 1012406-20 | 1012406-25 | 1012406-30 | |
| | 3.00 | 1012407-06 | 1012407-08 | 1012407-12 | 1012407-15 | 1012407-20 | 1012407-25 | 1012407-30 | |
| | 3.25 | 1012408-06 | 1012408-08 | 1012408-12 | 1012408-15 | 1012408-20 | 1012408-25 | 1012408-30 | |
| | 3.50 | 1012409-06 | 1012409-08 | 1012409-12 | 1012409-15 | 1012409-20 | 1012409-25 | 1012409-30 | |
| | 3.75 | 1012410-06 | 1012410-08 | 1012410-12 | 1012410-15 | 1012410-20 | 1012410-25 | 1012410-30 | |
| | 4.00 | 1012411-06 | 1012411-08 | 1012411-12 | 1012411-15 | 1012411-20 | 1012411-25 | 1012411-30 | |
| | 4.50 | - | - | 1012412-12 | 1012412-15 | - | - | - | |
| | 5.00 | - | - | 1012413-12 | 1012413-15 | - | - | - | |

Data on File at Abbott.

Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

COMPLIANCE CHART

| INFLATION PRESSURE (atm) | BALLOON DIAMETER (mm) | | | | | | | | | | | | |
|--------------------------|-----------------------|------|------|-------|------|------|------|------|------|------|------|------|------|
| | 1.20 | 1.50 | 2.00 | 2.25 | 2.50 | 2.75 | 3.00 | 3.25 | 3.50 | 3.75 | 4.00 | 4.50 | 5.00 |
| | MINI TREK™ | | | TREK™ | | | | | | | | | |
| 2 | 1.10 | 1.34 | 1.69 | 1.87 | 2.15 | 2.16 | 2.49 | 2.91 | 2.93 | 3.06 | 3.24 | 3.69 | 4.10 |
| 3 | 1.12 | 1.36 | 1.80 | 1.95 | 2.20 | 2.35 | 2.56 | 2.98 | 3.04 | 3.22 | 3.46 | 3.91 | 4.28 |
| 4 | 1.13 | 1.38 | 1.85 | 2.00 | 2.27 | 2.41 | 2.66 | 3.05 | 3.11 | 3.31 | 3.56 | 4.04 | 4.43 |
| 5 | 1.14 | 1.40 | 1.90 | 2.06 | 2.34 | 2.49 | 2.73 | 3.12 | 3.20 | 3.41 | 3.67 | 4.17 | 4.59 |
| 6 | 1.17 | 1.43 | 1.94 | 2.12 | 2.40 | 2.56 | 2.81 | 3.19 | 3.28 | 3.50 | 3.78 | 4.29 | 4.73 |
| 7 | 1.18 | 1.45 | 1.98 | 2.17 | 2.45 | 2.63 | 2.87 | 3.25 | 3.38 | 3.59 | 3.89 | 4.40 | 4.86 |
| 8 NP | 1.20 | 1.47 | 2.01 | 2.21 | 2.51 | 2.70 | 2.93 | 3.31 | 3.45 | 3.67 | 3.99 | 4.49 | 4.98 |
| 9 | 1.21 | 1.49 | 2.03 | 2.25 | 2.55 | 2.75 | 2.97 | 3.36 | 3.51 | 3.74 | 4.06 | 4.56 | 5.05 |
| 10 | 1.23 | 1.51 | 2.07 | 2.28 | 2.58 | 2.80 | 3.01 | 3.39 | 3.56 | 3.79 | 4.11 | 4.62 | 5.13 |
| 11 | 1.24 | 1.52 | 2.09 | 2.31 | 2.60 | 2.83 | 3.04 | 3.42 | 3.60 | 3.84 | 4.18 | 4.68 | 5.20 |
| 12 | 1.25 | 1.54 | 2.12 | 2.33 | 2.62 | 2.86 | 3.07 | 3.45 | 3.64 | 3.87 | 4.23 | 4.73 | 5.25 |
| 13 | 1.26 | 1.55 | 2.14 | 2.35 | 2.64 | 2.88 | 3.09 | 3.48 | 3.67 | 3.91 | 4.28 | 4.80 | 5.31 |
| 14 RBP | 1.27 | 1.56 | 2.17 | 2.37 | 2.67 | 2.91 | 3.11 | 3.50 | 3.70 | 3.94 | 4.32 | 4.86 | 5.36 |
| 15 | 1.28 | 1.57 | 2.19 | 2.39 | 2.69 | 2.92 | 3.14 | 3.53 | 3.73 | 3.99 | 4.38 | 4.92 | 5.41 |
| 16 | 1.29 | 1.59 | 2.21 | 2.41 | 2.71 | 2.96 | 3.16 | 3.56 | 3.75 | 4.02 | 4.44 | 4.98 | 5.50 |
| 17 | 1.30 | 1.60 | 2.24 | 2.43 | 2.73 | 2.99 | 3.18 | 3.58 | 3.78 | 4.05 | 4.49 | 5.04 | 5.56 |

Data on File at Abbott.

CAUTION: This product is intended for use by or under the direction of a physician. Prior to use, reference the Instructions for Use, inside the product carton (when available) or at vascular.eifu.abbott or at medical.abbott/manuals for more detailed information on Indications, Contraindications, Warnings, Precautions and Adverse Events. This material is intended for use with healthcare professionals only.

Information contained herein for DISTRIBUTION outside the U.S. ONLY. Check the regulatory status of the device in areas where CE marking is not the regulation in force.

Illustrations are artist's representations only and should not be considered as engineering drawings or photographs. Photos on file at Abbott.

Abbott International BVBA

Park Lane, Culliganlaan 2B, 1831 Diegem, Belgium, Tel: 32.2.714.14.11

™ Indicates a trademark of the Abbott Group of Companies.

‡ Indicates a third-party trademark, which is property of its respective owner.

www.cardiovascular.abbott

©2021 Abbott. All rights reserved. MAT-2111306 v1.0

