

Single Use Electrosurgical Snares SD-400

SnareMaster Plus

The Ultimate in Flexibility



SnareMaster Plus

Specific Cold Snare Polypectomy Design

Hexagonal Shape

- Allows for precise manipulation and control when capturing polyps, whilst maintaining an appropriate margin of normal tissue required for cold snare polypectomy



Thin Wire

- 0.3 mm thin wire with adequate stiffness facilitates clean, sharp, and cold cutting

Flexibility When You Need It


Cold and Hot Handle

- Equipped with an additional connector for hot polypectomy if required



Smoother Operation

- Special coating of the wire facilitates smooth operation
- Less force is required to cut the lesion sharply and quickly

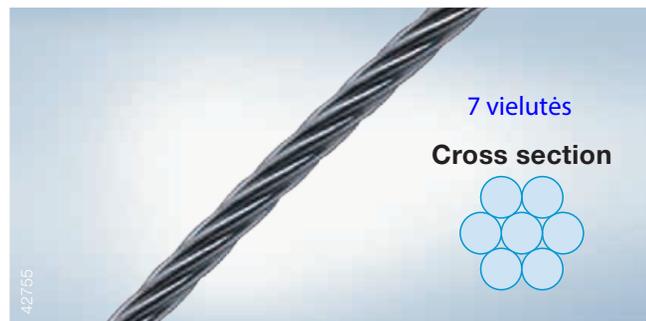


Durability

- Loop retains shape even after multiple cuts

Various sizes (10 mm, 15 mm)

- Various sizes enable cold and hot polypectomy of small lesions

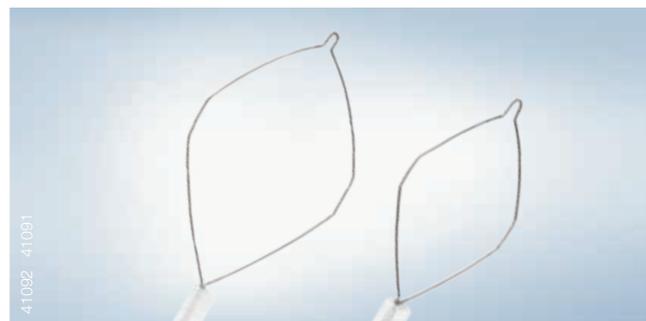

Specifications

Model	Loop design	Working length	Loop diameter	Min. Working Channel ø	Diameter of snare loop wire
SD-400U-10	Hexagon	2300 mm	10 mm	2.8 mm	0.3 mm
SD-400U-15	Hexagon	2300 mm	15 mm	2.8 mm	0.3 mm



SD-400U-10/15



E0428370EN · 11/19 · ABC

Olympus reserves the right of errors, modification and changes of the service and/or product offerings.



OLYMPUS EUROPA SE & CO. KG
 Postbox 10 49 08, 20034 Hamburg, Germany
 Wendenstrasse 14-18, 20097 Hamburg, Germany
 Phone: +49 40 23773-0, Fax: +49 40 233765
 www.olympus-europa.com