

ET CATHETERS

EC-PRO EMBRYO TRANSFER CATHETERS

Designed following the advice of some of the most reputable gynaecologists in the world. Complying with the high quality standards of Kitazato, EC-PRO catheters offer gynaecologists their improved visibility and superior control to maximise the implantation rate.

MAIN FEATURES

MALLEABLE GUIDE WITH MEMORY

Maintains its shape for a long period of time.

CATHETER WITH A SMOOTH TIP

For atraumatic insertion.

GUIDE WITH WIDENED ENTRANCE

Allows the catheter to be inserted easily.

SUPPORTED CATHETERS AVAILABLE WITHOUT THE GUIDE

Ideal for Double-Step transfer.

ERGONOMIC HANDLE

Allows safe and precise handling during the whole process.

GUIDE WITH LASER MARKINGS

Prevents toxicity after the use of ink.

APPROXIMATION MARKS ON THE CATHETER

For monitoring the insertion in the uterine cavity and facilitate the formation of the phases.

STYLETS AVAILABLE FOR THE WHOLE RANGE

Adds rigidity to the guide in the most difficult transfers.

EC-PRO ROUTINE

Shape the guide to comfortably and quickly access the uterine cavity and benefit from the precision provided by its ergonomic handle.

Supported catheter available, with and without the guide.






EC-PRO MASTER

35,36 p.d.

Ensures the correct placement of the embryos in the uterine cavity in the transfer and maximises the possibilities of implantation thanks to the highest visibility under ultrasound of the whole range.

The original design of this catheter adds the greatest visibility, from the handle to the tip, to the features of control, malleability and smoothness of the EC-PRO line.

Also available in a Supported version, guaranteeing easy insertion in the guide thanks to the cannula on the interior of the catheter.

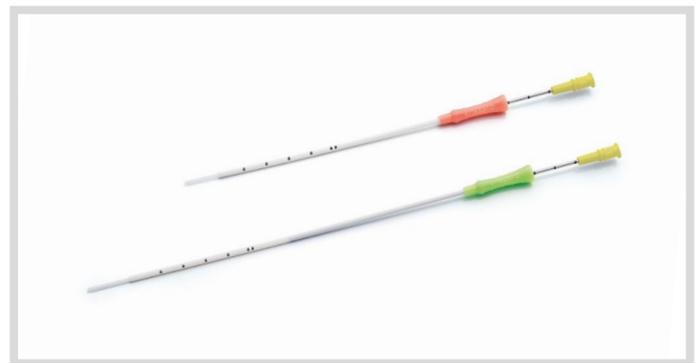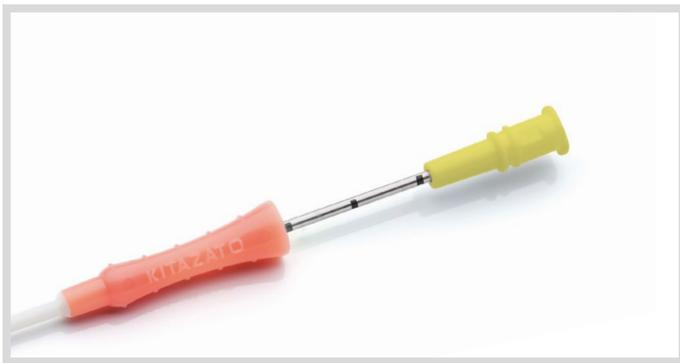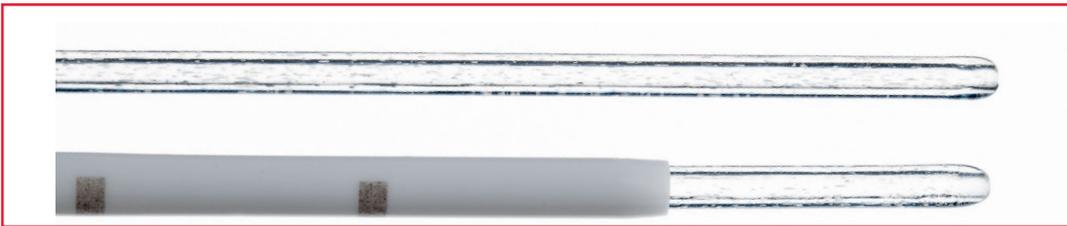


EC-PRO TRIAL

Closed tip catheter designed for examining the cervix before the transfer, and determining its difficulty. Avoids the entry of cervical mucus in the guide during channelling and guarantees an easy and clean insertion of the transfer catheter.



ET CATHETERS

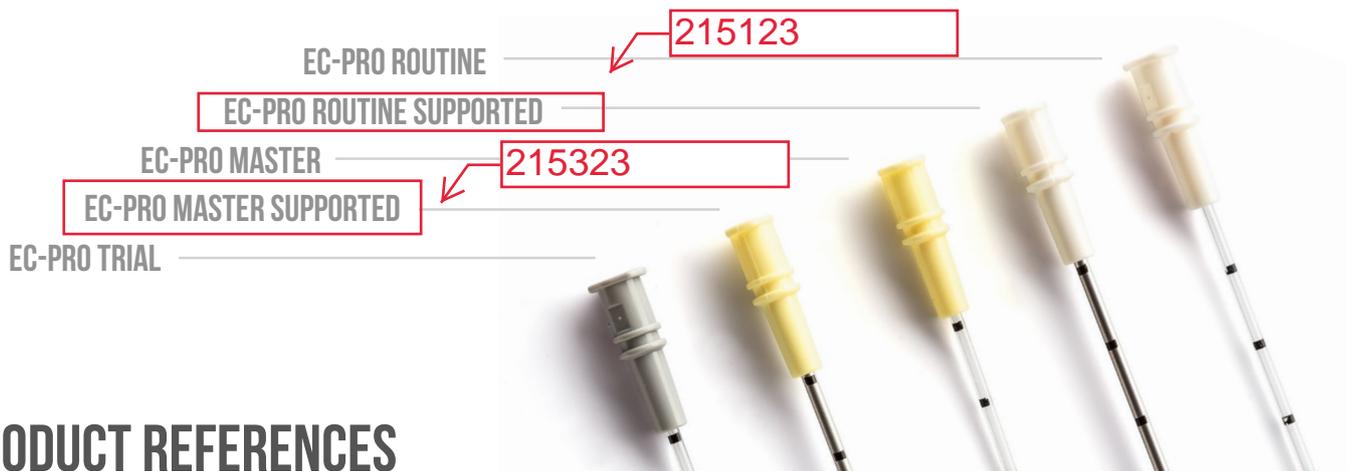

PRODUCT REFERENCES

Order number	Catheter Description	Length	Supported	High Visibility
EC-PRO				
215018	Routine	18cm	-	-
215023	Routine	23cm	-	-
215218	Master	18cm	-	x
215223	Master	23cm	-	x
215118	Routine Supported	18cm	x	-
215123	Routine Supported	23cm	x	-
215318	Master Supported	18cm	x	x
215323	Master Supported	23cm	x	x
205118	Routine Supported (catheter only)	18cm	x	-
205123	Routine Supported (catheter only)	23cm	x	-
205318	Master Supported (catheter only)	18cm	x	x
205323	Master Supported (catheter only)	23cm	x	x
275018	Trial	18cm	-	-
275023	Trial	23cm	-	-
294015	Stylet	15cm	-	-
294020	Stylet	20cm	-	-

Stiletas
naudojamas
kartu su 215323

*Kitazato Precurved ET Catheters are sold in boxes of 10 units.

**EC-PRO Diameter: 4,7Fr

QUALITY CONTROL

- MEA (Mouse Embryo Assay):
One cell assay \geq 80% after 96 hours
- Endotoxin: < 20EU/device
- Sterility (Bacteria, Fungi)
- Intracutaneous Reactivity Test
- Sensitization Test

COMPOSITION

(EC-PRO Routine ET Catheters)

Catheter: Polyurethane

Guide: Fluorocarbon Resin

Sterilization: Ethylene Oxide Gas Irradiation