

Celsite[®], Surecan[®], Cytocan

Access Port Systems, PICCs, Accessories
and Non-Coring Port Needles




Vascular Systems

Celsite[®], Surecan[®], Cytocan Access Port Systems, PICCs, Accessories and Non-Coring Port Needles

Product group	Product	Page
Access Port Systems	Implantation Sites	3
Access Port Systems for venous access	Celsite[®] PSU with titanium chamber and polysulphone body	4
	Celsite[®] Concept with titanium chamber, polysulphone body and 3 silicone suture areas for port fixation	5
	Celsite[®] Epoxy with titanium chamber and epoxy body	6
	Celsite[®] Discreet with titanium chamber and epoxy body	7
	Celsite[®] IMPLANTOFIX made of polysulphone	8
	Celsite[®] ECG for ECG assisted placement of the catheter	9
	Celsite[®] Valved Catheter	10
	Celsite[®] Double Port with titanium chamber and epoxy body	11
Access Port Systems for arterial access	Celsite[®] Arterial with polysulphone body	12
Access Port Systems for arterial access	Celsite[®] Anthron[®] Arterial	13
Access Port Systems for peritoneal access	Celsite[®] Peritoneal with titanium chamber and epoxy body	14
Access Port Systems for pleural access	Celsite[®] DRAINAPORT with titanium chamber and epoxy body	15
Access Port Systems for spinal or epidural access	Celsite[®] Spinal made of polysulphone with titanium filter	16
Celsite [®] , Surecan [®] , Cytocan	Characteristics, MRI, CECT	17
High Pressure Injections (CT)	Recommended flow rates for HP injections	18
Access Port Systems	Overview and type declaration	19
Celsite [®]	Accessories and separate kits	20-21
Celsite [®] PICC-Cel	Peripherally Inserted Central Catheters	22-23
Non-Coring port needles	Cytocan Surecan[®]	24-25
Non-Coring safety port needles	Surecan[®] Safety II	26-27

Access Port Systems

Implantation sites

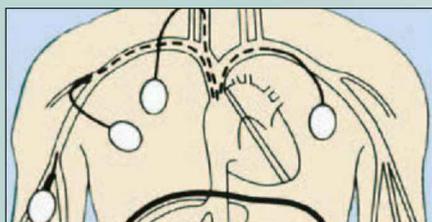

- ① **Venous access** for repeated intra-venous administration of, for example, chemotherapy, antibiotics and anti-viral drugs, total parenteral nutrition (TPN), blood sampling or transfusion




- ② **Arterial access** for intra-arterial administration of chemotherapy



- ③ **Epidural or intra-thecal access** for spinal administration of pain relieving drugs




- ④ **Peritoneal access** for loco-regional chemotherapy and (i. e. with Drainaport®) for hydration and drainage of malignant ascites

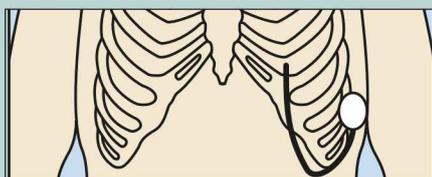

- ⑤ **Pleural access** for drainage of malignant pleural effusion (MPE)

Venous Access Port System

Celsite® PSU

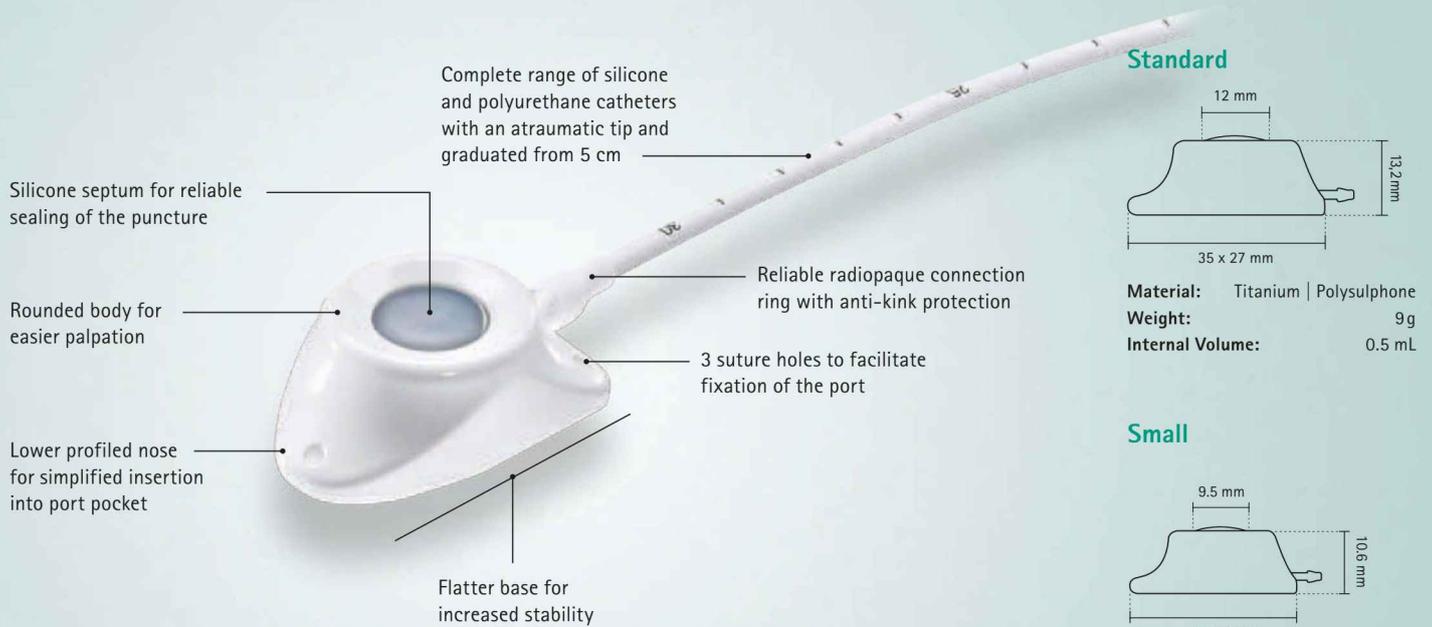

Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)**									Implantation technique	Type	Reference	Accessories see page 18/19
						Viscosity 5.8 mPa.s (cP)			Viscosity 11.4 mPa.s (cP)									
						19G	22G	22G	20G	19G	22G	20G	19G					
Standard																		
PUR	5 /1.7	1.1	900	22	10	2	5	6	1	3	5	Braunule, Seldinger	ST301C	4432096	②			
PUR	5 /1.7	1.1	370	22	10	2	5	6	1	3	5	OTW	ST3010TW	4433726	③			
PUR	6.5/2.1	1.4	800	28	11	2	5	7	1	4	6	Seldinger	ST301P	4430441	①			
PUR	6.5/2.1	1.4	800	28	11	2	5	7	1	4	6	Surgical cut-down	T301P	4430387	⑥			
Silicone***	6.5/2.2	1.0	800	24	10	2	6	7	1	4	5	Seldinger	ST301F	4430433	①			
Silicon	6.5/2.2	1.0	800	24	10	2	6	7	1	4	5	Surgical cut-down	T301F	4430000	⑥			
Silicone***	8.5/2.8	1.1	800	24	11	2	6	7	1	4	6	Seldinger	ST301	4430425	①			
Silicone	8.5/2.8	1.1	800	24	11	2	6	7	1	4	6	Surgical cut-down	T301	4430018	⑥			
PUR (high flow)***	8.5/2.8	1.6	800	39	12	2	6	8	1	4	7	Seldinger	ST301H	4432460	①			
PUR (high flow)	8.5/2.8	1.6	800	39	12	2	6	8	1	4	7	Surgical cut-down	T301H	4432452	⑥			
Silicone (high flow)	10 /3.2	1.6	800	38	12	2	6	9	1	4	6	Seldinger	ST301G	4433823	①			
Small																		
PUR	5 /1.7	1.1	900	19	10	2	5	7	1	3	5	Braunule, Seldinger	ST305C	4436962	②			
PUR	6.5/2.1	1.4	800	28	11	2	5	8	1	4	5	Seldinger	ST305P	4436946	①			
Silicone***	6.5/2.2	1.0	800	24	10	2	5	8	1	4	6	Seldinger	ST305	4433750	①			
Silicone	6.5/2.2	1.0	800	24	10	2	5	8	1	4	6	Surgical cut-down	T305	4436903	⑥			
Silicone	8.5/2.8	1.1	800	24	11	2	5	8	1	3	6	Seldinger	ST305L	4436920	①			
PUR (high flow)	8.5/2.8	1.6	800	39	12	2	6	9	1	4	6	Seldinger	ST305H	4433556	①			

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.


*** This version is also available as pre-connected system.

Venous Access Port System

Celsite® Concept




Standard

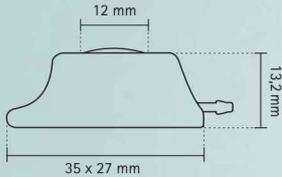

Material: Titanium | Polysulphone
Silicone
Weight: 8.6 g
Internal Volume: 0.5 mL

Small

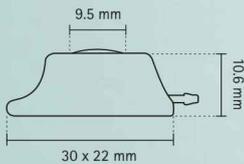

Material: Titanium | Polysulphone
Silicone
Weight: 4.6 g
Internal Volume: 0.25 mL




Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)**						Implantation technique	Type	Reference	Accessories see page 18/19	
						Viscosity 5.8 mPa.s (cP)			Viscosity 11.4 mPa.s (cP)							
						19G	22G	22G	20G	19G	22G					20G
Standard																
Silicone	6.5/2.2	1.0	800	24	11	2	6	7	1	4	5	Seldinger	ST501F	4437024	①	
Silicone	6.5/2.2	1.0	800	24	11	2	6	7	1	4	5	Surgical cut-down	T501F	4437021	⑥	
Silicone	8.5/2.8	1.1	800	24	11	2	6	7	1	4	6	Seldinger	ST501	4437022	①	
Silicone	8.5/2.8	1.1	800	24	11	2	6	7	1	4	6	Surgical cut-down	T501	4437020	⑥	
Small																
Silicone	6.5/2.2	1.0	800	24	11	2	5	8	1	4	6	Seldinger	ST505	4437027	①	
Silicone	8.5/2.8	1.1	800	24	11	2	5	8	1	3	6	Seldinger	ST505L	4437029	①	
PUR (high flow)	8.5/2.8	1.6	800	39	12	2	6	9	1	4	6	Seldinger	ST505H	4437028	①	

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.

Venous Access Port System

Celsite® Epoxy

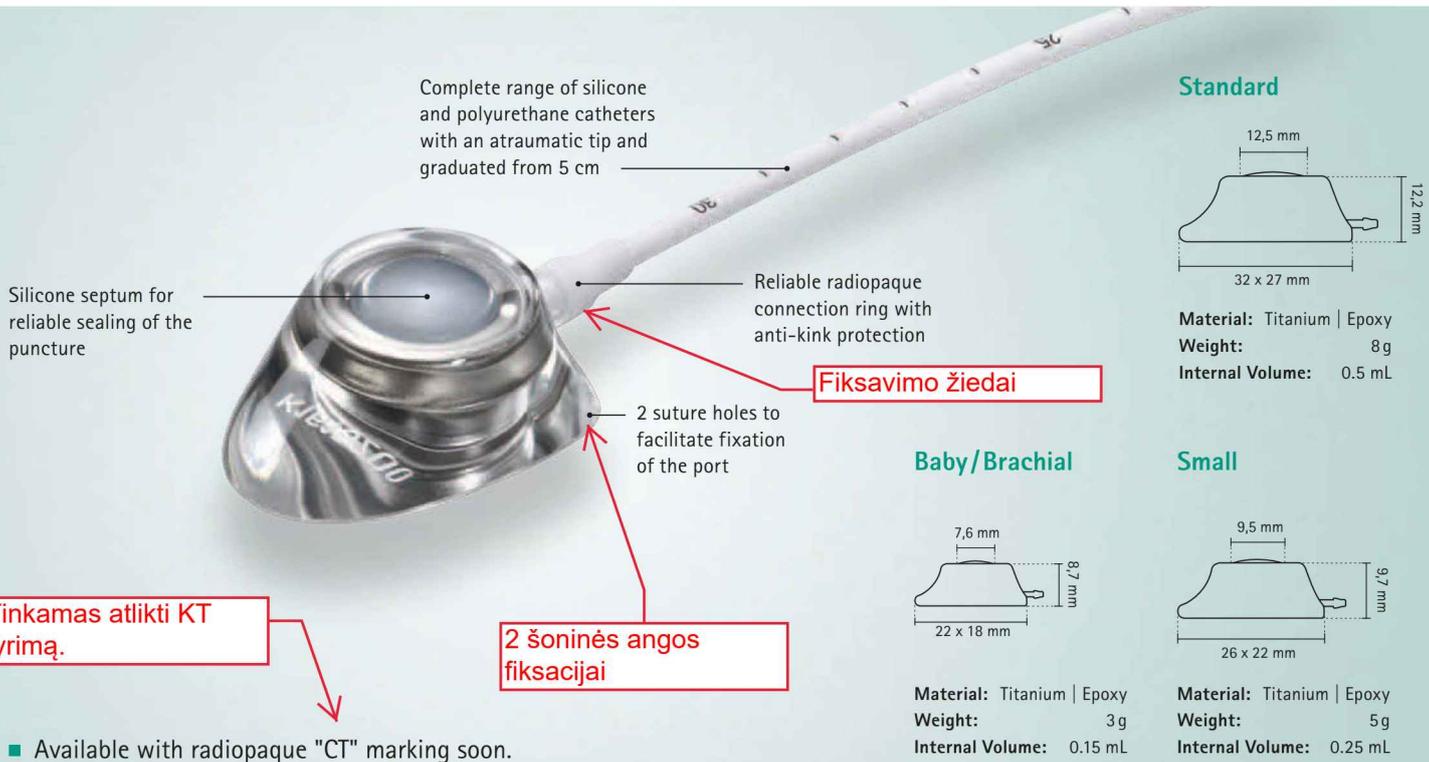

Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)	325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)**						Implantation technique	Type	Reference	Accessories see page 18/19
					Viscosity 5.8 mPa.s (cP)		Viscosity 11.4 mPa.s (cP)							
					19G	22G	22G	20G	19G	22G				
Standard														
PUR	5 /1.7	1.1	900	22 10	2	5	6	1	3	5	Braunule, Seldinger	ST201C	4432045	②
Silicone	6.5 /2.2	1.0	800	24 10	2	6	7	1	4	5	Surgical cut-down	T201F	4430034	⑥
Silicone	6.5 /2.2	1.0	800	24 10	2	6	7	1	4	5	Seldinger	ST201F	4430409	①
PUR	6.5 /2.1	1.4	800	28 11	2	5	7	1	4	6	Seldinger	ST201P	4430417	①
PUR (high flow)	8.5 /2.8	1.6	800	39 12	2	6	8	1	4	7	Seldinger	ST201H	4433149	①
Silicone	8.5 /2.8	1.1	800	24 11	2	6	7	1	4	6	Surgical cut-down	T201	4430026	⑥
Silicone	8.5 /2.8	1.1	800	24 11	2	6	7	1	4	6	Seldinger	ST201	4430395	①
Silicone (high flow)	10 /3.2	1.6	800	38 12	2	6	9	1	4	6	Seldinger	ST201G	4433807	①
Small														
Silicone***	6.5 /2.2	1.0	800	24 10	2	5	8	1	4	6	Seldinger	ST205	4430893	①
Silicone	6.5 /2.2	1.0	800	24 10	2	5	8	1	4	6	Surgical cut-down	T205	4430085	⑥
Baby/Brachial														
PUR	4.5 /1.5	0.8	800	14 8	2	4	-	1	3	-	Seldinger	Babyport®	4433742	④
PUR	5 /1.7	1.1	700	22 10	2	5	-	1	4	-	Seldinger, OTW	Brachial	4433734	⑩
Silicone	6 /2.0	1.2	600	27 12	2	5	-	1	4	-	Seldinger	Babyport® S	4433842	⑤

6

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.


*** This version is also available as pre-connected system.

Low Profile Venous Access Port System

Celsite® Discreet




Standard

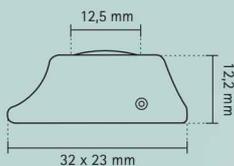

Material: Titanium | Epoxy
Weight: 7g
Internal Volume: 0.5 mL

Small

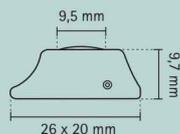

Material: Titanium | Epoxy
Weight: 4g
Internal Volume: 0.25 mL

- The unique low profile design with 90° connection provides a high level of discretion
- The surgical incision can be placed vertically to follow the subcutaneous traction lines
- Better cosmetic results due to a discreet scar
- 90° connection minimizes the risk of port rotation and therefore reduces the risk of catheter kinking
- Available with radiopaque "CT" marking soon.



Catheter	Exit can-nula	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)	325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)**						Implantation technique	Type	Reference	Accessories see page 18/19	
						Viscosity 5.8 mPa.s (cP)			Viscosity 11.4 mPa.s (cP)							
						19G	22G	22G	20G	19G	22G					20G
Standard																
Silicone	left	8.5/2.8	1.1	800	26	11	2	6	7	1	4	6	Seldinger	STL201L	4430144	⑦
Silicone	right	8.5/2.8	1.1	800	26	11	2	6	7	1	4	6	Seldinger	STR201L	4430145	⑦
PUR	left	8.5/2.8	1.6	800	39	12	2	6	8	1	4	7	Seldinger	STL201H	4440201	⑦
PUR	right	8.5/2.8	1.6	800	39	12	2	6	8	1	4	7	Seldinger	STR201H	4440202	⑦
Small																
Silicone	left	6.5/2.2	1.1	800	22	11	2	5	8	1	4	6	Seldinger	STL205F	4430146	⑦
Silicone	right	6.5/2.2	1.1	800	22	11	2	5	8	1	4	6	Seldinger	STR205F	4430147	⑦
PUR	left	6.5/2.1	1.4	800	28	11	2	5	8	1	4	5	Seldinger	STL205P	4440203	⑦
PUR	right	6.5/2.1	1.4	800	28	11	2	5	8	1	4	5	Seldinger	STR205P	4440204	⑦

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.



** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.

*** The kits for Brachial L and Brachial R contain Winged Surecan® in G22 x 15 mm (instead of G20 x 20mm).


Venous Access Port System

Celsite® IMPLANTOFIX

For repeated intra-venous administration of, for example, chemotherapy, antibiotics and anti-viral drugs, total parenteral nutrition (TPN), blood sampling or transfusion.

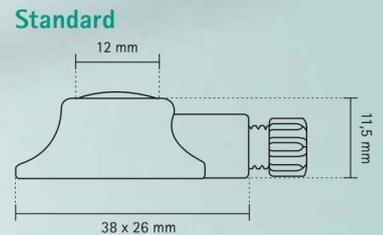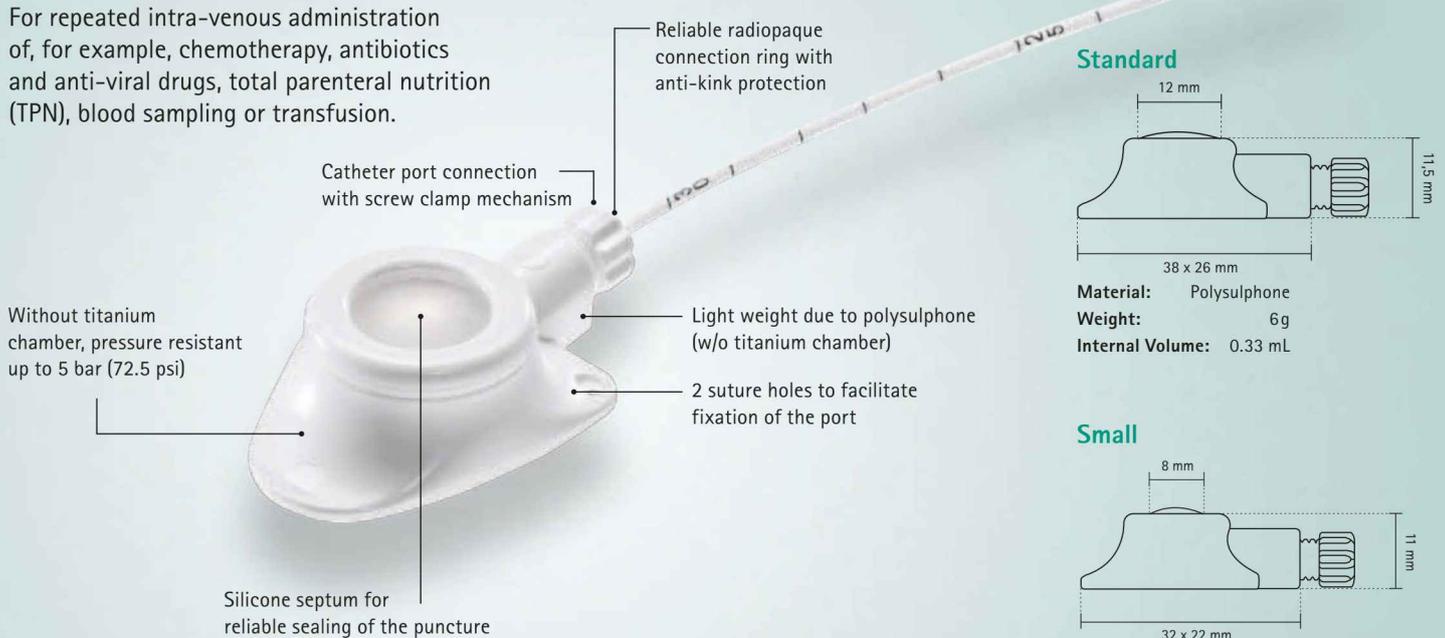

Material: Polysulphone
Weight: 6 g
Internal Volume: 0.33 mL

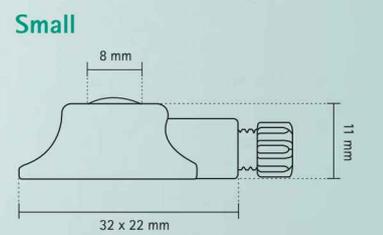

Material: Polysulphone
Weight: 4.3 g
Internal Volume: 0.15 mL

Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference	Accessories see page 18/19
				19G	22G				
Standard									
PUR	5 / 1.7	1.1	700	22	10	Surgical cut-down	IMPLANTOFIX	4430263	⑥
PUR	5 / 1.7	1.1	370	22	10	Seldinger, OTW	IMPLANTOFIX	4438604	⑬
PUR	5 / 1.7	1.1	700	22	10	Braunule	IMPLANTOFIX	4438620	⑪
Silicone	6 / 2.0	1.2	600	23	11	Seldinger	IMPLANTOFIX S	4438704	⑫
Small									
PUR	5 / 1.7	1.1	370	22	10	Seldinger, OTW	IMPLANTOFIX	4438647	⑬
PUR	5 / 1.7	1.1	700	22	10	Surgical cut-down	IMPLANTOFIX	4433521	⑥
PUR	5 / 1.7	1.1	700	22	10	Braunule	IMPLANTOFIX	4438663	⑪
Silicone	6 / 2.0	1.2	600	23	11	Seldinger	IMPLANTOFIX S	4438747	⑫

All IMPLANTOFIX products contain:

- 2 x Screw connectors
- 2 x Straight Surecan® 22G x 30 mm
- 1 x Spanner
- 1 x Vein lifter


Venous Access Port System

Celsite® ECG

Atlaiko iki 325 psi slėgi




Standard

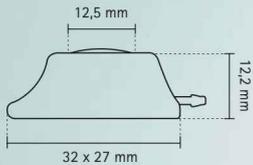

Material: Titanium | Epoxy
Weight: 8 g
Internal Volume: 0.5 mL

Small

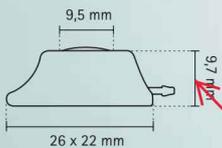

Material: Titanium | Epoxy
Weight: 5 g
Internal Volume: 0.25 mL

- The catheter can be placed via intra-atrial ECG detection. This allows accurate placement of the catheter tip into the superior vena cava without intraoperative fluoroscopy.
- Position of the catheter is defined by monitoring changes in the P-wave (see image below)
- No intraoperative x-ray control needed, which leads to time and cost savings
- No radiation exposure during intervention for theatre staff
- Specially marked guide-wire, catheter and single use ECG cable are included
- Universal adapter is required for connection with ECG monitor (B. Braun Certodyn Universal Adapter, Reference 4150228)
- Available with radiopaque "CT" marking soon.




aukštis 9,7mm, svoris 5g.



Localisation

- ① Maximal P-wave height is reached and maintained when the catheter enters into the right atrium. The amplitude of the P-wave will decrease if the catheter is advanced down into the inferior vena cava. After identifying the area where the P-wave begins to develop its maximal amplitude (which corresponds anatomically to the junction between superior vena cava and the right atrium) advance the catheter a further 2 cm.
- ② This is the final position of the catheter tip.

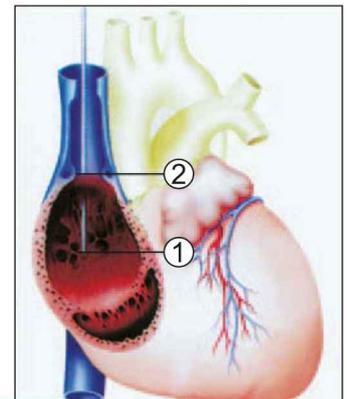

Skersmuo 6,5-8,5

ilgis 500 mm

Su titanu ir silikonine membrana

Rinkinio sudėtis yra 21 psl.

Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)	325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)**						Implantation technique	Type	Reference	series see page 18/19	
					Viscosity 5.8 mPa.s (cP)		Viscosity 11.4 mPa.s (cP)								
					19G	22G	22G	20G	19G	22G					20G
Standard															
Silicone	6.5 / 2.2	1.0	500	24	10	2	6	7	1	4	5	Seldinger (ECG)	ST201F ECG	4440140	⑨
Silicone	8.5 / 2.8	1.1	500	24	11	2	6	7	1	4	6	Seldinger (ECG)	ST201 ECG	4430140	⑨
Small															
Silicone	6.5 / 2.2	1.0	500	24	10	2	5	8	1	4	6	Seldinger (ECG)	ST205F ECG	4440111	⑨
Silicone	8.5 / 2.8	1.1	500	24	11	2	5	8	1	3	6	Seldinger (ECG)	ST205 ECG	4430111	⑨

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.

Venous Access Port Systems with valved catheter

Celsite® Valved Catheter



- For chemotherapy, administration of antibiotics, parenteral nutrition and blood sampling
- Allows easy infusion and aspiration
- Distal 3-way valve

Conditions:



Infusion




Aspiration



Closed


- Minimizes risk of catheter blockage due to thrombosis
- Anti-Reflex radiopaque silicone catheter

Standard

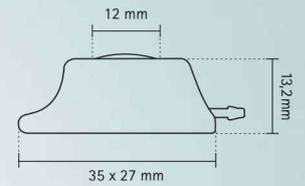

Material: Titanium | Polysulphone
Weight: 9 g
Internal Volume: 0.5 mL

Small

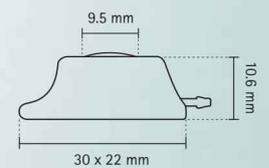

Material: Titanium | Polysulphone
Weight: 4.7 g
Internal Volume: 0.25 mL


Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference	Accessories see page 18/19
				19G	22G				
Standard									
Silicone	7.5 / 2.5	1.5	800	20	9	Seldinger	ST301V	4430092	⑦
Small									
Silicone	7.5 / 2.5	1.5	800	20	9	Seldinger	ST305V	4430095	⑦

Venous Access Port Systems

Celsite® Double Port




Standard

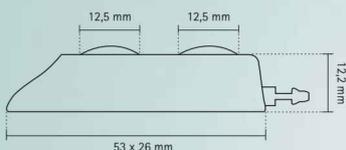

Material: Titanium | Epoxy
Weight: 14 g
Internal Volume: 0.5 mL x 2

Small

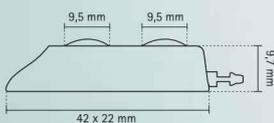

Material: Titanium | Epoxy
Weight: 7.5 g
Internal Volume: 0.25 mL x 2

- For simultaneous infusion of incompatible drugs
- For infusion with high flow rates by using both lumina
- Administration of continuous infusion and bolus injection
- Alternating puncture sites
- Profiled shape to be easily placed in a small pocket
- Small size facilitates implantation in paediatric and underweight patients
- Off-set silicone catheter tip ensures that no mixing of drugs occurs at the catheter tip
- Available with radiopaque "CT" marking soon.



Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		325 PSI Recommended maximum flow rates (mL/s) Contrast media at 37°C (325 psi = 22.4 bar)						Implantation technique	Type	Reference	Accessories see page 18/19	
						Viscosity 5.8 mPa.s (cP)			Viscosity 11.4 mPa.s (cP)							
						19G	22G	22G	20G	19G	22G					20G
Standard																
Silicone	10 / 3.2	1.2 x 2	800	24	10	2	5	8	1	4	6	Seldinger	ST401L	4430100	⑦	
Small																
Silicone	10 / 3.2	1.2 x 2	800	24	10	2	5	8	1	4	6	Seldinger	ST405L	4430101	⑦	

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

** With a catheter of 20 cm and Surecan® Safety II. For countries under CE mark only.

Venous Access Port Systems

Celsite® Arterial





- For loco-regional chemotherapy of liver tumours and hepatic artery infusion therapies
- The access port is implanted at the base of the ribs, while the catheter is introduced into the *arteria gastroduodenalis* so that the catheter tip is located in the *arteria hepatica*
- The radiopaque silicone catheter has three rings to facilitate immobilisation of the catheter in the artery


Accessories:

Every Access Port kit contains

- 2 Straight Surecan® needles 22G x 30 mm
- 1 vein lifter

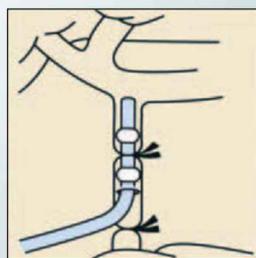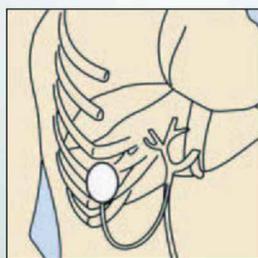


Celsite® Standard

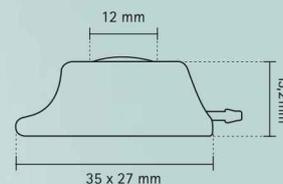

Material: Titanium | Polysulphone
Weight: 9g
Internal Volume: 0.5 mL

Celsite® IMPLANTOFIX Standard

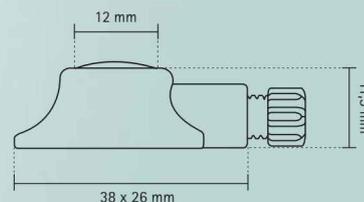


Material: Polysulphone
Weight: 6g
Internal Volume: 0.33 mL

Catheter	Access Port	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference
					19G	22G			
Standard									
Silicone	Celsite® (Titanium/ Polysulphone)	6.5 / 2.2	1.0	800	24	10	Surgical cut-down	T302	4430042
PUR	IMPLANTOFIX (Polysulphone)	5 / 1.7	1.1	700	22	10	Surgical cut-down	IMPLANTOFIX	4438817**

Arterial Access Port Systems


Celsite® Anthron® Arterial

Celsite® Small

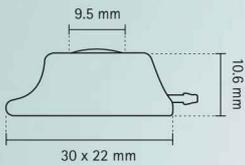

Material: Titanium | Polysulphone
Weight: 4.7 g
Internal Volume: 0.25 mL

Tapered catheter

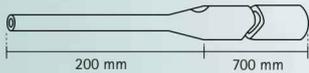


- The Anthron® catheter is specially designed for percutaneous access via the femoral or axillary artery, to the hepatic artery, for hepatic arterial infusion of chemotherapy (HAIC).
- Anthron® is a hydrophilic polyurethane catheter to which heparin is ionically bound. It is particularly indicated to help prevent catheter occlusion and catheter-related thrombosis.

Accessories:

Every Access Port kit contains

- 1 tunneling rod
- 1 vein lifter
- 2 Straight Surecan® needles 22G x 30 mm
- 1 Winged Surecan 20G x 20 mm

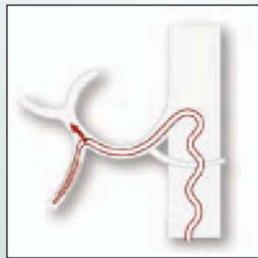

Catheter	OD (F)	OD (mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference
					19G	22G			
Small									
Tapered PUR Anthron®	Non tapered portion 5F Distal tip 2.7F	1.7/0.9	1.1/0.5	Total 900 Tapered 200	18	10	Percutaneous	R305-A5ST	4442465

* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

Peritoneal Access Port Systems

Celsite® Peritoneal



- For loco-regional chemotherapy of peritoneal metastases and ovarian cancer
- The access port is implanted at the base of the ribs and the catheter is placed at the required location inside the abdominal cavity
- The radiopaque silicone catheter with multiple perforations ensures optimal diffusion of infused drugs and reliable patency of the catheter


Accessories:

Every Access Port kit contains


- 2 Straight Surecan® needles 22G x 30 mm
- 1 vein lifter

The implantation accessories kit AP16F can be ordered separately (Reference 4430493; see page 21).

Standard

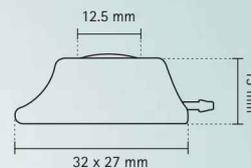

Material: Titanium | Epoxy
Weight: 10 g
Internal Volume: 0.5 mL

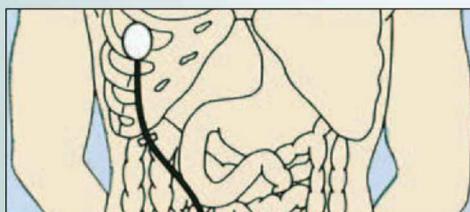


Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference
				19G	22G			
Standard								
Silicone	15 / 4.9	2.6	420	46	12	peritoneal	T203J	4430069


Peritoneal/Pleural Access Port System

Celsite® DRAINAPORT

Standard

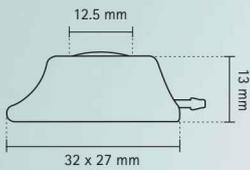

Material: Titanium | Epoxy
Weight: 10g
Internal Volume: 0.5 mL

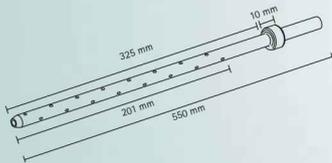

- For intra-peritoneal administration of chemotherapy, hydration, drainage of malignant ascites, or drainage of malignant pleural effusion
- Avoids repeated, painful puncture for drainage
- Improves quality of life and is an easy and effective solution for home care treatment
- Celsite® DRAINAPORT can be implanted percutaneously or by surgical cut-down technique
- **Catheter cuff** promotes tissue ingrowth to reduce infection risk and holds the catheter securely in place
- **Connection** is secured with the radiopaque titanium connection ring
- **Silicone septum** for reliable puncture and easy port location
- **Anatomic design** with delta shape profile, light weight and easy to suture
- Smooth, large and flexible **multiperforated silicone catheter** with 49 oval holes (Ø 1.1 x 1.6 mm) from the tip up to 20 cm to prevent blockage of the catheter and ensure optimal efficiency



Catheter	OD (F/mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference
				19G	22G			
Standard								
Silicone	15/4.9	2.6	550	46	12	peritoneal / pleural	T203J-1	4430169

Accessories:

Every Access Port kit contains

- 2 Straight Surecan® needles 22G x 33 mm

The implantation accessories kit AP16F can be ordered separately (Reference 4430493; see page 21).


* Gravity infusion of saline (0.9%) through a 22G respectively 19G needle from a height difference of 1 m and a catheter length of 40 cm.

Epidural/Intra-thecal Access Port System


Celsite® Spinal

- For spinal administration of pain relieving drugs
- The catheter is tunneled under the skin to the access port, which is implanted at the base of the ribs
- Light weight and comfortable
- Profiled shape design facilitates insertion
- Integrated 20 µm titanium filter prevents the passage of particles

Standard

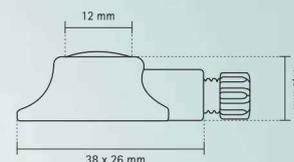


Material: Polysulphone | Titanium-Filter
Weight: 6g
Internal Volume: 0.33 mL

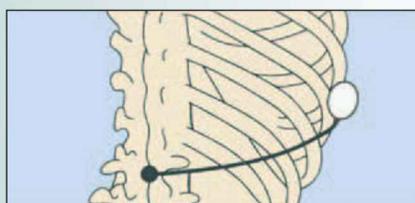

Catheter

Every Access Port kit contains 2 catheters:

- 1 multiperforated closed tip polyamide catheter (PA)
- 1 open tip polyurethane catheter (PUR) with a teflon-coated guide wire

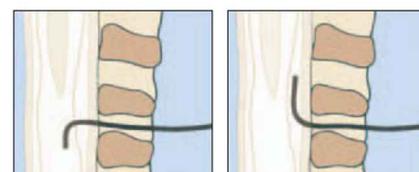

Catheter	OD	OD (mm)	ID (mm)	Length (mm)	Flow rate* (ml/min)		Implantation technique	Type	Reference
					19G	22G			
Standard									
PUR and PA	19G	1.05	0.6	1000	4	3	spinal/epidural	ST304-19	4430096
PUR and PA	20G	0.86	0.45	1000	1	1	spinal/epidural	ST304-20	4430097

Accessory kit:

- Screw connector (2x)
- Spanner
- Anti-kink device (2x)
- Tunneling rod
- Winged Surecan® needle 20G x 20 mm
- Omnifix syringe 10 mL
- Perican Tuohy needle 16G (ST304-19) or Tuohy needle 18G (ST304-20)
- Sterifix 0.2 µm filter
- Scalpel size 10 and 11
- Sterican needle 20G x 70 mm
- Perifix LOR syringe
- Straight Surecan® 22G x 30 mm (2x)

MR Compatibility and High Pressure Resistance

Celsite®



MR-Conditional

Suderintas su MRT

Non-clinical testing demonstrated that Celsite® Access Ports and Surecan®/Cytocan port needles (including Safety II and Ultrasite®) are MR Conditional. A patient with these devices can be scanned immediately after placement under the following conditions:

- Static magnetic field of 3-Tesla and 1.5-Tesla
- Maximum spatial gradient magnetic field of 710 Gauss/cm or less
- Maximum whole body averaged specific absorption rate (SAR) of 2.9 W/kg for 15 minutes of scanning

MR image quality may be compromised if the area of interest is in the exact same area or relatively close to the position of the devices. Therefore, optimization of MR imaging parameters to compensate for the presence of these devices may be necessary.

Please see instructions for use for general information and information on MRI-related heating.

Pressure Resistance



All venous Celsite® Access Ports with titanium chamber (except for valved catheters and Celsite® Implantofix) are pressure resistant up to 325 psi (22.4 bar).

Please see instructions for use for detailed device information regarding high pressure injection.

Be latekso

Material

All Celsite® Access Ports are latex-, PVC- and DEHP-free. All Surecan®/Cytocan needles are latex- and DEHP-free.

LATEX FREE

PVC FREE

DEHP FREE

Recommended maximum flow rates (mL/s)

Celsite® Access Port Systems

with Angled Surecan® needle and
Winged Surecan® needle without Y-site



Contrast media at 37°

Celsite® type		Viscosity 5.8 mPa.s (cP)			Viscosity 11.4 mPa.s (cP)		
		Needle size			Needle size		
		22G	20G	19G	22G	20G	19G
Celsite® Baby/ Brachial	Babyport® – Babyport® PC	2	4	-	1	3	-
	Brachial	2	4	-	1	3	-
	Brachial L – Brachial R	2	4	-	1	3	-
	Babyport® S	2	4	-	2	4	-
Double Port	ST405L	2	5	6	2	4	6
Celsite® Small	STL205P – STR205P	2	4	6	2	3	5
	ST305P	2	4	6	2	3	4
	ST305C	2	4	5	1	3	4
	T/ST305 – T/ST205 – ST505 – ST315 – ST215 – ST205F ECG	2	4	5	2	3	4
	STL205F – STR205F	2	4	5	2	3	4
	ST305L – ST505L – ST205ECG – ST315L	2	4	5	2	3	5
	ST305H – ST505H	2	5	7	2	4	6
Double Port	ST401L	2	5	7	2	4	6
Celsite® Standard	ST301C – ST501C – ST201C – ST3010TW	2	5	6	2	4	5
	T/ST301F – ST201F ECG – T/ST201F – T/ST501F – ST311F	2	5	6	2	4	6
	T/ST301P – ST201P	2	5	6	2	4	6
	T/ST301 – ST311 – T/ST201 – T/ST501 – ST201ECG	2	5	6	2	4	6
	STL201L – STR201L	2	5	6	2	4	6
	ST201H – T/ST301H – ST311H – ST501H	2	5	7	2	5	7
	STL201H – STR201H	2	5	7	2	5	7
	ST301G – ST201G – ST501G	2	5	8	2	5	7

Recommended maximum pressure (CT function) – 325 psi (22.4 bar)

Flow rates may vary depending on temperature of contrast media and length of the implanted catheter.

For flow rates of Surecan® Safety II, please refer to concerned Access Port product pages.

Overview and type declaration

Celsite®

Indication	Catheter	OD	Catheter material	Access Port type	Dead volume port	Dead volume catheter (mL/cm)	
VENOUS	Small catheters	5 F	Polyurethane	ST201C	0.50 mL	0.010 mL	
				ST301C, ST3010TW, ST501C			
				ST305C			
				4430263, 4438604, 4438620 (Implantofix)			
		6.5 F	Polyurethane	4438647, 4438663, 4433521 (Implantofix)		0.15 mL	0.015 mL
				ST201P, T301P, ST301P	0.50 mL		
		6 F	Silicone	ST305P, STL205P, STR205P		0.25 mL	0.011 mL
				4438704 (Implantofix)		0.33 mL	
		6.5 F	Silicone	4438747 (Implantofix)		0.15 mL	0.008 mL
				T201F, ST201F, T301F, ST301F, ST311F*, T501F, ST501F, ST201F ECG		0.50 mL	
Large and high flow catheters	8.5 F	Silicone	T205, ST205, ST215*, T305, ST305, ST315*, ST505		0.25 mL	0.010 mL	
			T201, ST201, T301, ST301, ST311*, T501, ST501, STL201L, STR201L		0.50 mL		
	8.5 F	Polyurethane	ST305L, ST505L		0.25 mL	0.020 mL	
			ST201H, T301H, ST301H, ST311H*, STL201H, STR201H		0.50 mL		
	10 F	Silicone	ST305H, ST505H		0.25 mL	0.020 mL	
			ST201G, ST301G		0.50 mL		
SPECIALITY VENOUS	Small catheters	4.5 F	Polyurethane	Babyport®	0.15 mL	0.005 mL	
		5 F	Polyurethane	Brachial	0.15 mL	0.010 mL	
		6 F	Silicone	Babyport® S	0.15 mL	0.011 mL	
		6.5 F	Silicone	STR205F, STL205F, ST205F ECG	0.25 mL	0.008 mL	
	Large and high flow catheters	8.5 F	Silicone	STR201L, STL201L, ST201 ECG		0.50 mL	0.010 mL
				ST205ECG		0.25 mL	
	Valved catheters	7.5 F	Silicone	ST301V		0.50 mL	0.018 mL
				ST305V		0.25 mL	
	Double port catheters	10 F	Silicone	ST401L		0.50 mL x 2	0.013 mL
				ST405L		0.25 mL x 2	
OTHER SPECIALITIES	Small arterial catheters	5 F	Polyurethane	4438817	0.33 mL	0.010 mL	
		6.5 F	Silicone	T302	0.50 mL	0.008 mL	
	Peritoneal catheters	15 F	Silicone	T203J, T203J-1		0.50 mL	0.053 mL
				19 G	Polyamide and polyurethane	ST304-19	
Spinal/epidural catheters	20 G	ST304-20					0.002 mL

* Pre-connected Access Port Systems

S

Accessories
S=Set

T

Reservoir
T=Titanium

R

Discreet
R=Right
L=Left

2

Body
2=Epoxy
2 suture holes
3=Polysulphone
3 suture holes
4=Double Port Epoxy
5=Polysulphone
with Silicone ears

0

Connection
0=with connection
ring
1=pre-connected

1

Indication
1=Venous (std)
2=Arterial
3=Peritoneal
& pleural
4=Spinal
5=Venous
(small)


F

Catheter
L=Silicone (large)
F=Silicone (small)
G=Silicone (high flow)
P=PUR
J=Silicone peritoneal/
pleural
H=PUR (high flow)
C=PUR Tecothane
V=Valve
OTW=Over the wire
ECG =ECG implantation
technique

Accessories

Venous accessories

Pieces	Implantation technique	Percutaneous					
		Seldinger		OTW	Seldinger/OTW	Seldinger/Braunule	Braunule
	Kit designation	Kit ①	Kit ②	Kit ③	Kit ③	Kit ②	Kit ①
2	Straight Surecan® needles	22 G x 30 mm	22 G x 30 mm				
1	Vein lifter	x	x	x	x	x	x
1	Puncture needle	18 G x 70 mm					
1	Splittocan needle					14 G x 80 mm	14 G x 80 mm
1	J guide wire with dispenser	0.035" x 50 cm	0.035" x 50 cm	0.035" x 70 cm	0.035" x 70 cm	0.035" x 50 cm	
1	Dilator			6F x 100 mm	6F x 100 mm		
1	Tear-away introducer	L 180/140 mm	L 180/140 mm			L 180/140 mm	
1	Tunnelling rod	x	x	x	x	x	x
1	Omnifix luer syringe	10 mL	10 mL				
1	Winged Surecan® needle	20 G x 20 mm		20 G x 20 mm		20 G x 20 mm	

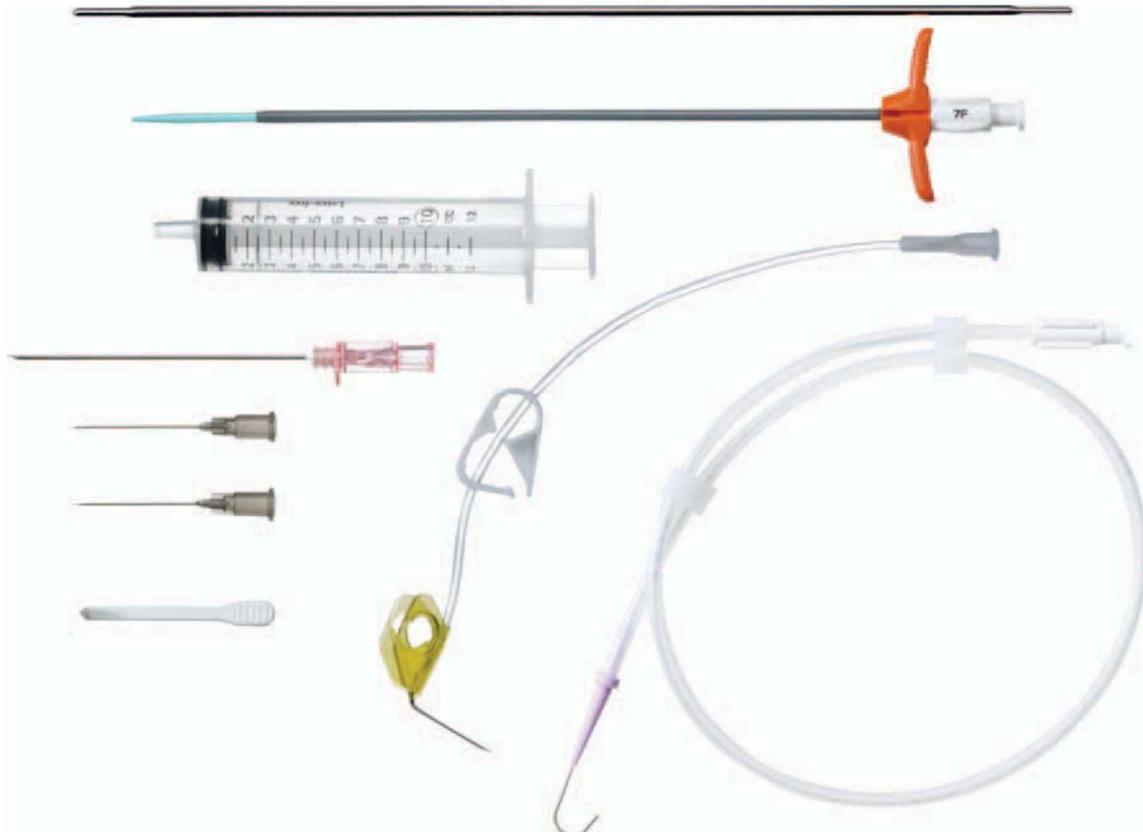

Accessories

Venous accessories

Pieces	Implantation technique	Surgical		Percutaneous			
		Cut-down	Seldinger/OTW	Seldinger			ECG
	Kit designation	Kit ⑥	Kit ⑩ (Brachial)	Kit ④ (Baby)	Kit ⑤ (Baby)	Kit ⑦	Kit ⑨
2	Straight Surecan® needles	22 G x 30 mm	22 G x 30 mm	22 G x 30 mm	22 G x 30 mm	22 G x 30 mm	22 G x 30 mm
1	Vein lifter	x	x	x	x	x	x
1	Puncture needle		18 G x 70 mm	20 G x 50 mm	18 G x 70 mm	18 G x 70 mm	18 G x 70 mm
1	Introcan needle			20 G x 32 mm			
1	J guide wire with dispenser		0.035" x 150 cm	0.025" x 50 cm	0.035" x 50 cm	0.035" x 50 cm	0.035" x 70 cm
1	ECG cable						x
1	Tear-away introducer/ Dilator		L 180/140 mm	L 80/50 mm	L 180/140 mm	180/140 mm	180/140 mm
1	Tunnelling rod		x	x	x	x	x
1	Omnifix luer syringe		10 mL	10 mL	10 mL	10 mL	10 mL
1	Winged Surecan® needle		22 G x 15 mm	22 G x 15 mm	22 G x 15 mm	20 G x 20 mm	20 G x 20 mm

Rinkinio sudėtis:

Separate accessory kits

	Reference	4430483	4430484	4430492	4430493
Pieces	Kit designation	AP 6F	AP 7F	AP 9F	AP 16F
1	Puncture needle	20 G x 50 mm Seldinger	18 G x 70 mm Seldinger	18 G x 70 mm Seldinger	18 G x 70 mm Seldinger
1	Introcan needle	20 G x 32 mm			
1	J guide wire with dispenser	0.025" x 50 cm	0.035" x 50 cm	0.035" x 50 cm	0.035" x 40 cm
1	Tear-away introducer/ Dilator	6F, short (80/50 mm)	7F x180/140 mm	9F x180/140 mm	16F with dilator 12F-14F
1	Tunnelling rod	x	x	x	x
1	Vein lifter	x	x	x	
1	Omnifix luer syringe	10 mL	10 mL	10 mL	10 mL
1	Winged Surecan® needle	22 G x 15 mm	20 G x 20 mm	20 G x 20 mm	19 G x 25 mm

Celsite® PICC-Cel


Celsite® PICC-Cel is a Peripherally Inserted Central venous polyurethane (PUR) Catheter, indicated for short to long-term drug infusion therapy (up to 3 months) such as antibiotherapy, blood sampling and transfusion, Total Parenteral Nutrition (TPN), antivirals or chemotherapy.

Celsite® PICC-Cel is also indicated for:

- Central Venous Pressure (CVP) monitoring
- High Pressure (HP) injection of contrast media (5 mL/sec, 300 psi max, 20.68 bar).

Catheter characteristics

- Single and double lumen catheters available:
Double lumen catheters allow simultaneous infusion of incompatible drugs
- PUR catheter marked every cm from 0 at the hub end helps monitoring of catheter progress and confirmation of complete catheter removal
- MR Safe (if using with Ultrasite® valve then the system is MR Conditional)
- Reverse taper design catheter from the hub for 10 cm reduces the risk of bleeding and increases resistance to kinking at skin exit
- Radiopaque

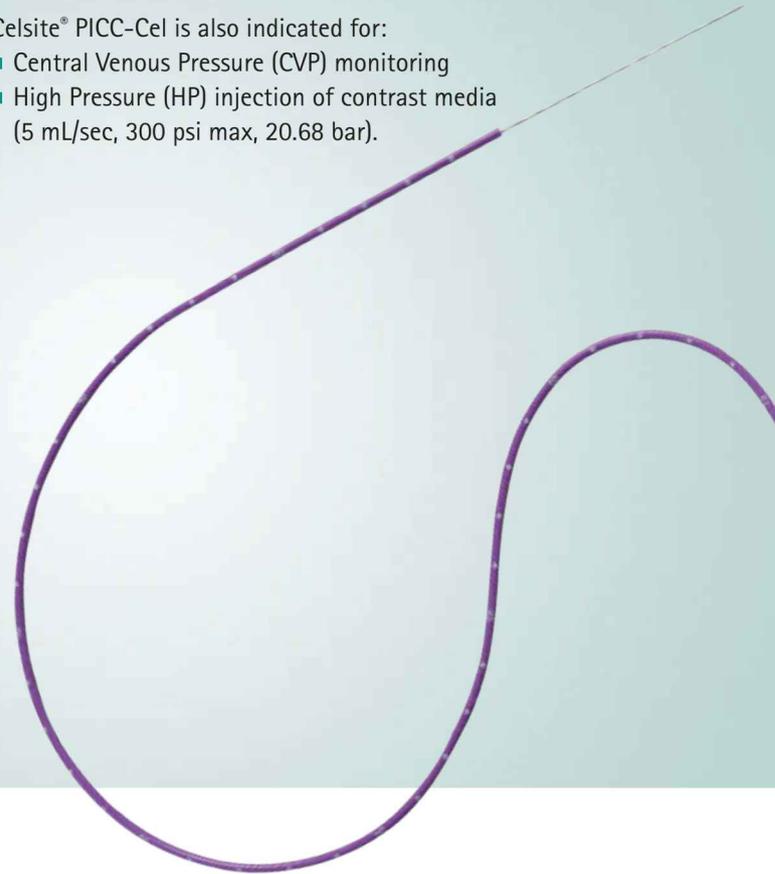

PICC-Cel range

Box of 5 units

N° lumen	Catheter/ peelable introducer	Catheter length	Maximum flow rate (300 psi/20 bar)	Reference
Single	4F	51 cm	5 mL/sec	04439001
Single	5F	61 cm	5 mL/sec	04439002
Double	5F	56 cm	5 mL/sec	04439003
Double	6F	61 cm	5 mL/sec	04439004

Implantation Accessories in the kit:



21G x 70 mm echogenic puncture needle for better visualisation under echo-guidance.



10 mL Omnifix® syringe



Scalpel (N° 11)



0.018" x 70 cm stainless steel guide-wire with 7cm floppy gold-plated tungsten tip, for introduction of the peelable introducer and to measure catheter length



Ultrasite® needle-free valve (1 for each line)



PICC-Cel Grip-Lok® Securement device



4F, 5F or 6F x 10 cm peelable introducer (color coded for easier identification)

Tape measure (66 cm/26") facilitates catheter measurement and trimming before implantation



- Luer lock hubs for reliable attachment, different color hubs allow for easier identification and management of therapies

- Pre-mounted stylet improves catheter advancement

- Side port for flushing during catheter placement


- Maximum pressure and flow rate marked on clamp

- Clear external tubing for easy confirmation of blood reflux, marked "Power Injectable" and with the gauge for CVP monitoring

- Catheter size clearly indicated on the fixation hub


PICC-Cel Grip-Lok® Securement device

Ref. 04439010 – Box of 50 units

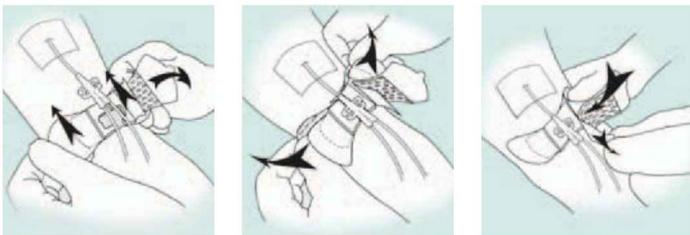

As suturing is not recommended* due to infection risk, one PICC-Cel Grip-Lok® is available in each Celsite® PICC-Cel box for safe and clean fixation of the PICC-Cel catheter at the skin exit.


PICC-Cel Grip-Lok securement devices can be ordered separately in boxes of 50 units.

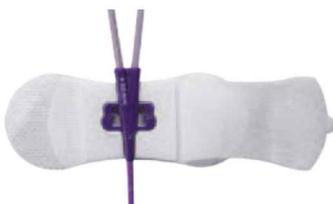

- The 2 holes on the hub are compatible with other securement devices



- Ultrasite® needle-free valve minimizes the risk of blood exposure and inadvertent air embolism

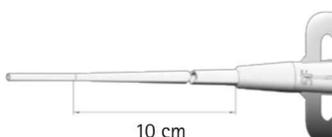

- 10 cm reverse taper catheter design reduces risk of bleeding



- Luer connector for flushing of the guide wire

* Infusion Nursing Standards of Practice for catheter stabilization J. Infus. Nursing 2006;29(IS):S1-S92

Non-Coring Needles for Access Ports

Surecan® | Cytocan

Angled Surecan® non-coring needle

- use for short-term infusions
- latex- and DEHP-free (hub)



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
19G	1.1	15	50	4438000
19G	1.1	20	50	4439430
19G	1.1	25	50	4439406
20G	0.9	15	50	4439929
20G	0.9	20	50	4439937
20G	0.9	25	50	4439945
20G	0.9	35	50	4434862
22G	0.7	15	50	4439813
22G	0.7	20	50	4439821
22G	0.7	25	50	4439830
22G	0.7	35	50	4434870

Straight Surecan® non-coring needle

- use for bolus injection or flushing of the Access Port
- latex- and DEHP-free (hub)



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
20G	0.9	40	100	4439953
20G	0.9	70	100	4439998
20G	0.9	90	100	4440000
22G	0.7	30	100	4439848
24G	0.55	25	100	4439414

Cytocan non-coring needle with fixation base

- use for long-term infusions
- flexible, transparent fixation base for reliable deployment
- latex- and DEHP-free (fixation base)
- extension tubing with clamp (PVC-free)
- tubing length cannula to connector: 250 +/- 10 mm



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
19G	1.1	15	25	4438035
19G	1.1	20	25	4438019
19G	1.1	25	25	4438027
20G	0.9	15	25	4439759
20G	0.9	20	25	4439767
20G	0.9	25	25	4439775
20G	0.9	40	25	4439777
22G	0.7	15	25	4439694
22G	0.7	20	25	4439635
22G	0.7	25	25	4439686

Non-Coring Winged Needles for Access Ports

Surecan®

Winged Surecan® non-coring needle

- High pressure resistant up to 325 psi (22.4 bar)
- use for long-term infusions
- flexible wings for relieved puncture and fixation
- latex- and DEHP-free (wings)
- extension tubing with clamp (PVC-free)
- tubing length cannula to connector: 200 +/- 10 mm



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
19G	1.1	15	15	4448286
19G	1.1	20	15	4448294
19G	1.1	25	15	4448308
20G	0.9	15	15	4448332
20G	0.9	20	15	4448340
20G	0.9	25	15	4448359
20G	0.9	30	15	4448367
22G	0.7	12	15	4448375
22G	0.7	15	15	4448383
22G	0.7	20	15	4448391
22G	0.7	25	15	4448405

Winged Surecan® LP non-coring needle

- use for long-term infusions
- flexible wings for relieved puncture and fixation
- latex- and DEHP-free (wings)
- extension tubing with clamp (PVC-free)
- tubing length cannula to connector: 200 +/- 10 mm



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
19G	1.1	15	15	3448286
19G	1.1	20	15	3448294
19G	1.1	25	15	3448308
20G	0.9	15	15	3448332
20G	0.9	20	15	3448340
20G	0.9	25	15	3448359
20G	0.9	30	15	3448367
22G	0.7	12	15	3448375
22G	0.7	15	15	3448383
22G	0.7	20	15	3448391
22G	0.7	25	15	3448405

Winged Surecan® non-coring needle with Y-site

- use for long-term infusions
- flexible wings for relieved puncture and fixation
- latex- and DEHP-free (wings)
- extension tubing with clamp (PVC-free)
- tubing length cannula to connector: 200 +/- 10 mm
- Y-site configuration



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
19G	1.1	20	15	4448430
19G	1.1	25	15	4448448
20G	0.9	15	15	4448472
20G	0.9	20	15	4448480
20G	0.9	25	15	4448499
22G	0.7	15	15	4448529
22G	0.7	20	15	4448537
22G	0.7	25	15	4448545
22G	0.7	30	15	4448553

Non-Coring Safety Needles for Access Ports

Surecan® Safety II

**DEHP
FREE**

**LATEX
FREE**

- Intuitive handling
- Minimizes risk of needle stick injuries
- High patient comfort

Your product advantages at a glance

- High pressure resistant up to 325 psi (22.4 bar)
- MR Conditional
- Ergonomic wings to provide a reliable and firm grip
- Low profile for a discreet appearance
- Transparent fixation base ensures an improved visibility of the puncture site and allows early identification of infections
- Non-absorbent foam pad for increased patient comfort
- Deployment of safety mechanism is easy and confirmed visually and audibly

Easy removal

- 1 Stabilise the needle base on the port
 - 2 Firmly pull the wings up until you hear a "Click"
- Safety mechanism is enabled



Surecan® Safety II non-coring safety needle

- tubing length cannula to connector:
190 +/- 10 mm



Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
G 19	1.1	15	20	4447000
G 19	1.1	20	20	4447001
G 19	1.1	25	20	4447002
G 19	1.1	32	20	4447003
G 19	1.1	38	20	4447004
G 20	0.9	15	20	4447005
G 20	0.9	20	20	4447006
G 20	0.9	25	20	4447007
G 20	0.9	32	20	4447008
G 20	0.9	38	20	4447009
G 22	0.7	15	20	4447010
G 22	0.7	20	20	4447011
G 22	0.7	25	20	4447012
G 22	0.7	32	20	4447013

Non-Coring Safety Needles for Access Ports

Surecan® Safety II

Surecan® Safety II non-coring safety needle with Ultrasite® and Y-site


- Y-site configuration
- tubing length Y-site to connector: 82 +/- 10 mm
- tubing length cannula to Y-site: 90 +/- 10 mm




Size	Cannula diameter (mm)	Cannula length (mm)	Sales unit-pcs.	Reference
G 19	1.1	15	20	4447028
G 19	1.1	20	20	4447029
G 19	1.1	25	20	4447030
G 19	1.1	32	20	4447031
G 19	1.1	38	20	4447032
G 20	0.9	15	20	4447033
G 20	0.9	20	20	4447034
G 20	0.9	25	20	4447035
G 20	0.9	32	20	4447036
G 22	0.7	15	20	4447038
G 22	0.7	20	20	4447039
G 22	0.7	25	20	4447040

Please note for high pressure injection through Surecan® Safety II needles with Y-site:

- Ensure that the Ultrasite® valve is firmly attached
- Close the clamp on the main line
- Inject through the Ultrasite® valve

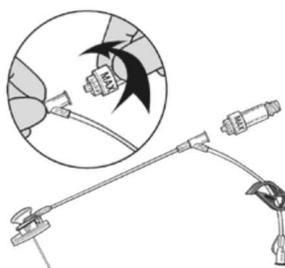

Surecan® Safety II with Ultrasite® Y



B. Braun Melsungen AG | Vascular Systems | Sieversufer 8 | 12359 Berlin | Germany
Phone +49 30 689897-0 | Fax +49 30 689897-30 | www.bbraun.com

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany
Phone +49 7461 95-0 | Fax +49 7461 95-26 00 | www.aesculap.com

Aesculap – a B. Braun company

Subject to technical changes. All rights reserved.
This brochure may only be used for the exclusive
purpose of obtaining information about our
products. Reproduction in any form partial or
otherwise is not permitted.

Brochure No. 6050179

1112/2.0/8