


ARGYLE™*

PERITONEAL DIALYSIS CATHETERS

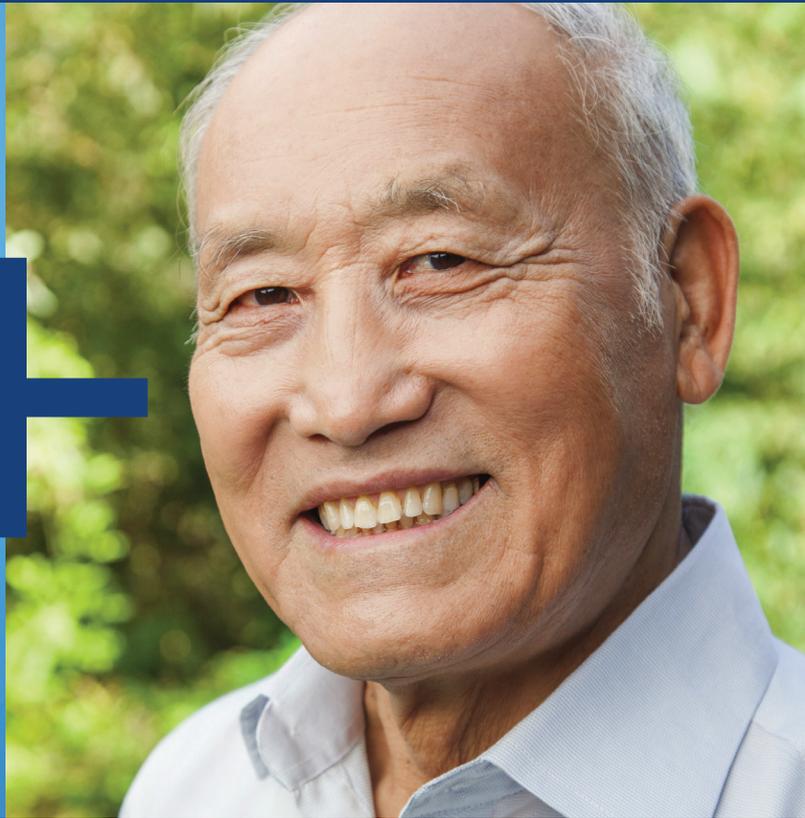

Medtronic
Further, Together

ARGYLE™* PERITONEAL DIALYSIS CATHETERS

Peritoneal dialysis catheters are intended to provide permanent access to the peritoneal cavity for patients to fill and drain their dialysate daily. Due to this daily, repetitive use, it's important to customize the catheter exit site for each patient to help facilitate easy patient and caregiver management. Optimized catheter placement and exit sites can help promote patient compliance for proper maintenance and may help reduce potential complications due to difficult to manage catheters. The Argyle™* peritoneal dialysis catheter line includes a variety of catheter configurations designed to enable patient customization to support long-term access preservation and patency.

TABLE OF CONTENTS

Catheter & Catheter Kit Ordering Guide	3-5
Tenckhoff Peritoneal Catheters	
Tenckhoff Peritoneal Catheters	6
1 Cuff	7
2 Cuffs.....	8
Oreopoulos-Zellerman Missouri Catheters	8
Swan Neck Peritoneal Catheters	
Swan Neck Tenckhoff Catheters	10
Swan Neck Curl Cath Catheters	11
Swan Neck Oreopoulos-Zellerman Catheters	11
Swan Neck Missouri Catheters	12
Moncrief-Popovich Swan Neck Curl Cath Catheters	12
Swan Neck Presternal Catheters	13-14
Curl Cath Peritoneal Catheters	
Curl Cath Catheter	15-16
1 Cuff	17
2 Cuffs	17
Pediatric Peritoneal Catheters	
Tenckhoff Catheters	19
Tenckhoff Curl Cath Catheters	19
Swan Neck Tenckhoff Catheters	20
Swan Neck Curl Cath Catheters	20
Swan Neck Missouri Catheters	21
Swan Neck Curl Cath Missouri Catheters	21
Accessories	22-25
Index	26-27

CATHETER & CATHETER KIT **ORDERING GUIDE**

CATHETER CONFIGURATION

The catheter products in this catalog can be ordered in either a catheter configuration or a catheter kit configuration.

Sterile Contents:

- Peritoneal Catheter (1)
 - Beta-Cap™ Adapter (1)
 - Beta-Cap™ Clamp (1)
-

When ordering a single catheter, you may also purchase the following accessories.

Single Catheter Accessories:

16Fr Pull-Apart Introducer Set with PTFE

8815544010	Includes:	5/case
	16Fr Pull-Apart Sheath†/Dilator (1)	
	18ga Introducer Needle (1)	
	12mL Syringe (1)	
	J/Straight 0.038 in Guide Wire (1)	
	†Made with PolyTetra FluoroEthylene (PTFE)	

Tunneling Stylet (Choose one tunneling stylet per catheter)

8817791001	Atraumatic Peritoneal Tunneling Stylet	5/case
8888415679	Faller Tunneling Stylet	1/case

Straightening Stylet, Non-sterile, Autoclavable (Choose one stylet per catheter)

8888415646	Insertion/Straightening Stylet, 44.5 cm	1/case
8888415661	Insertion/Straightening Stylet, 61 cm	1/case

Swan Neck Catheter Stencil (Choose one stencil if ordering the swan neck catheter)

8888412304	Swan Neck Tenckhoff Catheter Stencil, 2 Cuffs, Left/Right	1/case
8888412312	Swan Neck Missouri Catheter Stencil, 2 Cuffs, Left/Right	1/case

CATHETER KIT CONFIGURATION

Sterile Contents:

- Peritoneal Catheter (1)
- Beta-Cap™ Adapter, Cap and Clamp (1)
- 16Fr Pull-Apart Sheath/Dilator (1)
- 18ga Introducer Needle (1)
- J/Straight Guide Wire (1)
- 12mL Syringe (1)
- #11 Scalpel (1)
- Tunneling Stylet (1)
- 4 in x 4 in Gauze Sponges (6)

Peritoneal catheter kits are convenient and cost-efficient. Kits simplify inventory and save staff time when assembling supplies needed for catheter insertion.

Peritoneal catheter kits are recommended for percutaneous placement of the catheter.

The Peritoneal catheter kit provides an ideal alternative to catheter placement by laparotomy. Each kit features a 16Fr pull-apart introducer used to insert the catheter percutaneously using modified Seldinger technique. This technique uses a smaller incision — helping minimize the risk of postoperative dialysate leakage.

†Made with PolyTetra FluoroEthylene (PTFE)

When ordering a catheter kit, you may also purchase the following accessories.

Catheter Kit Accessories

Straightening Stylet (Choose one stylet per catheter)

8888415646	Insertion/Straightening Stylet, 44.5 cm	1/case
8888415661	Insertion/Straightening Stylet, 61 cm	1/case

Swan Neck Catheter Stencil (Choose one stencil if ordering the swan neck catheter)

8888412304	Swan Neck Tenckhoff Catheter Stencil, 2 Cuffs, Left/Right	1/case
8888412312	Swan Neck Missouri Catheter Stencil, 2 Cuffs, Left/Right	1/case

TENCKHOFF
**PERITONEAL
CATHETERS**

ARGYLE™* TENCKHOFF PERITONEAL CATHETERS

Recommended for peritoneal dialysis or intraperitoneal chemotherapy


Broad Selection — Choose from a variety of catheter lengths and cuff placements and types. Or choose the Argyle™** Universal Tenckhoff catheter, a unique catheter that can be trimmed to fit most patient sizes. The Tenckhoff catheter is offered in the following cuff configurations:

- One attached cuff – in either the subcutaneous or preperitoneal position.
- Two attached cuffs – one subcutaneous and one preperitoneal.
- Two disks with one or two cuffs – Oreopoulos-Zellerman, Missouri catheter style helps keep catheter in place and prevent omental wrapping.


Cost Efficient — For your convenience, Tenckhoff catheters are packaged one per carton, enabling you to stock multiple configurations at minimal cost. Catheters are also available in time-saving, cost-efficient kits. They contain the basic items needed for percutaneous insertion, featuring the pull-apart sheath/dilator and tunneling stylet.

Complete System — Each Tenckhoff catheter comes with an adapter, cap and clamp. The Ultem™** adapter forms a tight closure with standard CAPD transfer sets.


Tenckhoff Catheters — 1 Cuff

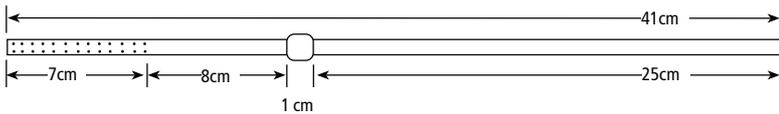

Catheter 8814843001 Tenckhoff Catheter, 1 Preperitoneal Cuff, 41 cm 1/case

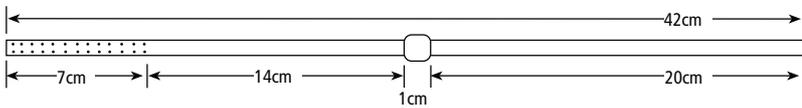

Catheter 8810889003 Tenckhoff Catheter, 1 Subcutaneous Cuff, 42 cm 1/case

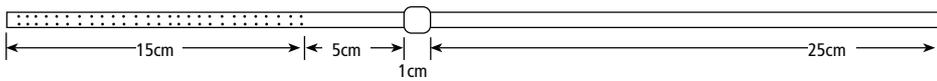


Catheter 8814843002 Tenckhoff Catheter, Universal, 1 Preperitoneal Cuff, 46 cm 1/case

Kit 8817278001 Tenckhoff Catheter, Universal, 1 Preperitoneal Cuff Kit, 46 cm 5/case


TENCKHOFF PERITONEAL CATHETERS

ARGYLE™* TENCKHOFF PERITONEAL CATHETERS (CONTINUED)

Tenckhoff Catheters — 2 Cuffs

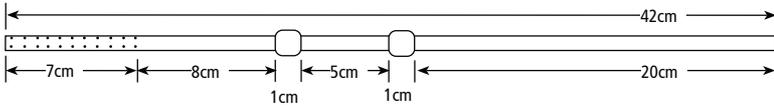


Catheter	8810888003	Tenckhoff Catheter, 2 Cuffs, 42 cm	1/case
Kit	8817278008	Tenckhoff Catheter, 2 Cuffs Kit, 42 cm	5/case

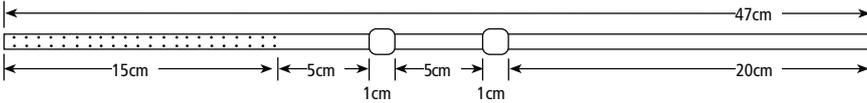

Catheter	8810888012	Tenckhoff Catheter, Universal, 2 Cuffs, 47 cm	1/case
----------	------------	---	--------

Oreopoulos-Zellerman Missouri Catheters

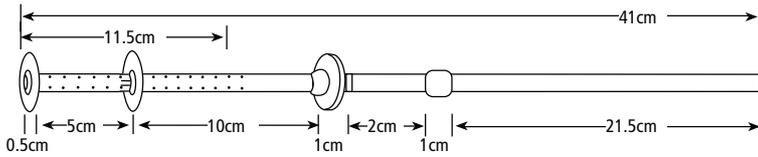

Catheter	8888411009	Oreopoulos-Zellerman Catheter, Missouri Catheter 1 Bead, 2 Disks, 2 Cuffs, 41 cm	1/case
----------	------------	---	--------



SWAN NECK PERITONEAL CATHETERS

SWAN NECK SPECIALTY PERITONEAL CATHETERS


Swan neck specialty peritoneal catheters give you the options and features you need to simplify access for peritoneal dialysis.

Broad Selection — Choose from a variety of catheter options: Tenckhoff, Curl Cath, Missouri, Oreopoulos-Zellerman and Presternal. Pediatric catheters are also available.

- The Swan Neck catheter employs a preformed 180 degree angle designed to keep both the internal and external portions of the catheter positioned downward. The angled catheter is designed to keep the internal portion from migrating out of position and to decrease stress on the exit site, which often occurs when a straight catheter is manually curved.
- The felt cuffs anchor the catheter subcutaneously and help establish a natural infection barrier through tissue ingrowth.
- A radiopaque stripe on the front of each catheter makes it easy to differentiate under X-ray. During insertion and postplacement, this stripe helps identify catheter twisting.
- Stencils are available for skin markings before surgery. Each stencil reflects the exact shape of the catheter and is critical to ensure proper catheter implantation.

Swan Neck Tenckhoff Catheters

The Swan Neck Tenckhoff catheter is a standard Tenckhoff with a preformed 180 degree angle. The angled catheter keeps the internal portion from migrating out of position and decreases stress on the exit site, which often occurs when a straight catheter is manually curved. Indicated for any patient requiring peritoneal dialysis or intraperitoneal chemotherapy.

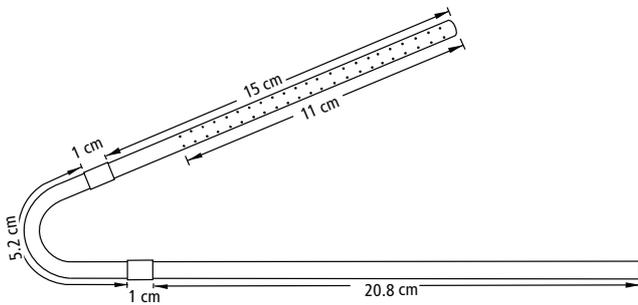


Catheter 8888412007 Swan Neck Tenckhoff Catheter, 2 Cuffs, Left, 43 cm 1/case

SWAN NECK PERITONEAL CATHETERS

Swan Neck Curl Cath Catheters

Swan Neck Curl Cath catheter employs the weight of the curl and the preformed bend to prevent internal catheter migration. The outside portion of the catheter is positioned downward, preventing cellular debris from collecting at the exit site.

Indicated for patients who have problems with catheter migration and inflow discomfort.

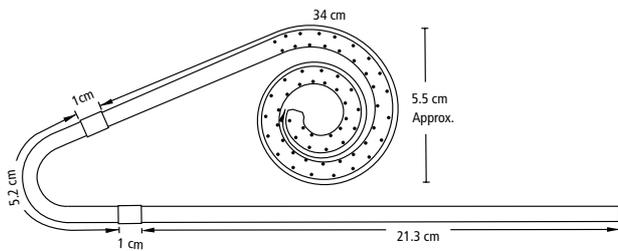


Catheter	8888413807	Swan Neck Curl Cath Catheter, 2 Cuffs, Left, 62.5 cm	1/case
Kit	8888413823	Swan Neck Curl Cath Catheter Kit, 2 Cuffs, Left, 62.5 cm	1/case
Catheter	8888413815	Swan Neck Curl Cath Catheter, 2 Cuffs, Right, 62.5 cm	1/case
Kit	8888413831	Swan Neck Curl Cath Catheter Kit, 2 Cuffs, Right, 62.5 cm	1/case

Swan Neck Oreopoulos-Zellerman Catheters

The Swan Neck Oreopoulos-Zellerman (also called Toronto Western) Catheter is available in Tenckhoff and Missouri configurations. These catheters are differentiated from other Swan Neck catheters by the two silicone disks added to the internal catheter segment. These disks are designed to prevent omental wrapping and to keep the catheter placed low in the peritoneum.

Indicated for patients who have repeated problems with catheter migration and omental occlusion.

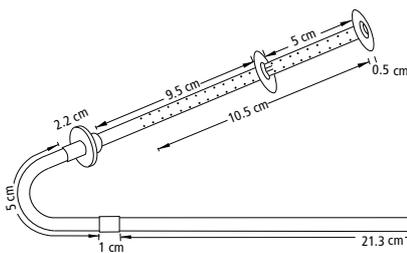


Catheter	8888413005	Swan Neck Oreopoulos-Zellerman Catheter, 2 Cuffs, Left, 44.5 cm	1/case
Catheter	8888413013	Swan Neck Oreopoulos-Zellerman Catheter, 2 Cuffs, Right, 44.5 cm	1/case

SWAN NECK SPECIALTY PERITONEAL CATHETERS (CONTINUED)


Swan Neck Missouri Catheters

The Swan Neck Missouri Catheter is the same basic design as the Swan Neck Tenckhoff and Curl Cath catheters with the addition of a silicone bead and a felt disk. The bead rests just inside the peritoneum to prevent dialysate leakage. The felt disk is placed just outside the peritoneum and is sutured to the rectus muscle. Fibroblasts grow into the fibers of the felt cuff to anchor the catheter. A second cuff in the subcutaneous tissue is designed to seal the exit site, helping to prevent bacteria from entering into the catheter tract.

Indicated for patients who experience dialysate leakage.

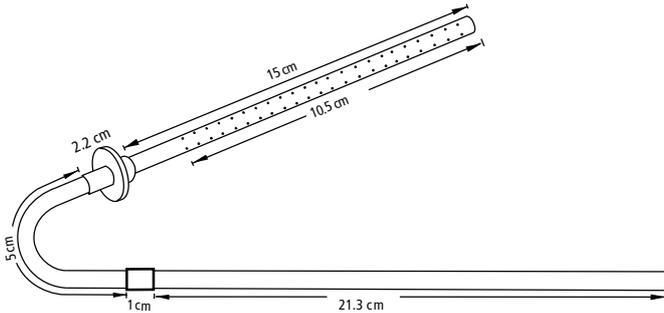


Catheter	8888412601	Swan Neck Missouri Catheter, 2 Cuffs, Left, 44.5 cm	1/case
Catheter	8888412619	Swan Neck Missouri Catheter, 2 Cuffs, Right, 44.5 cm	1/case

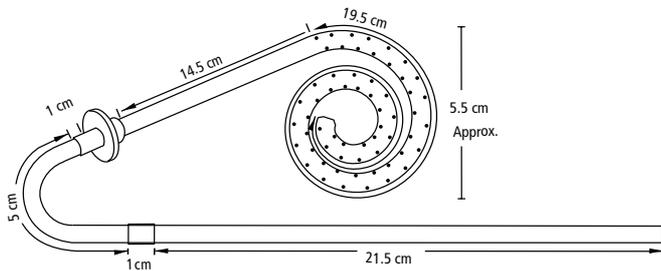

Catheter	8888413401	Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Left, 62.5 cm	1/case
Catheter	8888413419	Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Right, 62.5 cm	1/case

Moncrief-Popovich Swan Neck Curl Cath Catheters

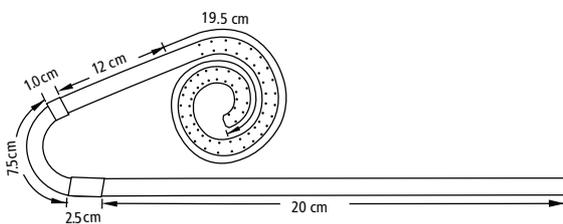

Catheter	8888414029	Moncrief-Popovich Swan Neck Curl Cath Catheter, 2 Cuffs, 62.5 cm	5/case
----------	------------	--	--------

SWAN NECK PERITONEAL CATHETERS


Swan Neck Presternal Catheters

The presteral (upper) segment of this catheter is tunneled into the subcutaneous tissue of the chest wall. The upper segment connects with the intraperitoneal (lower) segment by a double-barbed titanium connector. Both the presteral and intraperitoneal segments can be trimmed to custom fit the patient. The intraperitoneal segment is a variation of the Missouri catheter.


Indicated for:

- Patients with a more active lifestyle. As the exit site is located on the upper chest wall instead of at the belt line, this catheter is ideal for those patients who enjoy wading in a pool or soaking in a tub.
- Obese patients who may have difficulty manipulating a catheter at the belt line.
- Incontinent patients.

Assembled Presternal Catheter

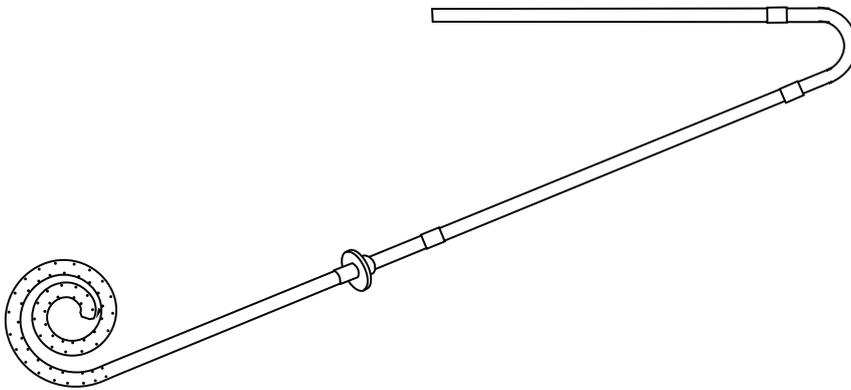


Presternal upper Segment 60.3 cm

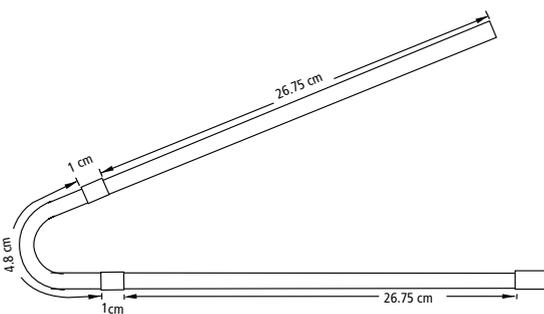

SWAN NECK SPECIALTY PERITONEAL CATHETERS (CONTINUED)

Swan Neck Presternal Catheters

Intraperitoneal Lower Segment 52 cm

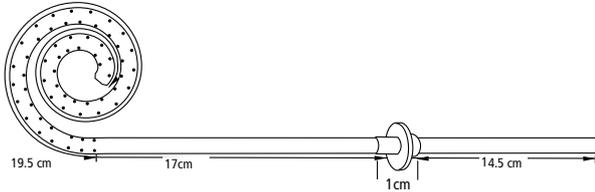

Catheter	8888414011	Swan Neck Presternal Catheter, 2 Cuffs, 112.3 cm (Swan Neck: 60.3 cm + 52 cm curl)	1/case
Catheter	8888111132	Swan Neck Pre-Sternal 2 Cuff, Curl Cath 1 Cuff, Tray	1/case
Catheter	8888121132	Presternal Peritoneal Dialysis Catheter Kit	5/case

When ordering the Swan Neck presteral catheter kit, also order the following procedural items:

8888412312	Swan Neck Missouri Catheter Stencil, 2 Cuffs, Left/Right	1/case
8888412304	Swan Neck Tenckhoff Catheter Stencil, 2 Cuffs, Left/Right	1/case
8888415679	Faller Tunneling Stylet	1/case
8888415661	Insertion/Straightening Stylet, 61 cm	1/case
8888415604	Titanium Two Part Adapter	1/case

Tunneler sold via Scanlaninternational.com.



SWAN NECK PERITONEAL CATHETERS

CURL CATH **CATHETERS**

CURL CATH CATHETERS

Recommended for peritoneal dialysis or intraperitoneal chemotherapy

The Curl Cath catheter, a modified Tenckhoff catheter, is suitable for acute or chronic peritoneal dialysis and intraperitoneal chemotherapy. It can be inserted either percutaneously or surgically.

The additional mass of the Curl Cath catheter keeps it low in the pelvic gutter, resulting in less catheter migration. The catheter's curled tubing and numerous inflow/outflow holes diffuse dialysate gently into the patient.

The Curl Cath catheter is made of translucent, medical-grade silicone rubber with a radiopaque stripe. It has a 4.9 mm OD and a 2.6 mm ID, and is available in adult or pediatric lengths. Unlike polyurethane catheters, the Curl Cath can be extended or repaired if necessary, eliminating the need for costly catheter replacement caused by external damage.

Felt cuffs anchor the catheter subcutaneously and help establish a natural infection barrier through tissue ingrowth. To meet the varying needs of patient and physician, two cuff configurations are offered:

- One attached cuff in the preperitoneal position
- Two attached cuffs, one subcutaneous and one preperitoneal

The Ultem™** adapter included in the packaging provides a tight closure with standard CAPD transfer sets. Titanium adapters are also available.

Curl Cath catheters are also available in convenient catheter kits, which provide the necessary components for inserting the catheter by modified Seldinger technique. Each Kit features a 16Fr Pull-Apart Introducer† and Peritoneal Catheter Tunneling Stylet.


†Made with PolyTetra FluoroEthylene (PTFE)

Product highlights include:


- Curled tubing with numerous inflow/outflow holes
- Can be inserted percutaneously or surgically
- Made of translucent, medical-grade silicone rubber
- Radiopaque stripe
- Can be extended or repaired if necessary
- Sterile

CURL CATH PERITONEAL CATHETERS

Curl Cath Catheters – 1 Cuff

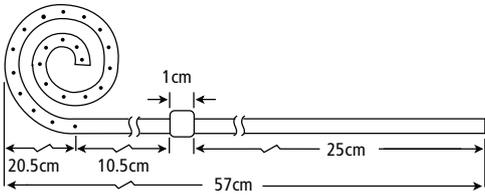

Catheter 8811313014 Curled Catheter, 1 Cuff, 57 cm 1/case

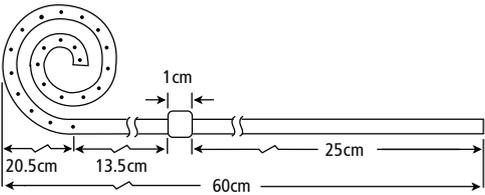


Catheter 8811313013 Curled Catheter, 1 Cuff, 60 cm 1/case

Catheter Kit 8817278010 Curled Catheter, 1 Cuff Kit, 60 cm 5/case


Catheter 8888411710 Curled Catheter, 1 Cuff, 62 cm 1/case

Curl Cath Catheters – 2 Cuffs

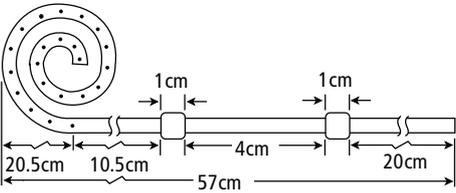

Catheter 8811313015 Curled Catheter, 2 Cuffs, 57 cm 1/case

Catheter Kit 8817278007 Curled Catheter, 2 Cuffs Kit, 57 cm 5/case

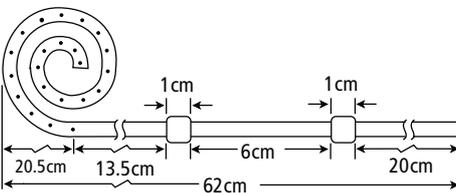

Catheter 8811313010 Curled Catheter, 2 Cuffs, 62 cm 1/case


Catheter Kit 8817278006 Curled Catheter, 2 Cuffs Kit, 62 cm 5/case

Catheter 8888411702 Curled Catheter, 2 Cuffs, Curled, 62 cm, Sterile 1/case

PEDIATRIC
**PERITONEAL
CATHETERS**


PEDIATRIC PERITONEAL CATHETERS

Pediatric Tenckhoff Catheters

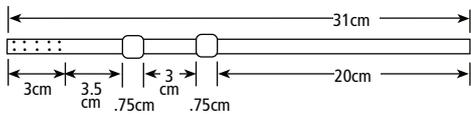


Catheter 8812329002 Pediatric Tenckhoff, 2 Cuffs, 31 cm **2.1** 1/case

Catheter 8888414227 Infant Tenckhoff, 2 Cuff, 31 cm, Sterile 1/case

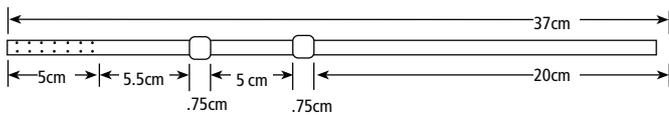


Catheter 8810890003 Pediatric Tenckhoff, 2 Cuffs, 37 cm **2.2** 1/case

Catheter 8888414201 Pediatric Tenckhoff, 2 Cuff, 30 cm, Sterile 1/case

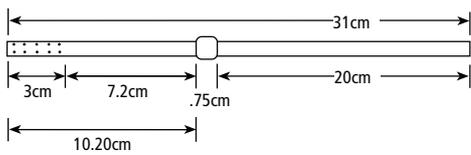


Catheter 8812329001 Pediatric Tenckhoff, 1 Subcutaneous Cuff, 31 cm 1/case

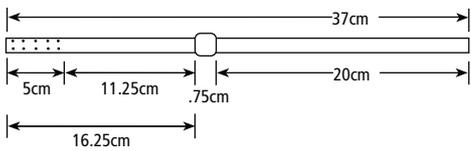

Catheter 8812321001 Pediatric Tenckhoff, 1 Subcutaneous Cuff, 37 cm 1/case

Pediatric Tenckhoff Curl Cath Catheter

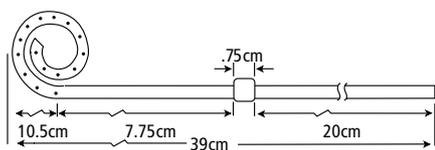


Catheter 8815677001 Pediatric Curl Cath Catheter, 1 Cuff, 39 cm 1/case

Catheter 8888417022 Pediatric Tenckhoff Curl Cath Catheter, 1 Cuff, 39 cm 1/case

PEDIATRIC PERITONEAL CATHETERS


PEDIATRIC PERITONEAL CATHETERS (CONTINUED)

Pediatric Swan Neck Tenckhoff Catheter

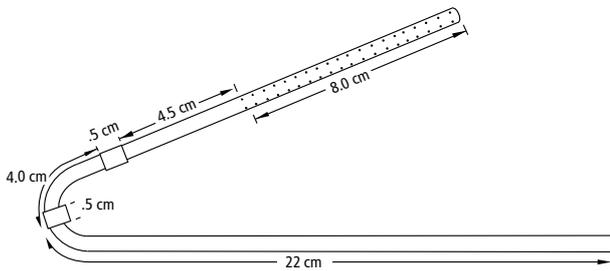


Catheter 8888410506 Pediatric Swan Neck Tenckhoff Catheter, 2 Cuffs, Left, 37.5 cm 1/case

Pediatric Swan Neck Curl Cath Catheters

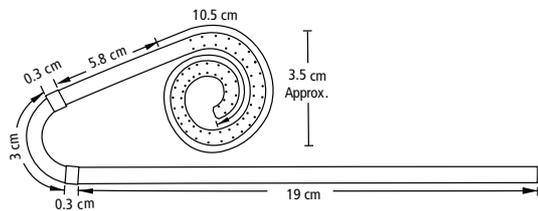


2.3

Catheter 8888413100 Infant Swan Neck Curl Cath Catheter, 2 Cuffs, 38.9 cm (Warady) 1/case

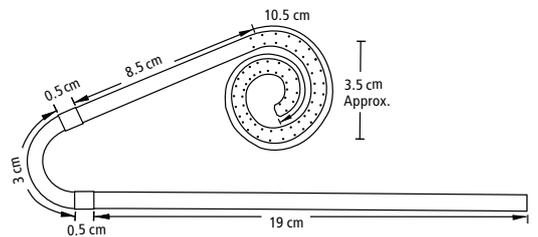


Catheter 8888414813 Pediatric Swan Neck Curl Cath Catheter, 2 Cuffs, Left, 42 cm 1/case

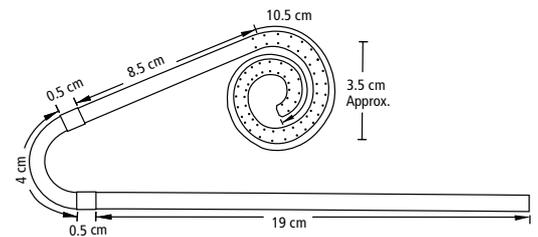

Catheter 8888413101 Pediatric Swan Neck Curl Cath Catheter, 2 Cuffs, 43 cm (Warady) 1/case


Catheter 8888414508 Pediatric, Tenckhoff, 1 Cuff, 34.25 cm 1/case

PEDIATRIC PERITONEAL CATHETERS

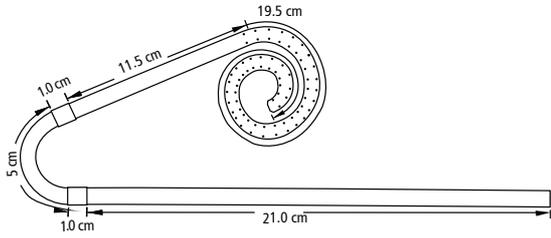


Catheter 8888413102 Adolescent Swan Neck Curl Cath Catheter, 2 Cuffs, 59 cm (Warady) 1/case

Pediatric Swan Neck Missouri Catheter

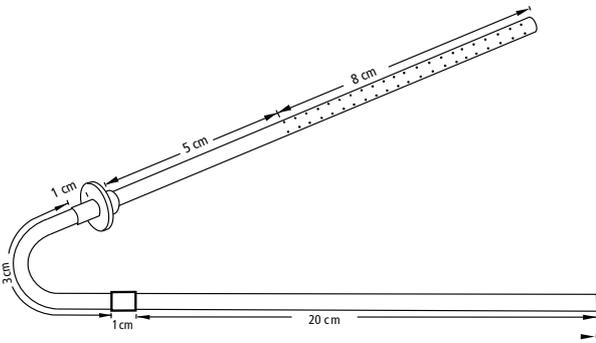

Catheter 8888412643 Pediatric Swan Neck Missouri Catheter, 2 Cuffs, Left, 38 cm 1/case

Pediatric Swan Neck Curl Cath Missouri Catheter

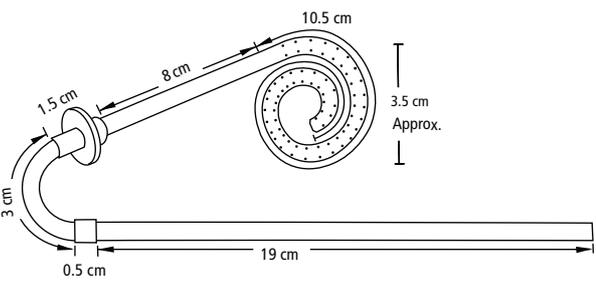

Catheter 8888415307 Pediatric Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Left, 42 cm 1/case

ACCESSORIES:
**PERITONEAL
CATHETERS**

ACCESSORIES — PERITONEAL CATHETERS

Peritoneal Catheter Straightening Stylet

Straightening stylets help keep the catheter straight for maneuvering during insertion into the peritoneum. Stylets come in a variety of sizes and can be used on all catheter types. Nonsterile. Autoclavable.

8888415646 Insertion/Straightening Stylet, 44.5 cm, for use with catheters 44.5 cm or shorter 1/case

8888415661 Insertion/Straightening Stylet, 61 cm, for use with catheters 61 cm or shorter 1/case

Atraumatic Peritoneal Tunneling Stylet

This atraumatic tunneling stylet penetrates tissue cleanly, creating a subcutaneous tunnel with minimal bleeding and tissue trauma. As the stylet dilates the tissue, the practitioner simultaneously pulls the catheter through the tunnel. Made of pliable stainless steel, it can be bent to create a customized, curved tunnel. The stylet is 5Fr X 23 cm (9"). Sterile. Single use only.

8817791001 Atraumatic Peritoneal Tunneling Stylet 5/case

Faller Tunneling Stylet

This stylet is an alternative to the atraumatic tunneling stylet. It has a cutting edge on the portion leading into the tunnel and a barbed portion on the opposite end, which connects to the catheter. Sterile.

8888415679 Faller Tunneling Stylet 1/case

Catheter Stencils

These stencils are used to provide a clear outline for the Swan Neck or Swan Neck Missouri catheters before implantation. Details include proper position on the abdominal wall and position of the silicone bead (Missouri) and also include two choices of exit sites so the physician can select the exit site most appropriate for each patient. Each stencil can be turned over to select left- or right-sided placement. Nonsterile.

8888412304 Swan Neck Tenckhoff Catheter Stencil, 2 Cuffs, Left/Right 1/case

8888412312 Swan Neck Missouri Catheter Stencil, 2 Cuffs, Left/Right 1/case

Titanium Adapter

An alternative to the Beta-Cap™ Adapter. Can be used on all Argyle™ peritoneal catheters. Sterile.

8888415604 Titanium Two Part Adapter 1/case

Scanlan™ Tunneler

Sold by Scanlaninternational.com.

Medical Grade Silicone Adhesive

This sterile adhesive is used to:

- Bond felt cuffs to the catheter
- Repair damaged catheters
- Extend the length of existing catheters

8810807001 Medical Grade Silicone Adhesive‡ 1/case

‡Must be ordered separately for any of the above uses.

ACCESSORIES — PERITONEAL CATHETERS

ACCESSORIES — PERITONEAL CATHETERS (CONTINUED)

Injection Sealing Caps

Used to close off the end of any Argyle™ Catheter. Universal luer fitting. Sterile. Latex-free.

8813791001 Injection Sealing Cap **50/case**

Beta-Cap™ Adapter

These adapters form a tight closure with standard CAPD transfer sets, or connect to a sealing cap. The Beta-Cap™ adapter is made of Ultem™ resin. Can be used with all Argyle™ Peritoneal Catheters. Sterile.

8814661001 Beta-Cap™ Adapter **25/case**

Beta-Cap™ Clamp

Used to clamp the peritoneal catheter tubing when not in use. Sterile.

8810805001 Beta-Cap™ Clamp **25/case**

Catheter Extender

Provides 6 inches of silicone tubing with a titanium connector. Not attached. Sterile.

8888415612 Adult Titanium Catheter Extender **1/case**

Peri-Patch™ Repair Kit

Use an Argyle™ Peritoneal Catheter Repair Kit to replace the external portion of a damaged catheter and eliminate the need for catheter replacement. All components are sterile. Kits are double packaged; the inner package can be placed directly on a sterile field.

NOTE: To repair an indwelling catheter, there must be at least 3 cm undamaged tubing at the exit site.

8888285001 Peri-Patch™ Repair Kit **1/case**

NOTE: Also order glue using code # **8810807001**

Sterile contents include:

- 1 15 cm Silicone Rubber Catheter Extension (2.6 mm ID, 4.9 mm OD) with Double-Barbed Connector
- 1 Beta-Cap™ Adapter, Cap and Clamp
- 1 Glue Mold

Introducer Set

Introducer sets are used when placing peritoneal catheters percutaneously.

8815544010 16Fr Pull-Apart Introducer Set **5/case**

Sterile contents include:

- 1 16Fr Pull-Apart Sheath/Dilator
- 1 18 ga Introducer Needle
- 1 12 mL Syringe
- 1 J/Straight 0.038" Guide Wire

†Made with PolyTetra FluoroEthylene (PTFE)

ACCESSORIES — PERITONEAL CATHETERS



INDEX

CODE	PRODUCT	PAGE	CODE	PRODUCT	PAGE
8810805001	Beta-Cap™ Clamp	17	8888412312	Swan Neck Missouri Catheter Stencil	2, 3, 10, 16
8810807001	Medical Grade Silicone Adhesive	16	8888412601	Swan Neck Missouri Catheter, 2 Cuffs, Left, 44.5 cm	8
8810888003	Tenckhoff Catheter, 2 Cuffs, 42 cm	5	8888412619	Swan Neck Missouri Catheter, 2 Cuffs, Right, 44.5 cm	8
8810888012	Tenckhoff Catheter, Universal, 2 Cuffs, 47 cm	5	8888412643	Pediatric Swan Neck Missouri Catheter, 2 Cuffs, Left, 38 cm	15
8810889003	Tenckhoff Catheter, 1 Subcutaneous Cuff, 42 cm	4	8888413005	Swan Neck Oreopoulos-Zellerman Catheter, 2 Cuffs, Left, 44.5 cm	7
8810890003	Pediatric Tenckhoff, 2 Cuffs, 37 cm	13	8888413013	Swan Neck Oreopoulos-Zellerman Catheter, 2 Cuffs, Right, 44.5 cm	7
8811313010	Curled Catheter, 2 Cuffs, 62 cm	12	8888413100	Infant Swan Neck Curl Cath Catheter, 2 Cuffs, 38.9 cm	14
8811313013	Curled Catheter, 1 Cuff, 60 cm	12	8888413101	Pediatric Swan Neck Curl Cath Catheter, 2 Cuffs, 43 cm	14
8811313014	Curled Catheter, 1 Cuff, 57 cm	12	8888413102	Adolescent Swan Neck Curl Cath Catheter, 2 Cuffs, 59 cm	14
8811313015	Curled Catheter, 2 Cuffs, 57 cm	12	8888413401	Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Left, 62.5 cm	8
8812321001	Pediatric Tenckhoff, 1 Subcutaneous Cuff, 37 cm	13	8888413419	Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Right, 62.5 cm	8
8812329001	Pediatric Tenckhoff, 1 Subcutaneous Cuff, 31 cm	13	8888413807	Swan Neck Curl Cath Catheter, 2 Cuffs, Left, 62.5 cm	7
8812329002	Pediatric Tenckhoff, 2 Cuffs, 31 cm	13	8888413815	Swan Neck Curl Cath Catheter, 2 Cuffs, Right, 62.5 cm	7
8813791001	Injection Sealing Cap	17	8888413823	Swan Neck Curl Cath Catheter Kit, 2 Cuffs, Left, 62.5 cm	7
8814661001	Beta-Cap™ Adapter	17	8888413831	Swan Neck Curl Cath Catheter Kit, 2 Cuffs, Right, 62.5 cm	7
8814843001	Tenckhoff Catheter, 1 Preperitoneal Cuff, 41 cm	4	8888414011	Swan Neck Presternal Catheter, 2 Cuffs, 112.3 cm	10
8814843002	Tenckhoff Catheter, Universal, 1 Preperitoneal Cuff, 46 cm	4	8888414029	Moncrief-Popovich Swan Neck Curl Cath Catheter, 2 Cuffs, 62.5 cm	8
8815544010	16Fr Pull-Apart Introducer Set	2, 17	8888414813	Pediatric Swan Neck Curl Cath Catheter, 2 Cuffs, Left, 42 cm	14
8815677001	Pediatric Curl Cath Catheter, 1 Cuff, 39 cm	13	8888415307	Pediatric Swan Neck Curl Cath Missouri Catheter, 2 Cuffs, Left, 42 cm	15
8817278001	Tenckhoff Catheter, Universal, 1 Preperitoneal Cuff Kit, 46 cm	4	8888415604	Titanium Two Part Adapter	10, 16
8817278006	Curled Catheter, 2 Cuffs Kit, 62 cm	12	8888415612	Adult Titanium Catheter Extender	17
8817278007	Curled Catheter, 2 Cuffs Kit, 57 cm	12	8888415646	Insertion/Straightening Stylet, 44.5 cm	2, 3, 16
8817278008	Tenckhoff Catheter, 2 Cuffs Kit, 42 cm	5	8888415661	Insertion/Straightening Stylet, 61 cm	2, 3, 10, 16
8817278010	Curled Catheter, 1 Cuff Kit, 60 cm	12	8888415679	Faller Tunneling Stylet	2, 10, 16
8817791001	Atraumatic Peritoneal Tunneling Stylet	2, 16			
8888111132	Swan Neck Pre-Sternal 2 Cuff, Curl Cath 1 Cuff, Tray 1/Case	14			
8888121132	Presternal Peritoneal Dialysis Catheter Kit 5/Case	14			
8888285001	Peri-Patch™ Repair Kit	17			
8888410506	Pediatric Swan Neck Tenckhoff Catheter, 2 Cuffs, Left, 37.5 cm	13			
8888411009	Oreopoulos-Zellerman Catheter, 1 Bead, 2 Disks, 2 Cuffs, 41 cm	5			
8888412007	Swan Neck Tenckhoff Catheter, 2 Cuffs, Left, 43 cm	6			
8888412304	Swan Neck Tenckhoff Catheter Stencil, 2 Cuffs, Left/Right	2, 3, 10, 16			



© 2020 Medtronic. All rights reserved. Medtronic, Medtronic logo and Further, Together are trademarks of Medtronic.™ Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company. Argyle™ is a trademark of Cardinal Health, Inc., used under license. 04/2020 – US180486(1) – [WF#1913136]

15 Hampshire Street
Mansfield, MA 02048

800.962.9888
508.261.8000

[medtronic.com/covidien](https://www.medtronic.com/covidien)

Medtronic