


LEGFLOW® XP - OTW

LEGFLOW® XP A STEP FORWARD IN THE TREATMENT OF PERIPHERAL LESIONS

▶ Safety and Efficacy clinically proven in real life patients:

- ▶ **ReFlow Study:** International multicenter registry in long (>20 cm) and complex (TASC C8D) lesions¹
- ▶ **Primary Patency:** 71.1% at 12 months and 65.6% at 24 months;
- ▶ **Freedom from TLR:** 79.9% at 12 months and 70.7% at 24 months;

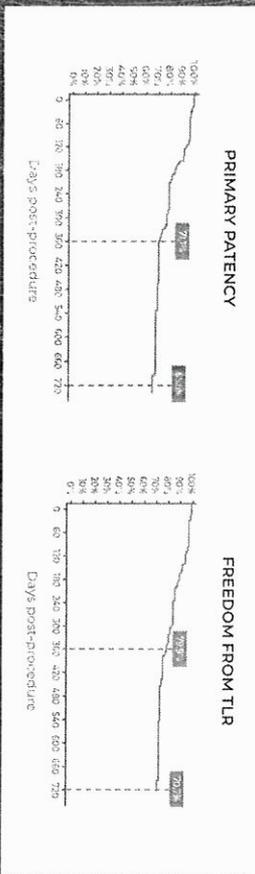

SAFEPAx® THE MOST ADVANCED COATING TECHNOLOGY AVAILABLE

A NEW DCB ERA PROVIDING:

- ▶ **No distal embolization** thanks to **No Flaking** effect for the treatment of the most fragile patients;
- ▶ **No Reflow** and **No Slow Flow** phenomena;
- ▶ **Lowest Wash-Off** effect among contemporaneous DCBs ensuring the optimal therapeutic treatment
- ▶ **Highest deliverability** for an unmatched navigability, even in the most complex and distal districts;
- ▶ **Safety and Efficacy** clinically proven on thousands of patients including PCTs with F.U. up to 5 years.

¹ Dr. Boersma data presented @ Charing Cross 2021
 Correlation data on file

ORDERING INFORMATIONS

LEGFLOW® XP OTW 0.014

3 p02.

UCL (cm)	Balloon Ø (mm)	Balloon Length (mm)											
		20	30	40	50	60	70	80	90	100	120	150	200
150	4.0	LXSS 3.0-20	LXSS 3.0-30	LXSS 3.0-40	LXSS 3.0-50	LXSS 3.0-60	LXSS 3.0-70	LXSS 3.0-80	LXSS 3.0-90	LXSS 3.0-100	LXSS 3.0-120	LXSS 3.0-150	LXSS 3.0-200
	5.0	LXSS 4.0-20	LXSS 4.0-30	LXSS 4.0-40	LXSS 4.0-50	LXSS 4.0-60	LXSS 4.0-70	LXSS 4.0-80	LXSS 4.0-90	LXSS 4.0-100	LXSS 4.0-120	LXSS 4.0-150	LXSS 4.0-200
	6.0	LXSS 5.0-20	LXSS 5.0-30	LXSS 5.0-40	LXSS 5.0-50	LXSS 5.0-60	LXSS 5.0-70	LXSS 5.0-80	LXSS 5.0-90	LXSS 5.0-100	LXSS 5.0-120	LXSS 5.0-150	LXSS 5.0-200
	7.0	LXSS 6.0-20	LXSS 6.0-30	LXSS 6.0-40	LXSS 6.0-50	LXSS 6.0-60	LXSS 6.0-70	LXSS 6.0-80	LXSS 6.0-90	LXSS 6.0-100	LXSS 6.0-120	LXSS 6.0-150	LXSS 6.0-200
	8.0	LXSS 7.0-20	LXSS 7.0-30	LXSS 7.0-40	LXSS 7.0-50	LXSS 7.0-60	LXSS 7.0-70	LXSS 7.0-80	LXSS 7.0-90	LXSS 7.0-100	LXSS 7.0-120	LXSS 7.0-150	LXSS 7.0-200
	9.0	LXSS 8.0-20	LXSS 8.0-30	LXSS 8.0-40	LXSS 8.0-50	LXSS 8.0-60	LXSS 8.0-70	LXSS 8.0-80	LXSS 8.0-90	LXSS 8.0-100	LXSS 8.0-120	LXSS 8.0-150	LXSS 8.0-200
	10.0	LXSS 9.0-20	LXSS 9.0-30	LXSS 9.0-40	LXSS 9.0-50	LXSS 9.0-60	LXSS 9.0-70	LXSS 9.0-80	LXSS 9.0-90	LXSS 9.0-100	LXSS 9.0-120	LXSS 9.0-150	LXSS 9.0-200
	12.0	LXSS 10.0-20	LXSS 10.0-30	LXSS 10.0-40	LXSS 10.0-50	LXSS 10.0-60	LXSS 10.0-70	LXSS 10.0-80	LXSS 10.0-90	LXSS 10.0-100	LXSS 10.0-120	LXSS 10.0-150	LXSS 10.0-200
	15.0	LXSS 12.0-20	LXSS 12.0-30	LXSS 12.0-40	LXSS 12.0-50	LXSS 12.0-60	LXSS 12.0-70	LXSS 12.0-80	LXSS 12.0-90	LXSS 12.0-100	LXSS 12.0-120	LXSS 12.0-150	LXSS 12.0-200
	2.0	LXSS 15.0-20	LXSS 15.0-30	LXSS 15.0-40	LXSS 15.0-50	LXSS 15.0-60	LXSS 15.0-70	LXSS 15.0-80	LXSS 15.0-90	LXSS 15.0-100	LXSS 15.0-120	LXSS 15.0-150	LXSS 15.0-200
120	2.0	LXSS 2.0-40	LXSS 2.0-60	LXSS 2.0-80	LXSS 2.0-100	LXSS 2.0-120	LXSS 2.0-150	LXSS 2.0-200					
	2.5	LXSS 2.5-40	LXSS 2.5-60	LXSS 2.5-80	LXSS 2.5-100	LXSS 2.5-120	LXSS 2.5-150	LXSS 2.5-200					
	3.0	LXSS 3.0-40	LXSS 3.0-60	LXSS 3.0-80	LXSS 3.0-100	LXSS 3.0-120	LXSS 3.0-150	LXSS 3.0-200					
	3.5	LXSS 3.5-40	LXSS 3.5-60	LXSS 3.5-80	LXSS 3.5-100	LXSS 3.5-120	LXSS 3.5-150	LXSS 3.5-200					
150	2.0	LXSS 2.0-40	LXSS 2.0-60	LXSS 2.0-80	LXSS 2.0-100	LXSS 2.0-120	LXSS 2.0-150	LXSS 2.0-200					
	2.5	LXSS 2.5-40	LXSS 2.5-60	LXSS 2.5-80	LXSS 2.5-100	LXSS 2.5-120	LXSS 2.5-150	LXSS 2.5-200					
	3.0	LXSS 3.0-40	LXSS 3.0-60	LXSS 3.0-80	LXSS 3.0-100	LXSS 3.0-120	LXSS 3.0-150	LXSS 3.0-200					
	3.5	LXSS 3.5-40	LXSS 3.5-60	LXSS 3.5-80	LXSS 3.5-100	LXSS 3.5-120	LXSS 3.5-150	LXSS 3.5-200					
4.0	LXSS 4.0-40	LXSS 4.0-60	LXSS 4.0-80	LXSS 4.0-100	LXSS 4.0-120	LXSS 4.0-150	LXSS 4.0-200						

LEGFLOW® XP

UCL (cm)	Balloon Ø (mm)	Balloon Length (mm)											
		20	30	40	50	60	70	80	90	100	120	150	200
15	2.0	LXSS 15.0-20	LXSS 15.0-30	LXSS 15.0-40	LXSS 15.0-50	LXSS 15.0-60	LXSS 15.0-70	LXSS 15.0-80	LXSS 15.0-90	LXSS 15.0-100	LXSS 15.0-120	LXSS 15.0-150	LXSS 15.0-200
	2.5	LXSS 15.0-20	LXSS 15.0-30	LXSS 15.0-40	LXSS 15.0-50	LXSS 15.0-60	LXSS 15.0-70	LXSS 15.0-80	LXSS 15.0-90	LXSS 15.0-100	LXSS 15.0-120	LXSS 15.0-150	LXSS 15.0-200
	3.0	LXSS 15.0-20	LXSS 15.0-30	LXSS 15.0-40	LXSS 15.0-50	LXSS 15.0-60	LXSS 15.0-70	LXSS 15.0-80	LXSS 15.0-90	LXSS 15.0-100	LXSS 15.0-120	LXSS 15.0-150	LXSS 15.0-200
	3.5	LXSS 15.0-20	LXSS 15.0-30	LXSS 15.0-40	LXSS 15.0-50	LXSS 15.0-60	LXSS 15.0-70	LXSS 15.0-80	LXSS 15.0-90	LXSS 15.0-100	LXSS 15.0-120	LXSS 15.0-150	LXSS 15.0-200

Contact your sales Representative for any inquiries and product availability
¹ UCL, usable catheter length

DATASHEET

LEGFLOW[®] XP Paclitaxel Releasing Peripheral Balloon Dilatation Catheter

General details: *3 pos.*

Manufacturer Company:	Cardionovum [®] GmbH
CE Mark released by:	CE 1434
Sterilization method:	Ethylene oxide
Usage:	Single-use
Storage recommendations:	0°C - 25°C away from light, humidity
Product Shelf Life:	2 years

Technical details:

Drug	
Coating Technology:	SAFEPAX [®]
Active Pharmaceutical Ingredient (API):	Paclitaxel
Excipient:	Shelloic Ammonium Salt
Dosage:	3 µg / mm ²

Catheter	
Catheter Architecture:	Bi-lumen OTW (Over The Wire)
Material:	Hub: Polycarbonate Shaft: PEBA Distal GW tubing: PEBA
Catheter Usable Lengths:	80 cm (LXS), 135 (LXM) and 150 cm (LXL)
Tip Length:	6 mm
Guidewire Compatibility:	0.014"
Introducer Sheath Compatibility:	<ul style="list-style-type: none"> • Diameters 3.0 – 5.0 [mm], length 20 – 200 [mm] – 6F • Diameter 6.0 [mm], length 20 – 80 [mm] – 6F • Diameter 6.0 – 7.0 [mm], length 100 – 200 [mm] – 6F • Diameters 7.0 – 10.0 [mm], length 20 – 80 [mm] – 6F • Diameter 8.0 [mm], length 100 – 200 [mm] – 8F • Diameter 12.0 [mm], length 20 – 80 [mm] – 8F

Balloon	
Radiopaque markers:	2 Pt - Ir
Rewrap Mechanism:	Ø3 & 4mm balloons: Three wing clockwise folded balloon with a rewrap (shape memory) Ø 5 to 12mm balloons: Six wing clockwise folded balloon with a rewrap (shape memory)
Nominal Pressure (NP):	<ul style="list-style-type: none"> • 8 bar for balloon diameters 3.00 – 9.00 with lengths 20 to 80 mm • 6 bar for balloon diameters 10.00 – 12.00 with lengths 20 to 80 mm • 6 bar for balloon diameters 4.00 – 7.00 with lengths 100 to 200 mm • 5 bar for balloon diameters 8.00 with lengths 100 to 200 mm
Rate Burst Pressure (RBP):	<ul style="list-style-type: none"> • 14 bar for balloon diameters 3.00 – 7.00 mm • 13 bar for balloon diameter 8.00 mm