

Over-The-Wire  
Design

104 pozivaja

Accordion CoAx™  
Stone Control Device

Simplify and Shorten  
*Endoscopic  
Lithotripsy*  
Procedures

Finally, a device to provide stone fragment  
localization and removal, **over a guidewire**,  
during endoscopic lithotripsy procedures...

**It's About Time!®**

PercSys.  
INNOVATIONS IN UROLOGY

# Accordion CoAx™ Stone Control Device

The Accordion CoAx Stone Control Device provides physicians stone fragment localization during endoscopic lithotripsy procedures:

- Prevents retrograde stone fragment migration and facilitates fragment removal during **ureteroscopic lithotripsy** procedures
- Prevents antegrade fragment migration into the ureter during **percutaneous nephrolithotomy (PCNL)**

## Ureteroscopic Lithotripsy

**Tracks** like a catheter


- Coaxial design allows **introduction over a 0.035" or 0.038" guidewire**
- Hydrophilic coating enables undeployed device to **track past impacted stones** like a low profile catheter
- Tapered tip **facilitates entry** into the ureteral orifice and past the stone

**Secures** like a basket

- Designed to **prevent retrograde migration** of the stone or small fragments
- **Withstands inadvertent contact** with holmium laser and intracorporeal lithotripter energy
- Linear actuation enables easy **release from intact stones** when required


**Clears** like nothing else

- **Conforms to and fills lumen** of the ureter
- Allows irrigant inflow to **rinse smaller lithotripsy fragments** distally
- **Sweeps larger fragments** into the bladder

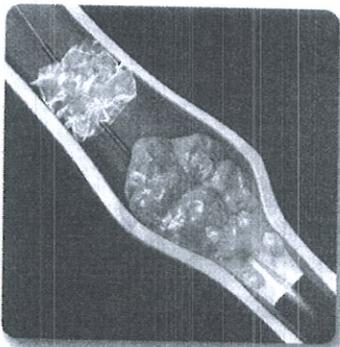


## Percutaneous Nephrolithotomy (PCNL)

- Can be **introduced over a 0.035" or 0.038" guidewire** using either an antegrade or retrograde approach
- Designed to **prevent migration of stone fragments** into the ureter
- **Conforms to and fills most structures** in the proximal urinary tract
- **Provides a stable platform** that is not easily displaced during the fragmentation process
- **Facilitates atraumatic entry** either through the skin or via the ureteral orifice
- Device advancement and film deployment **easily monitored using fluoroscopic guidance**

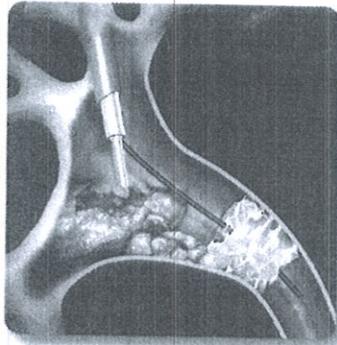

Using an **antegrade** approach, the low profile shaft allows introduction either alongside or inside the nephrostomy sheath

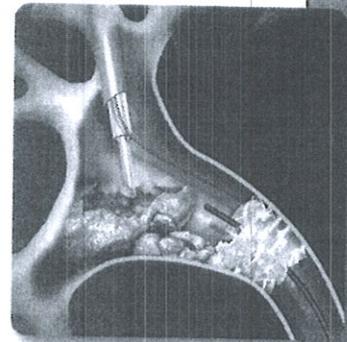

Using a **retrograde** approach, the guidewire can be removed and contrast media may be injected through the central lumen to conduct a retrograde pyelogram

### Ordering Information

Accordion CoAx Stone Control Device

Deployed OD: (nominal)	Length:	Guidewire Compatibility:	Order No:
10mm	75cm	0.035" or 0.038"	PA1305-01-10-05
15mm	75cm	0.035" or 0.038"	PA1305-01-15-05

PercSys  
3260 Hillview Avenue, Suite 100  
Palo Alto, CA 94304  
Toll Free: 1.877.987.6633  
Tel: 1.650.493.4200  
Fax: 1.650.493.4212

[www.percsys.com](http://www.percsys.com)

**PercSys**  
INNOVATIONS IN UROLOGY

© 2010 Percutaneous Systems, Inc. U.S. Patent # 7,462,183 B2; Additional U.S. and Foreign Patents Pending  
PercSys and It's About Time are registered trademarks and Accordion CoAx is a trademark of Percutaneous Systems, Inc.  
90-0949-01 (R 5/2010)