

# RADIFOCUS® Glidewire Advantage®

Guidewire

## PRODUCT SPECIFICATIONS

1. Diametras 0.014", 0.018", 0.035"  
 2. Ilgis 180 cm, 260 cm, 300 cm  
 3. PTFE danga (dangos ilgis 155, 235, 275 cm.)  
 4. 25 cm distalinė dalis, padengta hidrofiline danga.  
 5. Viršūnės smailėjimo ilgis - 5 cm.  
 6. Viršūnės modifikacijos: tiesus, kampu.

Product code	Shaft	Tip coil marker	Wire outer diameter	Overall length	Proximal spiral PTFE coating length	Distal hydrophilic coated length	Flexible tip length (taper)	Tip figure
RA*FA14181CM	Extra stiff	Yes	0.014" / 0.36 mm	180 cm	155 cm	25 cm	1 cm	Angled
RA*FA14301CM	Extra stiff	Yes	0.014" / 0.36 mm	300 cm	275 cm	25 cm	1 cm	Angled
RA*FS14301CM	Extra stiff	Yes	0.014" / 0.36 mm	300 cm	275 cm	25 cm	1 cm	Straight
RA*FA18181CM	Extra stiff	Yes	0.018" / 0.46 mm	180 cm	155 cm	25 cm	1 cm	Angled
RA*FA18301CM	Extra stiff	Yes	0.018" / 0.46 mm	300 cm	275 cm	25 cm	1 cm	Angled
RA*FS18301CM	Extra stiff	Yes	0.018" / 0.46 mm	300 cm	275 cm	25 cm	1 cm	Straight
RA*CA35185CM	Extra stiff	No	0.035" / 0.89 mm	180 cm	155 cm	25 cm	5 cm	Angled
RA*CA35265CM	Extra stiff	No	0.035" / 0.89 mm	260 cm	235 cm	25 cm	5 cm	Angled
RA*CS35185CM	Extra stiff	No	0.035" / 0.89 mm	180 cm	155 cm	25 cm	5 cm	Straight
RA*CS35265CM	Extra stiff	No	0.035" / 0.89 mm	260 cm	235 cm	25 cm	5 cm	Straight

Packed 1 per unit box.

## GLIDE WITH CONFIDENCE

- Radifocus® Guidewire M Standard type
- Radifocus® Guidewire M Stiff type
- Radifocus® Guidewire M Shapeable type
- Radifocus® Guidewire M Non-Vascular type
- Radifocus® Guidewire M with Radiopaque tip marker
- Radifocus® Guidewire GT with gold coil
- Radifocus® Guidewire Advantage Peripheral Guidewires
- Muso® Wire Stainless steel PTA Guidewire

poz 2. Specialus pravedėjas



RADIFOCUS®  
**Glidewire Advantage®**  
Guidewire

THIS IS YOUR  
**ADVANTAGE!**



**Terumo Corporation**  
44-1, 2-chome, Hatagaya  
Shibuya-ku, Tokyo 151-0072,  
Japan  
Tel.: +81 3 3374 8111  
Fax: +81 3 3374 8196

**Terumo Europe NV**  
Interleuvenlaan 40  
3001 Leuven, Belgium  
Tel.: +32 16 38 12 11  
Fax: +32 16 40 02 49

**Terumo Europe NV**  
Interventional Systems  
Liaison Office  
Immeuble Inside  
7-9 Rue Eugène et Armand Peugeot  
92500 Rueil Malmaison,  
France  
Tel.: +33 147 16 09 30  
Fax: +33 147 16 09 49

**Terumo Europe NV**  
Benelux Sales Division  
Belgium:  
Tel.: +32 16 39 25 80  
The Netherlands:  
Tel.: 0800 0220396

**Terumo Europe NV**  
Emerging Markets Division  
Tel.: +32 16 38 13 08

**Terumo Deutschland GmbH**  
Tel.: +49 6196 80 230

**Terumo Deutschland GmbH**  
Zweigniederlassung Spreitenbach  
Tel.: +41 56 419 10 10

**Terumo Deutschland GmbH**  
Zweigniederlassung Österreich  
Tel.: +43 2236 378020

**Terumo France S.A.S.**  
Tel.: +33 130 96 13 00

**Terumo Italia S.r.l.**  
Tel.: +39 0651 96 14 20

**Terumo Europe España SL**  
Tel.: +34 9021 01 298

**Terumo UK Ltd**  
Tel.: +44 1276 480440

**Terumo Russia LLC**  
Tel.: +7 495 988 4740

**Terumo Sweden AB**  
Tel.: +46 3174 85 880

**Terumo Corporation Dubai Branch**  
(Middle East and Africa)  
Tel.: +971 4 221 2220

**Terumo BCT Tıbbi Cihazlar Dağıtım ve Hizmetleri A.Ş.**  
Tel.: +90 216 590 0074



© Registered Trademark  
Published by Terumo Europe NV


www.terumo-europe.com

IS3886B03131TYV(1015)E



# HYBRID TECHNOLOGY FOR OUTSTANDING PERFORMANCE.

**Glidewire Advantage**, dedicated to all peripheral intervention levels, from sheath placement to ilio-femoro-popliteal and below-the-knee treatments. It combines maneuverability, crossability and device support to **reach and cross** the target lesion.

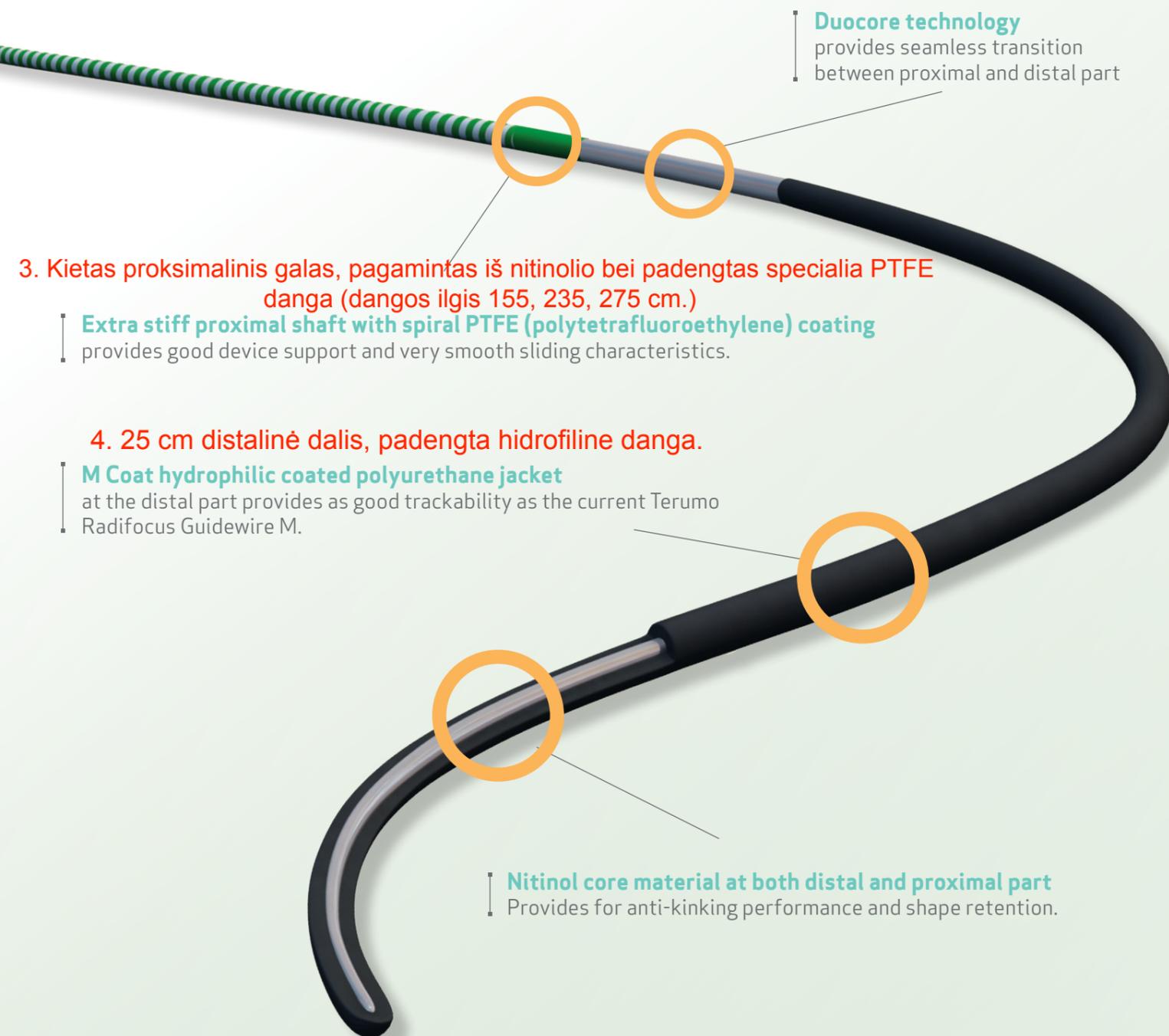

## 1<sup>ST</sup> ADVANTAGE: DURABILITY


Anti-kinking performance and shape retention capability due to Nitinol core material (both distal tip and proximal body).

## 2<sup>ND</sup> ADVANTAGE: CROSSABILITY

Crossing behaviour similar to a Terumo Radifocus Guidewire M - with the unmatched hydrophilic properties of M Coat at the distal 25 cm. **4. 25 cm distalinė dalis, padengta hidrofiline danga.**

## 3<sup>RD</sup> ADVANTAGE: STEERABILITY

Complete navigation and torque control due to the special treated Nitinol shaft-material.

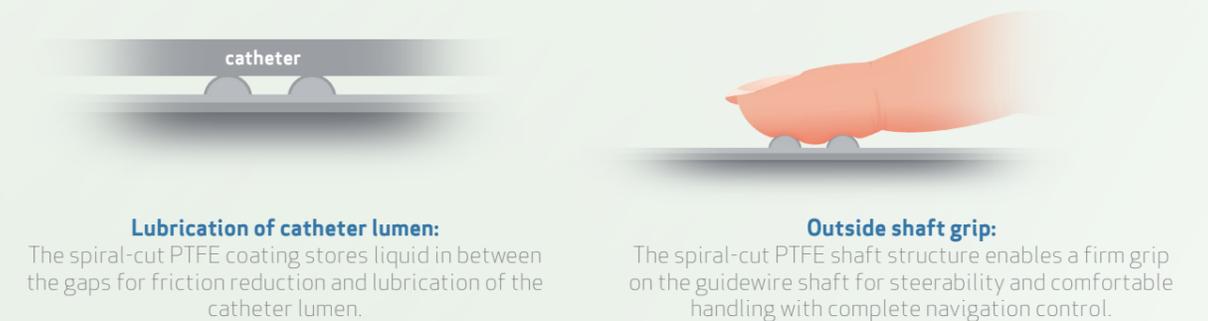

## 4<sup>TH</sup> ADVANTAGE: GUIDEWIRE SUPPORT

The extra stiff proximal shaft with the proprietary spiral-cut PTFE coating provides good device support and very smooth steerability even in complex anatomies. The proximal and distal shaft parts are connected seamlessly with unique DuoCore fusing technology to provide complete navigation control.

GLIDEWIRE ADVANTAGE OFFERS AN ALL LEVEL PERIPHERAL TREATMENT, FROM SHEATH PLACEMENT TO ILIACO-FEMORO-POPLITEAL AND BELOW-THE-KNEE TREATMENTS.

