


CLIMBER™

PTCA Guiding Catheter

Route to success

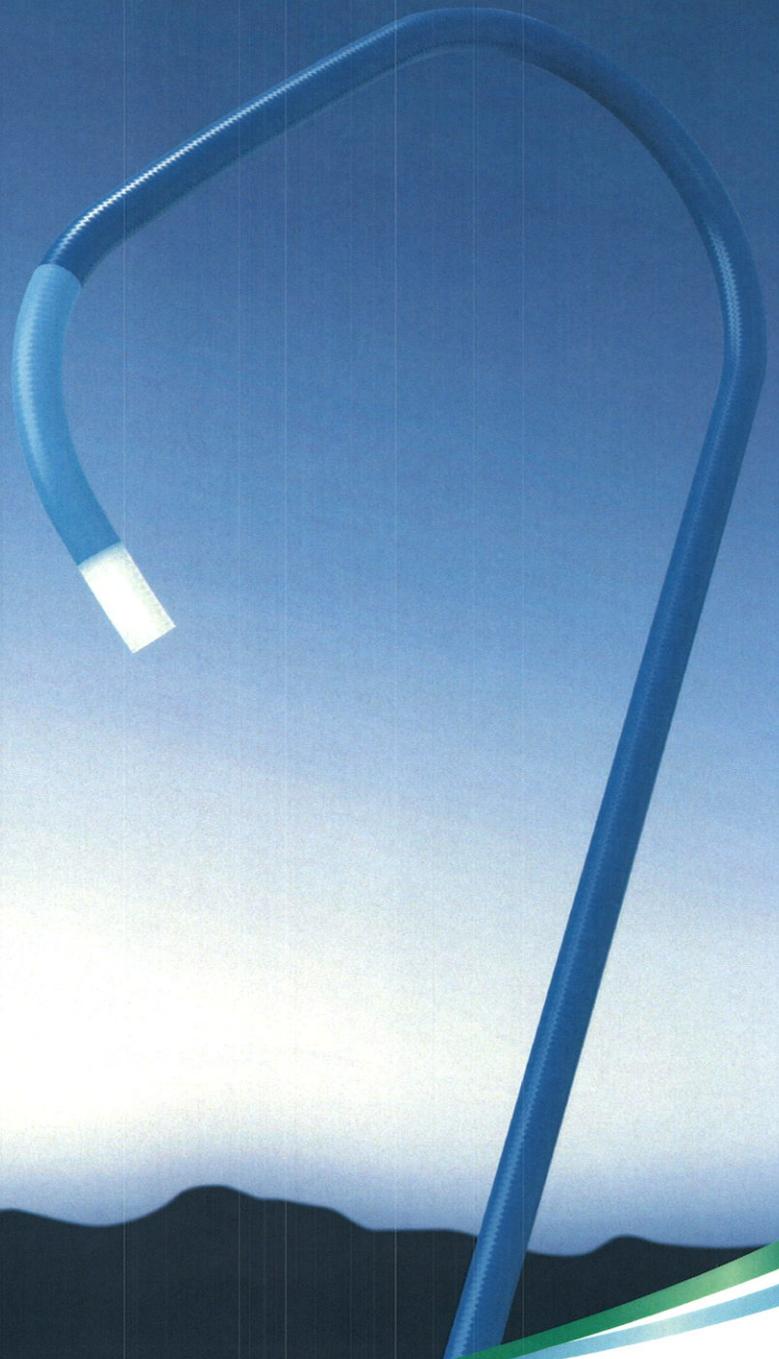

 **TERUMO®**

35,38 pirkimo dalys

Nuoroda internete:

<http://www.terumo-europe.com/en-emea/interventional-cardiology/coronary-intervention-products/ptca-guiding-catheter/climber%C2%AE-ptca-guiding-catheter>


New aspects in your daily PCI procedure

Hydrophilic coated new guiding catheter

- **Frictionless insertion** even in tortuous or calcified anatomy
- **Reduce risk of spasm** in transradial procedure
- **Provide tactile 1:1 torque control**
- Distal 7 cm and proximal 25 cm uncoated parts assure a **safe and comfortable procedure** while maintaining the benefit of **hydrophilic coating**



Hydrophilic coating

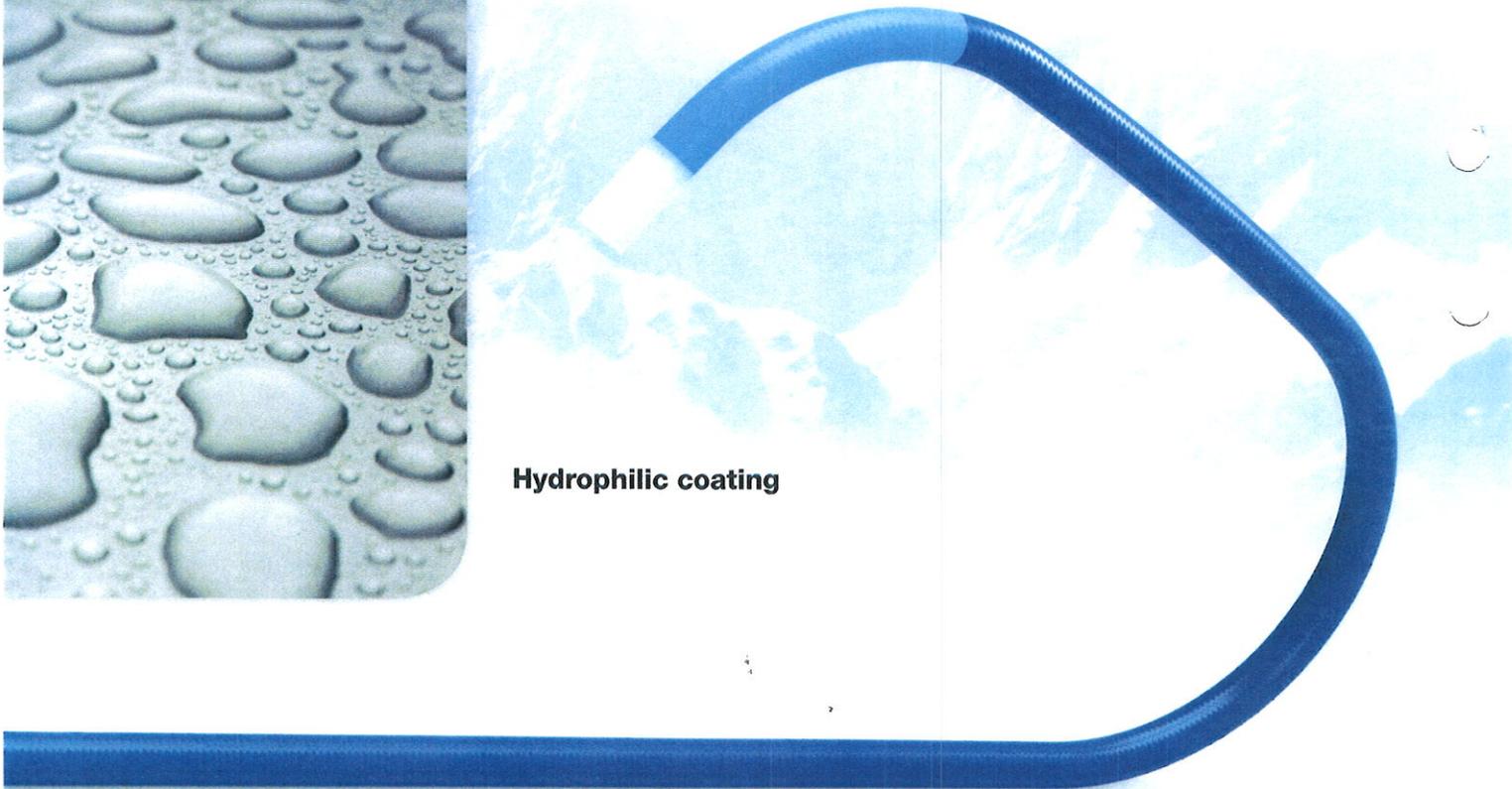

Ideal combination of largest inner lumen and passive backup support shaft




5Fr I.D.
0.058"/1.47 mm



6Fr I.D.
0.071"/1.80 mm



7Fr I.D.
0.082"/2.08 mm

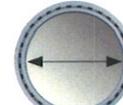

8Fr I.D.
0.091"/2.31 mm

- **Extra large 6 Fr inner lumen** allows to keep its system throughout procedure which contribute saving cost and contrast media
- **Flat braiding equipped inside shaft material** provides reliable backup support, pushability and excellent kink resistance
- **Polyamide outer shaft material** ensures excellent shape retention and minimal temperature softening of the catheter even in complex case treatment

Broad choices of curves maximizes your success

			5Fr / 1.76 mm	6Fr / 2.08 mm	7Fr / 2.37 mm	7Fr / 2.37 mm with Sidehole	8Fr / 2.68 mm	8Fr / 2.68 mm with Sidehole	
Standard curves	Amplatz Left	AL	0,75	GC-P5AL070N	GC-P6AL070N	GC-P7AL070N	GC-P7AL070A	GC-P8AL070N	GC-P8AL070A
		AL	1	GC-P5AL100N	GC-P6AL100N	GC-P7AL100N	GC-P7AL100A	GC-P8AL100N	GC-P8AL100A
		AL	1,5	GC-P5AL150N	GC-P6AL150N	GC-P7AL150N	GC-P7AL150A	GC-P8AL150N	GC-P8AL150A
		AL	2	GC-P5AL200N	GC-P6AL200N	GC-P7AL200N	GC-P7AL200A	GC-P8AL200N	GC-P8AL200A
		AL	3	GC-P5AL300N	GC-P6AL300N	GC-P7AL300N	GC-P7AL300A	GC-P8AL300N	GC-P8AL300A
	Amplatz Right	AR	1	GC-P5AR100N	GC-P6AR100N	GC-P7AR100N	GC-P7AR100A	GC-P8AR100N	GC-P8AR100A
		AR	2	GC-P5AR200N	GC-P6AR200N	GC-P7AR200N	GC-P7AR200A	GC-P8AR200N	GC-P8AR200A
	Judkins Left	JL	3	GC-P5JL300N	GC-P6JL300N	GC-P7JL300N	GC-P7JL300A	GC-P8JL300N	GC-P8JL300A
		JL	3,5	GC-P5JL350N	GC-P6JL350N	GC-P7JL350N	GC-P7JL350A	GC-P8JL350N	GC-P8JL350A
		JL	4	GC-P5JL400N	GC-P6JL400N	GC-P7JL400N	GC-P7JL400A	GC-P8JL400N	GC-P8JL400A
		JL	4,5	GC-P5JL450N	GC-P6JL450N	GC-P7JL450N	GC-P7JL450A	GC-P8JL450N	GC-P8JL450A
		JL	5	GC-P5JL500N	GC-P6JL500N	GC-P7JL500N	GC-P7JL500A	GC-P8JL500N	GC-P8JL500A
	Judkins Right	JR	3	GC-P5JR300N	GC-P6JR300N	GC-P7JR300N	GC-P7JR300A	GC-P8JR300N	GC-P8JR300A
		JR	3,5	GC-P5JR350N	GC-P6JR350N	GC-P7JR350N	GC-P7JR350A	GC-P8JR350N	GC-P8JR350A
		JR	4	GC-P5JR400N	GC-P6JR400N	GC-P7JR400N	GC-P7JR400A	GC-P8JR400N	GC-P8JR400A
		JR	4,5	GC-P5JR450N	GC-P6JR450N	GC-P7JR450N	GC-P7JR450A	GC-P8JR450N	GC-P8JR450A
		JR	5	GC-P5JR500N	GC-P6JR500N	GC-P7JR500N	GC-P7JR500A	GC-P8JR500N	GC-P8JR500A
		JR	6	GC-P5JR600N	GC-P6JR600N	GC-P7JR600N	GC-P7JR600A	GC-P8JR600N	GC-P8JR600A

CLIMBER™

		5Fr / 1.76 mm	6Fr / 2.08 mm	7Fr / 2.37 mm	7Fr / 2.37 mm with Sidehole	8Fr / 2.68 mm	8Fr / 2.68 mm with Sidehole	
Extra Backup curves	ExtraBackUp for left coronary	EBU 3	GC-P5EB300N	GC-P6EB300N	GC-P7EB300N	GC-P7EB300A	GC-P8EB300N	GC-P8EB300A
		EBU 3,25	GC-P5EB320N	GC-P6EB320N	GC-P7EB320N	GC-P7EB320A	GC-P8EB320N	GC-P8EB320A
		EBU 3,5	GC-P5EB350N	GC-P6EB350N	GC-P7EB350N	GC-P7EB350A	GC-P8EB350N	GC-P8EB350A
		EBU 3,75	GC-P5EB370N	GC-P6EB370N	GC-P7EB370N	GC-P7EB370A	GC-P8EB370N	GC-P8EB370A
		EBU 4	GC-P5EB400N	GC-P6EB400N	GC-P7EB400N	GC-P7EB400A	GC-P8EB400N	GC-P8EB400A
		EBU 4,25	GC-P5EB420N	GC-P6EB420N	GC-P7EB420N	GC-P7EB420A	GC-P8EB420N	GC-P8EB420A
		EBU 4,5	GC-P5EB450N	GC-P6EB450N	GC-P7EB450N	GC-P7EB450A	GC-P8EB450N	GC-P8EB450A
	Extra BackUp for Right coronary	EBU 4,75	GC-P5EB470N	GC-P6EB470N	GC-P7EB470N	GC-P7EB470A	GC-P8EB470N	GC-P8EB470A
		EBU 5	GC-P5EB500N	GC-P6EB500N	GC-P7EB500N	GC-P7EB500A	GC-P8EB500N	GC-P8EB500A
		XBRCA 3	GC-P5XR300N	GC-P6XR300N	GC-P7XR300N	GC-P7XR300A	GC-P8XR300N	
		XBRCA 3,5	GC-P5XR350N	GC-P6XR350N	GC-P7XR350N	GC-P7XR350A	GC-P8XR350N	
		XBRCA 3,75	GC-P5XR370N	GC-P6XR370N	GC-P7XR370N	GC-P7XR370A	GC-P8XR370N	
		XBRCA 4	GC-P5XR400N	GC-P6XR400N	GC-P7XR400N	GC-P7XR400A	GC-P8XR400N	
		XBRCA 4,25	GC-P5XR420N	GC-P6XR420N	GC-P7XR420N	GC-P7XR420A	GC-P8XR420N	
	Other curves	Multipurpose	XBRCA 4,5	GC-P5XR450N	GC-P6XR450N	GC-P7XR450N	GC-P7XR450A	GC-P8XR450N
MPA -			GC-P5MP000N	GC-P6MP000N	GC-P7MP000N	GC-P7MP000A		
HockeyStick		HS -	GC-P5HS000N	GC-P6HS000N	GC-P7HS000N	GC-P7HS000A		
		IMA	IM -	GC-P5IM000N	GC-P6IM000N	GC-P7IM000N	GC-P7IM000A	
Bypass		LCB -	GC-P5BP010N	GC-P6BP010N	GC-P7BP010N	GC-P7BP010A		
		RCB -	GC-P5BP020N	GC-P6BP020N	GC-P7BP020N	GC-P7BP020A		
Curves for Radial access		TIG	TIG 3	GC-P5TR300N	GC-P6TR300N			
			TIG 3,5	GC-P5TR350N	GC-P6TR350N			
		Tiger	TIG 4	GC-P5TR400N	GC-P6TR400N			
			TIG 4,5	GC-P5TR450N	GC-P6TR450N			
	TIG 5		GC-P5TR500N	GC-P6TR500N				

Specifications

Outer coating: Hydrophilic

Catheter length: 100 cm

Non coated at Distal 7 cm and Proximal 25 cm



Terumo Corporation
44-1, 2-chome, Hatagaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 38 12 11
Fax: +32 16 40 02 49

Terumo Europe Interventional Systems Liaison Office
Immeuble Krebs
8, rue Jean-Pierre Timbaud
CS 10720
78067 Saint Quentin
en Yvelines Cedex,
France
Tel.: +33 139 30 95 75
Fax: +33 139 30 95 70

European Offices

Terumo Europe N.V. Benelux Sales Division
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex, France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Italia S.r.l.
Via Simone Martini, 143/145
00142 Roma RM, Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo UK Ltd
Tamesis, The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH Zweigniederlassung Österreich
Liebermannstrasse F10 301
2345 Brunn am Gebirge, Österreich
Tel.: +43 2236 378 020
Fax: +43 2236 378 020 100

Terumo Deutschland GmbH Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid), Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Denmark Filial of Terumo Sweden AB
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda, Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

Terumo Russia LLC
BC Northern Tower
13 Floor, 10 Testovskaya Street
Moscow 123317, Russia
Tel.: + 7 495 988 4740
Fax: + 7 495 988 4739

Terumo Corporation Dubai Branch (Middle East and Africa)
Al Masraf Tower, 22nd Floor
P. O. Box 20291 Dubai
U.A.E.
Tel.: + 971 4 221 2220
Fax: + 971 4 221 3330

Addresses subject to change, please consult our Website.

© Registered Trademark
Published by Terumo Europe N.V.