

# Safer Power Morcellation

now a reality...

59 pirkimo dalis

## Laparoscopic Surgery

<http://www.safemorcellation.com/>

It is well known and accepted that the ability for the medical community to perform laparoscopic surgery offers tremendous benefits to the patient. In cases where a hysterectomy or myomectomy is indicated, the advantages to the patient include no large abdominal laparotomy incision, faster healing and recovery from the surgery, less postoperative pain, lower risk of surgical site infection and therefore a rapid return to normal activities.

Tissue morcellation, performed in combination with an isolation containment bag, is one of the most important advancement in laparoscopic surgery.

## Introducing Versator + MorSafe = A Winning Combination

The Versator tissue morcellator is unique in that it is a reusable product., offering you the ability to morcellate using a light weight, powerful disposable Handpiece while having the ability to reuse the rotor up to 10 times.

The MorSafe tissue morcellation bag has been designed and is specifically indicated for morcellation and as a result offers significant features and benefits that cannot be duplicated by other non-indicated bags being considered. Using MorSafe, you can now successfully prevent the dissemination of morcellated tissue around the peritoneal cavity.

The Versator tissue morcellator, when used in combination with the MorSafe tissue morcellation bag offers the latest and the best standard of patient care. And finally, for the first time, you can feel confident that you can morcellate large tissue quickly, efficiently and most important, safely.

## Morcellate Safely...



Versator and Morsafe are registered / trademarks of Veol Medical Technologies and the products are covered by various national and international approved and pending patents.

## TISSUE MORCELLATION BAG

Maišas vienkartinio naudojimo, skirtas saugiam morcellavimui

MorSafe™ is an innovative single-use disposable device that prevents dissemination of tissue in the peritoneal cavity, allowing for safe morcellation.

It is intended to be used during gynecological procedures such as laparoscopic myomectomy or laparoscopic hysterectomy.

Used in conjunction with the Versator® Tissue Morcellator, the system contains and removes the shredded benign tissue mass. It also prevents the CO<sub>2</sub> from entering the patient, potentially reducing CO<sub>2</sub> side effects.


The device has unique features to allow for quick deployment, insufflation, morcellation and spill-proof withdrawal of the bag.

MorSafe™, with its unique patented two port design, offers the surgeon superior visibility during the surgery compared to a single port approach. Designed to fit and take the shape of the abdomen, it has been constructed utilizing a special tear resistant material to prevent leakage.

It also contains a special ring in the bag opening to allow the surgeon ultimate control of the bag opening and easy access to the interior of the bag during surgery.

Su dviem praėjimais: vienas - laparoskopui, kitas - morceliatoriui.

Permatomas, atsparus įplėšimams ir įpjovimui

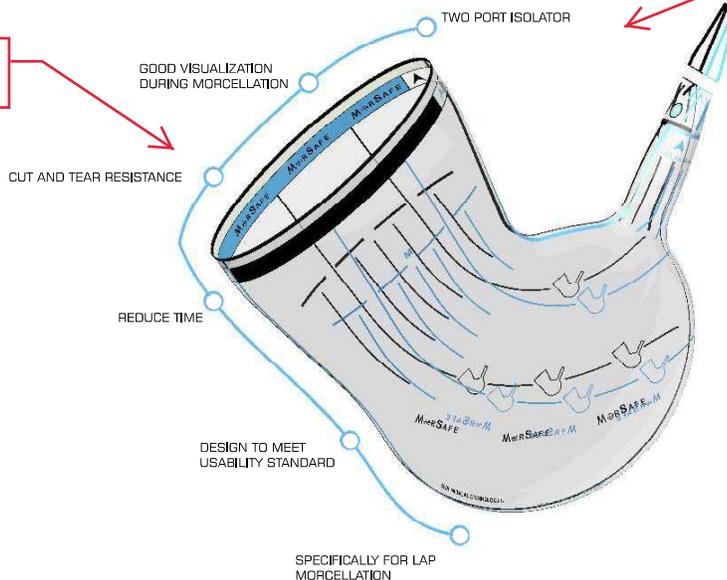

World's 1<sup>st</sup> product specifically designed for Safe Morcellation !

# SAFE TISSUE MORCELLATION SYSTEM

The safe tissue morcellation system comprises of VDU, Handpiece, Rotor and Tissue morcellation Bag. This unique package also includes accessories such as Obturator for Handpiece insertion and **Introducer for insertion of the Tissue Morcellation Bag.**

Komplekte maišo pravedējas

## PRODUCT SPECIFICATION

HANDPIECE	
Cutting Speed	Upto 2000 RPM
Cutting Diameter	15 mm
Cutting Modes	CUT, CUT PLUS & OFF
Weight	180 Grams
Core Guard poitions	Top, Right & Left
Working Length	160 mm
VERSATOR DRIVE UNIT - VDU	
Speed Range	250-2000 RPM
Power Rating	100-240VAC,350W,50-60Hz
Dimensions	345 x 215 x 110 mm
Weight	4.5 Kgs / 10 Lbs

TISSUE MORCELLATION BAG	
Number of Ports	Two
Sizes	Available in 3 sizes Small: 1600 cc Medium: 2100 cc Large: 2600 cc
ROTOR	
Rotor Length	2000 mm
Weight	460 Grams
No, of Reuse	10 cycles - autoclave
OBTURATOR	
Length	340 mm
Weight	26 Grams

Maišo tūris

## COMPETITIVE OVERVIEW

PARAMETERS	VERSATOR	GYNECARE MORCELLEX	KARLSTORZ ROTOCUT G1	LINA XCISE
Reusability	Partially Reusable	Single Use	Reusable	Single Use
Contained Morcellation	Morcellation Bag available in 3 sizes	Not available	Not available	Not available
Speed Range	250 - 2000 RPM	125 - 1000 RPM	500 - 1200 RPM	1000 RPM

For more details or to place an order, contact...



A-747, TTC Industrial Area, MIDC Pawane, Navi Mumbai, India - 400705  
Phone: +91 22 27634771, +91 9892332477 | www.veolmedtech.com  
info@safemorcellation.com

Copyright© Veol Medical Technologies 2017

Distributed by: