

# Safeflow

Premium Safety



The valve that makes infusion safety simple

# When infusion safety really matters

Safe, simple and cost-effective needle-free infusion system

## Safety first


Safeflow needle-free infusion system provides safe options for the preparation and administration of IV therapies. Safeflow's valve promotes clinical safety first and optimizes clinical infusion flow rates.

## Versatile and innovative system

Safeflow needle-free infusion system provides needle-free products for all your IV therapy needs hospital wide, including IV administration sets and extension lines.

Safeflow products feature easy access for Luer-Slip and Luer-Lock connections providing clinicians with needle-free access for faster, easier drug delivery, blood sampling and infusions.

Perfect for your clinical success!

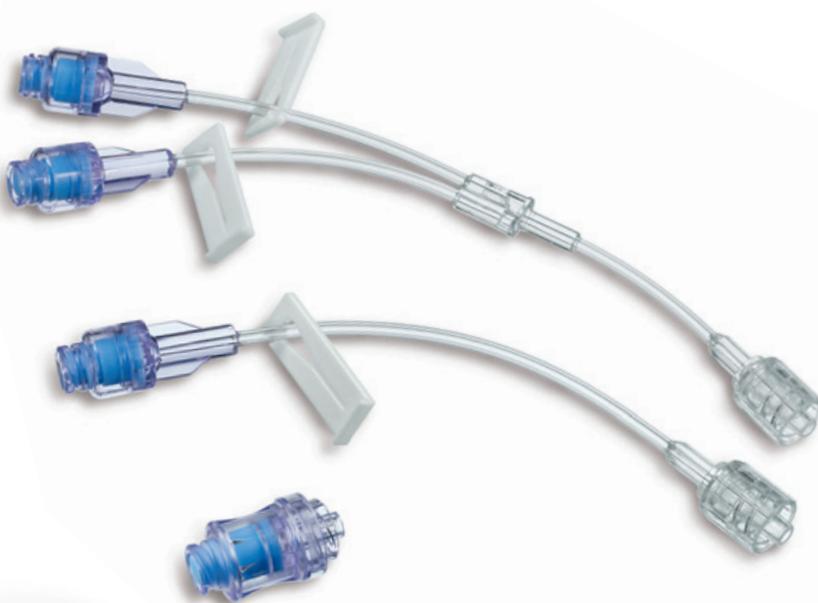

## Extension lines

Safeflow is always a closed system. The smallbore extension line allows manipulation away from the injection side, which ensures treatment of patient by infusion, transfusion, aspiration of blood

samples and intermittent injection. Tubings with smaller internal diameter are minimizing the contained volume and can be used for pediatric/neonate patients.




### PVC/DEHP-free Neutrapur® tubing

Particularly patient-friendly. No plasticizer leachable. Ecological disposal.



### Safeflow valve

Offers convenient and safe needle-free access for injection, aspiration or parallel infusions as a closed system.

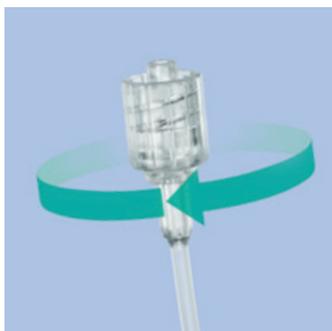

### Rotating Luer-Lock connector

Permits safer and more convenient connection without twisting the tube.




### Back check valve

Prevents backflow and therefore it offers higher safety during parallel pressure and gravity infusions.

## User benefits

- Swabbable
- Compatible with Luer-Slip and Luer-Lock
- Allows multiple access for direct aspiration or fluid injection
- Suitable for all IV Systems

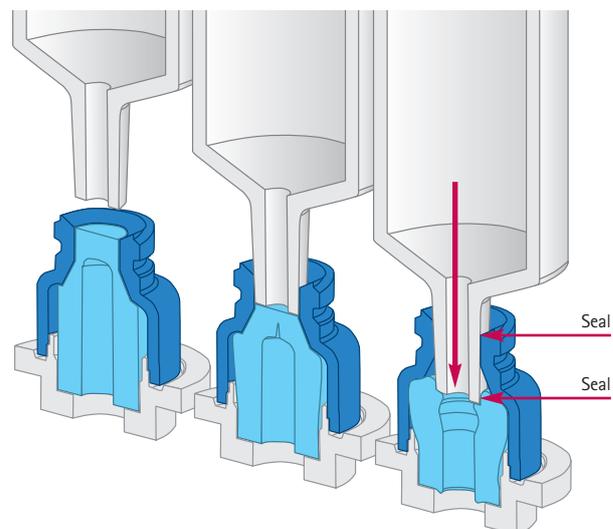

## Optimal flow rates

The Safeflow Luer activated valve features an innovative design that provides a straight fluid pathway resulting in uninterrupted fluid flow providing flow rates of 360 ml/min.

# Making infusion safety simple

## IV Accessories

Product Type	Description		Latex free	DEHP free	PVC free		Sales unit/pcs.	Code No. (REF)
	Safeflow valve		•	•	•		50	409100H
	<b>Extension Sets</b>	<b>Tubing inner Ø (mm)</b>	<b>Latex free</b>	<b>DEHP free tubing</b>	<b>PVC free tubing</b>	<b>Tubing length (cm)</b>	<b>Sales unit/pcs.</b>	<b>Code No. (REF)</b>
	Smallbore Extension Set with Safeflow valve	1.0	•	•	•	10	50	4097154
	Smallbore Y-Extension Set with two Safeflow valves	1.0	•	•	•	12	50	4097145
	Smallbore Y-Extension Set with two Safeflow valves and back check valve	1.0	•	•	•	12	50	4097148
	<b>IV Administration Sets</b>							
	Intrafix® SafeSet with Safeflow	3.0	•			210	100	4063004
	Intrafix® Primline Comfort with Safeflow	3.0	•			180	100	4062158
	Infusomat® Space Line with Safeflow	3.0	•	•	•	230	100	8700110SP
	Infusomat® Space Line, Dosifix®	3.0	•	•	•	270	40	8250245SP



Luer-Lock



Luer-Slip

Product Data of Safeflow	
Flow rate	360 ml/min
Latex-free	✓
Luer-Slip and Luer-Lock	✓
Complete system	✓
Priming volume	0.09 ml
Number of accesses	up to 200 accesses
Time in use	24 hours or as per hospital protocol

**B | BRAUN**  
SHARING EXPERTISE

B. Braun Melsungen AG

34209 Melsungen  
Germany  
Tel +49(0) 56 61 71-0  
www.bbraun.com