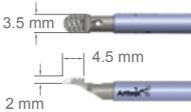
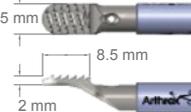


Shaver Blades and Burrs

Color bars reflect family color shown on the adapter

	Description	4 (Drill depth in mm)	6 (Drill depth in mm)
	<p>PowerPick™</p> <p>Initially designed to make quick work of microfracture procedures by using the shaver handpiece to power the Ø 1.5 mm trocar tip. The 30° and 45° angled shafts facilitate use in a variety of other applications, including ACL femoral tunnel location.</p>	<p>AR-8150PP-30 (30°)</p> <p>AR-8150PP-45 (45°)</p>	<p>AR-8150PX-45 (45°)</p>

	Description	3.5 (Width in mm)	4 (Width in mm)	5.5 (Width in mm)
	<p>PowerRasp™</p> <p>Perfect for arthroscopic bone debridement by reciprocating motion, creating a smooth and flat surface, eg, subacromial decompressions, notchplasties.</p>	<p>AR-8350PR</p> 	<p>AR-8400PR</p> 	<p>AR-8550PR</p> 

Retractable Hood Burrs


Color bars reflect family color shown on the adapter

	Description	4 (Ø in mm)	5 (Ø in mm)	5.5 (Ø in mm)
	<p>FlushCut™ Oval Burrs</p> <p>Allows the benefit of using the shaft as a cutting depth guide during procedures where a flat surface is desired, such as in an acromioplasty. By using the shaft as a depth stop and keeping the position parallel to the surface during use, the ideal acromioplasty can be completed quickly.</p>		<p>AR-8500RFOE</p>	<p>AR-8550RFOE</p>
	<p>Hip-Length Round Burrs</p> <p>Acetabular preparations can be made with the hood in place to protect the labrum and surrounding tissues when performing a hip labral repair procedure, and when necessary, the hood can be retracted to allow full burr head exposure to address other hip pathology, such as a CAM lesion.</p>	<p>AR-6400RRBE</p>	<p>AR-6500RRBE</p>	<p>AR-6550RRBE</p>


Shaver Blades and Burrs for Small Joint Handpiece AR-8330SJ

Color bars reflect family color shown on the adapter

	Description	1 (Ø in mm)	2 (Ø in mm)	3 (Ø in mm)	3.5 (Ø in mm)
	PowerRasp™ Resector, Small Hub Ideally suited for soft tissue, osteochondral and osteophyte resection, or bony site preparation.				AR-9350PR
	PowerPick™ Drill, Small Hub Excellent for micro-drilling subchondral bone defects.	AR-9100PP-00 (0°) AR-9100PP-45 (45°)			
	Torpedo™ Shaver Blade Tapered tip and scissor-like cutting action facilitate rapid, aggressive soft tissue resection.				AR-9350TD
	Dissector, Small Hub Sharp, teardrop-shaped outer window. The toothed inner cutting window aggressively attacks tissue, grabbing it and drawing it into the cutting area.			AR-9300DS	AR-9350DS
	Sabre, Small Hub Sharp, teardrop-shaped outer window. The smooth inner cutting window leaves crisp, clean edges for a more anatomic appearance following resection.		AR-9200SR	AR-9300SR	AR-9350SR
	Oval Burr, Small Hub Designed for rapid, aggressive bone resection in notchplasty, distal clavicle, and subacromial decompression procedures.			AR-9300BT	AR-9350BT
	Round Burr, Small Hub Ideally suited for soft tissue, osteochondral and osteophyte resection, or bony site preparation.			AR-9300BT	AR-9350BT

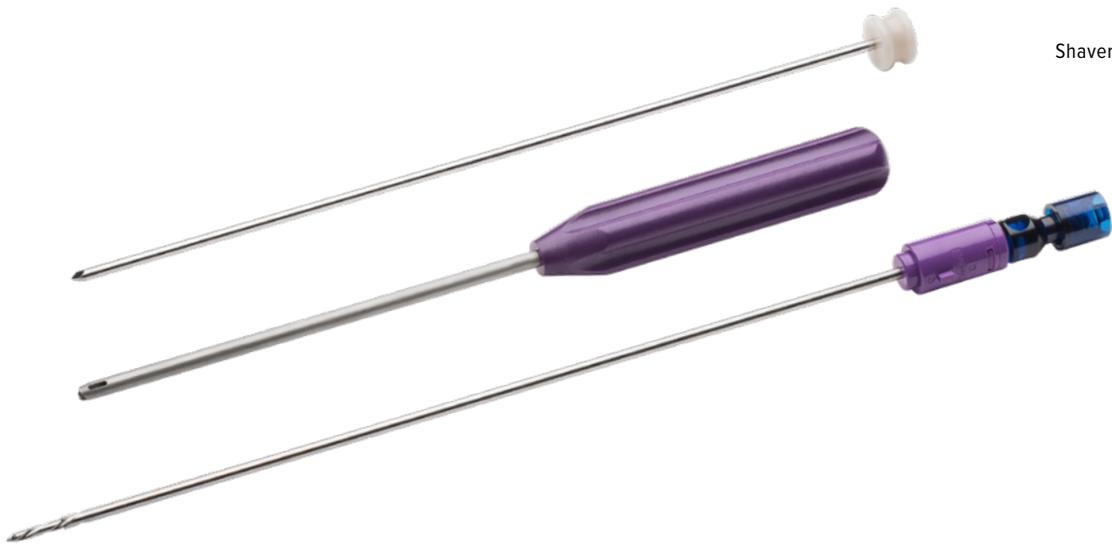

AR-1934DSSK
ShaverDrill disposables kit for
3 mm SutureTak anchor

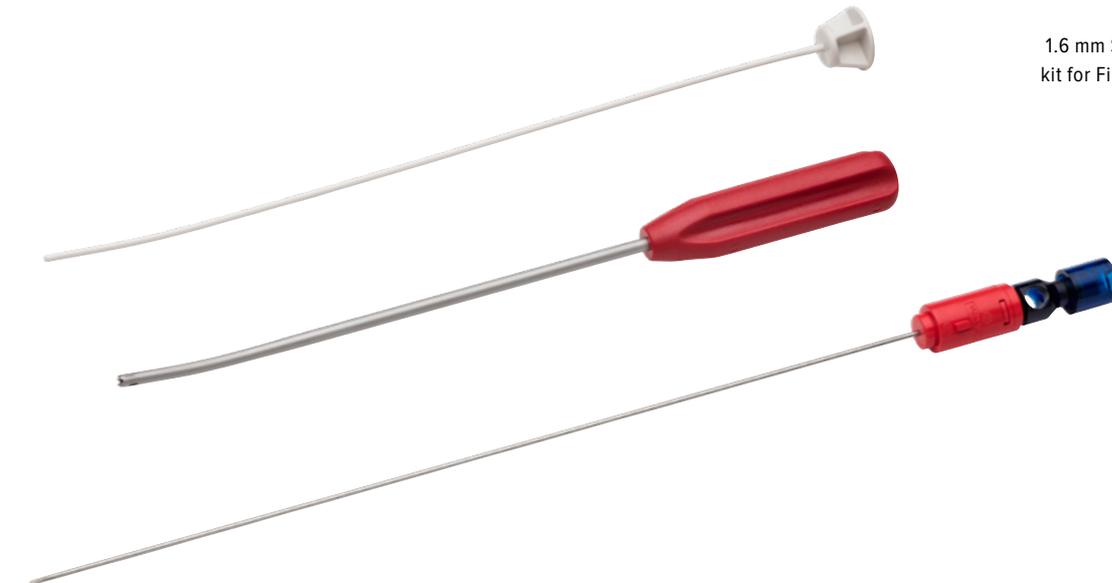

AR-3600NDSR-1
1.6 mm ShaverDrill disposables
kit for FiberTak soft anchor with
curved spear